

# WORLD INTELLECTUAL PROPERTY ORGANIZATION International Bureau



#### INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 6: C07H 21/04, C07K 14/00, 16/00, C12N

A1

(11) International Publication Number:

WO 98/06737

(43) International Publication Date:

19 February 1998 (19.02.98)

(21) International Application Number:

1/20, 15/00, 15/09, 15/63

PCT/US97/14593

(22) International Filing Date:

15 August 1997 (15.08.97)

(30) Priority Data:

60/024,050 08/706,408

16 August 1996 (16.08.96) 30 August 1996 (30.08.96)

US US

(60) Parent Applications or Grants

(63) Related by Continuation

US Filed on US Filed on

60/024,050 (CIP) 16 August 1996 (16.08.96) 08/706,408 (CIP)

30 August 1996 (30.08.96)

(71) Applicants (for all designated States except US): THE REGENTS OF THE UNIVERSITY OF CALIFORNIA [US/US]; 22nd floor, 300 Lakeside Drive, Oakland, CA 94612-3550 (US). STATE OF OREGON, acting by and through THE STATE BOARD OF HIGHER EDUCATION on behalf of THE UNIVERSITY OF OREGON [US/US]; Eugene, OR 97403 (US). AURORA BIOSCIENCES [US/US]; 11149 North Torrey Pines Road, La Jolla, CA 92037 (US).

(72) Inventors; and

(75) Inventors/Applicants (for US only): TSIEN, Roger, Y. [US/US]; 8535 Nottingham Place, La Jolla, CA 92037 (US). CUBITT, Andrew, B. [GB/GB]; 11083 Caminito Encanto, San Diego, CA 92131 (US). HEIM, Roger [CH/CH]; 530 Stratford Court #309B, Del Mar, CA 92014 (US). ORMO, Mats, F. [SE/SE]; Juristvägen 18, S-141 73 Huddinge (SE). REMINGTON, James, S. [US/US]; 2100 River Loop 1, Eugene, OR 97404 (US).

(74) Agent: HAILE, Lisa, A.; Fish & Richardson P.C., Suite 1400. 4225 Executive Square, La Jolla, CA 92037 (US).

(81) Designated States: AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, CA, CH, CN, CU, CZ, DE, DK, EE, ES, FI, GB, GE, GH, HU, IL, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, UA, UG, US, UZ, VN, ARIPO patent (GH, KE, LS, MW, SD, SZ, UG, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, ML, MR, NE, SN, TD, TG).

#### Published

With international search report.

Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.

(54) Title: LONG WAVELENGTH ENGINEERED FLUORESCENT PROTEINS

(57) Abstract

Engineered fluorescent proteins, nucleic acids encoding them and methods of use.

# FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AL	Albania	ES	Spain	LS	Lesotho	SI	Slovenia
AM	Armenia	FI	Finland	LT	Lithuania	SK	Slovakia
AT	Austria	FR	France	LU	Luxembourg	SN	Senegal
AU	Australia	GA	Gabon	LV	Latvia	SZ	Swaziland
AZ	Azerbaijan	GB	United Kingdom	MC	Monaco	TD	Chad
BA	Bosnia and Herzegovina	CE	Georgia	MD	Republic of Moldova	TG	Togo
BB	Barbados	GH	Ghana	MG	Madagascar	TJ	Tajikistan
BE	Belgium	GN	Guinea	MK	The former Yugoslav	TM	Turkmenistan
BF	Burkina Faso	GR	Greece		Republic of Macedonia	TR	Turkey
BG	Bulgaria	HU	Hungary	ML	Mali	TT	Trinidad and Tobago
BJ	Benin	1E	Ireland	MN	Mongolia	ŪA	Ukraine
BR	Brazil	IL	Israel	MR	Mauritania	UG	Uganda
BY	Belarus	IS	Iceland	MW	Melawi	US	United States of America
CA	Canada	ſΤ	Italy	MX	Mexico	UZ	Uzbekistan
CF	Central African Republic	JP	Japan	NB	Niger	VN	Vict Nam
CG	Congo	KE	Kenya	NL	Netherlands	YU	Yugoslavia
CH	Switzerland	KG	Kyrgyzstan	NO	Norway	zw	Zimbabwc
CI	Côte d'Ivoire	KP	Democratic People's	NZ	New Zealand		
CM	Cameroon		Republic of Korea	PL	Poland		
CN	China	KR	Republic of Korea	PT	Portugal		
CU	Cuba	KZ	Kazakstan	RO	Romania		
CZ	Czech Republic	LC	Saint Lucia	RU	Russian Federation		
DE	Germany	u	Liechtenstein	SD	Sudan		
DK	Denmark	LK	Sri Lanka	SE	Sweden		
EE	Estonia	LR	Liberia	SC	Singapore		

# LONG WAVELENGTH ENGINEERED FLUORESCENT PROTEINS

5

10

15

20

25

30

## **BACKGROUND OF THE INVENTION**

This application claims the benefit of the earlier filing date of a United States
60/024,050
provisional patent application serial number filed on August 16, 1996 entitled
"Long Wavelength Mutant Fluorescent Proteins" and patent application serial number
08/706,408filed on August 30, 1996 entitled "Long Wavelength Engineered Fluorescent
Proteins." both of which are herein incorporated by reference.

This invention was made in part with Government support under grant no. MCB 9418479 awarded by the National Science Foundation. The Government may have rights in this invention.

Fluorescent molecules are attractive as reporter molecules in many assay systems because of their high sensitivity and ease of quantification. Recently, fluorescent proteins have been the focus of much attention because they can be produced in vivo by biological systems, and can be used to trace intracellular events without the need to be introduced into the cell through microinjection or permeabilization. The green fluorescent protein of Aequorea victoria is particularly interesting as a fluorescent protein. A cDNA for the protein has been cloned. (D.C. Prasher et al., "Primary structure of the Aequorea victoria green-fluorescent protein," Gene (1992) 111:229-33.) Not only can the primary amino acid sequence of the protein be expressed from the cDNA, but the expressed protein can fluoresce. This indicates that the protein can undergo the cyclization and oxidation Aequorea green fluorescent protein believed to be necessary for fluorescence. ("GFP") is a stable, proteolysis-resistant single chain of 238 residues and has two absorption maxima at around 395 and 475 nm. The relative amplitudes of these two peaks is sensitive to environmental factors (W. W. Ward. Bioluminescence and Chemiluminescence (M. A. DeLuca and W. D. McElroy, eds) Academic Press pp. 235-242 (1981); W. W. Ward & S. H. Bokman Biochemistry 21:4535-4540 (1982); W. W. Ward et al. Photochem. Photobiol. 35:803-808 (1982)) and illumination history (A. B. Cubitt et al. Trends Biochem. Sci. 20:448-455 (1995)), presumably reflecting two or more ground states. Excitation at the primary absorption peak of 395 nm yields an emission maximum at 508 nm with a quantum yield of 0.72-0.85 (O. Shimomura and F.H. Johnson J. Cell. Comp. Physiol. 59:223 (1962);

J. G. Morin and J. W. Hastings, J. Cell. Physiol. 77:313 (1971); H. Morise et al. Biochemistry 13:2656 (1974); W. W. Ward Photochem. Photobiol. Reviews (Smith, K. C. ed.) 4:1 (1979); A. B. Cubitt et al. Trends Biochem. Sci. 20:448-455 (1995); D. C. Prasher Trends Genet. 11:320-323 (1995); M. Chalfie Photochem. Photobiol. 62:651-656 (1995); W. W. Ward. Bioluminescence and Chemiluminescence (M. A. DeLuca and W. D. McElroy, eds) Academic Press pp. 235-242 (1981); W. W. Ward & S. H. Bokman Biochemistry 21:4535-4540 (1982); W. W. Ward et al. Photochem. Photobiol. 35:803-808 (1982)). The fluorophore results from the autocatalytic cyclization of the polypeptide backbone between residues Ser<sup>65</sup> and Gly<sup>67</sup> and oxidation of the □-β bond of Tyr<sup>66</sup> (A. B.

5

10

15

20

30

Cubitt et al. Trends Biochem. Sci. 20:448-455 (1995); C. W. Cody et al. Biochemistry 32:1212-1218 (1993); R. Heim et al. Proc. Natl. Acad. Sci. USA 91:12501-12504 (1994)). Mutation of Ser<sup>65</sup> to Thr (S65T) simplifies the excitation spectrum to a single peak at 488 nm of enhanced amplitude (R. Heim et al. Nature 373:664-665 (1995)), which no longer gives signs of conformational isomers (A. B. Cubitt et al. Trends Biochem. Sci. 20:448-455 (1995)).

Fluorescent proteins have been used as markers of gene expression, tracers of cell lineage and as fusion tags to monitor protein localization within living cells. (M. Chalfie et al., "Green fluorescent protein as a marker for gene expression," Science 263:802-805; A.B. Cubitt et al., "Understanding, improving and using green fluorescent proteins," TIBS 20, November 1995, pp. 448-455. U.S. patent 5,491,084, M. Chalfie and D. Prasher. Furthermore, engineered versions of Aequorea green fluorescent protein have been identified that exhibit altered fluorescence characteristics, including altered excitation and emission maxima, as well as excitation and emission spectra of different shapes. (R. Heim et al., "Wavelength mutations and posttranslational autoxidation of green fluorescent protein," Proc. Natl. Acad. Sci. USA, (1994) 91:12501-04; R. Heim et al., "Improved green 25 fluorescence," Nature (1995) 373:663-665.) These properties add variety and utility to the arsenal of biologically based fluorescent indicators.

There is a need for engineered fluorescent proteins with varied fluorescent properties.

#### BRIEF DESCRIPTION OF THE DRAWINGS

Figs. 1A-1B. (A) Schematic drawing of the backbone of GFP produced by

Molscript (J.P. Kraulis, J. Appl. Cryst., 24:946 (1991)). The chromophore is shown as a ball and stick model. (B) Schematic drawing of the overall fold of GFP. Approximate residue numbers mark the beginning and ending of the secondary structure elements.

5

10

15

20

25

30

Figs. 2A-2C. (A) Stereo drawing of the chromophore and residues in the immediate vicinity. Carbon atoms are drawn as open circles, oxygen is filled and nitrogen is shaded. Solvent molecules are shown as isolated filled circles. (B) Portion of the final  $2F_o$ - $F_c$  electron density map contoured at  $1.0 \,\Box$ , showing the electron density surrounding the chromophore. (C) Schematic diagram showing the first and second spheres of coordination of the chromophore. Hydrogen bonds are shown as dashed lines and have the indicated lengths in Å. Inset: proposed structure of the carbinolamine intermediate that is presumably formed during generation of the chromophore.

Fig. 3 depicts the nucleotide sequence (SEQ ID NO:1) and deduced amino acid sequence (SEQ ID NO:2) of an Aequorea green fluorescent protein.

Fig. 4 depicts the nucleotide sequence (SEQ ID NO:3) and deduced amino acid sequence (SEQ ID NO:4) of the engineered *Aequorea*-related fluorescent protein S65G/S72A/T203Y utilizing preferred mammalian codons and optimal Kozak sequence.

Figs. 5-1 to 5-28 present the coordinates for the crystal structure of Aequorea-related green fluorescent protein S65T.

Fig. 6 shows the fluorescence excitation and emission spectra for engineered fluorescent proteins 20A and 10C (Table F). The vertical line at 528 nm compares the emission maxima of 10C, to the left of the line, and 20A, to the right of the line.

#### SUMMARY OF THE INVENTION

This invention provides functional engineered fluorescent proteins with varied fluorescence characteristics that can be easily distinguished from currently existing green and blue fluorescent proteins. Such engineered fluorescent proteins enable the simultaneous measurement of two or more processes within cells and can be used as fluorescence energy donors or acceptors when used to monitor protein-protein interactions through FRET. Longer wavelength engineered fluorescent proteins are particularly useful because photodynamic toxicity and auto-fluorescence of cells are significantly reduced at longer wavelengths. In particular, the introduction of the substitution T203X, wherein X is an aromatic amino acid, results in an increase in the excitation and emission wavelength

maxima of Aequorea-related fluorescent proteins.

In one aspect, this invention provides a nucleic acid molecule comprising a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of *Aequorea* green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid substitution located no more than about 0.5 nm from the chromophore of the engineered fluorescent protein, wherein the substitution alters the electronic environment of the chromophore, whereby the functional engineered fluorescent protein has a different fluorescent property than *Aequorea* green fluorescent protein.

10

15

20

25

5

In one aspect this invention provides a nucleic acid molecule comprising a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least a substitution at T203 and, in particular, T203X, wherein X is an aromatic amino acid selected from H, Y, W or F. said functional engineered fluorescent protein having a different fluorescent property than Aequorea green fluorescent protein. In one embodiment, the amino acid sequence further comprises a substitution at S65, wherein the substitution is selected from S65G, S65T, S65A, S65L, S65C, S65V and S65I. In another embodiment, the amino acid sequence differs by no more than the substitutions S65T/T203H; S65T/T203Y; \$72A/F64L/\$65G/T203Y; \$65G/V68L/Q69K/\$72A/T203Y; \$72A/\$65G/V68L/T203Y; S65G/S72A/T203Y; or S65G/S72A/T203W. In another embodiment, the amino acid sequence further comprises a substitution at Y66, wherein the substitution is selected from Y66H, Y66F, and Y66W. In another embodiment, the amino acid sequence further comprises a mutation from Table A. In another embodiment, the amino acid sequence further comprises a folding mutation. In another embodiment, the nucleotide sequence encoding the protein differs from the nucleotide sequence of SEQ ID NO:1 by the substitution of at least one codon by a preferred mammalian codon. In another embodiment, the nucleic acid molecule encodes a fusion protein wherein the fusion protein comprises a polypeptide of interest and the functional engineered fluorescent protein.

30

In another aspect, this invention provides a nucleic acid molecule comprising a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green

5

30

fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid substitution at L42, V61, T62, V68, Q69, Q94, N121, Y145, H148, V150, F165, I167, Q183, N185, L220, E222 (not E222G), or V224, said functional engineered fluorescent protein having a different fluorescent property than *Aequorea* green fluorescent protein. In one embodiment, amino acid substitution is:

L42X, wherein X is selected from C, F, H, W and Y, V61X, wherein X is selected from F, Y, H and C, T62X, wherein X is selected from A, V, F, S, D, N, Q, Y, H and C, V68X, wherein X is selected from F, Y and H, Q69X, wherein X is selected from K, R, E and G, 10 094X, wherein X is selected from D, E, H, K and N, N121X, wherein X is selected from F, H, W and Y, Y145X, wherein X is selected from W, C, F, L, E, H, K and Q, H148X, wherein X is selected from F, Y, N, K, Q and R, V150X, wherein X is selected from F, Y and H, 15 F165X, wherein X is selected from H, Q, W and Y, 1167X, wherein X is selected from F, Y and H, O183X, wherein X is selected from H, Y, E and K, N185X, wherein X is selected from D, E, H, K and Q, L220X, wherein X is selected from H, N, Q and T, 20 E222X, wherein X is selected from N and Q, or V224X, wherein X is selected from H, N, Q, T, F, W and Y.

In a further aspect, this invention provides an expression vector comprising expression control sequences operatively linked to any of the aforementioned nucleic acid molecules. In a further aspect, this invention provides a recombinant host cell comprising the aforementioned expression vector.

In another aspect, this invention provides a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid substitution located no more than about 0.5 nm from the chromophore of the engineered fluorescent protein, wherein the substitution alters the

electronic environment of the chromophore, whereby the functional engineered fluorescent protein has a different fluorescent property than Aequorea green fluorescent protein.

5

10

15

30

In another aspect, this invention provides a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least the amino acid substitution at T203, and in particular, T203X, wherein X is an aromatic amino acid selected from H, Y, W or F, said functional engineered fluorescent protein having a different fluorescent property than Aeguorea green fluorescent protein. In one embodiment, the amino acid sequence further comprises a substitution at S65, wherein the substitution is selected from S65G, S65T, S65A, S65L, S65C, S65V and S65I. In another embodiment, the amino acid sequence differs by no more than the substitutions S65T/T203H; S65T/T203Y; S72A/F64L/S65G/T203Y; S72A/S65G/V68L/T203Y; S65G/V68L/Q69K/S72A/T203Y; S65G/S72A/T203Y; or S65G/S72A/T203W. In another embodiment, the amino acid sequence further comprises a substitution at Y66, wherein the substitution is selected from Y66H, Y66F, and Y66W. In another embodiment, the amino acid sequence further comprises a folding mutation. In another embodiment, the engineered fluorescent protein is part of a fusion protein wherein the fusion protein comprises a polypeptide of interest and the functional engineered fluorescent protein.

In another aspect this invention provides a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid substitution at L42, V61, T62, V68, Q69, Q94, N121, Y145, H148, V150, F165, I167, Q183, N185, L220, E222, or V224, said functional engineered fluorescent protein having a different fluorescent property than Aequorea green fluorescent protein.

In another aspect, this invention provides a fluorescently labelled antibody comprising an antibody coupled to any of the aforementioned functional engineered fluorescent proteins. In one embodiment, the fluorescently labelled antibody is a fusion protein wherein the fusion protein comprises the antibody fused to the functional engineered fluorescent protein.

In another aspect, this invention provides a nucleic acid molecule comprising a nucleotide sequence encoding an antibody fused to a nucleotide sequence encoding a

functional engineered fluorescent protein of this invention.

5

10

15

20

25

30

In another aspect, this invention provides a fluorescently labelled nucleic acid probe comprising a nucleic acid probe coupled to a functional engineered fluorescent protein whose amino acid sequence of this invention. The fusion can be through a linker peptide.

PCT/US97/14593

In another aspect, this invention provides a method for determining whether a mixture contains a target comprising contacting the mixture with a fluorescently labelled probe comprising a probe and a functional engineered fluorescent protein of this invention; and determining whether the target has bound to the probe. In one embodiment, the target molecule is captured on a solid matrix.

In another aspect, this invention provides a method for engineering a functional engineered fluorescent protein having a fluorescent property different than Aequorea green fluorescent protein, comprising substituting an amino acid that is located no more than 0.5 nm from any atom in the chromophore of an Aequorea-related green fluorescent protein with another amino acid; whereby the substitution alters a fluorescent property of the protein. In one embodiment, the amino acid substitution alters the electronic environment of the chromophore.

In another aspect, this invention provides a method for engineering a functional engineered fluorescent protein having a different fluorescent property than Aequorea green fluorescent protein comprising substituting amino acids in a loop domain of an Aequorea-related green fluorescent protein with amino acids so as to create a consensus sequence for phosphorylation or for proteolysis.

In another aspect, this invention provides a method for producing fluorescence resonance energy transfer comprising providing a donor molecule comprising a functional engineered fluorescent protein this invention; providing an appropriate acceptor molecule for the fluorescent protein; and bringing the donor molecule and the acceptor molecule into sufficiently close contact to allow fluorescence resonance energy transfer.

In another aspect, this invention provides a method for producing fluorescence resonance energy transfer comprising providing an acceptor molecule comprising a functional engineered fluorescent protein of this invention; providing an appropriate donor molecule for the fluorescent protein; and bringing the donor molecule and the acceptor molecule into sufficiently close contact to allow fluorescence resonance energy

transfer. In one embodiment, the donor molecule is a engineered fluorescent protein whose amino acid sequence comprises the substitution T203I and the acceptor molecule is an engineered fluorescent protein whose amino acid sequence comprises the substitution T203X, wherein X is an aromatic amino acid selected from H, Y, W or F, said functional engineered fluorescent protein having a different fluorescent property than *Aequorea* green fluorescent protein.

In another aspect, this invention provides a crystal of a protein comprising a fluorescent protein with an amino acid sequence substantially identical to SEQ ID NO: 2, wherein said crystal diffracts with at least a 2.0 to 3.0 angstrom resolution.

In another embodiment, this invention provides computational method of designing a fluorescent protein comprising determining from a three dimensional model of a crystallized fluorescent protein comprising a fluorescent protein with a bound ligand, at least one interacting amino acid of the fluorescent protein that interacts with at least one first chemical moiety of the ligand, and selecting at least one chemical modification of the first chemical moiety to produce a second chemical moiety with a structure to either decrease or increase an interaction between the interacting amino acid and the second chemical moiety compared to the interaction between the interacting amino acid and the first chemical moiety.

In another embediment, this invention provides a computational method of modeling the three dimensional structure of a fluorescent protein comprising determining a three dimensional relationship between at least two atoms listed in the atomic coordinates of Figs. 5-1 to 5-28.

In another embodiment, this invention provides a device comprising a storage device and, stored in the device, at least 10 atomic coordinates selected from the atomic coordinates listed in Figs. 5-1 to 5-28. In one embodiment, the storage device is a computer readable device that stores code that receives as input the atomic coordinates. In another embodiment, the computer readable device is a floppy disk or a hard drive.

#### DETAILED DESCRIPTION OF THE INVENTION

#### 30 I. DEFINITIONS

5

10

15

20

25

Unless defined otherwise, all technical and scientific terms used herein have the same meaning as commonly understood by those of ordinary skill in the art to which this invention belongs. Although any methods and materials similar or equivalent to those described herein can be used in the practice or testing of the present invention, the preferred methods and materials are described. For purposes of the present invention, the following terms are defined below.

5

"Binding pair" refers to two moieties (e.g. chemical or biochemical) that have an affinity for one another. Examples of binding pairs include antigen/antibodies, lectin/avidin, target polynucleotide/probe oligonucleotide, antibody/anti-antibody, receptor/ligand, enzyme/ligand and the like. "One member of a binding pair" refers to one moiety of the pair, such as an antigen or ligand.

10

"Nucleic acid" refers to a deoxyribonucleotide or ribonucleotide polymer in either single- or double-stranded form, and, unless otherwise limited, encompasses known analogs of natural nucleotides that can function in a similar manner as naturally occurring nucleotides. It will be understood that when a nucleic acid molecule is represented by a DNA sequence, this also includes RNA molecules having the corresponding RNA sequence in which "U" replaces "T."

15

"Recombinant nucleic acid molecule" refers to a nucleic acid molecule which is not naturally occurring, and which comprises two nucleotide sequences which are not naturally joined together. Recombinant nucleic acid molecules are produced by artificial recombination, e.g., genetic engineering techniques or chemical synthesis.

20

25

Reference to a nucleotide sequence "encoding" a polypeptide means that the sequence, upon transcription and translation of mRNA, produces the polypeptide. This includes both the coding strand, whose nucleotide sequence is identical to mRNA and whose sequence is usually provided in the sequence listing, as well as its complementary strand, which is used as the template for transcription. As any person skilled in the art recognizes, this also includes all degenerate nucleotide sequences encoding the same amino acid sequence. Nucleotide sequences encoding a polypeptide include sequences containing introns.

30

"Expression control sequences" refers to nucleotide sequences that regulate the expression of a nucleotide sequence to which they are operatively linked. Expression control sequences are "operatively linked" to a nucleotide sequence when the expression control sequences control and regulate the transcription and, as appropriate, translation of the nucleotide sequence. Thus, expression control sequences can include appropriate

5

10

15

20

25

30

promoters, enhancers, transcription terminators, a start codon (i.e., ATG) in front of a protein-encoding gene, splicing signals for introns, maintenance of the correct reading frame of that gene to permit proper translation of the mRNA, and stop codons.

"Naturally-occurring" as used herein, as applied to an object, refers to the fact that an object can be found in nature. For example, a polypeptide or polynucleotide sequence that is present in an organism (including viruses) that can be isolated from a source in nature and which has not been intentionally modified by man in the laboratory is naturally-occurring.

"Operably linked" refers to a juxtaposition wherein the components so described are in a relationship permitting them to function in their intended manner. A control sequence "operably linked" to a coding sequence is ligated in such a way that expression of the coding sequence is achieved under conditions compatible with the control sequences, such as when the appropriate molecules (e.g., inducers and polymerases) are bound to the control or regulatory sequence(s).

"Control sequence" refers to polynucleotide sequences which are necessary to effect the expression of coding and non-coding sequences to which they are ligated. The nature of such control sequences differs depending upon the host organism; in prokaryotes, such control sequences generally include promoter, ribosomal binding site, and transcription termination sequence; in eukaryotes, generally, such control sequences include promoters and transcription termination sequence. The term "control sequences" is intended to include, at a minimum, components whose presence can influence expression, and can also include additional components whose presence is advantageous, for example, leader sequences and fusion partner sequences.

"Isolated polynucleotide" refers a polynucleotide of genomic, cDNA, or synthetic origin or some combination there of, which by virtue of its origin the "isolated polynucleotide" (1) is not associated with the cell in which the "isolated polynucleotide" is found in nature, or (2) is operably linked to a polynucleotide which it is not linked to in nature.

"Polynucleotide" refers to a polymeric form of nucleotides of at least 10 bases in length, either ribonucleotides or deoxynucleotides or a modified form of either type of nucleotide. The term includes single and double stranded forms of DNA.

The term "probe" refers to a substance that specifically binds to another substance (a "target"). Probes include, for example, antibodies, nucleic acids, receptors and

their ligands.

5

10

15

20

25

30

"Modulation" refers to the capacity to either enhance or inhibit a functional property of biological activity or process (e.g., enzyme activity or receptor binding); such enhancement or inhibition may be contingent on the occurrence of a specific event, such as activation of a signal transduction pathway, and/or may be manifest only in particular cell types.

The term "modulator" refers to a chemical (naturally occurring or non-naturally occurring), such as a synthetic molecule (e.g., nucleic acid, protein, non-peptide, or organic molecule), or an extract made from biological materials such as bacteria, plants, fungi, or animal (particularly mammalian) cells or tissues. Modulators can be evaluated for potential activity as inhibitors or activators (directly or indirectly) of a biological process or processes (e.g., agonist, partial antagonist, partial agonist, inverse agonist, antagonist, antineoplastic agents, cytotoxic agents, inhibitors of neoplastic transformation or cell proliferation, cell proliferation-promoting agents, and the like) by inclusion in screening assays described herein. The activity of a modulator may be known, unknown or partially known.

The term "test chemical" refers to a chemical to be tested by one or more screening method(s) of the invention as a putative modulator. A test chemical is usually not known to bind to the target of interest. The term "control test chemical" refers to a chemical known to bind to the target (e.g., a known agonist, antagonist, partial agonist or inverse agonist). Usually, various predetermined concentrations of test chemicals are used for screening, such as  $.01 \mu M$ ,  $.1 \mu M$ ,  $1.0 \mu M$ , and  $10.0 \mu M$ .

The term "target" refers to a biochemical entity involved a biological process.

Targets are typically proteins that play a useful role in the physiology or biology of an organism. A therapeutic chemical binds to target to alter or modulate its function. As used herein targets can include cell surface receptors, G-proteins, kinases, ion channels, phopholipases and other proteins mentioned herein.

The term "label" refers to a composition detectable by spectroscopic, photochemical, biochemical, immunochemical, or chemical means. For example, useful labels include <sup>32</sup>P, fluorescent dyes, fluorescent proteins, electron-dense reagents, enzymes (e.g., as commonly used in an ELISA), biotin, dioxigenin, or haptens and proteins for which antisera or monoclonal antibodies are available. For example, polypeptides of this invention can be made as detectible labels, by e.g., incorporating a them as into a polypeptide, and

used to label antibodies specifically reactive with the polypeptide. A label often generates a measurable signal, such as radioactivity, fluorescent light or enzyme activity, which can be used to quantitate the amount of bound label.

5

10

15

20

25

30

The term "nucleic acid probe" refers to a nucleic acid molecule that binds to a specific sequence or sub-sequence of another nucleic acid molecule. A probe is preferably a nucleic acid molecule that binds through complementary base pairing to the full sequence or to a sub-sequence of a target nucleic acid. It will be understood that probes may bind target sequences lacking complete complementarity with the probe sequence depending upon the stringency of the hybridization conditions. Probes are preferably directly labelled as with isotopes, chromophores, lumiphores, chromogens, fluorescent proteins, or indirectly labelled such as with biotin to which a streptavidin complex may later bind. By assaying for the presence or absence of the probe, one can detect the presence or absence of the select sequence or sub-sequence.

A "labeled nucleic acid probe" is a nucleic acid probe that is bound, either covalently, through a linker, or through ionic, van der Waals or hydrogen bonds to a label such that the presence of the probe may be detected by detecting the presence of the label bound to the probe.

The terms "polypeptide" and "protein" refers to a polymer of amino acid residues. The terms apply to amino acid polymers in which one or more amino acid residue is an artificial chemical analogue of a corresponding naturally occurring amino acid, as well as to naturally occurring amino acid polymers. The term "recombinant protein" refers to a protein that is produced by expression of a nucleotide sequence encoding the amino acid sequence of the protein from a recombinant DNA molecule.

The term "recombinant host cell" refers to a cell that comprises a recombinant nucleic acid molecule. Thus, for example, recombinant host cells can express genes that are not found within the native (non-recombinant) form of the cell.

The terms "isolated" "purified" or "biologically pure" refer to material which is substantially or essentially free from components which normally accompany it as found in its native state. Purity and homogeneity are typically determined using analytical chemistry techniques such as polyacrylamide gel electrophoresis or high performance liquid chromatography. A protein or nucleic acid molecule which is the predominant protein or nucleic acid species present in a preparation is substantially purified. Generally, an isolated

protein or nucleic acid molecule will comprise more than 80% of all macromolecular species present in the preparation. Preferably, the protein is purified to represent greater than 90% of all macromolecular species present. More preferably the protein is purified to greater than 95%, and most preferably the protein is purified to essential homogeneity, wherein other macromolecular species are not detected by conventional techniques.

5

10

15

20

25

30

The term "naturally-occurring" as applied to an object refers to the fact that an object can be found in nature. For example, a polypeptide or polynucleotide sequence that is present in an organism (including viruses) that can be isolated from a source in nature and which has not been intentionally modified by man in the laboratory is naturally-occurring.

The term "antibody" refers to a polypeptide substantially encoded by an immunoglobulin gene or immunoglobulin genes, or fragments thereof, which specifically bind and recognize an analyte (antigen). The recognized immunoglobulin genes include the kappa, lambda, alpha, gamma, delta, epsilon and mu constant region genes, as well as the myriad immunoglobulin variable region genes. Antibodies exist, e.g., as intact immunoglobulins or as a number of well characterized fragments produced by digestion with various peptidases. This includes, e.g., Fab' and F(ab)'2 fragments. The term "antibody," as used herein, also includes antibody fragments either produced by the modification of whole antibodies or those synthesized *de novo* using recombinant DNA methodologies.

The term "immunoassay" refers to an assay that utilizes an antibody to specifically bind an analyte. The immunoassay is characterized by the use of specific binding properties of a particular antibody to isolate, target, and/or quantify the analyte.

The term "identical" in the context of two nucleic acid or polypeptide sequences refers to the residues in the two sequences which are the same when aligned for maximum correspondence. When percentage of sequence identity is used in reference to proteins or peptides it is recognized that residue positions which are not identical often differ by conservative amino acid substitutions, where amino acids residues are substituted for other amino acid residues with similar chemical properties (e.g. charge or hydrophobicity) and therefore do not change the functional properties of the molecule. Where sequences differ in conservative substitutions, the percent sequence identity may be adjusted upwards to correct for the conservative nature of the substitution. Means for

making this adjustment are well known to those of skill in the art. Typically this involves scoring a conservative substitution as a partial rather than a full mismatch, thereby increasing the percentage sequence identity. Thus, for example, where an identical amino acid is given a score of 1 and a non-conservative substitution is given a score of zero, a conservative substitution is given a score between zero and 1. The scoring of conservative substitutions is calculated, e.g., according to known algorithm. See, e.g., Meyers and Miller, Computer Applic. Biol. Sci., 4: 11-17 (1988); Smith and Waterman (1981) Adv. Appl. Math. 2: 482; Needleman and Wunsch (1970) J. Mol. Biol. 48: 443; Pearson and Lipman (1988) Proc. Natl. Acad. Sci. USA 85: 2444; Higgins and Sharp (1988) Gene, 73: 237-244 and Higgins and Sharp (1989) CABIOS 5: 151-153; Corpet, et al. (1988) Nucleic Acids Research 16, 10881-90; Huang, et al. (1992) Computer Applications in the Biosciences 8, 155-65, and Pearson, et al. (1994) Methods in Molecular Biology 24, 307-31. Alignment is also often performed by inspection and manual alignment.

5

10

"Conservatively modified variations" of a particular nucleic acid sequence refers to those nucleic acids which encode identical or essentially identical amino acid 15 sequences, or where the nucleic acid does not encode an amino acid sequence, to essentially identical sequences. Because of the degeneracy of the genetic code, a large number of functionally identical nucleic acids encode any given polypeptide. For instance, the codons CGU, CGC, CGA, CGG, AGA, and AGG all encode the amino acid arginine. Thus, at every position where an arginine is specified by a codon, the codon can be altered to any of 20 the corresponding codons described without altering the encoded polypeptide. Such nucleic acid variations are "silent variations," which are one species of "conservatively modified variations." Every nucleic acid sequence herein which encodes a polypeptide also describes every possible silent variation. One of skill will recognize that each codon in a nucleic acid 25 (except AUG, which is ordinarily the only codon for methionine) can be modified to yield a functionally identical molecule by standard techniques. Accordingly, each "silent variation" of a nucleic acid which encodes a polypeptide is implicit in each described sequence. Furthermore, one of skill will recognize that individual substitutions, deletions or additions which alter, add or delete a single amino acid or a small percentage of amino acids (typically less than 5%, more typically less than 1%) in an encoded sequence are 30 "conservatively modified variations" where the alterations result in the substitution of an amino acid with a chemically similar amino acid. Conservative amino acid substitutions

PCT/US97/14593 WO 98/06737 15

providing functionally similar amino acids are well known in the art. The following six groups each contain amino acids that are conservative substitutions for one another:

- 1) Alanine (A), Serine (S), Threonine (T);
- 2) Aspartic acid (D), Glutamic acid (E);
- 3) Asparagine (N), Glutamine (Q);

5

15

20

25

30

- 4) Arginine (R), Lysine (K);
- 5) Isoleucine (I), Leucine (L), Methionine (M), Valine (V); and
- 6) Phenylalanine (F), Tyrosine (Y), Tryptophan (W).

The term "complementary" means that one nucleic acid molecule has the sequence of the binding partner of another nucleic acid molecule. Thus, the sequence 5'-10 ATGC-3' is complementary to the sequence 5'-GCAT-3'.

An amino acid sequence or a nucleotide sequence is "substantially identical" or "substantially similar" to a reference sequence if the amino acid sequence or nucleotide sequence has at least 80% sequence identity with the reference sequence over a given comparison window. Thus, substantially similar sequences include those having, for example, at least 85% sequence identity, at least 90% sequence identity, at least 95% sequence identity or at least 99% sequence identity. Two sequences that are identical to each other are, of course, also substantially identical.

A subject nucleotide sequence is "substantially complementary" to a reference nucleotide sequence if the complement of the subject nucleotide sequence is substantially identical to the reference nucleotide sequence.

The term "stringent conditions" refers to a temperature and ionic conditions used in nucleic acid hybridization. Stringent conditions are sequence dependent and are different under different environmental parameters. Generally, stringent conditions are selected to be about  $5\Box C$  to  $20\Box C$  lower than the thermal melting point  $(T_m)$  for the specific sequence at a defined ionic strength and pH. The T<sub>m</sub> is the temperature (under defined ionic strength and pH) at which 50% of the target sequence hybridizes to a perfectly matched probe.

The term "allelic variants" refers to polymorphic forms of a gene at a particular genetic locus, as well as cDNAs derived from mRNA transcripts of the genes and the polypeptides encoded by them.

The term "preferred mammalian codon" refers to the subset of codons from

PCT/US97/14593 WO 98/06737

among the set of codons encoding an amino acid that are most frequently used in proteins expressed in mammalian cells as chosen from the following list:

Amino Acid Preferred codons for high level mammalian expression

5		
	Gly	GGC,GGG
	Glu	GAG
	Asp	GAC
	Val	GUG,GUC
10	Ala	GCC,GCU
	Ser	AGC,UCC
	Lys	AAG
	Asn	AAC
	Met	AUG
15	Ile	AUC
	Thr	ACC
	Trp	UGG
	Cys	UGC
	Tyr	UAU,UAC
20	Leu	CUG
	Phe	UUC
	Arg	CGC,AGG,AGA
	Gln	CAG
	His	CAC
25	Pro	CCC

30

35

Fluorescent molecules are useful in fluorescence resonance energy transfer ("FRET"). FRET involves a donor molecule and an acceptor molecule. To optimize the efficiency and detectability of FRET between a donor and acceptor molecule, several factors need to be balanced. The emission spectrum of the donor should overlap as much as possible with the excitation spectrum of the acceptor to maximize the overlap integral. Also, the quantum yield of the donor moiety and the extinction coefficient of the acceptor should likewise be as high as possible to maximize R<sub>0</sub>, the distance at which energy transfer efficiency is 50%. However, the excitation spectra of the donor and acceptor should overlap as little as possible so that a wavelength region can be found at which the donor can be excited efficiently without directly exciting the acceptor. Fluorescence arising from direct excitation of the acceptor is difficult to distinguish from fluorescence arising from FRET. Similarly, the emission spectra of the donor and acceptor should overlap as little as possible so that the two emissions can be clearly distinguished. High fluorescence quantum yield of

the acceptor moiety is desirable if the emission from the acceptor is to be measured either as the sole readout or as part of an emission ratio. One factor to be considered in choosing the donor and acceptor pair is the efficiency of fluorescence resonance energy transfer between them. Preferably, the efficiency of FRET between the donor and acceptor is at least 10%, more preferably at least 50% and even more preferably at least 80%.

5

10

15

20

The term "fluorescent property" refers to the molar extinction coefficient at an appropriate excitation wavelength, the fluorescence quantum efficiency, the shape of the excitation spectrum or emission spectrum, the excitation wavelength maximum and emission wavelength maximum, the ratio of excitation amplitudes at two different wavelengths, the ratio of emission amplitudes at two different wavelengths, the excited state lifetime, or the fluorescence anisotropy. A measurable difference in any one of these properties between wild-type Aequorea GFP and the mutant form is useful. A measurable difference can be determined by determining the amount of any quantitative fluorescent property, e.g., the amount of fluorescence at a particular wavelength, or the integral of fluorescence over the emission spectrum. Determining ratios of excitation amplitude or emission amplitude at two different wavelengths ("excitation amplitude ratioing" and "emission amplitude ratioing", respectively) are particularly advantageous because the ratioing process provides an internal reference and cancels out variations in the absolute brightness of the excitation source, the sensitivity of the detector, and light scattering or quenching by the sample.

#### II. LONG WAVELENGTH ENGINEERED FLUORESCENT PROTEINS

### A. Fluorescent Proteins

5

10

15

20

25

30

As used herein, the term "fluorescent protein" refers to any protein capable of fluorescence when excited with appropriate electromagnetic radiation. This includes fluorescent proteins whose amino acid sequences are either naturally occurring or engineered (i.e., analogs or mutants). Many cnidarians use green fluorescent proteins ("GFPs") as energy-transfer acceptors in bioluminescence. A "green fluorescent protein," as used herein, is a protein that fluoresces green light. Similarly, "blue fluorescent proteins" fluoresce blue light and "red fluorescent proteins" fluoresce red light. GFPs have been isolated from the Pacific Northwest jellyfish, Aequorea victoria, the sea pansy, Renilla reniformis, and Phialidium gregarium. W.W. Ward et al., Photochem. Photobiol., 35:803-808 (1982); L.D. Levine et al., Comp. Biochem. Physiol., 72B:77-85 (1982).

A variety of Aequorea-related fluorescent proteins having useful excitation and emission spectra have been engineered by modifying the amino acid sequence of a naturally occurring GFP from Aequorea victoria. (D.C. Prasher et al., Gene, 111:229-233 (1992); R. Heim et al., Proc. Natl. Acad. Sci., USA, 91:12501-04 (1994); U.S. patent application 08/337,915, filed November 10, 1994; International application PCT/US95/14692, filed 11/10/95.)

As used herein, a fluorescent protein is an "Aequorea-related fluorescent protein" if any contiguous sequence of 150 amino acids of the fluorescent protein has at least 85% sequence identity with an amino acid sequence, either contiguous or non-contiguous, from the 238 amino-acid wild-type Aequorea green fluorescent protein of Fig. 3 (SEQ ID NO:2). More preferably, a fluorescent protein is an Aequorea-related fluorescent protein if any contiguous sequence of 200 amino acids of the fluorescent protein has at least 95% sequence identity with an amino acid sequence, either contiguous or non-contiguous, from the wild type Aequorea green fluorescent protein of Fig. 3 (SEQ ID NO:2). Similarly, the fluorescent protein may be related to Renilla or Phialidium wild-type fluorescent proteins using the same standards.

Aequorea-related fluorescent proteins include, for example and without limitation, wild-type (native) Aequorea victoria GFP (D.C. Prasher et al., "Primary structure of the Aequorea victoria green fluorescent protein," Gene, (1992) 111:229-33), whose nucleotide sequence (SEQ ID NO:1) and deduced amino acid sequence (SEQ ID NO:2) are presented in Fig. 3; allelic variants of this sequence, e.g., Q80R, which has the glutamine

residue at position 80 substituted with arginine (M. Chalfie et al., Science, (1994) 263:802-805); those engineered Aequorea-related fluorescent proteins described herein, e.g., in Table A or Table F, variants that include one or more folding mutations and fragments of these proteins that are fluorescent, such as Aequorea green fluorescent protein from which the two amino-terminal amino acids have been removed. Several of these contain different aromatic amino acids within the central chromophore and fluoresce at a distinctly shorter wavelength than wild type species. For example, engineered proteins P4 and P4-3 contain (in addition to other mutations) the substitution Y66H, whereas W2 and W7 contain (in addition to other mutations) Y66W. Other mutations both close to the chromophore region of the protein and remote from it in primary sequence may affect the spectral properties of GFP and are listed in the first part of the table below.

5

10

TABLE A

Clone	Mutation(s)	Excitation max (nm)	Emission max (nm)	Extinct. Coeff. (M <sup>-1</sup> cm <sup>-1</sup> )	Quantum yield
Wild type	None	395 (475)	508	21,000 (7,150)	0.77
P4	Y66H	383	447	13,500	0.21
P4-3	Y66H Y145F	381	445	14,000	0.38
W7	Y66W N146I M153T V163A N212K	433 (453)	475 (501)	18,000 (17,100)	0.67
W2	Y66W I123V Y145H H148R M153T V163A N212K	432 (453)	480	10,000 (9,600)	0.72
S65T	S65T	489	511	39,200	0.68
P4-1	S65T M153A	504 (396)	514	14,500 (8,600)	0.53

	K238E		
S65A	S65A	471	504
S65C	S65C	479	507
S65L	S65L	484	510
Y66F	Y66F	360	442
Y66W	Y66W	458	480

5

10

15

20

Additional mutations in *Aequorea*-related fluorescent proteins, referred to as "folding mutations," improve the ability of fluorescent proteins to fold at higher temperatures, and to be more fluorescent when expressed in mammalian cells, but have little or no effect on the peak wavelengths of excitation and emission. It should be noted that these may be combined with mutations that influence the spectral properties of GFP to produce proteins with altered spectral and folding properties. Folding mutations include: F64L, V68L, S72A, and also T44A, F99S, Y145F, N146I, M153T or A, V163A, I167T, S175G, S205T and N212K.

As used herein, the term "loop domain" refers to an amino acid sequence of an Aequorea-related fluorescent protein that connects the amino acids involved in the secondary structure of the eleven strands of the  $\Box$ -barrel or the central  $\Box$ -helix (residues 56-72) (see Fig. 1A and 1B).

As used herein, the "fluorescent protein moiety" of a fluorescent protein is that portion of the amino acid sequence of a fluorescent protein which, when the amino acid sequence of the fluorescent protein substrate is optimally aligned with the amino acid sequence of a naturally occurring fluorescent protein, lies between the amino terminal and carboxy terminal amino acids, inclusive, of the amino acid sequence of the naturally occurring fluorescent protein.

It has been found that fluorescent proteins can be genetically fused to other target proteins and used as markers to identify the location and amount of the target protein produced. Accordingly, this invention provides fusion proteins comprising a fluorescent protein moiety and additional amino acid sequences. Such sequences can be, for example, up to about 15, up to about 50, up to about 150 or up to about 1000 amino acids long. The

fusion proteins possess the ability to fluoresce when excited by electromagnetic radiation. In one embodiment, the fusion protein comprises a polyhistidine tag to aid in purification of the protein.

# B. Use Of The Crystal Structure Of Green Fluorescent Protein To Design Mutants Having Altered Fluorescent Characteristics

5

10

15

20

25

30

Using X-ray crystallography and computer processing, we have created a model of the crystal structure of Aequorea green fluorescent protein showing the relative location of the atoms in the molecule. This information is useful in identifying amino acids whose substitution alters fluorescent properties of the protein.

Fluorescent characteristics of Aequorea-related fluorescent proteins depend, in part, on the electronic environment of the chromophore. In general, amino acids that are within about 0.5 nm of the chromophore influence the electronic environment of the chromophore. Therefore, substitution of such amino acids can produce fluorescent proteins with altered fluorescent characteristics. In the excited state, electron density tends to shift from the phenolate towards the carbonyl end of the chromophore. Therefore, placement of increasing positive charge near the carbonyl end of the chromophore tends to decrease the energy of the excited state and cause a red-shift in the absorbance and emission wavelength maximum of the protein. Decreasing positive charge near the carbonyl end of the chromophore tends to have the opposte effect, causing a blue-shift in the protein's wavelengths.

Amino acids with charged (ionized D, E, K, and R), dipolar (H, N, Q, S, T, and uncharged D, E and K), and polarizable side groups (e.g., C, F, H, M, W and Y) are useful for altering the electronic environment of the chromophore, especially when they substitute an amino acid with an uncharged, nonpolar or non-polarizable side chain. In general, amino acids with polarizable side groups alter the electronic environment least, and; consequently, are expected to cause a comparatively smaller change in a fluorescent property. Amino acids with charged side groups alter the environment most, and, consequently, are expected to cause a comparatively larger change in a fluorescent property. However, amino acids with charged side groups are more likely to disrupt the structure of the protein and to prevent proper folding if buried next to the chromophore without any

additional solvation or salt bridging. Therefore charged amino acids are most likely to be tolerated and to give useful effects when they replace other charged or highly polar amino acids that are already solvated or involved in salt bridges. In certain cases, where substitution with a polarizable amino acid is chosen, the structure of the protein may make selection of a larger amino acid, e.g., W, less appropriate. Alternatively, positions occupied by amino acids with charged or polar side groups that are unfavorably oriented may be substituted with amino acids that have less charged or polar side groups. In another alternative, an amino acid whose side group has a dipole oriented in one direction in the protein can be substituted with an amino acid having a dipole oriented in a different

10 direction.

5

More particularly, Table B lists several amino acids located within about 0.5 nm from the chromophore whose substitution can result in altered fluorescent characteristics. The table indicates, underlined, preferred amino acid substitutions at the indicated location to alter a fluorescent characteristic of the protein. In order to introduce such substitutions, the table also provides codons for primers used in site-directed mutagenesis involving amplification. These primers have been selected to encode economically the preferred amino acids, but they encode other amino acids as well, as indicated, or even a stop codon, denoted by Z. In introducing substitutions using such degenerate primers the most efficient strategy is to screen the collection to identify mutants with the desired properties and then sequence their DNA to find out which of the possible substitutions is responsible. Codons are shown in double-stranded form with sense strand above, antisense strand below. In nucleic acid sequences, R=(A or g); Y=(C or T); M=(A or C); K=(g or T); S=(g or C); W=(A or T); H=(A, T, or C); B=(g, T, or C); V=(g, A, or C); D=(g, A, or T); N=(A, C, g, or T).

25

20

15

#### TABLE B

	Origina	ol position and presumed role	Change to	Codon	
30	L42	Aliphatic residue near C=N of chromophore	C <u>FH</u> LQR <u>WY</u> Z	5'YDS 3' 3'RHS 5'	
	V61	Aliphatic residue near central -CH= of chromophore FYHCL	R YDC	RHg	

	T62	Almost directly above center of chromophore bridge AVFS	КҮС	MRg
5			<u> D</u> ЕНК <u>И</u> Q	VAS BTS
			<u>FYHC</u> LR	YDC RHg
10	V68	Aliphatic residue near carbonyl and G67	<u>FYH</u> L	YWC RWg
	N121	Near C-N site of ring closure between T65 and G67 CFHLC	or <u>wy</u> z yds	RHS
15	Y145	Packs near tyrosine ring of chromophore	WCFL	TKS AMS
20			<u>DEHNKQ</u>	VAS BTS
	H148	H-bonds to phenolate oxygen	<u>FYN</u> I	wwc wwg
25			<u>KQR</u>	MRg KYC
	V150	Aliphatic residue near tyrosine ring of chromophore FYHL	YWC	RWg
30	F165	Packs near tyrosine ring	C <u>HO</u> R <u>WY</u> Z	YRS RYS
35	1167	Aliphatic residue near phenolate; I167T has effects	FYHL	YWC RWg
	T203	H-bonds to phenolic oxygen of chromophore	<u>FH</u> LQR <u>WY</u> Z	YDS RHS
40	E222	Protonation regulates ionization of chromophore	HKNQ	MAS KTS

Examples of amino acids with polar side groups that can be substituted with polarizable side groups include, for example, those in Table C.

24

#### TABLE C

	Origina	al position and presumed role	Change to	Codon
5	Q69	Terminates chain of H-bonding waters	KREG	RRg YYC
10	Q94	H-bonds to carbonyl terminus of chromophore	<u>DEHKN</u> Q	VAS BTS
10	Q183	Bridges Arg96 and center of chromophore bridge	<u>HY</u>	YAC RTG
15			<u>EK</u>	RAg YTC
	N185	Part of H-bond network near carbonyl of chromophore	DEHNKQ	VAS BTS

In another embodiment, an amino acid that is close to a second amino acid within about 0.5 nm of the chromophore can, upon substitution, alter the electronic properties of the second amino acid, in turn altering the electronic environment of the chromphore. Table D presents two such amino acids. The amino acids, L220 and V224, are close to E222 and oriented in the same direction in the  $\Box$  pleated sheet.

25

#### TABLE D

20	Origina	l position and presumed role	Change to	Codon
30	L220	Packs next to Glu222; to make GFP pH sensitive	<u>HKN</u> PQT	MMS KKS
35	V224	Packs next to Glu222; to make GFP pH sensitive	<u>нкирот</u>	MMS KKS
			C <u>FH</u> LQR <u>WY</u> Z	YDS RHS

One embodiment of the invention includes a nucleic acid molecule comprising a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of *Aequorea* green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least a substitution at Q69, wherein the functional engineered fluorescent protein has a different fluorescent property than *Aequorea* green fluorescent protein. Preferably, the substitution at Q69 is selected from the group of K, R, E and G. The Q69 substitution can be combined with other mutations to improve the properties of the protein, such as a functional mutation at S65.

5

10

15

20

One embodiment of the invention includes a nucleic acid molecule comprising a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of *Aequorea* green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least a substitution at E222, but not including E222G, wherein the functional engineered fluorescent protein has a different fluorescent property than *Aequorea* green fluorescent protein. Preferably, the substitution at E222 is selected from the group of N and Q. The E222 substitution can be combined with other mutations to improve the properties of the protein, such as a functional mutation at F64.

One embodiment of the invention includes a nucleic acid molecule comprising a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of *Aequorea* green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least a substitution at Y145, wherein the functional engineered fluorescent protein has a different fluorescent property than *Aequorea* green fluorescent protein.

Preferably, the substitution at Y145 is selected from the group of W, C, F, L, E, H, K and Q.

The Y145 substitution can be combined with other mutations to improve the properties of the protein, such as a Y66.

WO 98/06737 PCT/US97/14593

The invention also includes computer related embodiments, including computational methods of using the crystal coordinates for designing new fluorescent protein mutations and devices for storing the crystal data, including coordinates. For instance the invention includes a device comprising a storage device and, stored in the device, at least 10 atomic coordinates selected from the atomic coordinates listed in Figs. 5-1 to 5-28. More coordinates can be storage depending of the complexity of the calculations or the objective of using the coordinates (e.g. about 100, 1,000, or more coordinates). For example, larger numbers of coordinates will be desirable for more detailed representations of fluorescent protein structure. Typically, the storage device is a computer readable device that stores code that it receives as input the atomic coordinates. Although, other storage meand as known in the art are contemplated. The computer readable device can be a floppy disk or a hard drive.

5

10

5

10

15

20

25

30

#### C. Production Of Long Wavelength Engineered Fluorescent Proteins

Recombinant production of a fluorescent protein involves expressing a nucleic acid molecule having sequences that encode the protein.

In one embodiment, the nucleic acid encodes a fusion protein in which a single polypeptide includes the fluorescent protein moiety within a longer polypeptide. The longer polypeptide can include a second functional protein, such as FRET partner or a protein having a second function (e.g., an enzyme, antibody or other binding protein). Nucleic acids that encode fluorescent proteins are useful as starting materials.

The fluorescent proteins can be produced as fusion proteins by recombinant DNA technology. Recombinant production of fluorescent proteins involves expressing nucleic acids having sequences that encode the proteins. Nucleic acids encoding fluorescent proteins can be obtained by methods known in the art. Fluorescent proteins can be made by site-specific mutagenesis of other nucleic acids encoding fluorescent proteins, or by random mutagenesis caused by increasing the error rate of PCR of the original polynucleotide with 0.1 mM MnCl<sub>2</sub> and unbalanced nucleotide concentrations. See, e.g., U.S. patent application 08/337,915, filed November 10, 1994 or International application PCT/US95/14692, filed 11/10/95. The nucleic acid encoding a green fluorescent protein can be isolated by polymerase chain reaction of cDNA from A. victoria using primers based on the DNA sequence of A. victoria green fluorescent protein, as presented in Fig. 3. PCR methods are described in, for example, U.S. Pat. No. 4,683,195; Muilis et al. (1987) Cold Spring Harbor Symp. Quant. Biol. 51:263; and Erlich, ed., PCR Technology, (Stockton Press, NY, 1989).

The construction of expression vectors and the expression of genes in transfected cells involves the use of molecular cloning techniques also well known in the art. Sambrook et al., *Molecular Cloning -- A Laboratory Manual*, Cold Spring Harbor Laboratory, Cold Spring Harbor, NY, (1989) and *Current Protocols in Molecular Biology*, F.M. Ausubel et al., eds., (Current Protocols, a joint venture between Greene Publishing Associates, Inc. and John Wiley & Sons, Inc.). The expression vector can be adapted for function in prokaryotes or eukaryotes by inclusion of appropriate promoters, replication sequences, markers, etc.

Nucleic acids used to transfect cells with sequences coding for expression of the polypeptide of interest generally will be in the form of an expression vector including

expression control sequences operatively linked to a nucleotide sequence coding for expression of the polypeptide. As used, the term "nucleotide sequence coding for expression of" a polypeptide refers to a sequence that, upon transcription and translation of mRNA, produces the polypeptide. This can include sequences containing, e.g., introns.

Expression control sequences are operatively linked to a nucleic acid sequence when the expression control sequences control and regulate the transcription and, as appropriate, translation of the nucleic acid sequence. Thus, expression control sequences can include appropriate promoters, enhancers, transcription terminators, a start codon (i.e., ATG) in front of a protein-encoding gene, splicing signals for introns, maintenance of the correct reading frame of that gene to permit proper translation of the mRNA, and stop codons.

5

10

15

20

25

30

Methods which are well known to those skilled in the art can be used to construct expression vectors containing the fluorescent protein coding sequence and appropriate transcriptional/translational control signals. These methods include *in vitro* recombinant DNA techniques, synthetic techniques and *in vivo* recombination/genetic recombination. (See, for example, the techniques described in Maniatis, *et al.*, *Molecular Cloning A Laboratory Manual*, Cold Spring Harbor Laboratory, N.Y., 1989).

Transformation of a host cell with recombinant DNA may be carried out by conventional techniques as are well known to those skilled in the art. Where the host is prokaryotic, such as E. coli, competent cells which are capable of DNA uptake can be prepared from cells harvested after exponential growth phase and subsequently treated by the CaCl<sub>2</sub> method by procedures well known in the art. Alternatively, MgCl<sub>2</sub> or RbCl can be used. Transformation can also be performed after forming a protoplast of the host cell or by electroporation.

When the host is a eukaryote, such methods of transfection of DNA as calcium phosphate co-precipitates, conventional mechanical procedures such as microinjection, electroporation, insertion of a plasmid encased in liposomes, or virus vectors may be used. Eukaryotic cells can also be cotransfected with DNA sequences encoding the fusion polypeptide of the invention, and a second foreign DNA molecule encoding a selectable phenotype, such as the herpes simplex thymidine kinase gene. Another method is to use a eukaryotic viral vector, such as simian virus 40 (SV40) or bovine papilloma virus, to transiently infect or transform eukaryotic cells and express the protein. (Eukaryotic Viral

Vectors, Cold Spring Harbor Laboratory, Gluzman ed., 1982). Preferably, a eukaryotic host is utilized as the host cell as described herein.

Techniques for the isolation and purification of either microbially or eukaryotically expressed polypeptides of the invention may be by any conventional means such as, for example, preparative chromatographic separations and immunological separations such as those involving the use of monoclonal or polyclonal antibodies or antigen.

5

10

15

20

25

30

In one embodiment recombinant fluorescent proteins can be produced by expression of nucleic acid encoding for the protein in *E. coli. Aequorea*-related fluorescent proteins are best expressed by cells cultured between about 15 \(\text{C}\) and 30 \(\text{C}\) but higher temperatures (e.g. 37 \(\text{C}\)) are possible. After synthesis, these enzymes are stable at higher temperatures (e.g., 37 \(\text{C}\)) and can be used in assays at those temperatures.

A variety of host-expression vector systems may be utilized to express fluorescent protein coding sequence. These include but are not limited to microorganisms such as bacteria transformed with recombinant bacteriophage DNA, plasmid DNA or cosmid DNA expression vectors containing a fluorescent protein coding sequence; yeast transformed with recombinant yeast expression vectors containing the fluorescent protein coding sequence; plant cell systems infected with recombinant virus expression vectors (e.g., cauliflower mosaic virus, CaMV; tobacco mosaic virus, TMV) or transformed with recombinant plasmid expression vectors (e.g., Ti plasmid) containing a fluorescent protein coding sequence; insect cell systems infected with recombinant virus expression vectors (e.g., baculovirus) containing a fluorescent protein coding sequence; or animal cell systems infected with recombinant virus expression vectors (e.g., retroviruses, adenovirus, vaccinia virus) containing a fluorescent protein coding sequence, or transformed animal cell systems engineered for stable expression.

Depending on the host/vector system utilized, any of a number of suitable transcription and translation elements, including constitutive and inducible promoters, transcription enhancer elements, transcription terminators, etc. may be used in the expression vector (see, e.g., Bitter, et al., Methods in Enzymology 153:516-544, 1987). For example, when cloning in bacterial systems, inducible promoters such as pL of bacteriophage  $\Box$ , plac, ptrp, ptac (ptrp-lac hybrid promoter) and the like may be used. When cloning in mammalian cell systems, promoters derived from the genome of mammalian cells (e.g., metallothionein promoter) or from mammalian viruses (e.g., the

retrovirus long terminal repeat; the adenovirus late promoter; the vaccinia virus 7.5K promoter) may be used. Promoters produced by recombinant DNA or synthetic techniques may also be used to provide for transcription of the inserted fluorescent protein coding sequence.

5

10

15

20

25

30

In bacterial systems a number of expression vectors may be advantageously selected depending upon the use intended for the fluorescent protein expressed. For example, when large quantities of the fluorescent protein are to be produced, vectors which direct the expression of high levels of fusion protein products that are readily purified may be desirable. Those which are engineered to contain a cleavage site to aid in recovering fluorescent protein are preferred.

In yeast, a number of vectors containing constitutive or inducible promoters may be used. For a review see, Current Protocols in Molecular Biology, Vol. 2, Ed. Ausubel, et al., Greene Publish. Assoc. & Wiley Interscience, Ch. 13, 1988; Grant, et al., Expression and Secretion Vectors for Yeast, in Methods in Enzymology, Eds. Wu & Grossman, 31987, Acad. Press, N.Y., Vol. 153, pp.516-544, 1987; Glover, DNA Cloning, Vol. II, IRL Press, Wash., D.C., Ch. 3, 1986; and Bitter, Heterologous Gene Expression in Yeast, Methods in Enzymology, Eds. Berger & Kimmel, Acad. Press, N.Y., Vol. 152, pp. 673-684, 1987; and The Molecular Biology of the Yeast Saccharomyces, Eds. Strathern et al., Cold Spring Harbor Press, Vols. I and II, 1982. A constitutive yeast promoter such as ADH or LEU2 or an inducible promoter such as GAL may be used (Cloning in Yeast, Ch. 3, R. Rothstein In: DNA Cloning Vol.11, A Practical Approach, Ed. DM Glover, IRL Press, Wash., D.C., 1986). Alternatively, vectors may be used which promote integration of foreign DNA sequences into the yeast chromosome.

In cases where plant expression vectors are used, the expression of a fluorescent protein coding sequence may be driven by any of a number of promoters. For example, viral promoters such as the 35S RNA and 19S RNA promoters of CaMV (Brisson, et al., Nature 310:511-514, 1984), or the coat protein promoter to TMV (Takamatsu, et al., EMBO J. 6:307-311, 1987) may be used; alternatively, plant promoters such as the small subunit of RUBISCO (Coruzzi, et al., 1984, EMBO J. 3:1671-1680; Broglie, et al., Science 224:838-843, 1984); or heat shock promoters, e.g., soybean hsp17.5-E or hsp17.3-B (Gurley, et al., Mol. Cell. Biol. 6:559-565, 1986) may be used. These constructs can be introduced into plant cells using Ti plasmids, Ri plasmids, plant virus vectors, direct DNA transformation,

microinjection, electroporation, etc. For reviews of such techniques see, for example, Weissbach & Weissbach, Methods for Plant Molecular Biology, Academic Press, NY, Section VIII, pp. 421-463, 1988; and Grierson & Corey, Plant Molecular Biology, 2d Ed., Blackie, London, Ch. 7-9, 1988.

5

10

15

20

25

30

An alternative expression system which could be used to express fluorescent protein is an insect system. In one such system, Autographa californica nuclear polyhedrosis virus (AcNPV) is used as a vector to express foreign genes. The virus grows in Spodoptera frugiperda cells. The fluorescent protein coding sequence may be cloned into non-essential regions (for example, the polyhedrin gene) of the virus and placed under control of an AcNPV promoter (for example the polyhedrin promoter). Successful insertion of the fluorescent protein coding sequence will result in inactivation of the polyhedrin gene and production of non-occluded recombinant virus (i.e., virus lacking the proteinaceous coat coded for by the polyhedrin gene). These recombinant viruses are then used to infect Spodoptera frugiperda cells in which the inserted gene is expressed, see Smith, et al., J. Viol. 46:584, 1983; Smith, U.S. Patent No. 4,215,051.

Eukaryotic systems, and preferably mammalian expression systems, allow for proper post-translational modifications of expressed mammalian proteins to occur. Eukaryotic cells which possess the cellular machinery for proper processing of the primary transcript, glycosylation, phosphorylation, and, advantageously secretion of the gene product should be used as host cells for the expression of fluorescent protein. Such host cell lines may include but are not limited to CHO, VERO, BHK, HeLa, COS, MDCK, Jurkat, HEK-293, and WI38.

Mammalian cell systems which utilize recombinant viruses or viral elements to direct expression may be engineered. For example, when using adenovirus expression vectors, the fluorescent protein coding sequence may be ligated to an adenovirus transcription/translation control complex, e.g., the late promoter and tripartite leader sequence. This chimeric gene may then be inserted in the adenovirus genome by in vitro or in vivo recombination. Insertion in a non-essential region of the viral genome (e.g., region E1 or E3) will result in a recombinant virus that is viable and capable of expressing the fluorescent protein in infected hosts (e.g., see Logan & Shenk, Proc. Natl. Acad. Sci. USA, 81: 3655-3659, 1984). Alternatively, the vaccinia virus 7.5K promoter may be used. (e.g., see, Mackett, et al., Proc. Natl. Acad. Sci. USA, 79: 7415-7419, 1982; Mackett, et al., J.

Virol. 49: 857-864, 1984; Panicali, et al., Proc. Natl. Acad. Sci. USA 79: 4927-4931, 1982). Of particular interest are vectors based on bovine papilloma virus which have the ability to replicate as extrachromosomal elements (Sarver, et al., Mol. Cell. Biol. 1: 486, 1981). Shortly after entry of this DNA into mouse cells, the plasmid replicates to about 100 to 200 copies per cell. Transcription of the inserted cDNA does not require integration of the plasmid into the host's chromosome, thereby yielding a high level of expression. These vectors can be used for stable expression by including a selectable marker in the plasmid, such as the neo gene. Alternatively, the retroviral genome can be modified for use as a vector capable of introducing and directing the expression of the fluorescent protein gene in host cells (Cone & Mulligan, Proc. Natl. Acad. Sci. USA, 81:6349-6353, 1984). High level expression may also be achieved using inducible promoters, including, but not limited to, the metallothionine IIA promoter and heat shock promoters.

5

10

15

20

25

30

The invention can also include a localization sequence, such as a nuclear localization sequence, an endoplasmic reticulum localization sequence, a peroxisome localization sequence, a mitochondrial localization sequence, or a localized protein.

Localization sequences can be targeting sequences which are described, for example, in "Protein Targeting", chapter 35 of Stryer, L., Biochemistry (4th ed.). W.H. Freeman, 1995. The localization sequence can also be a localized protein. Some important localization sequences include those targeting the nucleus (KKKRK), mitochondrion (amino terminal MLRTSSLFTRRVQPSLFRNILRLQST-), endoplasmic reticulum (KDEL at C-terminus, assuming a signal sequence present at N-terminus), peroxisome (SKF at C-terminus), prenylation or insertion into plasma membrane (CaaX, CC, CXC, or CCXX at C-terminus), cytoplasmic side of plasma membrane (fusion to SNAP-25), or the Golgi apparatus (fusion to furin).

For long-term, high-yield production of recombinant proteins, stable expression is preferred. Rather than using expression vectors which contain viral origins of replication, host cells can be transformed with the fluorescent protein cDNA controlled by appropriate expression control elements (e.g., promoter, enhancer, sequences, transcription terminators, polyadenylation sites, etc.), and a selectable marker. The selectable marker in the recombinant plasmid confers resistance to the selection and allows cells to stably integrate the plasmid into their chromosomes and grow to form foci which in turn can be cloned and expanded into cell lines. For example, following the introduction of foreign DNA,

5

10

15

20

25

30

engineered cells may be allowed to grow for 1-2 days in an enriched media, and then are switched to a selective media. A number of selection systems may be used, including but not limited to the herpes simplex virus thymidine kinase (Wigler, et al., Cell, 11: 223, 1977), hypoxanthine-guanine phosphoribosyltransferase (Szybalska & Szybalski, Proc. Natl. Acad. Sci. USA, 48:2026, 1962), and adenine phosphoribosyltransferase (Lowy, et al., Cell, 22: 817, 1980) genes can be employed in tk', hgprt' or aprt' cells respectively. Also, antimetabolite resistance can be used as the basis of selection for dhfr, which confers resistance to methotrexate (Wigler, et al., Proc. Natl. Acad. Sci. USA, 77: 3567, 1980; O'Hare, et al., Proc. Natl. Acad. Sci. USA, 8: 1527, 1981); gpt, which confers resistance to mycophenolic acid (Mulligan & Berg, Proc. Natl. Acad. Sci. USA, 78: 2072, 1981; neo, which confers resistance to the aminoglycoside G-418 (Colberre-Garapin, et al., J. Mol. Biol., 150:1, 1981); and hygro, which confers resistance to hygromycin (Santerre, et al., Gene, 30: 147, 1984) genes. Recently, additional selectable genes have been described, namely trpB, which allows cells to utilize indole in place of tryptophan; hisD, which allows cells to utilize histinol in place of histidine (Hartman & Mulligan, Proc. Natl. Acad. Sci. USA, 85:8047, 1988); and ODC (ornithine decarboxylase) which confers resistance to the ornithine decarboxylase inhibitor, 2-(difluoromethyl)-DL-ornithine, DFMO (McConlogue L., In: Current Communications in Molecular Biology, Cold Spring Harbor Laboratory, ed., 1987).

DNA sequences encoding the fluorescence protein polypeptide of the invention can be expressed *in vitro* by DNA transfer into a suitable host celi. "Host cells" are cells in which a vector can be propagated and its DNA expressed. The term also includes any progeny of the subject host cell. It is understood that all progeny may not be identical to the parental cell since there may be mutations that occur during replication. However, such progeny are included when the term "host cell" is used. Methods of stable transfer, in other words when the foreign DNA is continuously maintained in the host, are known in the art.

The expression vector can be transfected into a host cell for expression of the recombinant nucleic acid. Host cells can be selected for high levels of expression in order to purify the fluorescent protein fusion protein. E. coli is useful for this purpose.

Alternatively, the host cell can be a prokaryotic or eukaryotic cell selected to study the activity of an enzyme produced by the cell. In this case, the linker peptide is selected to

5

10

15

20

25

30

include an amino acid sequence recognized by the protease. The cell can be, e.g., a cultured cell or a cell in vivo.

PCT/US97/14593

A primary advantage of fluorescent protein fusion proteins is that they are prepared by normal protein biosynthesis, thus completely avoiding organic synthesis and the requirement for customized unnatural amino acid analogs. The constructs can be expressed in *E. coli* in large scale for *in vitro* assays. Purification from bacteria is simplified when the sequences include polyhistidine tags for one-step purification by nickel-chelate chromatography. Alternatively, the substrates can be expressed directly in a desired host cell for assays *in situ*.

In another embodiment, the invention provides a transgenic non-human animal that expresses a nucleic acid sequence which encodes the fluorescent protein.

The "non-human animals" of the invention comprise any non-human animal having nucleic acid sequence which encodes a fluorescent protein. Such non-human animals include vertebrates such as rodents, non-human primates, sheep, dog, cow, pig, amphibians, and reptiles. Preferred non-human animals are selected from the rodent family including rat and mouse, most preferably mouse. The "transgenic non-human animals" of the invention are produced by introducing "transgenes" into the germline of the non-human animal. Embryonal target cells at various developmental stages can be used to introduce transgenes. Different methods are used depending on the stage of development of the embryonal target cell. The zygote is the best target for micro-injection. In the mouse, the male pronucleus reaches the size of approximately 20 micrometers in diameter which allows reproducible injection of 1-2 pl of DNA solution. The use of zygotes as a target for gene transfer has a major advantage in that in most cases the injected DNA will be incorporated into the host gene before the first cleavage (Brinster et al., Proc. Natl. Acad. Sci. USA 82:4438-4442, 1985). As a consequence, all cells of the transgenic non-human animal will carry the incorporated transgene. This will in general also be reflected in the efficient transmission of the transgene to offspring of the founder since 50% of the germ cells will harbor the transgene. Microinjection of zygotes is the preferred method for incorporating transgenes in practicing the invention.

The term "transgenic" is used to describe an animal which includes exogenous genetic material within all of its cells. A "transgenic" animal can be produced by cross-breeding two chimeric animals which include exogenous genetic material within cells used

10

15

20

in reproduction. Twenty-five percent of the resulting offspring will be transgenic *i.e.*, animals which include the exogenous genetic material within all of their cells in both alleles. 50% of the resulting animals will include the exogenous genetic material within one allele and 25% will include no exogenous genetic material.

Retroviral infection can also be used to introduce transgene into a non-human animal. The developing non-human embryo can be cultured in vitro to the blastocyst stage. During this time, the blastomeres can be targets for retro viral infection (Jaenich, R., Proc. Natl. Acad. Sci USA 73:1260-1264, 1976). Efficient infection of the blastomeres is obtained by enzymatic treatment to remove the zona pellucida (Hogan, et al. (1986) in Manipulating the Mouse Embryo, Cold Spring Harbor Laboratory Press, Cold Spring Harbor, N.Y.). The viral vector system used to introduce the transgene is typically a replication-defective retro virus carrying the transgene (Jahner, et al., Proc. Natl. Acad. Sci. USA 82:6927-6931, 1985; Van der Putten, et al., Proc. Natl. Acad. Sci USA 82:6148-6152, 1985). Transfection is easily and efficiently obtained by culturing the blastomeres on a monolayer of virus-producing cells (Van der Putten, supra; Stewart, et al., EMBO J. 6:383-388, 1987). Alternatively, infection can be performed at a later stage. Virus or virus-producing cells can be injected into the blastocoele (D. Jahner et al., Nature 298:623-628, 1982). Most of the founders will be mosaic for the transgene since incorporation occurs only in a subset of the cells which formed the transgenic nonhuman animal. Further, the founder may contain various retro viral insertions of the transgene at different positions in the genome which generally will segregate in the offspring. In addition, it is also possible to introduce transgenes into the germ line, albeit with low efficiency, by intrauterine retro viral infection of the midgestation embryo (D. Jahner et al., supra).

A third type of target cell for transgene introduction is the embryonal stem cell

(ES). ES cells are obtained from pre-implantation embryos cultured in vitro and fused with embryos (M. J. Evans et al. Nature 292:154-156, 1981; M.O. Bradley et al., Nature 309: 255-258, 1984; Gossler, et al., Proc. Natl. Acad. Sci USA 83: 9065-9069, 1986; and Robertson et al., Nature 322:445-448, 1986). Transgenes can be efficiently introduced into the ES cells by DNA transfection or by retro virus-mediated transduction. Such transformed ES cells can thereafter be combined with blastocysts from a nonhuman animal. The ES cells thereafter colonize the embryo and contribute to the germ line of the resulting chimeric animal. (For review see Jaenisch, R., Science 240: 1468-1474, 1988).

"Transformed" means a cell into which (or into an ancestor of which) has been introduced, by means of recombinant nucleic acid techniques, a heterologous nucleic acid molecule. "Heterologous" refers to a nucleic acid sequence that either originates from another species or is modified from either its original form or the form primarily expressed in the cell.

PCT/US97/14593

"Transgene" means any piece of DNA which is inserted by artifice into a cell. and becomes part of the genome of the organism (i.e., either stably integrated or as a stable extrachromosomal element) which develops from that cell. Such a transgene may include a gene which is partly or entirely heterologous (i.e., foreign) to the transgenic organism, or may represent a gene homologous to an endogenous gene of the organism. Included within this definition is a transgene created by the providing of an RNA sequence which is transcribed into DNA and then incorporated into the genome. The transgenes of the invention include DNA sequences which encode which encodes the fluorescent protein which may be expressed in a transgenic non-human animal. The term "transgenic" as used herein additionally includes any organism whose genome has been altered by in vitro manipulation of the early embryo or fertilized egg or by any transgenic technology to induce a specific gene knockout. The term "gene knockout" as used herein, refers to the targeted disruption of a gene in vivo with complete loss of function that has been achieved by any transgenic technology familiar to those in the art. In one embodiment, transgenic animals having gene knockouts are those in which the target gene has been rendered nonfunctional by an insertion targeted to the gene to be rendered non-functional by homologous recombination. As used herein, the term "transgenic" includes any transgenic technology familiar to those in the art which can produce an organism carrying an introduced transgene or one in which an endogenous gene has been rendered non-functional or "knocked out."

25

30

5

10

15

20

### III. USES OF ENGINEERED FLUORESCENT PROTEINS

The proteins of this invention are useful in any methods that employ fluorescent proteins.

The engineered fluorescent proteins of this invention are useful as fluorescent markers in the many ways fluorescent markers already are used. This includes, for example, coupling engineered fluorescent proteins to antibodies, nucleic acids or other receptors for use in detection assays, such as immunoassays or hybridization assays.

The engineered fluorescent proteins of this invention are useful to track the movement of proteins in cells. In this embodiment, a nucleic acid molecule encoding the fluorescent protein is fused to a nucleic acid molecule encoding the protein of interest in an expression vector. Upon expression inside the cell, the protein of interest can be localized based on fluorescence. In another version, two proteins of interest are fused with two engineered fluorescent proteins having different fluorescent characteristics.

5

10

15

20

25

30

The engineered fluorescent proteins of this invention are useful in systems to detect induction of transcription. In certain embodiments, a nucleotide sequence encoding the engineered fluorescent protein is fused to expression control sequences of interest and the expression vector is transfected into a cell. Induction of the promoter can be measured by detecting the expression and/or quantity of fluorescence. Such constructs can be used used to follow signaling pathways from receptor to promoter.

The engineered fluorescent proteins of this invention are useful in applications involving FRET. Such applications can detect events as a function of the movement of fluorescent donors and acceptor towards or away from each other. One or both of the donor/acceptor pair can be a fluorescent protein. A preferred donor and receptor pair for FRET based assays is a donor with a T203I mutation and an acceptor with the mutation T203X, wherein X is an aromatic amino acid-39, especially T203Y, T203W, or T203H. In a particularly useful pair the donor contains the following mutations: S72A, K79R, Y145F, M153A and T203I (with a excitation peak of 395 nm and an emission peak of 511 nm) and the acceptor contains the following mutations S65G, S72A, K79R, and T203Y. This particular pair provides a wide separation between the excitation and emission peaks of the donor and provides good overlap between the donor emission spectrum and the acceptor excitation spectrum. Other red-shifted mutants, such as those described herein, can also be used as the acceptor in such a pair.

In one aspect, FRET is used to detect the cleavage of a substrate having the donor and acceptor coupled to the substrate on opposite sides of the cleavage site. Upon cleavage of the substrate, the donor/acceptor pair physically separate, eliminating FRET. Assays involve contacting the substrate with a sample, and determining a qualitative or quantitative change in FRET. In one embodiment, the engineered fluorescent protein is used in a substrate for  $\Box$ -lactamase. Examples of such substrates are described in United States patent applications 08/407,544, filed March 20, 1995 and International Application

PCT/US96/04059, filed March 20, 1996. In another embodiment, an engineered fluorescent protein donor/acceptor pair are part of a fusion protein coupled by a peptide having a proteolytic cleavage site. Such tandem fluorescent proteins are described in United States patent application 08/594,575, filed January 31, 1996.

In another aspect, FRET is used to detect changes in potential across a membrane. A donor and acceptor are placed on opposite sides of a membrane such that one translates across the membrane in response to a voltage change. This creates a measurable FRET. Such a method is described in United States patent application 08/481,977, filed June 7, 1995 and International Application PCT/US96/09652, filed June 6, 1996.

5

10

15

20

The engineered protein of this invention are useful in the creation of fluorescent substrates for protein kinases. Such substrates incorporate an amino acid sequence recognizable by protein kinases. Upon phosphorylation, the engineered fluorescent protein undergoes a change in a fluorescent property. Such substrates are useful in detecting and measuring protein kinase activity in a sample of a cell, upon transfection and expression of the substrate. Preferably, the kinase recognition site is placed within about 20 amino acids of a terminus of the engineered fluorescent protein. The kinase recognition site also can be placed in a loop domain of the protein. (See, e.g. Figure 1B.) Methods for making fluorescent substrates for protein kinases are described in United States patent application 08/680,877, filed July 16, 1996.

A protease recognition site also can be introduced into a loop domain. Upon cleavage, fluorescent property changes in a measurable fashion.

The invention also includes a method of identifying a test chemical. Typically, the method includes contacting a test chemical a sample containing a biological entity labeled with a functional, engineered fluorescent protein or a polynucleotide encoding said functional, engineered fluorescent protein. By monitoring fluorescence (i.e. a fluorescent property) from the sample containing the functional engineered fluorescent protein it can be determined whether a test chemical is active. Controls can be included to insure the specificity of the signal. Such controls include measurements of a fluorescent property in the absence of the test chemical, in the presence of a chemical with an expected activity (e.g., a known modulator) or engineered controls (e.g., absence of engineered fluorescent protein, absence of engineered fluorescent protein polynucleotide or the absence of operably linkage of the engineered fluorescent protein).

The fluorescence in the presence of a test chemical can be greater or less than in the absence of said test chemical. For instance if the engineered fluorescent protein is used a reporter of gene expression, the test chemical may up or down regulate gene expression. For such types of screening, the polynucleotide encoding the functional, engineered fluorescent protein is operatively linked to a genomic polynucleotide or a re. Alternatively, the functional, engineered fluorescent protein is fused to second functional protein. This embodiment can be used to track localization of the second protein or to track protein-protein interactions using energy transfer.

20

25

30

15

5

10

### IV. PROCEDURES

Fluorescence in a sample is measured using a fluorimeter. In general, excitation radiation from an excitation source having a first wavelength, passes through excitation optics. The excitation optics cause the excitation radiation to excite the sample. In response, fluorescent proteins in the sample emit radiation which has a wavelength that is different from the excitation wavelength. Collection optics then collect the emission from the sample. The device can include a temperature controller to maintain the sample at a specific temperature while it is being scanned. According to one embodiment, a multi-axis translation stage moves a microtiter plate holding a plurality of samples in order to position different wells to be exposed. The multi-axis translation stage, temperature controller, auto-focusing feature, and electronics associated with imaging and data collection can be managed by an appropriately programmed digital computer. The computer also can

10

20

25

30

PCT/US97/14593 WO 98/06737 40

transform the data collected during the assay into another format for presentation. This process can be miniaturized and automated to enable screening many thousands of compounds.

Methods of performing assays on fluorescent materials are well known in the art and are described in, e.g., Lakowicz, J.R., Principles of Fluorescence Spectroscopy, New York: Plenum Press (1983); Herman, B., Resonance energy transfer microscopy, in: Fluorescence Microscopy of Living Cells in Culture, Part B, Methods in Cell Biology, vol. 30, ed. Taylor, D.L. & Wang, Y.-L., San Diego: Academic Press (1989), pp. 219-243; Turro, N.J., Modern Molecular Photochemistry, Menlo Park: Benjamin/Cummings Publishing Col, Inc. (1978), pp. 296-361.

The following examples are provided by way of illustration, not by way of limitation.

**EXAMPLES** 15

> As a step in understanding the properties of GFP, and to aid in the tailoring of GFPs with altered characteristics, we have determined the three dimensional structure at 1.9 Å resolution of the S65T mutant (R. Heim et al. Nature 373:664-665 (1995)) of A. victoria GFP. This mutant also contains the ubiquitous Q80R substitution, which accidentally occurred in the early distribution of the GFP cDNA and is not known to have any effect on the protein properties (M. Chalfie et al. Science 263:802-805 (1994)).

Histidine-tagged S65T GFP (R. Heim et al. Nature 373:664-665 (1995)) was overexpressed in JM109/pRSET<sub>B</sub> in 4 l YT broth plus ampicillin at 37□, 450 rpm and 5 1/min air flow. The temperature was reduced to 25  $\square$  at  $A_{595} = 0.3$ , followed by induction with 1mM isopropylthiogalactoside for 5h. Cell paste was stored at -80□ overnight, then was resuspended in 50 mM HEPES pH 7.9, 0.3 M NaCl, 5 mM 2-mercaptoethanol, 0.1 mM phenylmethyl-sulfonylfluoride (PMSF), passed once through a French press at 10,000 psi, then centrifuged at 20 K rpm for 45 min. The supernatant was applied to a Ni-NTA-agarose column (Oiagen), followed by a wash with 20 mM imidazole, then eluted with 100 mM imidazole. Green fractions were pooled and subjected to chymotryptic (Sigma) proteolysis (1:50 w/w) for 22 h at RT. After addition of 0.5 mM PMSF, the digest was reapplied to the

Ni column. N-terminal sequencing verified the presence of the correct N-terminal methionine. After dialysis against 20 mM HEPES, pH 7.5 and concentration to  $A_{490} = 20$ , rod-shaped crystals were obtained at RT in hanging drops containing 5  $\Box$ l protein and 5  $\Box$ l well solution, 22-26% PEG 4000 (Serva), 50 mM HEPES pH 8.0-8.5, 50 mM MgCl<sub>2</sub> and 10 mM 2-mercapto-ethanol within 5 days. Crystals were 0.05 mm across and up to 1.0 mm long. The space group is  $P2_12_12_1$  with a = 51.8, b = 62.8, c = 70.7 Å, Z=4. Two crystal forms of wild-type GFP, unrelated to the present form, have been described by M. A. Perrozo, K. B. Ward, R. B. Thompson, & W. W. Ward. J. Biol. Chem. 203, 7713-7716 (1988).

5

10

15

20

25

30

The structure of GFP was determined by multiple isomorphous replacement and anomalous scattering (Table E), solvent flattening, phase combination and crystallographic refinement. The most remarkable feature of the fold of GFP is an eleven stranded B-barrel wrapped around a single central helix (Fig. 1A and 1B), where each strand consists of approximately 9-13 residues. The barrel forms a nearly perfect cylinder 42 Å long and 24 Å in diameter. The N-terminal half of the polypeptide comprises three antiparallel strands, the central helix, and then 3 more anti-parallel strands, the latter of which (residues 118-123) is parallel to the N-terminal strand (residues 11-23). The polypeptide backbone then crosses the "bottom" of the molecule to form the second half of the barrel in a five-strand Greek Key motif. The top end of the cylinder is capped by three short, distorted helical segments, while one short, very distorted helical segment caps the bottom of the cylinder. The main-chain hydrogen bonding lacing the surface of the cylinder very likely accounts for the unusual stability of the protein towards denaturation and proteolysis. There are no large segments of the polypeptide that could be excised while preserving the intactness of the shell around the chromophore. Thus it would seem difficult to re-engineer GFP to reduce its molecular weight (J. Dopf & T.M. Horiagon Gene 173:39-43 (1996)) by a large percentage.

The p-hydroxybenzylideneimidazolidinone chromophore (C. W. Cody et al. Biochemistry 32:1212-1218 (1993)) is completely protected from bulk solvent and centrally located in the molecule. The total and presumably rigid encapsulation is probably responsible for the small Stokes' shift (i.e. wavelength difference between excitation and emission maxima), high quantum yield of fluorescence, inability of O<sub>2</sub> to quench the excited state (B.D. Nageswara Rao et al. Biophys. J. 32:630-632 (1980)), and resistance of the

chromophore to titration of the external pH (W. W. Ward. Bioluminescence and Chemiluminescence (M. A. DeLuca and W. D. McElroy, eds) Academic Press pp. 235-242 (1981); W. W. Ward & S. H. Bokman. Biochemistry 21:4535-4540 (1982); W. W. Ward et al. Photochem. Photobiol. 35:803-808 (1982)). It also allows one to rationalize why 5 fluorophore formation should be a spontaneous intramolecular process (R. Heim et al. Proc. Natl. Acad. Sci. USA 91:12501-12504 (1994)), as it is difficult to imagine how an enzyme could gain access to the substrate. The plane of the chromophore is roughly perpendicular (600) to the symmetry axis of the surrounding barrel. One side of the chromophore faces a surprisingly large cavity, that occupies a volume of approximately 135 Å<sup>3</sup> (B. Lee & F. M. 10 Richards. J. Mol. Biol. 55:379-400 (1971)). The atomic radii were those of Lee & Richards, calculated using the program MS with a probe radius of 1.4 Å. (M. L. Connolly, Science 221:709-713 (1983)). The cavity does not open out to bulk solvent. Four water molecules are located in the cavity, forming a chain of hydrogen bonds linking the buried side chains of Glu<sup>222</sup> and Gln<sup>69</sup>. Unless occupied, such a large cavity would be expected to de-stabilize the protein by several kcal/mol (S. J. Hubbard et al., Protein Engineering 7:613-626 (1994); 15 A. E. Eriksson et al. Science 255:178-183 (1992)). Part of the volume of the cavity might be the consequence of the compaction resulting from cyclization and dehydration reactions. The cavity might also temporarily accommodate the oxidant, most likely O<sub>2</sub> (A. B. Cubitt et al. Trends Biochem. Sci. 20:448-455 (1995); R. Heim et al. Proc. Natl. Acad. Sci. USA 91:12501-12504 (1994); S. Inouye & F.I. Tsuji. FEBS Lett. 351:211-214 (1994)), that 20 dehydrogenates the C-S bond of Tyr66. The chromophore, cavity, and side chains that contact the chromophore are shown in Figure 2A and a portion of the final electron density map in this vicinity in 2B.

The opposite side of the chromophore is packed against several aromatic and polar side chains. Of particular interest is the intricate network of polar interactions with the chromophore (Fig. 2C). His<sup>148</sup>, Thr<sup>203</sup> and Ser<sup>205</sup> form hydrogen bonds with the phenolic hydroxyl; Arg<sup>96</sup> and Gln<sup>94</sup> interact with the carbonyl of the imidazolidinone ring and Glu<sup>222</sup> forms a hydrogen bond with the side chain of Thr<sup>65</sup>. Additional polar interactions, such as hydrogen bonds to Arg<sup>96</sup> from the carbonyl of Thr<sup>62</sup>, and the side-chain carbonyl of Gln<sup>183</sup>, presumably stabilize the buried Arg<sup>96</sup> in its protonated form. In turn, this buried charge suggests that a partial negative charge resides on the carbonyl oxygen of the imidazolidinone ring of the deprotonated fluorophore, as has previously been suggested (W.

25

30

W. Ward. Bioluminescence and Chemiluminescence (M. A. DeLuca and W. D. McElroy, eds) Academic Press pp. 235-242 (1981); W. W. Ward & S. H. Bokman. Biochemistry 21:4535-4540 (1982); W. W. Ward et al. Photochem. Photobiol. 35:803-808 (1982)). Arg<sup>96</sup> is likely to be essential for the formation of the fluorophore, and may help catalyze the initial ring closure. Finally, Tyr<sup>145</sup> shows a typical stabilizing edge-face interaction with the benzyl ring. Trp<sup>57</sup>, the only tryptophan in GFP, is located 13 Å to 15 Å from the chromophore and the long axes of the two ring systems are nearly parallel. This indicates that efficient energy transfer to the latter should occur, and explains why no separate tryptophan emission is observable (D.C. Prasher et al. Gene 111:229-233 (1992). The two cysteines in GFP, Cys<sup>48</sup> and Cys<sup>70</sup>, are 24 Å apart, too distant to form a disulfide bridge. Cys<sup>70</sup> is buried, but Cys<sup>48</sup> should be relatively accessible to sulfhydryl-specific reagents. Such a reagent, 5,5'-dithiobis(2-nitrobenzoic acid), is reported to label GFP and quench its fluorescence (S. Inouye & F.I. Tsuji FEBS Lett. 351:211-214 (1994)). This effect was attributed to the necessity for a free sulfhydryl, but could also reflect specific quenching by the 5-thio-2-nitrobenzoate moiety that would be attached to Cys<sup>48</sup>.

5

10

15

20

25

30

Although the electron density map is for the most part consistent with the proposed structure of the chromophore (D.C. Prasher et al. Gene 111:229-233 (1992); C. W. Cody et al. Biochemistry 32:1212-1218 (1993)) in the cis [Z-] configuration, with no evidence for any substantial fraction of the opposite isomer around the chromophore double bond, difference features are found at >4  $\square$  in the final  $(F_0-F_c)$  electron density map that can be interpreted to represent either the intact, uncyclized polypeptide or a carbinolamine (inset to Fig. 2C). This suggests that a significant fraction, perhaps as much as 30% of the molecules in the crystal, have failed to undergo the final dehydration reaction. Confirmation of incomplete dehydration comes from electrospray mass spectrometry, which consistently shows that the average masses of both wild-type and S65T GFP (31,086±4 and 31,099.5±4 Da, respectively) are 6-7 Da higher than predicted (31,079 and 31,093 Da, respectively) for the fully matured proteins. Such a discrepancy could be explained by a 30-35% mole fraction of apoprotein or carbinolamine with 18 or 20 Da higher molecular weight The natural abundance of <sup>13</sup>C and <sup>2</sup>H and the finite resolution of the Hewlett-Packard 5989B electrospray mass spectrometer used to make these measurements do not permit the individual peaks to be resolved, but instead yields an average mass peak with a full width at half maximum of approximately 15 Da. The molecular weights shown include the His-tag,

which has the sequence MRGSHHHHHHH GMASMTGGQQM GRDLYDDDDK DPPAEF (SEQ ID NO:5). Mutants of GFP that increase the efficiency of fluorophore maturation might yield somewhat brighter preparations. In a model for the apoprotein, the Thr<sup>65</sup>-Tyr<sup>66</sup> peptide bond is approximately in the □-helical conformation, while the peptide of Tyr<sup>66</sup>-Gly<sup>67</sup> appears to be tipped almost perpendicular to the helix axis by its interaction with Arg<sup>96</sup>. This further supports the speculation that Arg<sup>96</sup> is important in generating the conformation required for cyclization, and possibly also for promoting the attack of Gly<sup>67</sup> on the carbonyl carbon of Thr<sup>65</sup> (A. B. Cubitt et al. *Trends Biochem. Sci.* 20:448-455 (1995)).

5

10

15

20

25

30

The results of previous random mutagenesis have implicated several amino acid side chains to have substantial effects on the spectra and the atomic model confirms that these residues are close to the chromophore. The mutations T203I and E222G have profound but opposite consequences on the absorption spectrum (T. Ehrig et al. FEBS Letters 367:163-166 (1995)). T203I (with wild-type Ser<sup>65</sup>) lacks the 475 nm absorbance peak usually attributed to the anionic chromophore and shows only the 395 nm peak thought to reflect the neutral chromophore (R. Heim et al. Proc. Natl. Acad. Sci. USA 91:12501-12504 (1994); T. Ehrig et al. FEBS Letters 367:163-166 (1995)). Indeed, Thr<sup>203</sup> is hydrogen-bonded to the phenolic oxygen of the chromophore, so replacement by Ile should hinder ionization of the phenolic oxygen. Mutation of Glu<sup>222</sup> to Gly (T. Ehrig et al. FEBS Letters 367:163-166 (1995)) has much the same spectroscopic effect as replacing Ser<sup>65</sup> by Gly, Ala, Cys, Val, or Thr, namely to suppress the 395 nm peak in favor of a peak at 470-490 nm (R. Heim et al. Nature 373:664-665 (1995); S. Delagrave et al. Bio/Technology 13:151-154 (1995)). Indeed Glu<sup>222</sup> and the remnant of Thr<sup>65</sup> are hydrogen-bonded to each other in the present structure, probably with the uncharged carboxyl of Glu<sup>222</sup> acting as donor to the side chain oxygen of Thr65. Mutations E222G, S65G, S65A, and S65V would all suppress such H-bonding. To explain why only wild-type protein has both excitation peaks, Ser<sup>65</sup>, unlike Thr<sup>65</sup>, may adopt a conformation in which its hydroxyl donates a hydrogen bond to and stabilizes Glu<sup>222</sup> as an anion, whose charge then inhibits ionization of the chromophore. The structure also explains why some mutations seem neutral. For example, Gln<sup>80</sup> is a surface residue far removed from the chromophore, which explains why its accidental and ubiquitous mutation to Arg seems to have no obvious intramolecular spectroscopic effect (M. Chalfie et al. Science 263:802-805 (1994)).

The development of GFP mutants with red-shifted excitation and emission

10

15

20

25

30

maxima is an interesting challenge in protein engineering (A. B. Cubitt et al. Trends Biochem. Sci. 20:448-455 (1995); R. Heim et al. Nature 373:664-665 (1995); S. Delagrave et al. Bio/Technology 13:151-154 (1995)). Such mutants would also be valuable for avoidance of cellular autofluorescence at short wavelengths, for simultaneous multicolor reporting of the activity of two or more cellular processes, and for exploitation of fluorescence resonance energy transfer as a signal of protein-protein interaction (R. Heim & R.Y. Tsien. Current Biol. 6:178-182 (1996)). Extensive attempts using random mutagenesis have shifted the emission maximum by at most 6 nm to longer wavelengths, to 514 nm (R. Heim & R.Y. Tsien. Current Biol. 6:178-182 (1996)); previously described "red-shifted" mutants merely suppressed the 395 nm excitation peak in favor of the 475 nm peak without any significant reddening of the 505 nm emission (S. Delagrave et al. Bio/Technology 13:151-154 (1995)). Because Thr<sup>203</sup> is revealed to be adjacent to the phenolic end of the chromophore, we mutated it to polar aromatic residues such as His, Tyr, and Trp in the hope that the additional polarizability of their 

systems would lower the energy of the excited state of the adjacent chromophore. All three substitutions did indeed shift the emission peak to greater than 520 nm (Table F). A particularly attractive mutation was T203Y/S65G/V68L/S72A, with excitation and emission peaks at 513 and 527 nm respectively. These wavelengths are sufficiently different from previous GFP mutants to be readily distinguishable by appropriate filter sets on a fluorescence microscope. The extinction coefficient, 36,500 M<sup>-1</sup>cm<sup>-1</sup>, and quantum yield, 0.63, are almost as high as those of S65T (R. Heim et al. Nature 373:664-665 (1995)).

Comparison of Aequorea GFP with other protein pigments is instructive.

Unfortunately, its closest characterized homolog, the GFP from the sea pansy Renilla reniformis (O. Shimomura and F.H. Johnson J. Cell. Comp. Physiol. 59:223 (1962); J. G. Morin and J. W. Hastings, J. Cell. Physiol. 77:313 (1971); H. Morise et al. Biochemistry 13:2656 (1974); W. W. Ward Photochem. Photobiol. Reviews (Smith, K. C. ed.) 4:1 (1979); W. W. Ward. Bioluminescence and Chemiluminescence (M. A. DeLuca and W. D. McElroy, eds) Academic Press pp. 235-242 (1981); W. W. Ward & S. H. Bokman Biochemistry 21:4535-4540 (1982); W. W. Ward et al. Photochem. Photobiol. 35:803-808 (1982)), has not been sequenced or cloned, though its chromophore is derived from the same FSYG sequence as in wild-type Aequorea GFP (R. M. San Pietro et al. Photochem. Photobiol. 57:63S (1993)). The closest analog for which a three dimensional structure is

available is the photoactive yellow protein (PYP, G. E. O. Borgstahl et al. Biochemistry 34:6278-6287 (1995)), a 14-kDa photoreceptor from halophilic bacteria. PYP in its native dark state absorbs maximally at 446 nm and transduces light with a quantum yield of 0.64. rather closely matching wild-type GFP's long wavelength absorbance maximum near 475 nm and fluorescence quantum yield of 0.72-0.85. The fundamental chromophore in both proteins is an anionic p-hydroxycinnamyl group, which is covalently attached to the protein via a thioester linkage in PYP and a heterocyclic iminolactam in GFP. Both proteins stabilize the negative charge on the chromophore with the help of buried cationic arginine and neutral glutamic acid groups, Arg52 and Glu46 in PYP and Arg96 and Glu222 in GFP, though in PYP the residues are close to the oxyphenyl ring whereas in GFP they are nearer the carbonyl end of the chromophore. However, PYP has an overall \(\sigma / \sigma \) fold with appropriate flexibility and signal transduction domains to enable it to mediate the cellular phototactic response, whereas GFP is a much more regular and rigid □-barrel to minimize parasitic dissipation of the excited state energy as thermal or conformational motions. GFP is an elegant example of how a visually appealing and extremely useful function, efficient fluorescence, can be spontaneously generated from a cohesive and economical protein structure.

### A. Summary Of GFP Structure Determination

5

10

15

20

25

30

Data were collected at room temperature in house using either Molecular Structure Corp. R-axis II or San Diego Multiwire Systems (SDMS) detectors (Cu Ki) and later at beamline X4A at the Brookhaven National Laboratory at the selenium absorption edge (iii = 0.979 Å) using image plates. Data were evaluated using the HKL package (Z. Otwinowski, in *Proceedings of the CCP4 Study Weekend: Data Collection and Processing*, L. Sawyer, N. Issacs, S. Bailey, Eds. (Science and Engineering Research Council (SERC), Daresbury Laboratory, Warrington, UK, (1991)), pp 56-62; W. Minor, XDISPLAYF (Purdue University, West Lafayette, IN, 1993)) or the SDMS software (A. J. Howard et al. *Meth. Enzymol.* 114:452-471 (1985)). Each data set was collected from a single crystal. Heavy atom soaks were 2 mM in mother liquor for 2 days. Initial electron density maps were based on three heavy atom derivatives using in-house data, then later were replaced with the synchrotron data. The EMTS difference Patterson map was solved by inspection, then used to calculate difference Fourier maps of the other derivatives. Lack of closure

10

15

20

25

in the Protein Data Bank (access code 1EMA).

refinement of the heavy atom parameters was performed using the Protein package (W. Steigemann, in Ph.D. Thesis (Technical University, Munich, 1974)). The MIR maps were much poorer than the overall figure of merit would suggest, and it was clear that the EMTS isomorphous differences dominated the phasing. The enhanced anomalous occupancy for the synchrotron data provided a partial solution to the problem. Note that the phasing power was reduced for the synchrotron data, but the figure of merit was unchanged. All experimental electron density maps were improved by solvent flattening using the program DM of the CCP4 (CCP4: A Suite of Programs for Protein Crystallography (SERC Daresbury Laboratory, Warrington WA4 4AD UK, 1979)) package assuming a solvent content of 38%. Phase combination was performed with PHASCO2 of the Protein package using a weight of 1.0 on the atomic model. Heavy atom parameters were subsequently improved by refinement against combined phases. Model building proceeded with FRODO and O (T. A. Jones et al. Acta. Crystallogr. Sect. A 47:110 (1991); T. A. Jones, in Computational Crystallography D. Sayre, Ed. (Oxford University Press, Oxford, 1982) pp. 303-317) and crystallographic refinement was performed with the TNT package (D. E. Tronrud et al. Acta Cryst. A 43:489-503 (1987)). Bond lengths and angles for the chromophore were estimated using CHEM3D (Cambridge Scientific Computing). Final refinement and model building was performed against the X4A selenomethione data set, using (2F<sub>o</sub>-F<sub>c</sub>) electron density maps. The data beyond 1.9 Å resolution have not been used at this stage. The final model contains residues 2-229 as the terminal residues are not visible in the electron density map, and the side chains of several disordered surface residues have been omitted. Density is weak for residues 156-158 and coordinates for these residues are unreliable. This disordering is consistent with previous analyses showing that residues 1 and 233-238 are dispensible but that further truncations may prevent fluorescence (J. Dopf & T.M. Horiagon. Gene 173:39-43 (1996)). The atomic model has been deposited

Table E

# **Diffraction Data Statistics**

Crystal	Resoluti on (Å)	Total obs	Unique obs	Compl. (%) <sup>3</sup>	Compl. (shell) <sup>b</sup>	Rmerge (%)°	Riso (%) <sup>d</sup>
R-axix II							
Native	2.0	51907	13582	80	69	4.1	5.8
EMTS*	2.6	17727	6787	87	87	5.7	20.6
SeMet	2.3	44975	10292	92	88	10.2	9.3
Multiwire							
HGI4-Se	3.0	15380	4332	84	79	7.2	28.8
X4a					•		
SeMet	1.8	126078	19503	80	55	9.3	9.4
EMTS	2.3	57812	9204	82	66	7.2	26.3

# **Phasing Statistics**

Derivative	Resolution (Å)	Number of sites	Phasing power <sup>f</sup>	Phasing Power(shell)	FOM <sup>8</sup>	FOM (shell)
In House						
<b>EMTS</b>	3.0	2	2.08	2.08	0.77	.072
SeMet	3.0	4	1.66	1.28	•	-
HGI4-Se	3.0	9	1.77	1.90	-	-
X4a						
EMTS	3.0	2	1.36	1.26	0.77	.072
SeMet	3.0	4	1.31	1.08	• .	

# **Atomic Model Statistics**

	Protein atoms	1790
5	Solvent atoms 9-	1
	Resol. range (Å)	20-1.9
	Number of reflections (F > 0) 17	7676
	Completeness	84.
	R. factor <sup>(h)</sup>	0.175
10	Mean B-value (Ų)	24.1
	Deviations from ideality	
	Bond lengths (Å)	0.014
	Bond angles (□)	1.9
	Restrained B-values (Å <sup>2</sup> )	4.3
15	Ramachandran outliers	0

Notes:

WO 98/06737 PCT/US97/14593

- (a) Completeness is the ratio of observed reflections to theoretically possible expressed as a percentage.
- (b) Shell indicates the highest resolution shell, typically 0.1-0.4 Å wide.
- (c) Rmerge =  $\Box$  |I <I>| /  $\Box$  I, where <I> is the mean of individual observations of intensities I.
- (d) Riso =  $\square |I_{DER} I_{NAT}| / \square I_{NAT}$

5

- (e) Derivatives were EMTS=ethymercurithiosalicylate (residues modified Cys<sup>48</sup> and Cys<sup>70</sup>), SeMet=selenomethionine substituted protein (Met<sup>1</sup> and Met<sup>233</sup> could not be located); HgI<sub>4</sub>-SeMet = double derivative HgI<sub>4</sub> on SeMet background.
- 10 (f) Phasing power =  $\langle F_H \rangle / \langle E \rangle$  where  $\langle F_H \rangle = r.m.s.$  heavy atom scattering and  $\langle E \rangle = lack$  of closure.
  - (g) FOM, mean figure of merit
  - (h) Standard crystallographic R-factor,  $R = \Box ||F_{obs}| |F_{catc}|| / \Box |F_{obs}|$

# B. Spectral properties of Thr<sup>203</sup> ("T203") mutants compared to S65T The mutations F64L, V68L and S72A improve the folding of GFP at 37 (B. P. Cormack et al. Gene 173:33 (1996)) but do not significantly shift the emission spectra.

**TABLE F** 

Clone	Mutations	Excitation max.(nm)	Extinction coefficient (10 <sup>3</sup> M <sup>-1</sup> cm <sup>-1</sup> )	Emission max.(nm)
S65T	S65T	489	39.2	511
5B	T203H/S65T	512	19.4	524
6C	T203Y/S65T	513	14.5	525
10B	T203Y/F64L/S65G/S72A	513	30.8	525
10C	T203Y/F65G/V68L/S72A	513	36.5	527
11	T203W/S65G/S72A	502	33.0	512

WO 98/0673	<b>7</b> 51			PCT/US97/14593
12H	T203Y/S65G/S72A	513	36.5	527
20A	T203Y/S65G/V68L/Q69K/S72A	515	46.0	527

5

The present invention provides novel long wavelength engineered fluorescent proteins. While specific examples have been provided, the above description is illustrative and not restrictive. Many variations of the invention will become apparent to those skilled in the art upon review of this specification. The scope of the invention should, therefore, be determined not with reference to the above description, but instead should be determined with reference to the appended claims along with their full scope of equivalents.

All publications and patent documents cited in this application are 10 incorporated by reference in their entirety for all purposes to the same extent as if each individual publication or patent document were so individually denoted.

#### SEQUENCE LISTING

- (1) GENERAL INFORMATION:
  - (i) APPLICANT: The Regents of the University of California et al.
  - (ii) TITLE OF INVENTION: LONG WAVELENGTH MUTANT FLUORESCENT PROTEINS
  - (iii) NUMBER OF SEQUENCES: 4
  - (iv) CORRESPONDENCE ADDRESS:
    - (A) ADDRESSEE: Fish & Richardson P.C.
    - (B) STREET: 4225 Executive Square, Suite 1400
    - (C) CITY: La Jolla
    - (D) STATE: CA
    - (E) COUNTRY: USA
    - (F) ZIP: 92037
  - (v) COMPUTER READABLE FORM:
    - (A) MEDIUM TYPE: Floppy disk
    - (B) COMPUTER: IBM PC compatible
    - (C) OPERATING SYSTEM: PC-DOS/MS-DOS
    - (D) SOFTWARE: PatentIn Release #1.0, Version #1.30
  - (vii) PRIOR APPLICATION DATA:
    - (A) APPLICATION NUMBER: US 60/024,050
    - (B) FILING DATE: 16-AUG-1996
    - (C) CLASSIFICATION:
  - (vii) PRIOR APPLICATION DATA:
    - (A) APPLICATION NUMBER: US 08/706,408
    - (B) FILING DATE: 30-AUG-1996
    - (C) CLASSIFICATION:
  - (viii) ATTORNEY/AGENT INFORMATION:

    - (A) NAME: Haile, Lisa A.
      (B) REGISTRATION NUMBER: 38,347
    - (C) REFERENCE/DOCKET NUMBER: 07257/056WO1
    - (ix) TELECOMMUNICATION INFORMATION:
      - (A) TELEPHONE: 619/678-5070
      - (B) TELEFAX: 619/678-5099
- (2) INFORMATION FOR SEQ ID NO:1:
  - (i) SEQUENCE CHARACTERISTICS:
    - (A) LENGTH: 716 base pairs
    - (B) TYPE: nucleic acid
    - (C) STRANDEDNESS: single (D) TOPOLOGY: linear
  - (ix) FEATURE:

    - (A) NAME/KEY: CDS (B) LOCATION: 1..714

	(xi)	SEC	UENC	E DE	SCRI	PTIC	N: 9	EQ 1	D NC	):1:						
ATG Met 1	AGT Ser	AAA Lys	GGA Gly	GAA Glu 5	GAA Glu	CTT Leu	TTC Phe	ACT Thr	GCA Ala 10	GTT Val	GTC Val	CCA Pro	ATT Ile	CTT Leu 15	GTT Val	48
GAA Glu	TTA Leu	GAT Asp	GGT Gly 20	GAT Asp	GTT Val	AAT Asn	GGG Gly	CAC His 25	AAA Lys	TTT Phe	TCT Ser	GTC Val	AGT Ser 30	GGA Gly	GAG Glu	96
GGT Gly	GAA Glu	GGT Gly 35	GAT Asp	GTA Val	ACA Thr	TAC Tyr	GGA Gly 40	AAA Lys	CTT Leu	ACC Thr	CTT Leu	AAA Lys 45	TTT Phe	ATT Ile	TGC Cys	144
ACT Thr	ACT Thr 50	GGA Gly	AAA Lys	CTA Leu	CCT Pro	GTT Val 55	CCA Pro	TGG Trp	CCA Pro	ACA Thr	CTT Leu 60	GTC Val	ACT Thr	ACT Thr	TTC Phe	192
TCT Ser 65	TAT Tyr	GGT Gly	GTT Val	CAA Gln	TGC Cys 70	TTT Phe	TCA Ser	AGA Arg	TAC Tyr	CCA Pro 75	GAT Asp	CAT His	ATG Met	AAA Lys	CGG Arg 80	240
CAT His	GAC Asp	TTT Phe	TTC Phe	AAG Lys 85	AGT Ser	GCC Ala	ATG Met	CCC Pro	GAA Glu 90	GGT Gly	TAT Tyr	GTA Val	CAG Gln	CAA Gln 95	AGA Arg	288
ACT Thr	ATA Ile	TTT Phe	TTC Phe 100	AAA Lys	GAT Asp	GAC Asp	GGG Gly	AAC Asn 105	TAC Tyr	AAG Lys	ACA Thr	CGT Arg	GCT Ala 110	GAA Glu	GTC Val	336
AAG Lys	TTT Phe	GAA Glu 115	GGT Gly	GAT Asp	ACC Thr	CTT Leu	GTT Val 120	AAT Asn	AGA Arg	ATC Ile	GAG Glu	TTA Leu 125	AAA Lys	GGT Gly	ATT Ile	384
GAT Asp	TTT Phe 130	AAA Lys	GAA Glu	GAT Asp	GGA Gly	AAC Asn 135	ATT Ile	CTT Leu	GGA Gly	CAT His	AAA Lys 140	TTG Leu	GAA Glu	TAC Tyr	AAC Asn	432
TAT Tyr 145	AAC Asn	TCA Ser	CAC His	AAT Asn	GTA Val 150	TAC Tyr	ATC Ile	ATG Met	GCA Ala	GAC Asp 155	AAA Lys	CAA Gln	AAG Lys	AAT Asn	GGA Gly 160	480
ATC Ile	AAA Lys	GTT Val	AAC Asn	TTC Phe 165	AAA Lys	ATT Ile	AGA Arg	CAC His	AAC Asn 170	ATT Ile	GAA Glu	GAT Asp	GGA Gly	AGC Ser 175	GTT Val	528
CAA Gln	CTA Leu	GCA Ala	GAC Asp 180	TAT Tyr	TAT Tyr	CAA Gln	CAA Gln	AAT Asn 185	Thr	CCA Pro	ATT	CTC Leu	GAT Asp 190	GGC Gly	CCT Pro	576
GTC Val	CTT Leu	TTA Leu 195	Pro	GAC Asp	AAC Asn	CAT His	TAC Tyr 200	Leu	TCC Ser	ACA Thr	CAA Gln	TCT Ser 205	GCC Ala	CTT Leu	TCG Ser	624
AAA Lys	GAT Asp 210	Pro	AAC Asn	GAA Glu	AAG Lys	AGA Arg 215	Asp	CAC	ATG Met	GTC Val	CTT Leu 220	Leu	GAG Glu	TTT Phe	GTA Val	672
ACA Thr 225	Ala	GCT Ala	GGG Gly	ATT Ile	ACA Thr 230	CAT His	GGC Gly	ATG Met	GAT Asp	GAA Glu 235	Leu	TAC	AAA Lys			714
TA																716

## (2) INFORMATION FOR SEQ ID NO:2:

- (i) SEQUENCE CHARACTERISTICS:
  - (A) LENGTH: 238 amino acids
    (B) TYPE: amino acid

- (D) TOPOLOGY: linear
- (ii) MOLECULE TYPE: protein
- (xi) SEQUENCE DESCRIPTION: SEQ ID NO:2:

Met Ser Lys Gly Glu Glu Leu Phe Thr Ala Val Val Pro Ile Leu Val

Glu Leu Asp Gly Asp Val Asn Gly His Lys Phe Ser Val Ser Gly Glu 20 25 30

Gly Glu Gly Asp Val Thr Tyr Gly Lys Leu Thr Leu Lys Phe Ile Cys 35 40 45

Thr Thr Gly Lys Leu Pro Val Pro Trp Pro Thr Leu Val Thr Thr Phe 50 60

Ser Tyr Gly Val Gln Cys Phe Ser Arg Tyr Pro Asp His Met Lys Arg 65 70 75 80

His Asp Phe Phe Lys Ser Ala Met Pro Glu Gly Tyr Val Gln Gln Arg 85 90 95

Thr Ile Phe Phe Lys Asp Asp Gly Asn Tyr Lys Thr Arg Ala Glu Val

Lys Phe Glu Gly Asp Thr Leu Val Asn Arg Ile Glu Leu Lys Gly Ile 115 120 125

Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly His Lys Leu Glu Tyr Asn 130 140

Tyr Asn Ser His Asn Val Tyr Ile Met Ala Asp Lys Gln Lys Asn Gly 145 155 160

Ile Lys Val Asn Phe Lys Ile Arg His Asn Ile Glu Asp Gly Ser Val

Gln Leu Ala Asp Tyr Tyr Gln Gln Asn Thr Pro Ile Leu Asp Gly Pro 180 185 190

Val Leu Leu Pro Asp Asn His Tyr Leu Ser Thr Gln Ser Ala Leu Ser 195 200 205

Lys Asp Pro Asn Glu Lys Arg Asp His Met Val Leu Leu Glu Phe Val 210 220

Thr Ala Ala Gly Ile Thr His Gly Met Asp Glu Leu Tyr Lys

- (2) INFORMATION FOR SEQ ID NO:3:
  - (i) SEQUENCE CHARACTERISTICS:
    - (A) LENGTH: 720 base pairs
    - (B) TYPE: nucleic acid
    - (C) STRANDEDNESS: single
    - (D) TOPOLOGY: linear
  - (ix) FEATURE:
    - (A) NAME/KEY: CDS
    - (B) LOCATION: 1..720

	(xi)	SEC	UENC	E DE	SCRI	PTIC	N: 5	SEQ I	מ סו	: 3 :						
ATG Met	GTG Val 240	AGC Ser	AAG Lys	GGC Gly	GAG Glu	GAG Glu 245	CTG Leu	TTC Phe	ACC Thr	GGG Gly	GTG Val 250	GTG Val	CCC Pro	ATC Ile	CTG Leu	48
GTC Val 255	GAG Glu	CTG Leu	GAC Asp	GGC Gly	GAC Asp 260	GTA Val	AAC Asn	GGC Gly	CAC His	AAG Lys 265	TTC Phe	AGC Ser	GTG Val	TCC Ser	GGC Gly 270	96
GAG Glu	GGC Gly	GAG Glu	GGC Gly	GAT Asp 275	GCC Ala	ACC Thr	TAC Tyr	GGC Gly	AAG Lys 280	CTG Leu	ACC Thr	CTG Leu	AAG Lys	TTC Phe 285	ATC Ile	144
TGC Cys	ACC Thr	ACC Thr	GGC Gly 290	AAG Lys	CTG Leu	CCC Pro	GTG Val	CCC Pro 295	TGG Trp	CCC Pro	ACC Thr	CTC Leu	GTG Val 300	ACC Thr	ACC Thr	192
TTC Phe	GGC Gly	TAC Tyr 305	GGC Gly	GTG Val	CAG Gln	TGC Cys	TTC Phe 310	GCC Ala	CGC Arg	TAC Tyr	CCC Pro	GAC Asp 315	CAC His	ATG Met	AAG Lys	240
CAG Gln	CAG Gln 320	GAC Asp	TTC Phe	TTC Phe	AAG Lys	TCC Ser 325	GCC Ala	ATG Met	CCC Pro	GAA Glu	GGC Gly 330	TAC Tyr	GTC Val	CAG Gln	GAG Glu	288
CGC Arg 335	ACC Thr	ATC Ile	TTC Phe	TTC Phe	AAG Lys 340	GAC Asp	GAC Asp	GGC Gly	AAC Asn	TAC Tyr 345	AAG Lys	ACC Thr	CGC Arg	GCC Ala	GAG Glu 350	336
GTG Val	AAG Lys	TTC Phe	GAG Glu	GGC Gly 355	GAC Asp	ACC Thr	CTG Leu	GTG Val	AAC Asn 360	CGC Arg	ATC Ile	GAG Glu	CTG Leu	AAG Lys 365	GGC Gly	384
ATC Ile	GAC Asp	TTC Phe	AAG Lys 370	GAC Asp	GAC Asp	GGC Gly	AAC Asn	ATC Ile 375	CTG Leu	GGG Gly	CAC His	AAG Lys	CTG Leu 380	GAG Glu	TAC Tyr	432
AAC Asn	TAC Tyr	AAC Asn 385	AGC Ser	CAC His	AAC Asn	GTC Val	TAT Tyr 390	ATC Ile	ATG Met	GCC Ala	GAC Asp	AAG Lys 395	CAG Gln	AAG Lys	AAC Asn	480
GGC Gly	ATC Ile 400	AAG Lys	GTG Val	AAC Asn	TTC Phe	AAG Lys 405	ATC Ile	CGC Arg	CAC His	AAC Asn	ATC Ile 410	GAG Glu	GAC Asp	GGC Gly	AGC Ser	528
GTG Val 415	Gln	CCC Pro	GCC Ala	GAC Asp	CAC His 420	TAC Tyr	CAG Gln	CAG Gln	AAC Asn	ACC Thr 425	Pro	ATC Ile	GGC Gly	GAC Asp	GGC Gly 430	576
CCC Pro	GTG Val	CTG Leu	CTG Leu	CCC Pro 435	GAC Asp	AAC Asn	CAC His	TAC	CTG Leu 440	AGC Ser	TAC Tyr	CAG Gln	TCC Ser	GCC Ala 445	CTG Leu	624
AGC Ser	AAA Lys	GAC Asp	CCC Pro 450	AAC Asn	GAG Glu	AAG Lys	CGC Arg	GAT Asp 455	His	ATG Met	GTC Val	CTG Leu	CTG Leu 460	GAG Glu	TTC Phe	672
GTG Val	ACC Thr	GCC Ala 465	GCC Ala	GGG Gly	ATC Ile	ACT Thr	CAC His 470	Gly	ATG Met	GAC Asp	GAG Glu	CTG Leu 475	TAC Tyr	AAG Lys	TAA *	720

## (2) INFORMATION FOR SEQ ID NO:4:

- (i) SEQUENCE CHARACTERISTICS:

  (A) LENGTH: 240 amino acids
  (B) TYPE: amino acid
  (D) TOPOLOGY: linear

WO 98/06737 PCT/US97/14593

- (ii) MOLECULE TYPE: protein
- (xi) SEQUENCE DESCRIPTION: SEQ ID NO:4:

Met Val Ser Lys Gly Glu Glu Leu Phe Thr Gly Val Val Pro Ile Leu
1 10 15

56

Val Glu Leu Asp Gly Asp Val Asn Gly His Lys Phe Ser Val Ser Gly 20 25 30

Glu Gly Glu Gly Asp Ala Thr Tyr Gly Lys Leu Thr Leu Lys Phe Ile 35 40 45

Cys Thr Thr Gly Lys Leu Pro Val Pro Trp Pro Thr Leu Val Thr Thr 50 55 60

Phe Gly Tyr Gly Val Gln Cys Phe Ala Arg Tyr Pro Asp His Met Lys 65 70 75 80

Gln Gln Asp Phe Phe Lys Ser Ala Met Pro Glu Gly Tyr Val Gln Glu 85 90 95

Arg Thr Ile Phe Phe Lys Asp Asp Gly Asn Tyr Lys Thr Arg Ala Glu 100 105 110

Val Lys Phe Glu Gly Asp Thr Leu Val Asn Arg Ile Glu Leu Lys Gly 115 120 125

Ile Asp Phe Lys Asp Asp Gly Asn Ile Leu Gly His Lys Leu Glu Tyr 130 140

Asn Tyr Asn Ser His Asn Val Tyr Ile Met Ala Asp Lys Gln Lys Asn 145 150 155 160

Gly Ile Lys Val Asn Phe Lys Ile Arg His Asn Ile Glu Asp Gly Ser 165 170 175

Val Gln Pro Ala Asp His Tyr Gln Gln Asn Thr Pro Ile Gly Asp Gly
180 185 190

Pro Val Leu Leu Pro Asp Asm His Tyr Leu Ser Tyr Gln Ser Ala Leu 195 200 205

Ser Lys Asp Pro Asn Glu Lys Arg Asp His Met Val Leu Leu Glu Phe 210 215 220

Val Thr Ala Ala Gly Ile Thr His Gly Met Asp Glu Leu Tyr Lys \*
225 230 240

## WHAT IS CLAIMED IS:

1. A nucleic acid molecule comprising a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of *Aequorea* green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least the substitution T203X, wherein X is an aromatic amino acid selected from H, Y, W or F, said functional engineered fluorescent protein having a different fluorescent property than *Aequorea* green fluorescent protein.

2. The nucleic acid molecule of claim 1 wherein the amino acid sequence further comprises a substitution at S65, wherein the substitution is selected from S65G, S65T, S65A, S65L, S65C, S65V and S65I.

3. The nucleic acid molecule of claim 1 wherein the amino acid sequence differs by no more than the substitutions S65T/T203H; S65T/T203Y; S72A/F64L/S65G/T203Y; S72A/S65G/V68L/T203Y; S65G/V68L/Q69K/S72A/T203Y; S65G/S72A/T203Y; or S65G/S72A/T203W.

 4. The nucleic acid molecule of claim 1 or 2 wherein the amino acid sequence further comprises a substitution at Y66, wherein the substitution is selected from Y66H, Y66F, and Y66W.

1 2

5. The nucleic acid molecule of claim 1 or 2 wherein the amino acid sequence further comprises a mutation from Table A.

1 2

6. The nucleic acid molecule of claim 1 or 2 wherein the amino acid sequence further comprises a folding mutation.

7.

1

WO 98/06737 PCT/US97/14593 58

The nucleic acid molecule of any of claims 1-3 wherein the

2	nucleotide sequence encoding the protein differs from the nucleotide sequence of SEQ ID
3	NO:1 by the substitution of at least one codon by a preferred mammalian codon.
1	8. The nucleic acid molecule of any of claims 1-3 encoding a fusion
2	protein wherein the fusion protein comprises a polypeptide of interest and the functional
3	engineered fluorescent protein.
1	9. An expression vector comprising expression control sequences
2	operatively linked to a nucleic acid molecule comprising a nucleotide sequence encoding a
3	functional engineered fluorescent protein whose amino acid sequence is substantially
4	identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2)
5	and which differs from SEQ ID NO:2 by at least the amino acid substitution T203X,
6	wherein X is an aromatic amino acid selected from H, Y, W or F, said functional engineered
7	fluorescent protein having a different fluorescent property than Aequorea green fluorescent
8	protein.
1	10. The expression vector of claim 9 wherein the amino acid sequence
2	further comprises a substitution at S65, wherein the substitution is selected from S65G,
3	S65T, S65A, S65L, S65C, S65V and S65I.
1	11. The expression vector of claim 9 wherein the amino acid sequence
2	differs by no more than the substitutions S65T/T203H; S65T/T203Y;
3	S72A/F64L/S65G/T203Y; S72A/S65G/V68L/T203Y; S65G/V68L/Q69K/S72A/T203Y,
4	S65G/S72A/T203Y; or S65G/S72A/T203W.
1	12. The expression vector of claim 10 or 11 wherein the amino acid
2	sequence further comprises a substitution at Y66, wherein the substitution is selected from
3	Y66H, Y66F, and Y66W.

The expression vector of claim 10 or 11 wherein the amino acid 13. 1 sequence further comprises a mutation from Table A. 2 The expression vector of claim 9 or 10 wherein the amino acid 14. 3 sequence further comprises a folding mutation. 4 The expression vector of any of claims 9-11 wherein the nucleotide 15. 1 sequence encoding the protein differs from the nucleotide sequence of SEQ ID NO:1 by the 2 substitution of at least one codon by a preferred mammalian codon. 3 The expression vector of any of claims 9-11 encoding a fusion 16. 1 protein wherein the fusion protein comprises a polypeptide of interest and the functional 2 engineered fluorescent protein. 3 A recombinant host cell comprising an expression vector that 17. 1 comprises expression control sequences operatively linked to a nucleic acid molecule 2 comprising a nucleotide sequence encoding a functional engineered fluorescent protein 3 whose amino acid sequence is substantially identical to the amino acid sequence of 4 Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 5 by at least the amino acid substitution T203X, wherein X is an aromatic amino acid selected 6 from H, Y, W or F, said functional engineered fluorescent protein having a different 7 fluorescent property than Aequorea green fluorescent protein. 8 The recombinant host cell of claim 17 wherein the amino acid 18. 1 sequence further comprises a substitution at S65, wherein the substitution is selected from 2

S65G, S65T, S65A, S65L, S65C, S65V and S65I.

3

1	19.	The recombinant host cell of claim 17 wherein the amino acid
2	sequence differs by n	o more than the substitutions S65T/T203H; S65T/T203Y;
3	S72A/F64L/S65G/T2	03Y; S72A/S65G/V68L/T203Y; S65G/V68L/Q69K/S72A/T203Y;
4	S65G/S72A/T203Y;	or S65G/S72A/T203W.
1	20.	The recombinant host cell of claim 17 or 18 wherein the amino acid
2	sequence further com	prises a substitution at Y66, wherein the substitution is selected from
3	Y66H, Y66F, and Y6	66W.
1	21.	The recombinant host cell of claim 17 or 18 wherein the amino acid
2	sequence further con	nprises a mutation from Table A.
		•
1	22.	The recombinant host cell of claim 17 or 18 wherein the amino acid
2	sequence further con	nprises a folding mutation.
1	23.	The recombinant host cell of any of claims 17-19 wherein the
2	nucleotide sequence	encoding the protein differs from the nucleotide sequence of SEQ ${\rm I\!D}$
3	NO:1 by the substitu	ation of at least one codon by a preferred mammalian codon.
1	24.	The recombinant host cell of any of claims 17-19 encoding a fusion
2	protein wherein the	fusion protein comprises a polypeptide of interest and the functional
3	engineered fluoresc	ent protein.
1	25.	The recombinant host cell of any of claims 17-19 which is a
2	prokaryotic cell.	·
	-	
1	26.	The recombinant host cell of any of claims 17-19 which is a
2	eukarvotic cell	

PCT/US97/14593 WO 98/06737 61

1	27. A functional engineered fluorescent protein whose amino acid
2	sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent
3	protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least the amino acid
4	substitution T203X, wherein X is an aromatic amino acid selected from H, Y, W or F, said
5	functional engineered fluorescent protein having a different fluorescent property than
6	Aequorea green fluorescent protein.
1	28. The protein of claim 27 wherein the amino acid sequence further
2	comprises a substitution at S65, wherein the substitution is selected from S65G, S65T,
3	S65A, S65L, S65C, S65V and S65I.
1	29. The protein of claim 27 wherein the amino acid sequence differs by
2	no more than the substitutions S65T/T203H; S65T/T203Y; S72A/F64L/S65G/T203Y;
3	S72A/S65G/V68L/T203Y; S65G/V68L/Q69K/S72A/T203Y; S65G/S72A/T203Y; or
4	S65G/S72A/T203W.
1	30. The protein of claim 27 or 28 wherein the amino acid sequence
2	further comprises a substitution at Y66, wherein the substitution is selected from Y66H,
3	Y66F, and Y66W.
1	31. The protein of claim 27 or 28 wherein the amino acid sequence
2	further comprises a folding mutation.
1	32. The protein of any of claims 27-29 which is a fusion protein wherein
2	the fusion protein comprises a polypeptide of interest and the functional engineered
3	fluorescent protein.

l	33. A fluorescently labelled antibody comprising an antibody coupled to
2	a functional engineered fluorescent protein whose amino acid sequence is substantially
3	identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2)
4	and which differs from SEQ ID NO:2 by at least the amino acid substitution T203X,
5	wherein X is an aromatic amino acid selected from H, Y, W or F, said functional engineered
6	fluorescent protein having a different fluorescent property than Aequorea green fluorescent
7	protein.
1	34. The fluorescently labelled antibody of claim 33 wherein the amino
2	acid sequence further comprises a substitution at S65, wherein the substitution is selected
3	from S65G, S65T, S65A, S65L, S65C, S65V and S65I.
1	35. The fluorescently labelled antibody of claim 33 wherein the amino
2	acid sequence differs by no more than the substitutions S65T/T203H; S65T/T203Y;
3	S72A/F64L/S65G/T203Y; S72A/S65G/V68L/T203Y; S65G/V68L/Q69K/S72A/T203Y;
4	S65G/S72A/T203Y; or S65G/S72A/T203W.
1	36. The fluorescently labelled antibody of claim 33 or 34 wherein the
2	amino acid sequence further comprises a substitution at Y66, wherein the substitution is
3	selected from Y66H, Y66F, and Y66W.
1	37. The fluorescently labelled antibody of any of claims 33-35 which is a
2	fusion protein wherein the fusion protein comprises the antibody fused to the functional
3	engineered fluorescent protein.

1	38. A nucleic acid molecule comprising a nucleotide sequence encoding
2	an antibody fused to a nucleotide sequence encoding a functional engineered fluorescent
3	protein whose amino acid sequence is substantially identical to the amino acid sequence of
4	Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2
5	by at least the amino acid substitution T203X, wherein X is an aromatic amino acid selected
6	from H, Y, W or F, said functional engineered fluorescent protein having a different
7	fluorescent property than Aequorea green fluorescent protein.
1	39. The nucleic acid molecule of claim 38 wherein the amino acid
2	sequence further comprises a substitution at S65, wherein the substitution is selected from
3	S65G, S65T, S65A, S65L, S65C, S65V and S65I.
1	40. The nucleic acid molecule of claim 38 wherein the amino acid
2	sequence differs by no more than the substitutions S65T/T203H; S65T/T203Y;
3	S72A/F64L/S65G/T203Y; S72A/S65G/V68L/T203Y; S65G/V68L/Q69K/S72A/T203Y;
4	S65G/S72A/T203Y; or S65G/S72A/T203W.
•	
1	41. The nucleic acid molecule of claim 38 or 39 wherein the amino acid
2	sequence further comprises a substitution at Y66, wherein the substitution is selected from
3	Y66H, Y66F, and Y66W.
1	42. A fluorescently labelled nucleic acid probe comprising a nucleic acid
2	probe coupled to a functional engineered fluorescent protein whose amino acid sequence is
3	substantially identical to the amino acid sequence of Aequorea green fluorescent protein
4	(SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least the amino acid substitution
5	T203X, wherein X is an aromatic amino acid selected from H, Y, W or F, said functional
6	engineered fluorescent protein having a different fluorescent property than Aequorea green
7	fluorescent protein.

WO 98/06737 64 PCT/US97/14593

l	43. The fluorescently labelled nucleic acid probe of claim 42 wherein the
2	amino acid sequence further comprises a substitution at S65, wherein the substitution is
3	selected from S65G, S65T, S65A, S65L, S65C, S65V and S65I.
1	44. The fluorescently labelled nucleic acid probe of claim 42 wherein the
2	amino acid sequence differs by no more than the substitutions S65T/T203H; S65T/T203Y;
3	S72A/F64L/S65G/T203Y; S72A/S65G/V68L/T203Y; S65G/V68L/Q69K/S72A/T203Y;
4	S65G/S72A/T203Y; or S65G/S72A/T203W.
1	45. The nucleic acid molecule of claim 42 or 43 wherein the amino acid
2	sequence further comprises a substitution at Y66, wherein the substitution is selected from
3	Y66H, Y66F, and Y66W.
4	
1	46. A nucleic acid molecule comprising a nucleotide sequence encoding
2	a functional engineered fluorescent protein whose amino acid sequence is substantially
3	identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2)
4	and which differs from SEQ ID NO:2 by at least an amino acid substitution at L42, V61,
5	T62, V68, Q69, Q94, N121, Y145, H148, V150, F165, I167, Q183, N185, L220, E222 (not
6	E222G), or V224, said functional engineered fluorescent protein having a different
7	fluorescent property than Aequorea green fluorescent protein.
1	The nucleic acid molecule of claim 46 wherein the amino acid
2	substitution is:
3	L42X, wherein X is selected from C, F, H, W and Y,
4	V61X, wherein X is selected from F, Y, H and C,
5 ·	T62X, wherein X is selected from A, V, F, S, D, N, Q, Y, H and C,
6	V68X, wherein X is selected from F, Y and H,
7	Q69X, wherein X is selected from K, R, E and G,
8	Q94X, wherein X is selected from D, E, H, K and N,

9	NIZIX, wherein X is selected from F, A, w and T,
10	Y145X, wherein X is selected from W, C, F, L, E, H, K and Q,
11	H148X, wherein X is selected from F, Y, N, K, Q and R,
12	V150X, wherein X is selected from F, Y and H,
13	F165X, wherein X is selected from H, Q, W and Y,
14	1167X, wherein X is selected from F, Y and H,
15	Q183X, wherein X is selected from H, Y, E and K,
16	N185X, wherein X is selected from D, E, H, K and Q,
17	L220X, wherein X is selected from H, N, Q and T,
18	E222X, wherein X is selected from N and Q or
19	V224X, wherein X is selected from H, N, Q, T, F, W and Y.
20	
	48. An expression vector comprising expression control sequences
1 2	operatively linked to a nucleic acid molecule of comprising a nucleotide sequence encoding
3	a functional engineered fluorescent protein whose amino acid sequence is substantially
4	identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2)
5	and which differs from SEQ ID NO:2 by at least an amino acid substitution at L42, V61,
6	T62, V68, Q69, Q94, N121, Y145, H148, V150, F165, I167, Q183, N185, L220, E222 (not
7	E222G), or V224, said functional engineered fluorescent protein having a different
8	fluorescent property than Aequorea green fluorescent protein.
1	49. The expression vector of claim 48 wherein the amino acid
2	substitution is:
3	L42X, wherein X is selected from C, F, H, W and Y,
4 .	V61X, wherein X is selected from F, Y, H and C,
5	T62X, wherein X is selected from A, V, F, S, D, N, Q, Y, H and C,
6	V68X, wherein X is selected from F, Y and H,
7	Q69X, wherein X is selected from K, R, E and G,
8	Q94X, wherein X is selected from D, E, H, K and N,

9	N121X, wherein X is selected from F, H, W and Y,
LO	Y145X, wherein X is selected from W, C, F, L, E, H, K and Q,
11	H148X, wherein X is selected from F, Y, N, K, Q and R,
12	V150X, wherein X is selected from F, Y and H,
13	F165X, wherein X is selected from H, Q, W and Y,
14	1167X, wherein X is selected from F, Y and H,
15	Q183X, wherein X is selected from H, Y, E and K,
16	N185X, wherein X is selected from D, E, H, K and Q,
17	L220X, wherein X is selected from H, N, Q and T,
18	E222X, wherein X is selected from N and Q or
19	V224X, wherein X is selected from H, N, Q, T, F, W and Y.
1	50. A recombinant host cell comprising an expression vector that
2	comprises expression control sequences operatively linked to a nucleic acid molecule
3	comprising a nucleotide sequence encoding a functional engineered fluorescent protein
4	whose amino acid sequence is substantially identical to the amino acid sequence of
5	Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2
6	by at least an amino acid substitution at L42, V61, T62, V68, Q69, Q94, N121, Y145,
7	H148, V150, F165, I167, Q183, N185, L220, E222 (not E222G), or V224, said functional
8	engineered fluorescent protein having a different fluorescent property than Aequorea green
9	fluorescent protein.
1	51. The recombinant host cell of claim 50 wherein the amino acid
2	substitution is:
3	L42X, wherein X is selected from C, F, H, W and Y,
4	V61X, wherein X is selected from F, Y, H and C,
5	T62X, wherein X is selected from A, V, F, S, D, N, Q, Y, H and C,
6	V68X, wherein X is selected from F, Y and H,
7	O69X, wherein X is selected from K, R, E and G,

8	Q94X, wherein X is selected from D, E, H, K and N,
9	N121X, wherein X is selected from F, H, W and Y,
10	Y145X, wherein X is selected from W, C, F, L, E, H, K and Q,
11	H148X, wherein X is selected from F, Y, N, K, Q and R,
12	V150X, wherein X is selected from F, Y and H,
13	F165X, wherein X is selected from H, Q, W and Y,
14	I167X, wherein X is selected from F, Y and H,
15	Q183X, wherein X is selected from H, Y, E and K,
16	N185X, wherein X is selected from D, E, H, K and Q,
17	L220X, wherein X is selected from H, N, Q and T,
18	E222X, wherein X is selected from N and Q or
19	V224X, wherein X is selected from H, N, Q, T, F, W and Y.
20	
1	52. A functional engineered fluorescent protein whose amino acid
2	sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent
3	protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid
4	substitution at L42, V61, T62, V68, Q69, Q94, N121, Y145, H148, V150, F165, I167,
5	Q183, N185, L220, E222 (E222G), or V224, said functional engineered fluorescent protein
6	having a different fluorescent property than Aequorea green fluorescent protein.
1	53. The functional engineered fluorescent protein of claim 52 wherein the
2	amino acid substitution is:
3	L42X, wherein X is selected from C, F, H, W and Y,
4	V61X, wherein X is selected from F, Y, H and C,
5 .	T62X, wherein X is selected from A, V, F, S, D, N, Q, Y, H and C,
6	V68X, wherein X is selected from F, Y and H,
7	Q69X, whèrein X is selected from K, R, E and G,
8	Q94X, wherein X is selected from D, E, H, K and N,

9	N121X, wherein X is selected from F, H, W and Y,
10	Y145X, wherein X is selected from W, C, F, L, E, H, K and Q,
11	H148X, wherein X is selected from F, Y, N, K, Q and R,
12	V150X, wherein X is selected from F, Y and H,
13	F165X, wherein X is selected from H, Q, W and Y,
14	I167X, wherein X is selected from F, Y and H,
15	Q183X, wherein X is selected from H, Y, E and K,
16	N185X, wherein X is selected from D, E, H, K and Q,
17	L220X, wherein X is selected from H, N, Q and T,
18	E222X, wherein X is selected from N and Q or
19	V224X, wherein X is selected from H, N, Q, T, F, W and Y.
1	54. A fluorescently labelled antibody comprising an antibody coupled to
2	a functional engineered fluorescent protein whose amino acid sequence is substantially
3	identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2)
4	and which differs from SEQ ID NO:2 by at least an amino acid substitution at L42, V61,
5	T62, V68, Q69, Q94, N121, Y145, H148, V150, F165, I167, Q183, N185, L220, E222 (no
6	E222G), or V224, said functional engineered fluorescent protein having a different
7	fluorescent property than Aequorea green fluorescent protein.
1	55. The antibody of claim 54 wherein the amino acid substitution is:
2	L42X, wherein X is selected from C, F, H, W and Y,
3	V61X, wherein X is selected from F, Y, H and C,
4	T62X, wherein X is selected from A, V, F, S, D, N, Q, Y, H and C,
5	V68X, wherein X is selected from F, Y and H,
6	Q69X, wherein X is selected from K, R, E and G,
7	Q94X, wherein X is selected from D, E, H, K and N,
8	N121X, wherein X is selected from F, H, W and Y,

9	Y145X, wherein X is selected from W, C, F, L, E, H, K and Q,
LO	H148X, wherein X is selected from F, Y, N, K, Q and R,
11	V150X, wherein X is selected from F, Y and H,
12	F165X, wherein X is selected from H, Q, W and Y,
13	I167X, wherein X is selected from F, Y and H,
14	Q183X, wherein X is selected from H, Y, E and K,
15	N185X, wherein X is selected from D, E, H, K and Q,
16	L220X, wherein X is selected from H, N, Q and T,
17	E222X, wherein X is selected from N and Q or
18	V224X, wherein X is selected from H, N, Q, T, F, W and Y.
1	56. A nucleic acid molecule comprising a nucleotide sequence encoding
2	an antibody fused to a nucleotide sequence encoding a functional engineered fluorescent
3	protein whose amino acid sequence is substantially identical to the amino acid sequence of
4	Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2
5	by at least an amino acid substitution at L42, V61, T62, V68, Q69, Q94, N121, Y145,
6	H148, V150, F165, I167, Q183, N185, L220, E222 (not E222G), or V224, said functional
7	engineered fluorescent protein having a different fluorescent property than Aequorea green
8	fluorescent protein.
1	57. The nucleic acid molecule of claim 56 wherein the amino acid
2	substitution is:
3	L42X, wherein X is selected from C, F, H, W and Y,
4	V61X, wherein X is selected from F, Y, H and C,
5	T62X, wherein X is selected from A, V, F, S, D, N, Q, Y, H and C,
6	V68X, wherein X is selected from F, Y and H,
7	Q69X, wherein X is selected from K, R, E and G,
8	Q94X, wherein X is selected from D, E, H, K and N,
9	N121X, wherein X is selected from F, H, W and Y,

10	Y145X, wherein X is selected from W, C, F, L, E, H, K and Q,
11	H148X, wherein X is selected from F, Y, N, K, Q and R,
12	V150X, wherein X is selected from F, Y and H,
13	F165X, wherein X is selected from H, Q, W and Y,
14	1167X, wherein X is selected from F, Y and H,
15	Q183X, wherein X is selected from H, Y, E and K,
16	N185X, wherein X is selected from D, E, H, K and Q,
17	L220X, wherein X is selected from H, N, Q and T,
18	E222X, wherein X is selected from N and Q or
19	V224X, wherein X is selected from H, N, Q, T, F, W and Y.
1	58. A fluorescently labelled nucleic acid probe comprising a nucleic acid
2	probe coupled to a functional engineered fluorescent protein whose amino acid sequence is
3	substantially identical to the amino acid sequence of Aequorea green fluorescent protein
4	(SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid substitution
5	at L42, V61, T62, V68, Q69, Q94, N121, Y145, H148, V150, F165, I167, Q183, N185,
6	L220, E222 (E222G), or V224, said functional engineered fluorescent protein having a
7	different fluorescent property than Aequorea green fluorescent protein.
1	59. The probe of claim 58 wherein the amino acid substitution is:
2	L42X, wherein X is selected from C, F, H, W and Y,
3	V61X, wherein X is selected from F, Y, H and C,
4	T62X, wherein X is selected from A, V, F, S, D, N, Q, Y, H and C,
5	V68X, wherein X is selected from F, Y and H,
6	Q69X, wherein X is selected from K, R, E and G,
7	Q94X, wherein X is selected from D, E, H, K and N,
8	N121X, wherein X is selected from F, H, W and Y,
9	Y145X, wherein X is selected from W, C, F, L, E, H, K and Q,

10	H148X, wherein X is selected from F, Y, N, K, Q and R,
11	V150X, wherein X is selected from F, Y and H,
12	F165X, wherein X is selected from H, Q, W and Y,
13	1167X, wherein X is selected from F, Y and H,
14	Q183X, wherein X is selected from H, Y, E and K,
15	N185X, wherein X is selected from D, E, H, K and Q,
16	L220X, wherein X is selected from H, N, Q and T,
17	E222X, wherein X is selected from N and Q or
18	V224X, wherein X is selected from H, N, Q, T, F, W and Y.
1	60. A method for determining whether a mixture contains a target
2	comprising:
3	contacting the mixture with a fluorescently labelled probe comprising
4	a probe and a functional engineered fluorescent protein of claim 27 or claim 52; and
5	determining whether the target has bound to the probe.
1	61. The method of any of claim 60 the target is bound to a solid matrix.
1	
2	62. A method for engineering a functional engineered fluorescent protein
3	having a fluorescent property different than Aequorea green fluorescent protein, comprising
4	substituting an amino acid that is located no more than 0.5 nm from any atom in the
5	chromophore of an Aequorea-related green fluorescent protein with another amino acid;
6	whereby the substitution alters a fluorescent property of the protein.
1	63. The method of claim 62 wherein the amino acid substitution alters the
2	electronic environment of the chromophore.
3	

1	64.	A method for engineering a functional engineered fluorescent protein
2	having a different flu	orescent property than Aequorea green fluorescent protein comprising
3	substituting amino ac	cids in a loop domain of an Aequorea-related green fluorescent protein
4	with amino acids so	as to create a consensus sequence for phosphorylation or for
5	proteolysis.	
1	65.	A method for producing fluorescence resonance energy transfer
2	comprising:	
3		providing a donor molecule comprising a functional engineered
4	fluorescent protein o	f claim 27 or claim 52;
5		providing an appropriate acceptor molecule for the fluorescent
6	protein; and	
7		bringing the donor molecule and the acceptor molecule into
8	sufficiently close co	ntact to allow fluorescence resonance energy transfer.
1	66.	A method for producing fluorescence resonance energy transfer
2	comprising:	
3		providing an acceptor molecule comprising a functional engineered
4	fluorescent protein o	of claim 27 or claim 52;
5		providing an appropriate donor molecule for the fluorescent protein;
6	and	
7		bringing the donor molecule and the acceptor molecule into
8	sufficiently close co	entact to allow fluorescence resonance energy transfer.
1	67.	The method of claim 66 wherein the donor molecule is a engineered
2	fluorescent protein	whose amino acid sequence comprises the substitution T203I and the
3	acceptor molecule i	s a nutant fluorescent protein whose amino acid sequence comprises the
4	substitution T203X	, wherein X is an aromatic amino acid selected from H, Y, W or F, said
5	functional engineer	ed fluorescent protein having a different fluorescent property than
6	Aequorea green flu	orescent protein.

A nucleic acid molecule comprising a nucleotide sequence encoding 68. a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid substitution located no more than about 0.5 nm from the chromophore of the engineered fluorescent protein, wherein the substitution alters the electronic environment of the chromophore, whereby the functional engineered fluorescent protein has a different fluorescent property than Aequorea green fluorescent protein.

operatively linked to a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid substitution located no more than about 0.5 nm from the chromophore of the engineered fluorescent protein, wherein the substitution alters the electronic environment of the chromophore, whereby the functional engineered fluorescent protein has a different fluorescent property than Aequorea green fluorescent protein.

- 70. A functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least an amino acid substitution located no more than about 0.5 nm from the chromophore of the engineered fluorescent protein, wherein the substitution alters the electronic environment of the chromophore, whereby the functional engineered fluorescent protein has a different fluorescent property than Aequorea green fluorescent protein.
- 71. A crystal of a protein comprising a fluorescent protein with an amino acid sequence substantially identical to SEQ ID NO: 2, wherein said crystal diffracts with at least a 2.0 to 3.0 angstrom resolution.

1	72. The crystal of claim 71, wherein the fluorescent protein has at least
2	200 amino acids, a completeness value of at least 80% and has a crystal stability within
3	0.5% of its unit cell dimensions.
1	73. The crystal of claim 71, wherein the amino acid sequence comprises a
2	substitution at S65, wherein the substitution is selected from S65G, S65T, S65A, S65L,
3	S65C, S65V and S65I.
1	74. The crystal of claim 71, wherein said crystal has the following unit
2	cell dimensions in angstroms: a = 51.8, b= 62.8 and c= 70.7 with a space group of P2 2 2
	and an $\square$ angle of 90.00 $\square$ , a $\square$ angle of 90.00 $\square$ and a $\square$ angle of 90.00 $\square$ and the crystal has
3	a diffraction limit where 90% or greater of the potential reflections can be used to determine
4	the coordinates of the atoms.
5	the coordinates of the atoms.
:	75. A computational method of designing a fluoresent protein
2	comprising:
3	determining from a three dimensional model of a crystallized
4	fluorescent protein comprising a fluorescent protein with a bound ligand, at least one
5	interacting amino acid of the fluorescent protein that interacts with at least one first
6	chemical moiety of the ligand, and
7	selecting at least one chemical modification of the first chemical
8	moiety to produce a second chemical moiety with a structure to either decrease or increase
9	an interaction between the interacting amino acid and the second chemical moiety compared
10	to the interaction between the interacting amino acid and the first chemical moiety.
1	76. The computational method of claim 75, further comprising generating
2	the three dimensional model of the crystallized protein comprising a fluorescent protein
3	with an amino acid sequence substantially identical to SEQ ID NO:2.

1	77. The computational method of claim 75, wherein the selecting selects
2	the first chemical moiety that interacts with at least one of the amino acids listed in Figs. 5-1
3	to 5-28.
1	78. The computational method of claim 75, wherein the chemical
2	modification enhances hydrogen bonding interaction, charge interaction, hydrophobic
3	interaction, Van Der Waals interaction or dipole interaction between the second chemical
4	moiety and the interacting amino acid compared to the first chemical moiety and the
5	interacting amino acid.
1	79. A computational method of modeling the three dimensional structure
2	of a fluorescent protein comprising determining a three dimensional relationship between at
3	least two atoms listed in the atomic coordinates of Figs. 5-1 to 5-28.
1	80. The computational method of claim 79, wherein the determining
2	comprises determining the three dimensional structure of a fluorescent protein with an
3	amino acid sequence at least 80% identical to SEQ ID NO:2.
4	
1	81. The computational method of claim 79, wherein the determining
2	comprises determining the three dimensional structure of a fluorescent protein with an
3	amino acid sequence at least 95% identical to SEQ ID NO:2.
1	82. The computational method of claim 79, wherein the determining
2	comprises determining the three dimensional relationship of at least 1500 atoms listed in
3	Figs. 5-1 to 5-28.
1	83. A device comprising a storage device and, stored in the device, at
2	least 10 atomic coordinates selected from the atomic coordinates listed in Figs. 5-1 to 5-28.

- 1 84. The device of claim 83, wherein the storage device is a computer
- 2 readable device that stores code that receives as input the atomic coordinates.

The device of claim 84, wherein computer readable device is a floppy disk or a hard drive.

- 86. A nucleic acid molecule comprising a nucleotide sequence encoding a functional engineered fluorescent protein whose amino acid sequence is substantially identical to the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and which differs from SEQ ID NO:2 by at least a substitution at Q69, wherein said functional engineered fluorescent protein has a different fluorescent property than Aequorea green fluorescent protein.
- 9 87. The nucleic acid molecule of claim 86, wherein said substitution at Q69 is selected from the group of K, R, E and G.
- 11 88. The nucleic acid molecule of claim 86, wherein said amino acid sequence further comprises a function mutation at S65.
- 13 89. A nucleic acid molecule comprising a nucleotide sequence encoding a functional
  14 engineered fluorescent protein whose amino acid sequence is substantially identical to
  15 the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and
  16 which differs from SEQ ID NO:2 by at least a substitution at E222, but not including
  17 E222G, wherein said functional engineered fluorescent protein has a different
  18 fluorescent property than Aequorea green fluorescent protein.
- 19 90. The nucleic acid molecule of claim 89, wherein said substitution at E222 is selected from the group of N and Q.
- 21 91. The nucleic acid molecule of claim 89, wherein said amino acid sequence further comprises a function mutation at F64.
- 23 92. A nucleic acid molecule comprising a nucleotide sequence encoding a functional
  24 engineered fluorescent protein whose amino acid sequence is substantially identical to
  25 the amino acid sequence of Aequorea green fluorescent protein (SEQ ID NO:2) and
  26 which differs from SEQ ID NO:2 by at least a substitution at Y145, wherein said
  27 functional engineered fluorescent protein has a different fluorescent property than
  28 Aequorea green fluorescent protein.
- 29 93. The nucleic acid molecule of claim 92, wherein said substitution at Y145 is selected from the group of W, C, F, L, E, H, K and Q.

31 94. The nucleic acid molecule of claim 92, wherein said amino acid sequence further comprises a function mutation at Y66.

- 33 95. A method of identifying a test chemical, comprising:
- 34 contacting a test chemical a sample containing a biological entity labeled with a
- functional, engineered fluorescent protein or a polynucleotide encoding said functional,
- 36 engineered fluorescent protein, and
- detecting fluorescence of said functional engineered fluorescent protein.
- The method of claim 95, wherein said fluorescence in the presence of a test chemical is greater than in the absence of said test chemical.
- The method of claim 96, wherein said polynucleotide encoding said functional,
- engineered fluorescent protein is operatively linked to a genomic polynucleotide.
- 98. The method of claim 95, wherein said functional, engineered fluorescent protein is fused to second functional protein.
- The method of claim 96, wherein said polynucleotide encoding said functional, engineered fluorescent protein is operatively linked to a response element.
- The method of claim 96, wherein said polynucleotide encoding said functional,
- engineered fluorescent protein is operatively linked to a response element in a
- 48 mammalian cell.

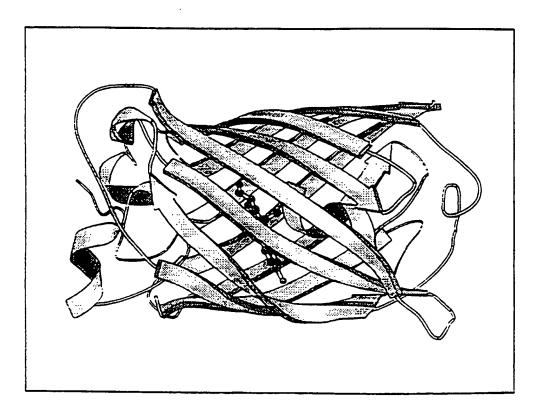


Figure 1a

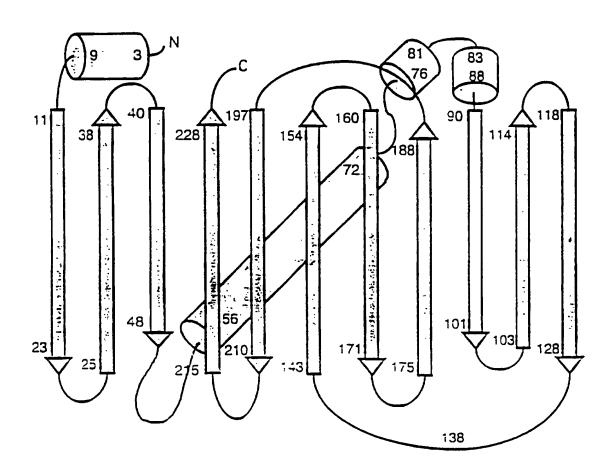


Figure 1b

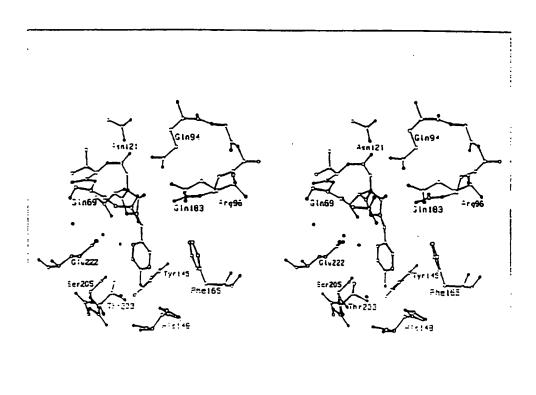


Figure 2a

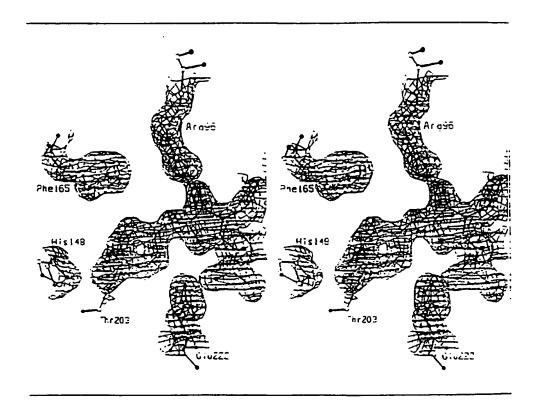


Figure 2b

#### (xi) SEQUENCE DESCRIPTION:

				`~'	, ,														
	%0: ! %0: 1		ATG Het 1	AGT Ser	ida Lys	GGA Gly	GAA Glu 5	Glu	ren C11	TTC Phe	ACT Thr	564 617 10	Val	GTC Val	EEA Pro	ATT ile	CTT Leu 15	Val	<b>-8</b>
																	Gly	GAS Glu	75
						ASD					Lys					Phe		TSE Cys	:44
				ACT Thr 50	Gly					Pro					Val			TTC Phe	172
				Tyr														223 P1A 08	240
				GAC ASD				Ser					Gly						225
				ATA Ile			LYS					lyr							336
				Phe		Gly					Asn								384
			Asp	177 Phe 130	LYS	Glu	ÀSD	Gly	Asn 135	He	Leu	Gly	His	Lys 143	Leu	Glu	Tyr	Asn	-32
			Tyr 145		Ser	His	ASR	vai 150	Tyr	ile	net	Ala	155	LYS	ผเก	LYS	Asn	:50	483
			Ιle	LYS	Val	ASR	Phe 165	LYS	::e	Arg	His	170	lle	Clu	ASO	Gly	Ser 175	val	528
			Gln	Leu	Ala	ASD 180	HIS	Iyr	C:n	SIn	AS:1	The	Pro	11e	Gly	450 190	Cly	Pro	576
			Val	Leu	Leu 195	Pro	ASD	Asn	H:5	1 yr 200	Leu	Ser	ihr	Gln	Ser 2 <b>05</b>	Ala	Leu	Ser	524
		(	Lys	GAT ASD 210	Pro	Asn	Giu	LYS	Arg 215	CZA	nis	Het	Val	Leu 220	Leu	Glu	Phe		57Z 717
			7hr 225	GCT Ala	Ala	Gly	He	Thr 230	# : S	S:Y	ret	CZA	61u 235	Leu	iye	Lys	16		

Figure 3

### 7/36

TODDY 3655, STOR numanized codon usage. With an additional amino acid after the start met to provide optimal kozak sequence

; 18 2" 36 45 54 ATS STS ASS AND GOS GAG GAG GTS TTS ASS STS GTS GTS GES ATS STS GTS GAG
Met Val Ser Lys Gly Glu Glu Leu Pne Tnr Gly Val Val Pro Ile Leu Val Glu
£3 77 81 90 99 108
TTS GAR SEE GAR GTA AAC GGE CAR AAG TTE AGE GTO TEE GGE GAG GGE GAG GGE
Leu Asp Gly Asp Val Asp Gly Ris Lys Phe Ser Val Ser Gly Glu Gly Glu Gly
117 126 135 144 153 162 SAT SEE ACC TAC GOC AAG CTG ACC CTG AAG TTC ATC TSC ACC ACC GCC AAG CTG
Asp Ala Thr Tyr Gly Lys Leu Thr Leu Lys Phe Ile Cys Thr Thr Gly Lys Leu
171 - 180 189 198 207 216
CCC GTG CCC TGG CCC ACC CTC GTG ACC ACC TTC GGC TAC GGC GTG CAG TGC TTC
Pro Val Pro Trp Pro Thr Leu Val Thr Thr Phe Gly Tyr Gly Val Gln Cye Phe
225 234 243 252 261 270 ccc csc csc csc csc csc csc csc csc cs
Ala Arg Tym Pro Asp His Met Lys Glm His Asp Pne Phe Lys Ser Ala Met Pro
279 298 297 306 315 324
GAA GGC TAC GTC CAG GAG CGC ACC ATC TTC TTC AAG GAC GAC GGC AAC TAC AAG
Glu Gly Tyr Val Gin Glu Arg The Ile Phe Phe Lys Asp Asp Gly Asn Tyr Lys
144 144 194
113 142 151 160 169 178 ACC CTC CTC GAG CTC AAC CTC AAC CTC AAC CTC CTC CTC CTC
Thr Arg Ala Glu Val Lys Phe Glu Gly Asp Thr Leu Val Ash Arg Ile Glu Leu
102 196 405 414 423 432
ANG GOT ATT GAC TITE ANG GAG GAC GOT AND ATT ETS GGS CAT ANG CITS GAG TAC
Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Els Lys Leu Glu Tyr
Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Els Lys Leu Glu Tyr
AAG GSC ATC GAC TTC AAG GAC GAC GSC AAC ATC CTS GGS CAC AAG CTG GAG TAC Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Eis Lys Leu Glu Tyr  441 450 459 468 477 486 AAC TAC AGC GAC GAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC
ANG GOT ATT GAC TITE ANG GAG GAC GGC AND ATT STO GGG CAC ANG CTG GAG TAC Lys Gly lie Asp Phe Lys Glu Asp Gly Ash Ile Leu Gly Eis Lys Leu Glu Tyr 441 450 459 468 477 486 AND TAC AND GGC CAC AND GTC TAT ATC ATG GGC GAC ANG CAG ANG AND GGC ATC
AAG GOT ATT GAC TITC AAG GAG GAC GGC AAC ATC STO GGO CAC AAG CTG GAG TAC Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Eis Lys Leu Glu Tyr  441 450 459 468 477 486 AAC TAC AAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC ASN Tyr Asn Ser His Asn Val Tyr Ile Met Ala Asp Lys Gln Lys Asn Gly Ile 455 504 513 522 531 540
AAG GOT ATT GAC TTC AAG GAG GAC GGC AAC ATC CTG GGG CAC AAG CTG GAG TAC Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Eis Lys Leu Glu Tyr  441 450 459 468 477 486 AAC TAC AAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC ASN Tyr Ash Ser His Ash Val Tyr Ile Met Ala Asp Lys Gln Lys Ash Gly Ile  485 504 513 522 531 540 AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAG CTC GCC
AAG GSC ATC GAC TTC AAG GAC GAC GSC AAC ATC CTS GGS CAC AAG CTG GAG TAC Lys Gly lie Asp Phe Lys Glu Asp Gly Asn lie Leu Gly Eis Lys Leu Glu Tyr  441 450 459 468 477 486 AAC TAC AAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC ASN Tyr Asn Ser His Asn Val Tyr lie Met Ala Asp Lys Gln Lys Asn Gly lie  495 504 513 522 531 540 AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAG CTC GCC
AAG GSC ATC GAC TTC AAG GAC GAC GSC AAC ATC CTS GGS CAC AAG CTG GAG TAC Lys Gly lie Asp Phe Lys Glu Asp Gly Asn fle Leu Gly Eis Lys Leu Glu Tyr  441 450 459 468 477 486 AAC TAC AAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC ASN Tyr Asn Ser His Asn Val Tyr fle Met Ala Asp Lys Gln Lys Asn Gly Ile  495 504 513 522 531 540 AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAC CTC GCC Lys Val Asn Phe Lys fle Arg His Asn Ile Glu Asp Gly Ser Val Gln Leu Ala  549 558 567 576 585 594
AAG GOT ATT GAC TTC AAG GAG GAC GGC AAC ATC CTG GGG CAC AAG CTG GAG TAC  Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Eis Lys Leu Glu Tyt  441
AAG GOO ATT GAC TTO AAG GAC GAC GGC AAC ATC CTG GGG CAC AAG CTG GAG TAC Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Eis Lys Leu Glu Tyr  441 450 459 468 477 486 AAC TAC AAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC ASN Tyr Asn Ser His Asn Val Tyr Ile Met Ala Asp Lys Gln Lys Asn Gly Ile  455 504 513 522 531 540 AAG GTG AAC TTT AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAG CTC GCC Lys Val Asn Phe Lys Ile Arg His Asn Ile Glu Asp Gly Ser Val Gln Leu Ala  549 558 567 576 585 594 GAC CAC TAC CAG CAG CAC AAC ACC GGC GAC GGC GTG CTG CTG CCC GAC
Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Els Lys Leu Glu Tyr  441 450 459 468 477 486  AAC TAC AAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC  Asn Tyr Asn Ser His Asn Val Tyr lie Met Ala Asp Lys Gln Lys Asn Gly Ile  485 504 513 522 531 540  AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAG CTC GCC  Lys Val Asn Phe Lys Ile Arg His Asn Ile Glu Asp Gly Ser Val Gln Leu Ala  549 558 567 576 585 594  GAC CAC TAC CAG CAG AAC ACC CGC ATC GGC GAC GGC GTG CTG CTG CCC GAC  Asp His Tyr Gln Gln Asn Tor Pro Ile Gly Asp Gly Pro Val Leu Leu Pro Asp
Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Els Lys Leu Glu Tyr  441 450 459 468 477 486  AAC TAC AAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC  Asn Tyr Asn Ser His Asn Val Tyr lie Met Ala Asp Lys Gln Lys Asn Gly Ile  495 504 513 522 531 540  AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAG CTC GCC  Lys Val Asn Phe Lys Ile Arg His Asn Ile Glu Asp Gly Ser Val Gln Leu Ala  549 558 567 576 585 594  GAC CAC TAC CAG CAG AAC ACC CCC ATC GGC GAC GGC CCC GTG CTG CCC  ASP His Tyr Gln Gln Asn Tnr Pro ILe Gly Asp Gly Pro Val Leu Leu Pro Asp  603 612 621 630 639 648
Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Els Lys Leu Glu Tyr  441 450 459 468 477 486  AAC TAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC  ASN Tyr Asn Ser His Asn Val Tyr lie Met Ala Asp Lys Gln Lys Asn Gly Ile  495 504 513 522 531 540  AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAG CTC GCC  Lys Val Asn Phe Lys Ile Arg His Asn Ile Glu Asp Gly Ser Val Gin Leu Ala  549 558 567 576 585 594  GAC CAC TAC CAG CAG AAC ACC CCC ATC GGC GAC GGC CCC GTG CTG CCC  ASP His Tyr Gln Gln Asn Tnr Pro ILe Gly Asp Gly Pro Val Leu Leu Pro Asp  603 612 621 630 639 648  AAC CAC TAC CTG AGC TAC CAG TAC GGC CTG AAA GAC CCC AAC GAG CGC
Lys Gly lie Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Els Lys Leu Glu Tyr  441 450 459 468 477 486  AAC TAC AAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC  Asn Tyr Asn Ser His Asn Val Tyr lie Met Ala Asp Lys Gln Lys Asn Gly Ile  495 504 513 522 531 540  AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAG CTC GCC  Lys Val Asn Phe Lys Ile Arg His Asn Ile Glu Asp Gly Ser Val Gln Leu Ala  549 558 567 576 585 594  GAC CAC TAC CAG CAG AAC ACC CCC ATC GGC GAC GGC CCC GTG CTG CCC  ASP His Tyr Gln Gln Asn Tnr Pro ILe Gly Asp Gly Pro Val Leu Leu Pro Asp  603 612 621 630 639 648
AAG GSC ATC GAC TTC AAG GAC GAC GSC AAC ATC CTS GGS CAC AAG CTG GAG TAC  Lys Gly lie Asp Phe Lys Glu Asp Gly Asn fle Leu Gly Eis Lys Leu Glu Tyr  441
AAG GSC ATC GAC TTC AAG GAC GAC GSC AAC ATC CTS GGS CAC AAG CTG GAG TAC  Lys Gly Ile Asp Phe Lys Glu Asp Gly Asn Ile Leu Gly Els Lys Leu Glu Tyr  441
AAG GSC ATC GAC TTC AAG GAC GAC GSC AAC ATC CTS GGS CAC AAG CTG GAG TAC  Lys Gly lie Asp Phe Lys Glu Asp Gly Asn fle Leu Gly Eis Lys Leu Glu Tyr  441
LYS GLY 11e ASP Phe LYS GLU ASP GLY ASR ILE LEU GLY ELS LYS LEU GLU TYT  441 450 459 468 477 486  ARC THE ARC AGE CAC ARC GTE TAT ATC ATG GGC GAC AAG CAG AAG AGC GGC ATC  ASR TYT ASD SET HIS ASR VAL TYT ILE MET ALA ASP LYS GLR LYS ASR GLY ILE  495 504 513 522 531 540  AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CTG CTC  LYS VAL ASR PHE LYS ILE ATG HIS ASR ILE GLU ASP GLY SET VAL GLR LEU ALA  549 558 567 576 585 594  GAC CAC TAC CAG CAG AAC ACC CCC ATC GGC GAC GGC GGC GTG CTG CTG CTG  ASP HIS TYT GLR GLR ASR TRT PTO ILE GLY ASP GLY PTO VAL LEU LEU PTO ASP  603 612 621 630 639 648  AAC CAC TAC CTG AGC TAC CAG TCC GGC CTG AGC AAC GAC GAC GAC  ASR HIS TYT LEU SET TYT GLR SET ALA LEU SET LYS ASP PTO ASR GLU LYS ATG  GAT CAC ATG GTC CTG CTG GAG TTC GTG ACC GGC ATC CAC GAG ATC  ASP HIS TYT LEU SET TYT GLR SET ALA LEU SET LYS ASP PTO ASR GLU LYS ATG  GAT CAC ATG GTC CTG CTG CTG GAG TTC GTG ACC GGC ATC CAC GGC ATC  ASP HIS RET VAL LEU LEU GLU PHE VAL TAT ALA ALA GLY ILE THT HIS GLY MET
ANG GOT ATT GAC TIT ANG GAG GAC GGC AAC ATC CTG GGG CAC RAG CTG GAG TAC LYS GLY ILE ASP Phe Lys Glu ASP GLY ASR ILE Leu GLY ELS LYS Leu GLU TYT  441 450 459 468 477 486  AAC TAC AAC AGC CAC AAC GTC TAT ATC ATG GCC GAC AAG CAG AAG AAC GGC ATC  ASR TYT ASD SET HIS ASR VAL TYT ILE MET ALA ASP LYS GLR LYS ASR GLY ILE  495 504 513 522 531 540  AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CAG CTC GCC  LYS VAL ASR Phe Lys ILE ATG HIS ASR ILE GLU ASP GLY SET VAL GLR LEU ALA  549 558 567 576 585 594  GAC CAC TAC CAG CAG AAC ACC CCC ATC GGC GAC GGC CCC GTG CTG CTG GCC  ASP HIS TYT GLR GLR ASR TAT PRO ILE GLY ASP GLY PRO VAL LEU LEU PRO ASP  AAC CAC TAC CTG AGC TAC CAG TCC GGC CTG AGC AAC GAC GGC  ASR HIS TYT LEU SET TYT GLR SET ALA LEU SET LYS ASP PTO ASR GLU LYS ATG  GAT CAC ATG GTC CTG CTG GAG TTC GTG ACC GCC GGG ATC CAC GAC ATC  ASP HIS TYT LEU SET TYT GLR SET ALA LEU SET LYS ASP PTO ASR GLU LYS ATG  GAT CAC ATG GTC CTG CTG GAG TTC GTG GCC GGG ATC ACC CAC GGC ATC  ASP HIS MET VAL LEU LEU GLU PHE VAL TAT ALA ALA GLY ILE THY HIS GLY MET
LYS GLY 11e ASP Phe LYS GLU ASP GLY ASR ILE LEU GLY ELS LYS LEU GLU TYT  441 450 459 468 477 486  ARC THE ARC AGE CAC ARC GTE TAT ATC ATG GGC GAC AAG CAG AAG AGC GGC ATC  ASR TYT ASD SET HIS ASR VAL TYT ILE MET ALA ASP LYS GLR LYS ASR GLY ILE  495 504 513 522 531 540  AAG GTG AAC TTC AAG ATC CGC CAC AAC ATC GAG GAC GGC AGC GTG CTG CTC  LYS VAL ASR PHE LYS ILE ATG HIS ASR ILE GLU ASP GLY SET VAL GLR LEU ALA  549 558 567 576 585 594  GAC CAC TAC CAG CAG AAC ACC CCC ATC GGC GAC GGC GGC GTG CTG CTG CTG  ASP HIS TYT GLR GLR ASR TRT PTO ILE GLY ASP GLY PTO VAL LEU LEU PTO ASP  603 612 621 630 639 648  AAC CAC TAC CTG AGC TAC CAG TCC GGC CTG AGC AAC GAC GAC GAC  ASR HIS TYT LEU SET TYT GLR SET ALA LEU SET LYS ASP PTO ASR GLU LYS ATG  GAT CAC ATG GTC CTG CTG GAG TTC GTG ACC GGC ATC CAC GAG ATC  ASP HIS TYT LEU SET TYT GLR SET ALA LEU SET LYS ASP PTO ASR GLU LYS ATG  GAT CAC ATG GTC CTG CTG CTG GAG TTC GTG ACC GGC ATC CAC GGC ATC  ASP HIS RET VAL LEU LEU GLU PHE VAL TAT ALA ALA GLY ILE THT HIS GLY MET

Figure 4

FIG 5-1

CRYST1 ORIGX1 ORIGX2 ORIGX3 SCALE1 SCALE2 SCALE3 ATOM	1 2	1.000 0.000 0.000 0.019 0.000 N CA	000 000 317 000 000 SER SER	0.000000 1.000000 0.000000 0.000000 0.015912 0.000000 2	0.000000 0.000000 1.000000 0.000000 0.014151 28.888 27.638		90.00 0.00000 0.00000 0.00000 0.00000 0.00000 52.301 52.516 51.644	1.00 85.05 1.00 80.05 1.00 85.35
ATOM ATOM	3 4	С О	SER SER	2 2	26.499 26.606	9.656	50.915	1.00 84.55
ATOM	5	СВ	SER	2	27.783	11.635	52.378	1.00 70.97
ATOM	6	OG	SER	2 3	27.690	12.033	51.012 51.731	1.00 44.08
atom Atom	7 8	N CA	LYS	3	25.418 24:141	10.191	51.036	1.00 87.15
ATOM	9	c	LYS	3	24.214	10.266	49.497	1.00 76.86
MOTA	10	0	LYS	3	24.107	9.258	48.774	1.00 78.27 1.00 89.44
ATOM	11	CB	LYS	3 3	23.127 21.768	11.240	51.521 51.949	1.00 89.44
ATOM ATOM	12 13	CD CD	LYS	3	20.681	11.781	51.987	1.00 76.53
ATOM	14	CE	LYS	3	20.711	12.655	53.243	1.00 68.55
MOTA	15	NZ	LYS	3	20.816	14.103	52.953	1.00 46.24
ATOM	16	N	GLY	4 4	24.318 24.297	11.495	49.015 47.605	1.00 45.97
MOTA MOTA	17 18	CA C	GLY	4	25.425	11.205	46.796	1.00 31.90
MOTA	19	ō	GLY	4	25.234	10.923	45.619	1.00 33.63
ATOM	20	N	GLU	5 5	26.606 27.821	11.082	47.420 46.726	1.00 32.54
atom atom	21 22	CA C	GLU	5 5	27.521	9.590	45.616	1.00 28.43
ATOM	23	ō	GLU	5	27.850	9.803	÷4.444	1.00 26.12
ATOM	24	CB	CLU	5	28.873	10.053	47.718	1.00 33.53
MOTA	25 26	CG CD	GLU	5 5	30.337 31.311	10.461 9.584	47.425 48.170	1.00 90.82
atom Atom	27		GLU	5	31.508	9.677	49.381	1.00 74.80
ATOM	28	OE2	GLU	5	31.839	8.653	47.403	1.00100.00
ATOM	29	N	GLU	6	26.883	9.499 7.410	46.017 45.150	1.00 23.57
MOTA MOTA	30 31	CA C	GLU	6 6	26.479 25.561	7.837	43.979	1.00 31.13
ATOM	32	ō	GLU	6	25.479	7.142	42.955	1.00 30.96
MOTA	33	CB	GLU	6	25.780	5.330	45.992	1.00 35.64
ATOM	34	CG N	GLU LEU	6 7	25.260 24.864	6.893 3.966	47.338 44.138	1.00 22.25
ATOM ATOM	35 36	CA	LEU	7	23.954	9.456	-3.089	1.00 21.61
ATOM	37	C	LEU	7	24.693	10.061	41.917	1.00 15.90
ATOM	38	0	LEU	7	24.152	10.250	40.836 43.665	1.00 18.38
ATOM ATOM	39 40	CS	LEU	7 7	23.050 21.672	10.548	44.098	1.00 22.41
MOTA	41	CD1		7	21.597	8.536	44.074	1.00 31.64
MOTA	42	CD2		7	21.332	10.591	45.485	1.00 33.14
ATOM	43	N	SHE	8 9	25.944 26.740	10.407	42.157 41.159	1.00 20.75
. ATOM ATOM	44 45	CA C	PHE	8	27.818	10.333	40.427	1.00 30.59
ATOM	46	ŏ,	PHE	9	28.590	10.856	39.600	1.00 30.05
MOTA	47	CB	PHE	9	27.309	12.375	41.820	1.00 15.95
ATOM ATOM	48 49	CG	PHE	3 3	26.222 25.672	13.355	42.163 43.447	1.00 13.29
MOTA	50	CD2		š	25.725	14.227	+1.129	1.00 13.12
ATOM	51	CE 1	PHE	8	24.661	14.290	43.772	1.00 15.14
ATOM .	52	CEZ		3 3	24.712	15.137	41.499 42.794	1.00 13.19 1.00 5.69
ATOM ATOM	53 54	) 11	PHE		24.192 27.798	15.170	42.699	1.00 27.35
ATOM	55	CA	THR	ż	28.704	8.122	40.175	1.00 34.93
MOTA	56	С	THR		28.709	998	33.636	1.00 45.11
ATOM	57	0	THR	ě	29.642	7.452	13.062	1.00 50.55
MOTA:	58 59	C3	THR THR		28.447 19.629	6.798 6.330	45.892 41.827	1.00 40.40
ATOM ATOM	50				17.801	5.779	19.959	1.00 39.76

ATOM 61 N GLY 10 27.690 8.510 37.956 1.00 30.53 ATOM 63 C GLY 10 27.691 8.458 16.507 1.00 23.21 ATOM 63 C GLY 10 27.144 9.746 35.914 1.00 16.55 ATOM 64 O GLY 10 27.011 10.729 36.617 1.00 25.70 ATOM 65 N VAL 11 26.299 10.863 31.971 1.00 22.28 ATOM 66 CA VAL 11 26.299 10.863 31.971 1.00 22.28 ATOM 66 CA VAL 11 26.299 10.863 31.971 1.00 22.28 ATOM 67 C VAL 11 22.4758 11.020 34.455 1.00 20.43 ATOM 68 O VAL 11 22.972 10.062 34.455 1.00 20.43 ATOM 69 CB VAL 11 22.972 10.062 34.455 1.00 20.43 ATOM 70 CGI VAL 11 22.972 10.062 31.734 1.00 31.75 ATOM 71 CGI VAL 11 27.490 10.048 32.015 1.00 31.85 ATOM 72 N VAL 12 22.407 13.624 34.516 1.00 31.75 ATOM 73 N VAL 12 22.407 13.624 34.516 1.00 13.85 ATOM 74 C VAL 12 22.407 13.624 34.516 1.00 14.37 ATOM 75 C VAL 12 22.407 13.624 34.516 1.00 14.37 ATOM 76 CB VAL 12 22.407 13.624 34.516 1.00 14.37 ATOM 77 CGI VAL 12 23.000 12.561 35.433 1.00 12.88 ATOM 77 CGI VAL 12 23.007 14.639 34.179 1.00 14.72 ATOM 78 CGI VAL 12 23.007 14.639 34.179 1.00 14.72 ATOM 78 CGI VAL 12 23.007 14.639 34.179 1.00 15.01 ATOM 79 N PRO 13 20.118 11.936 37.728 1.00 15.01 ATOM 80 CA PRO 13 20.118 11.936 37.728 1.00 16.55 ATOM 80 C PRO 13 20.118 11.936 37.728 1.00 16.57 ATOM 80 C PRO 13 20.118 11.936 31.235 1.00 16.75 ATOM 80 C PRO 13 19.737 15.468 35.337 1.00 15.14 ATOM 80 C PRO 13 19.737 15.468 35.337 1.00 15.14 ATOM 80 C PRO 13 19.737 15.468 35.337 1.00 15.14 ATOM 80 C PRO 13 19.737 15.468 35.337 1.00 15.14 ATOM 80 C PRO 13 19.737 15.468 35.337 1.00 15.14 ATOM 80 C PRO 13 19.737 15.468 35.337 1.00 15.14 ATOM 80 C PRO 13 19.737 15.468 35.337 1.00 15.14 ATOM 80 C PRO 13 19.737 15.468 35.337 1.00 15.14 ATOM 80 C PRO 13 19.737 15.468 33.337 1.00 15.14 ATOM 80 C PRO 13 19.043 12.422 33.550 1.00 16.69 ATOM 90 C ELU 15 16.694 2.199 33.499 1.00 1.00 1.04 ATOM 90 C ELU 15 16.695 30.909 33.499 1.00 11.04 ATOM 90 C ELU 15 16.919 30.909 33.499 1.00 11.04 ATOM 90 C ELU 15 16.919 30.909 35.909 1.00 11.04 ATOM 90 C ELU 15 16.919 30.909 35.909 1.00 13.69 ATOM 100 C VAL 16 16.919 2.1472 34.689 1.00 11.65 ATOM 100 C VAL 16									
NOM   63   C	MOTA	61	N	GLY	10	27.690	8.510	37.956	1.00 30.53
ATOM 63 C GLY 10 27.144 9.746 35.914 1.00 16.55 ATOM 65 N VAL 11 26.835 9.719 34.629 1.00 16.39 ATOM 66 CA VAL 11 26.209 10.863 33.971 1.00 22.28 ATOM 66 CA VAL 11 26.209 10.863 33.971 1.00 22.28 ATOM 67 C VAL 11 26.209 10.863 33.971 1.00 22.28 ATOM 69 CB VAL 11 23.972 10.062 34.455 1.00 20.43 ATOM 69 CB VAL 11 23.972 10.062 34.455 1.00 20.43 ATOM 70 CGI VAL 11 23.972 10.062 34.455 1.00 20.43 ATOM 71 CG2 VAL 11 27.480 10.048 32.015 1.00 31.75 ATOM 71 CG2 VAL 11 27.480 10.048 32.015 1.00 31.75 ATOM 72 N VAL 12 24.417 12.227 34.931 1.00 22.12 ATOM 73 CA VAL 12 22.007 12.561 34.455 1.00 13.85 ATOM 74 C VAL 12 22.007 14.639 34.179 1.00 13.42 ATOM 75 O VAL 12 22.007 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 22.007 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 22.007 13.624 34.516 1.00 14.37 ATOM 76 CB VAL 12 22.007 13.639 1.00 15.01 ATOM 77 CGI VAL 12 22.007 13.639 1.00 15.01 ATOM 78 CG2 VAL 12 23.007 13.639 1.00 15.01 ATOM 78 CG2 VAL 12 23.008 13.362 37.422 1.00 17.57 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 19.927 15.468 35.337 1.00 15.14 ATOM 85 CD PRO 13 19.927 15.468 35.337 1.00 15.14 ATOM 86 C ILE 14 20.196 16.766 33.557 1.00 17.36 ATOM 88 C ILE 14 20.196 16.766 33.557 1.00 17.36 ATOM 88 C ILE 14 20.196 16.766 33.557 1.00 17.36 ATOM 89 C ILE 14 22.201 18.169 35.015 1.00 17.69 ATOM 90 CG ILE 14 22.201 18.199 33.4425 1.00 1.00 15.41 ATOM 90 CG ILE 14 22.201 18.995 34.342 1.00 16.54 ATOM 90 CG ILE 14 22.201 18.995 34.342 1.00 15.14 ATOM 90 CG ILE 14 20.763 20.997 35.091 1.00 15.01 ATOM 90 CG ILE 14 20.763 20.997 35.094 1.00 15.43 ATOM 90 CG ILE 14 20.763 20.997 34.421 1.00 1.56 ATOM 90 CG ILE 14 20.763 20.997 34.421 1.00 1.56 ATOM 90 CG ILE 14 20.763 20.997 34.421 1.00 1.56 ATOM 90 CG ILE 14 20.763 20.907 35.991 1.00 10.00 ATOM 90 CG ILE 14 20.763 20.907 34.421 1.00 1.55 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 16.57 ATOM 102 CG ILE 15 15.008 20.098 31.913 1.00 16.57 ATOM 103 CA VAL 16									1.00 23.21
ATOM 64 O GLY 10 27.011 10.729 36.617 1.00 25.70 ATOM 65 N VAL 11 26.815 9.719 34.629 1.00 16.39 ATOM 66 CA VAL 11 26.825 9.719 34.629 1.00 16.39 ATOM 67 C VAL 11 24.758 11.020 34.479 1.00 22.28 ATOM 69 CB VAL 11 22.972 10.062 34.456 1.00 20.43 ATOM 69 CB VAL 11 22.972 10.062 34.456 1.00 30.87 ATOM 70 CG1 VAL 11 26.173 10.664 32.467 1.00 30.87 ATOM 71 CG2 VAL 11 22.17 10.664 32.467 1.00 30.87 ATOM 72 N VAL 12 22.417 12.227 34.931 1.00 31.75 ATOM 73 CA VAL 12 22.407 13.624 34.516 1.00 31.75 ATOM 75 O VAL 12 22.407 13.624 34.516 1.00 14.37 ATOM 76 CB VAL 12 22.000 13.662 37.422 1.00 14.53 ATOM 77 CG1 VAL 12 23.007 14.639 34.179 1.00 14.73 ATOM 78 CG2 VAL 12 23.200 13.672 36.839 1.00 15.01 ATOM 79 N PRO 13 21.180 13.382 34.066 1.00 14.73 ATOM 79 N PRO 13 21.180 13.382 34.066 1.00 14.73 ATOM 80 CA PRO 13 20.493 14.382 33.255 1.00 10.75 ATOM 81 C PRO 13 20.418 11.936 37.728 1.00 15.51 ATOM 82 O PRO 13 20.418 12.32.70 13.673 35.503 1.00 15.41 ATOM 83 CB PRO 13 19.727 15.468 35.337 1.00 15.41 ATOM 85 CD PRO 13 19.225 13.707 32.745 1.00 17.36 ATOM 85 CD PRO 13 19.225 13.707 32.745 1.00 15.41 ATOM 86 C PRO 13 19.235 13.707 32.745 1.00 15.41 ATOM 87 C PRO 13 19.031 12.422 33.550 1.00 15.41 ATOM 89 O ILE 14 18.768 18.769 35.439 1.00 15.01 ATOM 89 O ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 89 O ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 20.196 16.766 33.557 1.00 12.93 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 10.00 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.05 ATOM 99 CG ILE 14 18.724 18.878 32.399 1.00 11.05 ATOM 100 CD1 ILE 15 15.406 29.406 33.4355 1.00 15.57 ATOM 101 CO UNI ILE 15 15.406 29.406 33.4355 1.									1.00 16.55
ATOM 65 N VAL 11 26.835 9.719 34.629 1.00 16.39 ATOM 66 CA VAL 11 26.29 10.863 33.971 1.00 22.60 ATOM 68 C VAL 11 26.758 11.020 34.479 1.00 29.60 ATOM 69 CB VAL 11 23.972 10.062 34.456 1.00 29.60 ATOM 70 CG1 VAL 11 25.912 11.980 31.734 1.00 31.85 ATOM 71 CG2 VAL 11 25.912 11.980 31.734 1.00 31.85 ATOM 72 N VAL 12 22.407 12.227 34.931 1.00 30.87 ATOM 73 CA VAL 12 23.060 12.561 35.433 1.00 12.88 ATOM 74 C VAL 12 23.060 12.561 35.433 1.00 12.88 ATOM 75 O VAL 12 23.007 14.639 34.156 1.00 14.37 ATOM 76 CB VAL 12 23.007 14.639 34.179 1.00 13.74 ATOM 77 CG1 VAL 12 23.007 14.639 34.179 1.00 13.75 ATOM 78 CG2 VAL 12 23.007 14.639 34.179 1.00 13.74 ATOM 78 CG2 VAL 12 23.007 14.639 34.179 1.00 15.01 ATOM 78 CG2 VAL 12 23.007 14.639 34.179 1.00 15.01 ATOM 78 CG2 VAL 12 23.007 14.639 34.179 1.00 15.01 ATOM 78 CG2 VAL 12 23.007 14.639 34.179 1.00 15.01 ATOM 80 CA PRO 13 20.493 14.382 34.066 1.00 14.55 ATOM 80 CA PRO 13 20.493 14.382 34.066 1.00 14.55 ATOM 81 C PRO 13 20.493 14.382 34.066 1.00 16.55 ATOM 82 O PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 84 CG PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 85 CD PRO 13 19.225 13.707 32.745 1.00 17.62 ATOM 86 N LEE 14 19.893 17.991 34.266 1.00 14.91 ATOM 88 C LEE 14 19.893 17.991 34.266 1.00 14.91 ATOM 88 C LEE 14 19.893 17.991 34.266 1.00 12.08 ATOM 99 CG LEE 14 18.768 19.769 35.00 10.00 15.64 ATOM 90 CS LLE 14 19.893 17.991 34.266 1.00 12.08 ATOM 90 CS LLE 14 22.271 18.109 15.00 15.00 10.00 ATOM 90 CS LLE 14 19.893 17.991 34.266 1.00 12.08 ATOM 90 CG LEU 15 16.811 20.136 23.955 1.00 10.156 ATOM 90 CG LEU 15 16.811 20.136 23.955 1.00 10.158 ATOM 90 CG LEU 15 15.462 19.450 34.255 1.00 17.82 ATOM 90 CG LEU 15 16.811 20.136 3.955 1.00 10.00 10.00 ATOM 90 CG LEU 15 16.812 20.773 13.479 1.00 16.57 ATOM 100 CD1 LEU 15 15.567 39.90 34.401 1.00 15.37 ATOM 101 CC LEU 15 15.602 20.908 13.913 1.00 10.23 ATOM 102 C VAL 16 18.226 26.075 34.403 1.00 15.73 ATOM 103 CA VAL 16 18.226 26.075 34.403 1.00 15.73 ATOM 104 C CLU 17 14.301 27.529 25.317 3									1.00 25.70
ATOM 66 CA VAL 11 26.209 10.863 39.971 1.00 22.28 ATOM 67 C VAL 11 24.758 11.023 34.479 1.00 29.64 ATOM 68 O VAL 11 23.972 10.062 34.456 1.00 20.43 ATOM 69 CB VAL 11 26.173 10.664 32.467 1.00 30.87 ATOM 70 CG1 VAL 11 25.912 11.980 31.734 1.00 31.75 ATOM 71 CG2 VAL 11 27.480 10.048 32.015 1.00 31.85 ATOM 72 N VAL 12 24.417 12.227 34.931 1.00 20.12 ATOM 73 CA VAL 12 23.080 12.561 35.433 1.00 20.12 ATOM 75 O VAL 12 23.080 12.561 35.433 1.00 20.12 ATOM 75 O VAL 12 23.007 13.624 34.516 1.00 14.72 ATOM 76 CB VAL 12 23.007 13.624 34.516 1.00 14.72 ATOM 77 CG1 VAL 12 22.007 13.624 34.516 1.00 14.72 ATOM 77 CG1 VAL 12 22.007 13.623 34.799 1.00 15.40 ATOM 78 CG2 VAL 12 23.070 13.662 37.742 1.00 17.57 ATOM 78 CG2 VAL 12 23.070 13.662 37.742 1.00 17.57 ATOM 79 N PRO 13 20.493 14.382 34.066 1.00 14.72 ATOM 80 CA PRO 13 20.493 14.382 34.066 1.00 14.72 ATOM 81 C PRO 13 20.493 14.382 34.066 1.00 14.72 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.7.65 ATOM 82 C PRO 13 19.225 13.707 12.745 1.00 17.36 ATOM 85 CD PRO 13 19.225 13.707 22.745 1.00 17.36 ATOM 85 CD PRO 13 19.225 13.707 22.745 1.00 17.36 ATOM 86 C PRO 13 19.235 12.195 34.340 1.00 15.41 ATOM 86 C PRO 13 19.043 12.422 33.550 1.00 17.36 ATOM 86 C PRO 13 19.043 12.422 33.550 1.00 17.36 ATOM 86 C PRO 13 19.043 12.422 33.550 1.00 17.36 ATOM 86 C PRO 13 19.043 12.422 33.550 1.00 17.36 ATOM 87 C ALE 14 19.893 17.991 34.266 1.00 12.93 ATOM 87 C ALE 14 19.893 17.991 34.266 1.00 12.93 ATOM 89 C LEE 14 19.893 17.991 34.266 1.00 12.93 ATOM 90 CS LE 14 19.893 17.991 34.266 1.00 12.93 ATOM 90 CS LE 14 12.109 18.905 34.325 1.00 16.54 ATOM 90 CS LE 14 12.109 18.905 34.325 1.00 16.54 ATOM 90 CS LE 14 12.109 18.905 34.325 1.00 16.54 ATOM 90 CS LE 14 12.109 18.905 34.325 1.00 16.54 ATOM 90 CS LE 14 12.109 18.905 34.325 1.00 16.54 ATOM 90 CS LE 14 12.109 18.905 34.325 1.00 16.54 ATOM 90 CS LE 14 12.109 18.905 34.325 1.00 16.54 ATOM 90 CS LE 15 16.816 20.277 34.806 1.00 12.38 ATOM 90 CS LEU 15 16.816 20.277 34.806 1.00 12.58 ATOM 90 CS LEU 15 16.816 20.277 34.806 1.00 12.58 ATOM 90 CS LEU									. 00 16 39
ATOM 67 C VAL 11 24.758 11.020 34.479 1.00 29.60 ATOM 68 C VAL 11 22.972 10.062 34.456 1.00 20.48 ATOM 70 CG1 VAL 11 25.972 10.062 34.456 1.00 20.48 ATOM 71 CG2 VAL 11 25.972 10.063 34.456 1.00 20.18 ATOM 72 N VAL 12 24.417 12.227 34.931 1.00 31.85 ATOM 73 CA VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 74 C VAL 12 22.407 13.624 34.516 1.00 14.37 ATOM 75 O VAL 12 23.007 14.639 34.179 1.00 13.47 ATOM 76 CB VAL 12 23.007 14.639 34.179 1.00 13.74 ATOM 77 CG1 VAL 12 22.000 13.662 37.422 1.00 17.3 ATOM 78 CG2 VAL 12 23.200 13.662 37.422 1.00 17.3 ATOM 78 CG2 VAL 12 23.000 13.662 37.422 1.00 17.3 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 79 N PRO 13 21.180 13.382 34.066 1.00 14.73 ATOM 80 CA PRO 13 20.493 14.382 34.066 1.00 14.73 ATOM 81 C PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 85 CD PRO 13 20.315 12.195 34.340 1.00 17.66 ATOM 85 CD PRO 13 19.023 12.195 34.340 1.00 17.66 ATOM 86 K CLE 14 19.893 17.991 34.266 1.00 14.73 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 88 C LE 14 18.766 13.760 33.596 1.00 12.93 ATOM 89 O ILE 14 18.766 13.760 33.596 1.00 12.93 ATOM 89 O ILE 14 18.766 13.760 33.596 1.00 12.93 ATOM 91 CC1 ILE 14 22.271 18.169 33.955 1.00 16.15 ATOM 93 CD1 ILE 14 22.271 18.169 33.955 1.00 16.16 ATOM 93 CD1 ILE 14 22.191 18.199 19.007 35.015 1.00 16.15 ATOM 94 CD1 LEU 15 15.693 20.490 33.679 1.00 12.04 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 16.16 ATOM 96 C LEU 15 16.812 20.493 34.365 1.00 16.15 ATOM 97 C CLEU 15 16.811 20.136 33.955 1.00 16.15 ATOM 98 C LEU 15 17.080 21.599 35.015 1.00 16.15 ATOM 99 C CLEU 15 16.811 20.136 33.955 1.00 16.16 ATOM 90 C S LEU 15 16.811 20.136 33.955 1.00 16.16 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 12.25 ATOM 100 CD1 LEU 15 15.00 20.98 31.913 1.00 49.50 ATOM 100 CD1 LEU 15 15.00 20.98 31.913 1.00 49.50 ATOM 101 CD2 LEU 15 15.00 20.98 31.913 1.00 10.16 ATOM 101 CD2 LEU 15 15.00 20.985 31.919 1.00 10.56 ATOM 107 CG1 VAL 16 18.255 20.985 31.919 1									
ATOM 68 O VAL 11 23.972 10.062 34.456 1.00 20.43 ATOM 69 CB VAL 11 26.173 10.664 12.467 1.00 30.87 ATOM 70 CG1 VAL 11 25.912 11.980 31.734 1.00 31.75 ATOM 71 CG2 VAL 11 27.460 10.048 32.015 1.00 31.87 ATOM 72 N VAL 12 24.417 12.227 34.931 1.00 20.12 ATOM 73 CA VAL 12 22.060 12.561 35.433 1.00 12.88 ATOM 74 C VAL 12 22.067 13.624 34.516 1.00 14.37 ATOM 75 O VAL 12 23.007 13.624 34.516 1.00 14.37 ATOM 75 O VAL 12 23.007 13.624 34.516 1.00 14.37 ATOM 76 CB VAL 12 23.007 13.624 34.516 1.00 14.37 ATOM 77 CG1 VAL 12 22.007 13.602 37.7422 1.00 17.57 ATOM 77 CG1 VAL 12 22.007 13.603 37.728 1.00 16.53 ATOM 78 CG2 VAL 12 23.071 13.603 37.728 1.00 16.53 ATOM 79 N PRO 13 21.180 13.382 34.066 1.00 14.72 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.75 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.75 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.765 ATOM 83 CB PRO 13 20.116 15.589 34.141 1.00 7.65 ATOM 84 CC PRO 13 19.205 13.707 12.745 1.00 17.36 ATOM 85 CD PRO 13 19.205 13.707 12.745 1.00 17.36 ATOM 85 CD PRO 13 19.205 13.707 12.745 1.00 17.36 ATOM 85 CD PRO 13 19.205 13.707 12.745 1.00 17.36 ATOM 85 CD PRO 13 19.205 13.707 12.745 1.00 17.36 ATOM 86 CI LE 14 19.893 17.991 34.266 1.00 12.93 ATOM 86 CI LE 14 19.893 17.991 34.266 1.00 12.93 ATOM 89 C ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 89 C ILE 14 18.768 13.760 33.556 1.00 10.60 1.00 1.00 1.00 1.00 1.00 1.									
ATOM 69 CB VAL 11 26.173 10.664 32.467 1.00 30.87 ATOM 70 CG1 VAL 11 27.480 10.048 31.734 1.00 31.75 ATOM 71 CG2 VAL 11 27.480 10.048 32.015 1.00 33.85 ATOM 72 N VAL 12 24.417 12.227 34.931 1.00 20.18 ATOM 73 CA VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 74 C VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 75 O VAL 12 23.070 13.624 34.516 1.00 14.37 ATOM 76 CB VAL 12 23.070 13.624 34.516 1.00 16.37 ATOM 76 CB VAL 12 23.070 13.662 37.422 1.00 17.5 ATOM 77 CG1 VAL 12 23.070 13.077 36.839 1.00 15.01 ATOM 78 CG2 VAL 12 23.270 13.077 36.839 1.00 15.01 ATOM 79 N PRO 13 21.180 13.382 34.066 1.00 14.75 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 16.55 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 82 O PRO 13 19.797 15.468 15.337 1.00 15.14 ATOM 85 CD PRO 13 19.255 13.070 32.745 1.00 17.64 ATOM 85 CD PRO 13 19.255 13.070 32.745 1.00 17.64 ATOM 85 CD PRO 13 19.255 13.070 32.745 1.00 17.64 ATOM 86 K ILE 14 19.893 17.991 34.266 1.00 14.74 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 88 C ILE 14 18.768 13.769 33.550 1.00 15.41 ATOM 89 O ILE 14 18.768 13.769 33.550 1.00 12.93 ATOM 89 O ILE 14 18.768 13.769 33.550 1.00 12.93 ATOM 90 CS ILE 14 18.768 13.769 33.550 1.00 12.93 ATOM 90 CS ILE 14 18.768 13.769 33.550 1.00 12.93 ATOM 90 CS ILE 14 18.768 13.769 33.550 1.00 12.93 ATOM 90 CS ILE 14 22.271 18.169 35.015 1.00 16.08 ATOM 90 CS ILE 14 20.195 16.766 33.557 1.00 16.04 ATOM 90 CS ILE 14 20.195 16.760 33.556 1.00 12.93 ATOM 90 CS ILE 14 20.196 18.769 33.595 1.00 16.04 ATOM 90 CS ILE 14 20.196 18.769 33.595 1.00 16.05 ATOM 90 CS ILE 14 20.196 18.769 33.595 1.00 16.05 ATOM 90 CS ILE 14 20.196 18.769 33.595 1.00 16.05 ATOM 90 CS ILE 14 20.793 18.169 35.015 1.00 16.05 ATOM 90 CS ILE 14 20.793 18.169 30.00 33.955 1.00 16.05 ATOM 90 CS ILE 14 20.793 18.169 30.00 33.955 1.00 16.05 ATOM 90 CS ILE 15 16.691 20.196 33.999 1.00 10.05 ATOM 90 CS ILE 15 1.00 16.694 30.00 33.999 1.00 10.05 ATOM 90 CS ILE 15 1.00 80.80 11.90 11.90 11.90 11.90 11.90 11.90 11.90 11.90 11.90 11.90									1.00 29.60
ATOM 70 CG1 VAL 11 25.912 11.980 31.734 1.00 31.75 ATOM 71 CG2 VAL 11 27.480 10.048 32.051 1.00 33.85 ATOM 72 N VAL 12 24.417 12.227 34.931 1.00 20.12 ATOM 73 CA VAL 12 22.080 12.561 35.403 1.00 12.88 ATOM 74 C VAL 12 22.080 12.561 35.403 1.00 12.88 ATOM 75 O VAL 12 22.007 13.624 34.516 1.00 14.37 ATOM 75 O VAL 12 22.007 13.624 34.516 1.00 13.42 ATOM 77 CG1 VAL 12 23.007 14.639 34.179 1.00 13.42 ATOM 77 CG1 VAL 12 23.007 13.639 1.00 15.01 ATOM 77 CG1 VAL 12 23.007 13.639 1.00 15.01 ATOM 78 CG2 VAL 12 23.001 13.662 37.422 1.00 17.57 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 15.51 ATOM 79 N PRO 13 20.180 13.382 34.066 1.00 14.72 ATOM 80 CA PRO 13 20.493 14.182 33.265 1.00 10.76 ATOM 81 C PRO 13 20.186 13.382 34.066 1.00 14.72 ATOM 80 CA PRO 13 20.186 13.382 34.066 1.00 10.76 ATOM 82 O PRO 13 19.797 15.468 35.337 12.00 15.14 ATOM 85 CD PRO 13 19.255 13.707 32.745 1.00 17.36 ATOM 84 CG PRO 13 19.255 13.707 32.745 1.00 17.36 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 86 CT LE 14 20.196 16.766 33.557 1.00 19.69 ATOM 86 CT LE 14 19.893 17.91 34.266 1.00 12.93 ATOM 86 CT LE 14 18.768 13.750 33.596 1.00 12.93 ATOM 88 C LE 14 19.893 17.91 34.266 1.00 12.93 ATOM 89 O ILE 14 20.196 16.766 33.557 1.00 15.41 ATOM 80 CT LE 14 18.768 13.750 33.596 1.00 12.08 ATOM 89 O ILE 14 22.271 18.189 35.015 1.00 16.54 ATOM 90 CS ILE 14 20.783 20.207 35.084 1.00 11.04 ATOM 90 CS ILE 14 20.783 20.207 35.084 1.00 11.05 ATOM 90 CS ILE 14 20.783 20.207 35.084 1.00 11.05 ATOM 90 CD ILE 14 20.783 20.207 35.084 1.00 11.05 ATOM 90 CD ILE 14 20.783 20.207 35.084 1.00 11.05 ATOM 90 CD ILE 14 20.783 20.207 35.084 1.00 11.05 ATOM 90 CD ILE 14 20.783 20.207 35.084 1.00 11.05 ATOM 90 CD ILE 15 16.811 20.136 33.955 1.00 16.54 ATOM 90 CD ILE 15 16.811 20.136 33.955 1.00 16.54 ATOM 90 CD ILE 14 20.783 20.207 35.084 1.00 11.05 ATOM 90 CD ILE 15 15.086 20.207 35.084 1.00 11.05 ATOM 90 CD ILE 15 15.085 20.207 35.084 1.00 11.56 ATOM 90 CD ILE 15 15.085 20.207 35.084 1.00 11.85 ATOM 90 CD ILE 15 15.085 20.207 35.901 1.00 10.56 ATOM 90 CD ILE 15 1									00 20.43
ATOM 71 CG2 VAL 11 27.480 10.048 32.015 1.00 33.85 ATOM 72 N VAL 12 24.471 12.227 34.931 1.00 20.12 ATOM 73 CA VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 74 C VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 75 O VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 75 O VAL 12 23.070 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 23.070 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 23.070 13.662 37.422 1.00 17.57 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 15.51 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 15.51 ATOM 78 N PRO 13 21.180 13.382 34.066 1.00 14.72 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 20.161 15.589 34.141 1.00 7.65 ATOM 82 O PRO 13 20.161 15.589 34.141 1.00 7.65 ATOM 82 O PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 84 CG PRO 13 19.043 12.422 33.550 1.00 19.54 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.54 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 19.54 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 19.54 ATOM 87 CA TLE 14 19.893 17.991 34.266 1.00 12.08 ATOM 89 O LE 14 18.768 13.760 33.596 1.00 12.08 ATOM 89 O LE 14 18.768 13.760 33.596 1.00 12.08 ATOM 91 CG1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 92 CG2 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 93 CD1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 93 CD1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.54 ATOM 95 CA LEU 15 17.090 18.930 34.265 1.00 16.54 ATOM 95 CA LEU 15 17.090 18.930 34.265 1.00 16.54 ATOM 96 C LEU 15 17.090 18.000 33.696 1.00 12.08 ATOM 97 C LEU 15 17.090 18.000 33.696 1.00 12.08 ATOM 95 CA LEU 15 17.090 18.000 33.696 1.00 12.08 ATOM 96 C LEU 15 17.090 18.000 33.696 1.00 12.08 ATOM 96 C LEU 15 17.090 18.000 33.696 1.00 12.05 ATOM 96 C LEU 15 17.090 18.000 33.695 1.00 16.54 ATOM 96 C LEU 15 17.090 18.000 33.695 1.00 16.54 ATOM 100 CD1 LEU 15 17.090 18.000 33.695 1.00 16.54 ATOM 100 CD1 LEU 15 17.090 18.000 33.695 1.00 10.05 ATOM 100 CD1 LEU 15 17.090 18.000 33.695 1.00 10.05 ATOM 100 CD1 LEU 15 17.090 28.575 33.995 1.00 10.23 AT	MOTA	69	CB	VAL					
ATOM 73 CA VAL 12 24.417 12.227 34.931 1.00 20.12 88 ATOM 73 CA VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 74 C VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 75 CO VAL 12 23.007 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 23.007 14.639 34.179 1.00 15.01 ATOM 76 CB VAL 12 23.270 13.077 36.839 1.00 15.01 ATOM 77 CGI VAL 12 23.200 13.662 37.422 1.00 17.57 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 79 N PRO 13 20.188 13.382 34.066 1.00 14.72 ATOM 80 CA PRO 13 20.186 13.382 34.066 1.00 14.72 ATOM 81 C PRO 13 20.116 15.589 34.141 1.00 7.65 ATOM 82 C PRO 13 20.116 15.589 34.141 1.00 7.65 ATOM 82 C PRO 13 19.225 13.707 32.745 1.00 17.36 ATOM 82 C PRO 13 19.245 13.707 32.745 1.00 17.36 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 85 CD PRO 13 20.315 12.195 34.340 1.00 15.41 ATOM 86 N LEE 14 19.893 17.91 34.266 1.00 12.93 ATOM 86 N LEE 14 19.893 17.91 34.266 1.00 12.93 ATOM 88 C LEE 14 18.768 13.570 33.596 1.00 12.93 ATOM 80 C CB LEE 14 18.768 18.760 33.597 1.00 12.93 ATOM 80 C CB LEE 14 18.768 18.760 33.595 1.00 16.54 ATOM 90 CB LEE 14 18.768 18.760 33.595 1.00 16.54 ATOM 90 CB LEE 14 18.768 18.760 33.595 1.00 16.54 ATOM 90 CB LEE 14 22.109 18.993 37.91 34.266 1.00 12.93 ATOM 90 CB LEE 14 22.109 18.993 37.91 34.256 1.00 12.93 ATOM 90 CB LEE 14 22.109 18.905 34.325 1.00 16.54 ATOM 90 CB LEE 14 22.109 18.905 34.325 1.00 16.54 ATOM 90 CB LEE 14 22.109 18.905 34.325 1.00 16.54 ATOM 90 CB LEE 14 22.109 18.905 34.325 1.00 16.54 ATOM 90 CB LEE 14 22.109 18.905 34.325 1.00 16.54 ATOM 90 CB LEE 14 22.404 33.605 34.325 1.00 16.54 ATOM 90 CB LEE 14 22.404 33.605 34.935 1.00 16.54 ATOM 90 CB LEE 14 22.404 33.605 34.935 1.00 16.55 ATOM 90 CB LEE 15 16.811 20.136 33.955 1.00 16.54 ATOM 90 CB LEE 15 16.811 20.136 33.955 1.00 16.54 ATOM 90 CB LEE 15 16.811 20.136 33.955 1.00 16.55 ATOM 90 CB LEE 15 16.811 20.136 33.955 1.00 16.55 ATOM 90 CB LEE 15 16.812 20.136 33.955 1.00 16.55 ATOM 90 CB LEE 15 16.812 20.136 33.955 1.00 16.55 ATOM 90 CB LEE 15 16.812 20.136 33.955 1.00 10.05 ATOM 90 CB LEE 15 16.812 20.137 34.406	ATOM	70	CG1	VAL	11	25.912	11.980	31.734	1.00 31.75
ATOM 73 CA VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 74 C VAL 12 22.407 13.624 34.516 1.00 14.57 ATOM 75 O VAL 12 23.007 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 23.007 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 23.007 13.662 37.422 1.00 17.57 ATOM 77 CG1 VAL 12 22.000 13.662 37.422 1.00 17.57 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 85 CD PRO 13 19.043 12.222 33.550 1.00 19.69 ATOM 85 CD PRO 13 19.043 12.222 33.550 1.00 19.69 ATOM 85 CD PRO 13 19.043 12.222 33.550 1.00 19.69 ATOM 86 N ILE 14 20.196 16.766 33.575 1.00 15.49 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.08 ATOM 89 O ILE 14 18.768 18.760 33.596 1.00 12.08 ATOM 89 O ILE 14 18.768 18.760 33.596 1.00 12.08 ATOM 90 CS ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 90 CS ILE 14 22.271 18.169 35.015 1.00 16.58 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 99 CD ILE 14 20.183 20.305 34.325 1.00 16.54 ATOM 99 CD ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 99 CD ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 99 CD ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 99 CD ILE 15 17.899 19.307 34.421 1.00 13.85 ATOM 99 CD ILE 15 17.899 19.307 34.421 1.00 13.85 ATOM 99 CD ILE 15 17.899 19.307 34.421 1.00 13.85 ATOM 99 CD ILE 15 17.899 19.307 34.421 1.00 13.85 ATOM 99 CD ILEU 15 17.800 21.599 35.901 1.00 40.50 ATOM 99 CD ILEU 15 17.800 21.599 35.901 1.00 40.50 ATOM 99 CD ILEU 15 17.800 21.599 35.901 1.00 40.50 ATOM 99 CD ILEU 15 17.800 21.599 35.901 1.00 40.50 ATOM 100 CDI LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 17.000 27.27 ATOM 104 CD VAL 16 18.226 26.075 34.497 1.00 10.56 ATOM 105 CD VAL 16 18.226 26.075 34.497 1.00 10.57 ATOM 106 CB VAL 16 18.226 26.075 34.497 1.00 10.51 ATOM 1	MOTA	71	CG2	VAL	11	27.480	10.048		
ATOM 73 CA VAL 12 23.080 12.561 35.433 1.00 12.88 ATOM 74 C VAL 12 22.407 13.624 34.179 1.00 14.37 ATOM 75 O VAL 12 23.007 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 23.007 14.639 34.179 1.00 15.01 ATOM 77 CG1 VAL 12 22.000 13.662 37.422 1.00 17.57 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 82 O PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.25 13.707 32.745 1.00 17.36 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 17.36 ATOM 85 CD PRO 13 20.315 12.195 34.340 1.00 15.41 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 17.36 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.08 ATOM 89 O ILE 14 18.768 13.760 33.596 1.00 12.08 ATOM 89 O ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 90 CS ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 90 CS ILE 14 22.271 18.169 35.015 1.00 16.58 ATOM 90 CS ILE 14 22.271 18.169 35.015 1.00 16.58 ATOM 90 CS ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 90 CS ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 99 CD ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 99 CD ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 99 CD ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 99 CD ILE 15 17.899 19.307 34.421 1.00 13.85 ATOM 99 CD ILE 15 16.915 21.474 34.685 1.00 16.15 ATOM 99 CD ILE 15 17.899 19.307 34.421 1.00 13.85 ATOM 99 CD ILE 15 16.915 21.474 34.685 1.00 14.82 ATOM 99 CD ILEU 15 15.008 20.998 31.913 1.00 49.22 ATOM 100 CDI LEU 15 15.008 20.998 31.913 1.00 49.22 ATOM 101 CD2 LU 15 17.080 21.599 35.901 1.00 10.56 ATOM 99 CD ILEU 15 17.080 21.599 35.901 1.00 10.56 ATOM 103 CA VAL 16 16.964 33.905 34.479 1.00 10.56 ATOM 104 C VAL 16 16.964 33.905 34.479 1.00 10.57 ATOM 106 CD UAL 16 18.226 20.975 34.691 1.00 10.57 ATOM 107 CGI VAL 16 18.226 20.975 34.691 1.00 10.57 ATOM 108 CG2 VAL 16 19.520 23.965 34.691		72	N	VAL	12	24.417	12.227	34.931	1.00 20.12
ATOM 74 C VAL 12 22.407 13.624 34.516 1.00 14.37 ATOM 75 O VAL 12 23.007 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 23.270 13.077 36.839 1.00 15.01 ATOM 77 CG1 VAL 12 22.000 13.663 37.422 1.00 17.57 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 17.57 ATOM 79 N PRO 13 20.180 13.382 34.066 1.00 14.72 ATOM 80 CA PRO 13 20.493 14.182 33.265 1.00 10.76 ATOM 81 C PRO 13 20.166 15.589 34.141 1.00 7.65 ATOM 82 O PRO 13 19.797 15.468 35.337 1.00 15.76 ATOM 82 O PRO 13 19.797 15.468 35.337 1.00 15.41 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.41 ATOM 84 CG PRO 13 19.255 13.707 32.745 1.00 17.36 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 19.69 ATOM 86 N ILE 14 18.768 13.760 33.596 1.00 12.93 ATOM 88 C ILE 14 18.768 13.760 33.596 1.00 12.93 ATOM 88 C ILE 14 18.768 13.760 33.596 1.00 12.93 ATOM 89 O ILE 14 18.768 13.750 33.596 1.00 12.63 ATOM 91 CG1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 93 CD1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 93 CD1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.55 ATOM 94 N LEU 15 17.899 19.007 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 16.54 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 13.85 ATOM 99 CG LEU 15 16.815 21.474 34.685 1.00 13.85 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 10.06 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 10.56 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 10.56 ATOM 90 CG LEU 15 15.462 19.450 34.285 1.00 10.56 ATOM 90 CG LEU 15 15.462 19.450 34.285 1.00 10.56 ATOM 100 CD1 LEU 15 15.279 20.440 33.679 1.00 10.25 ATOM 100 CD1 LEU 15 15.279 20.440 33.679 1.00 10.25 ATOM 100 CD1 LEU 15 15.008 20.098 31.919 1.00 10.26 ATOM 100 CD1 LEU 15 15.008 20.098 31.919 1.00 10.26 ATOM 100 CD1 LEU 15 15.008 20.098 31.919 1.00 10.26 ATOM 100 CD1 LEU 15 15.008 20.098 31.919 1.00 10.26 ATOM 100 CD1 LEU 15 15.008 20.098 31.919 1.00 10.26 ATOM 100 CD1 LEU 15 15.008 20.098 31.919 1.00 10.26 ATOM 100 CD1 LEU 15 15.008 20.098 31.919 1.00 10.26 ATOM 100 CD1 LEU 15 15.008 20.098 31		73	CA	VAL	12	23.080	12.561	35.433	1.00 12.88
ATOM 75 O VAL 12 23.007 14.639 34.179 1.00 13.42 ATOM 76 CB VAL 12 23.270 13.077 36.68.99 1.00 15.01 ATOM 77 CG1 VAL 12 23.270 13.662 37.422 1.00 17.57 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 79 N PRO 13 21.180 13.382 34.066 1.00 14.77 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 14.77 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 14.77 ATOM 82 O PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 86 R ILE 14 20.196 16.766 33.557 1.00 14.91 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.08 ATOM 88 C ILE 14 18.768 18.760 33.596 1.00 12.08 ATOM 89 O ILE 14 18.768 18.760 33.596 1.00 12.08 ATOM 90 CS ILE 14 20.783 20.207 35.084 1.00 16.54 ATOM 91 CG1 ILE 14 20.783 20.207 35.084 1.00 16.54 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 16.15 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 16.08 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 16.08 ATOM 97 O LEU 15 16.811 20.136 33.955 1.00 16.15 ATOM 98 C LEU 15 16.811 20.136 33.955 1.00 16.15 ATOM 99 CG LEU 15 16.811 20.136 33.955 1.00 16.08 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 16.15 ATOM 90 CG LEU 15 16.811 20.136 33.955 1.00 16.08 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 16.15 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 16.54 ATOM 90 CG LEU 15 15.008 20.098 31.913 1.00 40.50 ATOM 100 CD1 LEU 15 15.008 20.098 31.913 1.00 40.50 ATOM 100 CD1 LEU 15 15.008 20.098 31.913 1.00 40.50 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 10.25 ATOM 102 CG LEU 17 13.904 26.144 34.870 1.00 16.54 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.25 ATOM 104 C VAL 16 18.873 24.668 34.098 1.00 11.55 ATOM 105 CG LU 17 13.904 26.147 34.685 1.00 11.55 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 11.53 ATOM 108 CG LU 17 13.904 26.147 33.169 1.00 10.55 ATOM 109 CG LEU 15 15.008 20.098 31.917 1.00 10.55 ATOM 101 CG2 UAL 16 18.273 24.668 34.098 1.00 11.55 ATOM 102 CG LU 17 13.904 26.440 33.679			•			22.407	13.624	34.516	1.00 14.37
ATOM 76 CB VAL 12 23.270 13.077 36.839 1.00 15.01 ATOM 77 CG1 VAL 12 22.000 13.662 37.422 1.00 17.57 ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 79 N PRO 13 21.180 13.382 34.066 1.00 14.75 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 82 CA PRO 13 20.116 15.589 34.141 1.00 7.65 ATOM 82 CA PRO 13 19.797 15.468 15.337 1.00 15.14 ATOM 83 CB PRO 13 19.251 13.707 32.745 1.00 17.36 ATOM 84 CG PRO 13 19.043 12.422 33.550 1.00 17.36 ATOM 85 CD PRO 13 20.315 12.195 34.340 1.00 15.41 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 15.41 ATOM 86 N ILE 14 18.768 18.760 33.557 1.00 12.93 ATOM 88 C ILE 14 18.768 18.760 33.557 1.00 12.93 ATOM 88 C ILE 14 18.768 18.760 33.559 1.00 11.04 ATOM 90 CB ILE 14 18.768 18.760 33.559 1.00 11.04 ATOM 91 CG1 ILE 14 22.271 18.169 25.015 1.00 11.54 ATOM 91 CG1 ILE 14 22.271 18.169 25.015 1.00 16.54 ATOM 93 CD1 ILE 14 22.271 18.169 25.015 1.00 16.54 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.55 ATOM 94 N LEU 15 77.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.915 21.474 34.685 1.00 13.85 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 13.85 ATOM 99 CG LEU 15 16.811 20.136 33.955 1.00 16.56 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 13.85 ATOM 99 CG LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 15.462 39.450 34.225 1.00 16.56 ATOM 99 CG LEU 15 15.462 39.450 34.225 1.00 16.56 ATOM 99 CG LEU 15 15.462 39.450 34.225 1.00 16.56 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 99 CG LEU 15 15.462 39.450 34.225 1.00 14.82 ATOM 100 CD1 LEU 15 17.080 21.599 35.901 1.00 10.56 ATOM 100 CD1 LEU 15 17.080 21.599 35.901 1.00 10.56 ATOM 100 CD1 LEU 15 17.080 21.599 35.901 1.00 10.56 ATOM 100 CD1 LEU 15 18.296 30.490 31.919 1.00 10.56 ATOM 100 CD1 LEU 15 18.296 30.890 31.919 1.00 10.56 ATOM 100 CD1 LEU 15 18.296 30.890 31.919 1.00 10.56 ATOM 100 CD1 LEU 15 18.296 30.890 31.479				VAL			14.639	34.179	
ATOH 78 CG1 VAL 12 22.000 13.662 37.422 1.00 17.57 ATOH 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOH 79 N PRO 13 21.180 13.382 34.066 1.00 14.72 ATOH 80 CA PRO 13 20.493 14.382 34.066 1.00 14.72 ATOH 80 CA PRO 13 20.493 14.382 33.265 1.00 10.765 ATOH 82 O PRO 13 19.797 15.468 35.337 1.00 15.14 ATOH 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOH 83 CB PRO 13 19.225 13.707 32.745 1.00 15.14 ATOH 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOH 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOH 86 N ILE 14 20.196 16.766 33.557 1.00 14.91 ATOH 87 CA 1LE 14 19.893 17.991 34.266 1.00 12.08 ATOH 88 C ILE 14 18.768 18.760 33.555 1.00 16.491 ATOH 88 C ILE 14 18.768 18.760 33.556 1.00 12.08 ATOH 90 CS ILE 14 22.109 18.905 34.325 1.00 16.08 ATOH 90 CS ILE 14 22.109 18.905 34.325 1.00 16.08 ATOH 90 CS ILE 14 22.109 18.905 34.325 1.00 16.08 ATOH 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.04 ATOH 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.64 ATOH 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.65 ATOH 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOH 96 C LEU 15 16.811 20.136 33.955 1.00 16.15 ATOH 97 C LEU 15 16.811 20.136 33.955 1.00 16.08 ATOH 99 CG LEU 15 16.811 20.136 33.955 1.00 14.04 ATOH 99 CG LEU 15 16.811 20.136 33.955 1.00 14.02 ATOH 99 CG LEU 15 16.811 20.136 33.955 1.00 14.02 ATOH 99 CG LEU 15 16.811 20.136 33.955 1.00 14.02 ATOH 99 CG LEU 15 16.811 20.136 33.955 1.00 14.02 ATOH 99 CG LEU 15 15.462 19.450 34.225 1.00 21.25 ATOH 99 CG LEU 15 15.462 19.450 34.225 1.00 21.25 ATOH 99 CG LEU 15 15.462 19.450 34.225 1.00 16.08 ATOH 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOH 102 CD LEU 15 15.008 20.098 31.913 1.00 40.50 ATOH 90 CG LEU 15 15.008 20.098 31.913 1.00 40.50 ATOH 90 CG LEU 15 15.008 20.098 31.913 1.00 10.56 ATOH 101 CD2 LEU 15 15.008 20.098 31.913 1.00 40.50 ATOH 102 CD LEU 15 15.008 20.098 31.913 1.00 40.50 ATOH 103 CA VAL 16 16.865 22.556 33.919 1.00 40.50 ATOH 104 C VAL 16 15.347 24.748 22.904 1.00 16.72 ATOH 105 C VAL 16 18.226 26.075 34.40 33.679 1.00 19.38 ATOH 107 CG LVAL 16 18.226 26.07									1.00 15.01
ATOM 78 CG2 VAL 12 23.781 11.936 37.728 1.00 16.55 ATOM 79 N PRO 13 21.180 13.382 34.066 1.00 14.72 ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 20.116 15.589 34.141 1.00 7.56 ATOM 82 CO PRO 13 19.797 15.468 35.37 1.00 15.14 ATOM 83 CB PRO 13 19.255 13.707 32.745 1.00 17.36 ATOM 84 CG PRO 13 19.043 12.422 33.550 1.00 19.641 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.641 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 15.44 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 12.93 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 89 C ILE 14 18.768 18.760 33.596 1.00 12.08 ATOM 89 O ILE 14 18.768 18.760 33.596 1.00 12.08 ATOM 90 CS ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 91 CG1 ILE 14 20.271 18.169 25.015 1.00 16.54 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 11.56 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 16.15 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 17.080 21.250 35.901 1.00 10.00 ATOM 99 CG LEU 15 13.279 20.440 33.679 1.00 10.56 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 13.379 20.440 33.679 1.00 46.97 ATOM 105 CVAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 107 CG1 VAL 16 18.266 27.474 35.571 1.00 10.53 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 107 CG1 VAL 16 18.266 27.474 35.571 1.00 10.53 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 107 CG1 VAL 16 18.266 26.075 34.468 1.00 14.24 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 107 CG1 VAL 16 18.274 24.748 32.904 1.00 16.57 ATOM 107 CG1 VAL 16 18.266 26.075 34.468 3.00 9.00 12.58 ATOM 107 CG1 VAL 16									
ATOM 79 N PRO 13 20.493 14.382 34.066 1.00 14.72 ATOM 81 C PRO 13 20.493 14.382 34.066 1.00 10.76 ATOM 81 C PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 82 O PRO 13 19.797 15.468 35.337 1.00 17.36 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 17.36 ATOM 84 CG PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 15.41 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 88 C ILE 14 18.768 18.760 33.557 1.00 14.91 ATOM 89 C ILE 14 18.768 18.760 33.595 1.00 10.293 ATOM 89 O ILE 14 18.768 18.760 33.595 1.00 11.04 ATOM 90 CS ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 90 CS ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 91 CGI ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CDI ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 16.54 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 16.54 ATOM 97 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 16.812 20.399 31.911 1.00 40.50 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 16.54 ATOM 99 CG LEU 15 15.462 19.450 34.255 1.00 16.54 ATOM 99 CG LEU 15 15.462 19.450 34.255 1.00 16.52 ATOM 99 CG LEU 15 15.462 19.450 34.255 1.00 10.00 ATOM 99 CG LEU 15 15.462 19.450 34.255 1.00 10.23 ATOM 100 CD1 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 101 CD2 LEU 15 15.367 24.744 33.699 1.00 10.23 ATOM 104 C VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 105 O VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 108 CG2 VAL 16 18.273 24.768 34.691 1.00 15.74 ATOM 110 CG LU 17 14.331 27.524 36.755 1.00 13.61 ATOM 110 CG LU 17 14.331 27.524 36.755 1.00 13.61 ATOM 111 C GLU 17 14.331 27.524 36.755 1.00 13.61 ATOM 112 C GLU 17 12.250 23.965 34.474 35.571 1.00 13.61 ATOM 114 CG GLU 17 12.250 23.965 34.499 1.00 12.58 ATOM 115 CD GLU 17 12.250 23.766 32.									
ATOM 80 CA PRO 13 20.493 14.382 33.265 1.00 10.76 ATOM 81 C PRO 13 20.116 15.589 34.141 1.00 7.65 ATOM 82 O PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 84 CG PRO 13 19.255 13.707 32.745 1.00 17.36 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 14.91 ATOM 86 N ILE 14 18.831 17.991 34.266 1.00 12.93 ATOM 87 CA ILE 14 18.768 13.760 33.596 1.00 12.08 ATOM 89 O ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 90 CB ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 91 CGI ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 92 CGI ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 93 CDI ILE 14 23.642 18.836 34.738 1.00 11.56 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.915 21.474 34.685 1.00 13.55 ATOM 97 O LEU 15 16.915 21.474 34.685 1.00 13.62 ATOM 99 CG LEU 15 16.915 21.474 34.685 1.00 13.62 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 14.82 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 12.25 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 12.25 ATOM 100 CDI LEU 15 15.279 20.440 33.679 1.00 40.50 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 10.00 ATOM 100 CDI LEU 15 15.008 20.098 31.913 1.00 40.50 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 40.50 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 104 C VAL 16 15.347 24.778 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.778 34.063 1.00 9.47 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.25 ATOM 107 CGI VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.500 25.402 35.344 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.1478 34.407 1.00 15.37 ATOM 101 CA GLU 17 14.086 27.474 35.571 1.00 14.91 ATOM 110 CA GLU 17 14.086 27.474 35.571 1.00 14.91 ATOM 111 C GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 112 C GLU 17 11.602 25.557 32.990 1.00 38.75 ATOM 114 CG GLU 17 12.200 23.766 32.127 1.00 14.91 ATOM 115 CD GLU 17 12.200 23.766 32.127 1.00 14.91 ATOM 110 CA GLU 17 12.200 23.766 32.127 1.00 14.91 ATOM 120 C LEU 18 12.956 30.998 33.									
ATOM 82 C PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.225 13.707 32.745 1.00 17.36 ATOM 84 CG PRO 13 19.225 13.707 32.745 1.00 17.36 ATOM 85 CD PRO 13 19.225 13.707 32.745 1.00 17.36 ATOM 86 CB PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 85 CD PRO 13 20.315 12.195 34.340 1.00 15.41 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 14.91 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 88 C ILE 14 18.768 18.760 33.595 1.00 12.08 ATOM 89 O ILE 14 18.768 18.760 33.595 1.00 12.08 ATOM 89 O ILE 14 18.768 18.760 33.595 1.00 12.08 ATOM 90 CB ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 91 CG1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 92 CG2 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 93 CD1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 11.56 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 O LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 33.679 1.00 40.50 ATOM 99 CB LEU 15 15.088 20.098 31.913 1.00 49.22 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 CD VAL 16 16.885 22.556 33.919 1.00 10.23 ATOM 103 CA VAL 16 16.885 22.556 33.919 1.00 10.23 ATOM 104 C VAL 16 16.885 22.556 33.919 1.00 10.23 ATOM 105 CD VAL 16 18.226 26.075 34.681 1.00 12.58 ATOM 107 CG1 VAL 16 18.226 26.075 34.681 1.00 12.58 ATOM 108 CG2 VAL 16 15.3147 24.748 32.904 1.00 15.47 ATOM 105 CD VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 18.250 23.9945 34.692 1.00 14.24 ATOM 110 CA GLU 17 14.086 27.474 35.571 1.00 12.38 ATOM 115 CD GLU 17 14.086 27.474 35.571 1.00 15.37 ATOM 116 OE1 CLU 17 14.08									
ATOM 82 O PRO 13 19.797 15.468 35.337 1.00 15.14 ATOM 83 CB PRO 13 19.225 13.707 32.745 1.00 17.36 ATOM 84 CG PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 85 CD PRO 13 20.315 12.195 34.340 1.00 15.41 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 14.91 ATOM 87 CA LE 14 19.893 17.991 34.266 1.00 12.93 ATOM 88 C ILE 14 18.768 18.760 33.556 1.00 12.93 ATOM 89 O ILE 14 18.768 18.760 33.556 1.00 12.08 ATOM 90 CB ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 90 CB ILE 14 22.1109 18.905 34.325 1.00 16.54 ATOM 91 CGI ILE 14 22.711 18.169 35.015 1.00 18.08 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CDI ILE 14 23.642 18.836 34.738 1.00 11.56 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 2.62 ATOM 97 O LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 10.00 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 10.00 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 10.00 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 2.62 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 101 CD2 LEU 15 15.262 19.450 34.285 1.00 2.02 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 40.50 ATOM 103 CA VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 104 C VAL 16 15.347 24.748 32.904 1.00 10.23 ATOM 105 O VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 106 CB VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 107 CGI VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 13.904 26.144 34.870 1.00 14.74 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 16.72 ATOM 101 CA GLU 17 13.904 26.144 34.870 1.00 14.53 ATOM 110 CA GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 110 CA GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 110 CA GLU 17 12.650 25.577 32.950 1.00 83.46 ATOM 110 CA GLU 17 12.650 25.577 32.950 1.00 83.46 ATOM 110 CA GLU 17 12.220 23.766 83.100 1.00 14.51 ATOM 110 CA GLU 17 12.220 23.766 83 1.00 14.53 ATOM 120 C LEU 18 12.952 30.855 71 34.805 1.00 14.53 ATOM 120 C LEU 18 12.952 30.855 71 3									
ATOM 84 CG PRO 13 19.225 13.707 32.745 1.00 17.36 ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 85 CD PRO 13 20.315 12.195 34.340 1.00 15.41 ATOM 86 % ILE 14 20.196 16.766 33.557 1.00 14.91 ATOM 86 % ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 88 C ILE 14 18.768 13.760 33.596 1.00 12.08 ATOM 89 O ILE 14 18.768 13.760 33.596 1.00 12.08 ATOM 89 O ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 90 CS ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 91 CG1 ILE 14 22.771 18.169 35.015 1.00 16.54 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 C LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 21.25 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 99 CG LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 10.56 ATOM 106 CB VAL 16 18.226 26.075 34.691 1.00 12.25 ATOM 106 CB VAL 16 18.226 26.075 34.691 1.00 12.25 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 15.74 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 15.37 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 11.36 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 11.36 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB CLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB CLU 17 14.331 27.524 36.765 1.00 15.37 ATOM 114 CC CLU 17 14.331 27.524 36.765 1.00 15.37 ATOM 115 CD CLU 17 14.331 27.524 36.765 1.00 15.37 ATOM 114 CC CLU 17 14.331 27.524 36.765 1.00 15.37 ATOM 114 CC CLU 17 14.331 27.52									
ATOM 85 CD PRO 13 19.043 12.422 33.550 1.00 19.69 ATOM 85 CD PRO 13 20.315 12.195 34.340 1.00 15.41 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 14.91 ATOM 87 CA LLE 14 19.893 17.991 34.266 1.00 12.08 ATOM 89 O ILE 14 18.768 13.769 13.00 11.04 ATOM 90 CB ILE 14 18.768 13.769 13.00 11.04 ATOM 90 CB ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 91 CG1 ILE 14 22.271 18.169 25.015 1.00 18.08 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 16.54 ATOM 97 O LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 15.442 19.541 33.199 1.00 40.50 ATOM 99 CG LEU 15 15.442 19.541 33.199 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 104 C VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 105 O VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 106 CB VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 109 N GLU 17 13.904 26.644 34.870 1.00 15.74 ATOM 100 CA CLU 17 14.331 27.524 36.765 1.00 13.67 ATOM 101 CA CLU 17 14.331 27.524 36.765 1.00 14.24 ATOM 109 N GLU 17 13.904 26.144 34.870 1.00 15.74 ATOM 101 CA CLU 17 14.331 27.524 36.765 1.00 13.67 ATOM 114 CG CLU 17 14.331 27.524 36.765 1.00 13.67 ATOM 115 CD GLU 17 14.331 27.524 36.765 1.00 13.67 ATOM 110 CA LEU 15 13.890 28.571 34.4805 1.00 15.74 ATOM 110 CA LEU 18 13.990 28.571 34.4805 1.00 15.74 ATOM 110 CA LEU 18 13.990 28.571 34.805 1.00 13.61 ATOM 110 CA LEU 18 13.990 28.571 34.805 1.00 13.61 ATOM 110 CA LEU 18 13.990 28.571 34.805 1.00 13.61 ATOM 110 CA LEU 18 13.990 28.571 34.805 1.00 13.61 ATOM 110 CA LEU 18 13.990 28.571 34.805 1.00 14.24 ATOM 110 CA LEU 18 13.990 28.571 34.805 1.00 14.15 ATOM 110 CA LEU 18 13.990 28.571 34.805 1.00 14.91 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 14.91 ATOM 121 CA LEU 18 13.990 28.571 34.8	ATOM								1.00 15.14
ATOM 85 CD PRO 13 20.315 12.195 34.34000 15.41 ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 14.91 ATOM 86 N ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 88 C ILE 14 18.768 18.760 33.596 1.00 12.93 ATOM 89 O ILE 14 18.768 18.760 33.596 1.00 12.08 ATOM 90 CB ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 91 CG1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 91 CG1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 11.56 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 16.25 ATOM 99 CB LEU 15 16.915 21.474 34.685 1.00 10.00 ATOM 99 CB LEU 15 15.462 19.450 34.225 1.00 21.25 ATOM 99 CB LEU 15 14.412 19.541 33.199 1.00 40.50 ATOM 99 CB LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 90 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 40.50 ATOM 103 CA VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.896 23.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 15.347 24.748 32.904 1.00 10.23 ATOM 104 C VAL 16 15.347 24.748 32.904 1.00 10.52 ATOM 105 CB VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 106 CB VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 110 CA CLU 17 13.904 26.144 34.870 1.00 15.37 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 14.43 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 14.34 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 25.344 1.00 14.53 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 15.37 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 14.31 ATOM 110 CA CLU 17 14.331 27.524 36.765 1.00 14.31 ATOM 110 CA CLU 18 13.990 28.571 34.695 1.00 14.32 ATOM 110 CA CLU 18 13.9	MOTA	83							
ATOM 86 N ILE 14 20.196 16.766 33.557 1.00 14.91 ATOM 87 CA ILE 14 19.893 17.991 34.266 1.00 12.93 ATOM 88 C ILE 14 18.768 13.760 33.596 1.00 12.08 ATOM 89 O ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 90 CB ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 91 CG1 ILE 14 22.271 18.169 35.015 1.00 18.08 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 O LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 CG LEU 15 15.462 19.450 34.225 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 21.25 ATOM 99 CG LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.25 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.25 ATOM 105 O VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 106 CB VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 107 CG1 VAL 16 15.2147 24.688 34.098 1.00 10.25 ATOM 108 CG2 VAL 16 15.226 26.075 34.691 1.00 12.25 ATOM 108 CG2 VAL 16 15.202 23.945 34.628 1.00 12.25 ATOM 109 N GLU 17 13.904 26.144 34.870 1.00 15.37 ATOM 101 CG CGU 17 14.331 27.524 36.765 1.00 13.61 ATOM 111 C GGU 17 14.331 27.524 36.765 1.00 13.61 ATOM 112 O GLU 17 12.436 24.778 32.990 1.00 15.37 ATOM 113 CB GLU 17 12.436 24.573 33.105 1.00 15.37 ATOM 114 CG GLU 17 12.436 24.573 33.105 1.00 15.37 ATOM 115 CD GLU 17 12.436 24.573 33.105 1.00 15.37 ATOM 116 OE1 CLU 18 13.990 28.571 34.805 1.00 17.92 ATOM 117 OE2 GLU 17 11.865 24.573 33.105 1.00 19.27 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 15.37 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.92 ATOM 124 CD LEU 13 15.553 32.049 35.579 1.00 14.31 ATOM 122 CB LEU 13 15.553 32.049 35.579 1.00 14.31 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 14.31	MOTA	84	CG	PRO					
ATOM 88 C ILE 14 18.768 17.791 34.286 1.00 12.08 ATOM 89 O ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 90 CB ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 91 CG1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 97 O LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CB LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 104 C VAL 16 15.347 24.748 32.904 1.00 10.54 ATOM 105 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 108 CG2 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 109 N GLU 17 13.904 26.144 34.870 1.00 12.85 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 12.85 ATOM 110 CA GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 13.61 ATOM 114 CG GLU 17 14.086 27.474 35.571 1.00 13.61 ATOM 115 CD GLU 17 14.086 27.474 35.571 1.00 15.74 ATOM 116 OE1 GLU 17 12.436 24.773 34.063 1.00 14.24 ATOM 117 OE2 GLU 17 12.260 25.402 35.344 1.00 14.15 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 14.93 ATOM 120 C LEU 18 12.553 30.978 33.917 1.00 14.31 ATOM 121 O LEU 18 13.990 28.571 34.805 1.00 14.93 ATOM 122 CG LEU 18 12.553 30.978 33.917 1.00 14.31 ATOM 122 CG LEU 18 15.533 32.049 35.579 1.00 14.93 ATOM 123 CG LEU 18 15.533 32.049 35.579 1.00 14.93	MOTA	85	CD	PRO	13	20.315			
ATOM 88 C ILE 14 18.768 17.791 34.286 1.00 12.08 ATOM 89 O ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 90 CB ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 91 CG1 ILE 14 22.271 18.169 35.015 1.00 16.54 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 97 O LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CB LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 104 C VAL 16 15.347 24.748 32.904 1.00 10.54 ATOM 105 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 108 CG2 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 109 N GLU 17 13.904 26.144 34.870 1.00 12.85 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 12.85 ATOM 110 CA GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 13.61 ATOM 114 CG GLU 17 14.086 27.474 35.571 1.00 13.61 ATOM 115 CD GLU 17 14.086 27.474 35.571 1.00 15.74 ATOM 116 OE1 GLU 17 12.436 24.773 34.063 1.00 14.24 ATOM 117 OE2 GLU 17 12.260 25.402 35.344 1.00 14.15 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 14.93 ATOM 120 C LEU 18 12.553 30.978 33.917 1.00 14.31 ATOM 121 O LEU 18 13.990 28.571 34.805 1.00 14.93 ATOM 122 CG LEU 18 12.553 30.978 33.917 1.00 14.31 ATOM 122 CG LEU 18 15.533 32.049 35.579 1.00 14.93 ATOM 123 CG LEU 18 15.533 32.049 35.579 1.00 14.93	MOTA	86	11		14				1.00 14.91
ATOM 89 O ILE 14 18.724 18.878 32.399 1.00 11.04 ATOM 90 CB ILE 14 21.109 18.905 34.325 1.00 16.54 ATOM 91 CG1 ILE 14 22.271 18.169 25.015 1.00 16.54 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 O LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CB LEU 15 14.412 19.541 33.199 1.00 40.50 ATOM 100 CD1 LEU 15 15.20.440 33.679 1.00 46.57 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.23 ATOM 103 CA VAL 16 16.885 22.556 33.919 1.00 10.23 ATOM 104 C VAL 16 15.347 24.748 32.904 1.00 10.25 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.225 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 15.344 1.00 15.37 ATOM 115 CD GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 116 OE1 GLU 17 12.650 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.650 25.557 32.950 1.00 83.46 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 15.37 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 15.37 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 15.61 ATOM 120 C LEU 18 14.116 29.914 35.401 1.00 15.61 ATOM 120 C LEU 18 14.116 29.914 35.401 1.00 15.61 ATOM 121 C LEU 18 14.116 29.914 35.401 1.00 15.61 ATOM 122 CB LEU 13 15.682 30.085 35.057 1.00 14.91 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 14.91 ATOM 122 CB LEU 13 15.682 30.03 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 30.03 35.005 1.00 19.27	MOTA	87	CA	ILE	14	19.893	17.991	34.266	1.00 12.93
ATOM 90 CB ILE 14 21.109 18.905 54.325 1.00 16.54 ATOM 91 CG1 ILE 14 20.783 20.207 35.084 1.00 18.08 ATOM 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 24.685 1.00 3.62 ATOM 97 O LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 99 CB LEU 15 15.462 19.450 34.2E5 1.00 21.25 ATOM 99 CG LEU 15 14.412 19.541 33.199 1.00 40.50 ATOM 99 CG LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.734 24.748 22.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 109 N GLU 17 13.904 26.144 34.870 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 14.43 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.031 27.524 36.765 1.00 15.37 ATOM 114 CG GUU 17 14.031 27.524 36.765 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 110 CE GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 112 C GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 113 CB GLU 17 12.456 24.178 34.447 1.00 15.37 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 121 C LEU 18 13.990 28.571 34.488 1.00 14.91 ATOM 122 CB LEU 18 13.5533 32.049 35.579 1.00 14.31 ATOM 123 CG LEU 18 15.5533 32.049 35.579 1.00 14.31 ATOM 124 CD1 LEU 18 15.5682 33.033 34.4488 1.00 12.30	ATOM	88	С	ILE	14	18.768	13.760	33.596	1.00 12.08
ATOM 90 CB ILE 14 21.109 18.905 54.325 1.00 16.54 ATOM 91 CG1 ILE 14 20.783 20.207 35.084 1.00 18.08 ATOM 93 CD1 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 24.685 1.00 3.62 ATOM 97 O LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 99 CB LEU 15 15.462 19.450 34.2E5 1.00 21.25 ATOM 99 CG LEU 15 14.412 19.541 33.199 1.00 40.50 ATOM 99 CG LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.734 24.748 22.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 109 N GLU 17 13.904 26.144 34.870 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 14.43 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.031 27.524 36.765 1.00 15.37 ATOM 114 CG GUU 17 14.031 27.524 36.765 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 110 CE GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 112 C GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 113 CB GLU 17 12.456 24.178 34.447 1.00 15.37 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 121 C LEU 18 13.990 28.571 34.488 1.00 14.91 ATOM 122 CB LEU 18 13.5533 32.049 35.579 1.00 14.31 ATOM 123 CG LEU 18 15.5533 32.049 35.579 1.00 14.31 ATOM 124 CD1 LEU 18 15.5682 33.033 34.4488 1.00 12.30	ATOM	89	0	ILE	14	18.724	18.878	32.399	1.00 11.04
ATOM 91 CG1 ILE 14 22.271 18.169 25.015 1.00 16.08 ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 O LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 O LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 24.25 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 40.50 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.56 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 22.904 1.00 16.72 ATOM 106 CB VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.24 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.031 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 14.331 27.524 36.765 1.00 15.37 ATOM 114 CG GLU 17 14.331 27.524 36.765 1.00 15.37 ATOM 115 CD GLU 17 14.331 27.524 36.765 1.00 15.37 ATOM 116 OE1 GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 117 OE2 GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 118 N LEU 18 13.900 28.571 34.807 1.00 18.75 ATOM 118 N LEU 18 13.990 28.571 34.695 1.00 19.87 ATOM 119 CA LEU 18 13.990 28.571 34.695 1.00 19.87 ATOM 120 C LEU 13 15.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 12 18 14.116 29.914 35.401 1.00 16.61 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 14.91 ATOM 123 CG LEU 13 15.528 30.938 32.917 1.00 14.31 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.5682 33.033 34.438 1.00 12.30		90	CB	ILE	14	21.109	18.905	34.325	1.00 16.54
ATOM 92 CG2 ILE 14 20.783 20.207 35.084 1.00 11.56 ATOM 93 CD1 ILE 14 23.642 18.836 34.738 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 O LEU 15 16.812 21.474 34.685 1.00 3.62 ATOM 99 CG LEU 15 15.462 19.450 35.901 1.00 10.00 ATOM 99 CG LEU 15 15.462 19.450 34.225 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 33.199 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.56 ATOM 104 C VAL 16 15.347 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.728 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 108 CG2 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.24 ATOM 100 CA CLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 110 CA CLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.436 24.178 34.447 1.00 14.15 ATOM 115 CD GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 116 0E1 GLU 17 12.436 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 14.91 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 14.91 ATOM 119 CA LEU 18 12.265 30.855 35.057 1.00 14.91 ATOM 120 C LEU 18 12.265 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.23					14		18.169	35.015	
ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 16.15 ATOM 94 N LEU 15 17.899 19.307 34.421 1.00 13.85 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 97 O LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 99 CB LEU 15 15.462 19.450 34.225 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.541 33.199 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.56 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.25 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.25 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 100 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 110 CA GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 115 CD GLU 17 12.650 25.402 35.344 1.00 15.37 ATOM 115 CD GLU 17 12.650 25.402 35.344 1.00 15.37 ATOM 116 OE1 GLU 17 12.650 25.402 35.344 1.00 15.37 ATOM 117 OE2 GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 49.50 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 14.91 ATOM 110 CA LEU 18 13.990 28.571 34.805 1.00 14.31 ATOM 111 C GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 112 C GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.260 30.855 35.055 1.00 18.75 ATOM 115 CD GLU 17 12.260 30.855 35.055 1.00 14.91 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 14.91 ATOM 121 C GLU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 123 CG LEU 18 15.533 32.049 35.579 1.00 14.93								35.084	
ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 97 C LEU 15 16.812 21.474 24.685 1.00 21.62 ATOM 99 CB LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 98 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 CB VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 108 CG2 VAL 16 18.273 24.668 34.098 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 13.904 26.144 34.870 1.00 12.58 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 115 CD GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 116 OE1 GLU 17 12.650 25.402 15.344 1.00 14.53 ATOM 115 CD GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 116 OE1 GLU 17 12.650 25.402 15.344 1.00 15.37 ATOM 116 OE1 GLU 17 12.650 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.650 25.557 32.950 1.00 83.46 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 121 CG LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 122 CB LEU 18 13.990 28.571 34.401 1.00 16.61 ATOM 122 CB LEU 18 13.990 28.571 34.401 1.00 16.61 ATOM 122 CB LEU 18 13.590 28.571 34.401 1.00 16.61 ATOM 122 CB LEU 18 13.590 28.571 34.409 1.00 11.356 ATOM 122 CB LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 122 CB LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.5682 33.033 34.4488 1.00 12.338									
ATOM 95 CA LEU 15 16.811 20.136 33.955 1.00 14.82 ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 97 O LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 98 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 22.904 1.00 16.72 ATOM 106 CB VAL 16 15.347 24.748 22.904 1.00 16.72 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 108 CG2 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 15.344 1.00 14.15 ATOM 115 CD GLU 17 12.650 25.402 15.344 1.00 14.15 ATOM 115 CD GLU 17 12.650 25.402 15.344 1.00 14.15 ATOM 115 CD GLU 17 12.650 25.402 15.344 1.00 14.15 ATOM 115 CD GLU 17 12.650 25.402 15.344 1.00 14.15 ATOM 116 OE1 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 14.31 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 14.31 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 120 C LEU 18 12.962 30.855 35.007 1.00 14.31 ATOM 120 C LEU 18 12.962 30.855 35.007 1.00 14.31 ATOM 121 C LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 14.33 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.5632 33.033 34.488 1.00 12.33									
ATOM 96 C LEU 15 16.915 21.474 34.685 1.00 3.62 ATOM 97 O LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 98 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 14.412 19.541 33.199 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.53 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.285 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 116 0E1 GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 117 OE2 GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 49.50 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 49.50 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.22 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 C LEU 18 13.990 38.571 34.805 1.00 14.24 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 14.91 ATOM 122 CB LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 122 CB LEU 13 15.682 33.033 34.438 1.00 12.38									. 00 14 82
ATOM 99 C LEU 15 17.080 21.509 35.901 1.00 10.00 ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 14.412 19.541 33.199 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.57 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.262 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 18.252 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.036 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 15.37 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.60 25.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 49.50 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 121 O LEU 18 13.990 28.571 34.805 1.00 14.91 ATOM 122 CB LEU 18 12.585 30.978 33.917 1.00 14.91 ATOM 123 CG LEU 18 12.585 30.630 35.005 1.00 14.91 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.5682 33.033 34.438 1.00 12.30									
ATOM 99 CB LEU 15 15.462 19.450 34.285 1.00 21.25 ATOM 99 CG LEU 15 14.412 19.541 33.199 1.00 40.50 ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.4688 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 108 CG2 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 18 15.682 33.033 34.438 1.00 21.40									
ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 40.50 ATOM 101 CD2 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.85 ATOM 108 CG2 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.24 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 12.585 30.855 35.057 1.00 14.91 ATOM 121 C LEU 18 12.585 30.855 35.057 1.00 14.31 ATOM 122 CB LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 122 CB LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.1033 34.438 1.00 21.43									
ATOM 100 CD1 LEU 15 13.279 20.440 33.679 1.00 46.97 ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.74 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 12 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.533 32.049 35.579 1.00 14.31 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.582 33.033 34.438 1.00 12.39									
ATOM 101 CD2 LEU 15 15.008 20.098 31.913 1.00 49.22 ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 121 O LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 122 CB LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 121 C LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 14.91 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 14.31 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39									
ATOM 102 N VAL 16 16.885 22.556 33.919 1.00 10.56 ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.24 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 15.74 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 120 C LEU 18 14.116 29.914 35.401 1.00 14.31 ATOM 120 C LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.533 32.049 35.579 1.00 14.31 ATOM 122 CB LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 19.30									
ATOM 103 CA VAL 16 16.964 23.905 34.479 1.00 10.23 ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.24 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 115 CD GLU 17 12.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 83.46 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 120 C LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 122 CB LEU 18 15.533 32.049 35.579 1.00 14.31 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39									
ATOM 104 C VAL 16 15.716 24.727 34.063 1.00 9.47 ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.20 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.32 ATOM 120 C LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 121 O LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 122 CB LEU 18 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39									
ATOM 105 O VAL 16 15.347 24.748 32.904 1.00 16.72 ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 115 CD GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 0E1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.22 ATOM 120 C LEU 18 13.990 28.571 34.805 1.00 17.22 ATOM 121 O LEU 18 13.990 28.571 34.805 1.00 14.91 ATOM 122 CB LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39	ATOM								
ATOM 106 CB VAL 16 18.273 24.668 34.098 1.00 12.85 ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39	ATOM								
ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 12 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39	MOTA	105	0	VAL	16			32.904	
ATOM 107 CG1 VAL 16 18.226 26.075 34.691 1.00 12.58 ATOM 108 CG2 VAL 16 19.520 23.945 34.628 1.00 14.24 ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 115 CD GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 121 O LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39	MOTA	106	CB	VAL	16	18.273	24.668	34.098	
ATOM 109 N GLU 17 15.059 25.317 35.060 1.00 14.43 ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39	ATOM	107	CG1	VAL	16	18.226	26.075	34.691	1.00 12.58
ATOM 110 CA GLU 17 13.904 26.144 34.870 1.00 13.61 ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39	MOTA	108	CG2	VAL	16	19.520	23.945	34.628	
ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39	MOTA	109	11	GLU	17	15.059	25.317	35.060	1.00 14.43
ATOM 111 C GLU 17 14.086 27.474 35.571 1.00 9.38 ATOM 112 O GLU 17 14.331 27.524 36.765 1.00 15.74 ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 117 OE2 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 15.682 33.033 34.438 1.00 12.39	ATOM	110	CA	GLU	17	13.904	26.144	34.870	1.00 13.61
ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.22 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 125 CD2 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38		111	С	GLU	17	14.086	27.474	35.571	1.00 9.38
ATOM 113 CB GLU 17 12.650 25.402 35.344 1.00 14.15 ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38			0	GLU	17	14.331		36.765	
ATOM 114 CG GLU 17 12.436 24.178 34.447 1.00 15.37 ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.22 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38			СВ		<u>:</u> 7	12.650	25.402		1.00 14.15
ATOM 115 CD GLU 17 11.865 24.573 33.105 1.00 49.50 ATOM 116 OE1 GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.22 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 124 CD1 LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38									1.00 15.37
ATOM 116 OE1'GLU 17 11.160 25.557 32.950 1.00 83.46 ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38									. 00 49 50
ATOM 117 OE2 GLU 17 12.220 23.766 32.127 1.00 38.75 ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38									
ATOM 118 N LEU 18 13.990 28.571 34.805 1.00 17.82 ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38									. 00 39 75
ATOM 119 CA LEU 18 14.116 29.914 35.401 1.00 16.61 ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38									
ATOM 120 C LEU 18 12.962 30.855 35.057 1.00 14.91 ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38				7 = 11					. 00 1/.52
ATOM 121 O LEU 18 12.585 30.978 33.917 1.00 14.31 ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38									
ATOM 122 CB LEU 13 15.426 30.630 35.005 1.00 13.56 ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38					- 3	12.962			14.91
ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38				LEU	- 3				
ATOM 123 CG LEU 13 15.533 32.049 35.579 1.00 19.27 ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38 ATOM 126 N ASP 19 12.480 31.551 36.082 1.00 17.88 ATOM 127 CA ASP 19 11.476 32.577 35.940 1.00 19.57					.3				
ATOM 124 CD1 LEU 13 16.740 32.182 36.489 1.00 21.40 ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 12.38 ATOM 126 N ASP 19 12.480 31.551 36.082 1.00 17.88 ATOM 127 CA ASP 19 11.476 32.577 35.940 1.00 19.57	MOTA				. 3				1.00 19.27
ATOM 125 CD2 LEU 13 15.682 33.033 34.438 1.00 13.38 ATOM 126 N ASP 19 12.480 31.551 36.082 1.00 17.88 ATOM 127 CA ASP 19 11.476 32.577 35.940 1.00 19.57	MOTA				:,3				1.00 21.40
ATOM 126 N ASP 19 12.480 31.551 36.082 1.00 17.58 ATOM 127 CA ASP 19 11.476 32.577 35.940 1.00 19.57	MOTA		CD2	LEU	:3	15.682	33.033	34.438	1.00 13.39
ATOM 127 CA ASP 19 11.476 32.577 35.940 1.00 19.57			:1		19	12.430	31.551	36.082	
			CA	ASP	19				

ATOM	128	С	ASP	19	12.098	33.896	36.360	1.00 11.65
MOTA	129	0	ASP	19	12.486	34.044	37.493	1.00 16.52
MOTA	130	CB	ASP	19	10.234	32.305	36.847	1.00 24.92
MOTA	131	CG	ASP	19	9.305	31.262	36.282	1.00 38.46
MOTA	132		ASP	19	9.572	30.587	36.989	1.00 61.49
MOTA	133		ASP	19 20	9.337 12.178	31.189 34.863	34.949 35.471	1.00 22.44
MOTA	134 135	N CA	GLY GLY	20	12.178	36.101	35.908	1.00 19.52
MOTA MOTA	136	C	GLY	20	12.048	37.385	35.538	1.00 19.35
ATOM	137	Ö	GLY	20	11.240	37.443	34.628	1.00 18.22
ATOM	138	N	ASP	21	12.401	38.407	36.286	1.00 13.19
ATOM	139	CA	ASP	21	11.908	39.737	36.112	1.00 15.36
ATOM	140	С	ASP	21	13.039	40.683	36.424	1.00 12.77
MOTA	141	0	ASP	21	13.517	40.742	37.569	1.00 15.18
ATOM	142	CB	ASP	21	10.701	40.036	37.040	1.00 22.26
ATOM	143	CG	ASP	21	10.230	41.491 42.407	37.022 36.557	1.00 30.80
ATOM	144 145	OD1	ASP ASP	21 21	10.878 9.062	41.658	37.604	1.00 45.92
MOTA MOTA	146	N	VAL	22	13.464	41.393	35.397	1.00 19.66
ATOM	147	CA	VAL	22	14.524	42.388	35.542	1.00 25.10
MOTA	148	c	VAL	22	14.010	43.780	35.154	1.00 18.25
MOTA	149	0	VAL	22	13.769	44.062	33.955	1.00 15.10
MOTA	150	CB	VAL	22	15.803	42.012	34.750	1.00 26.57
MOTA	151	CG1	VAL	22	16.861	43.127	34.896	1.00 24.27
MOTA	152	CG2	VAL	22	16.365	40.710	35.297	1.00 22.98
ATOM	153	N	ASN	23	13.823	44.641 45.993	36.166 35.908	1.00 25.32
MOTA	154 155	C <b>A</b>	ASN ASN	23 23	13.319 11.987	45.958	35.142	1.00 32.31
atoh atom	156	Ö	ASN	23	11.774	46.730	34.187	1.00 30.47
ATOM	157	CB	ASN	23	14.344	46.831	35.096	1.00 31.26
MOTA	158	CG	ASN	23	15.374	47.607	35.938	1.00 24.72
MOTA	159	OD1	ASN	23	15.795	47.183	37.024	1.00 27.22
MOTA	160	ND2	ASN	23	15.829	48.723	35.389	1.00 41.15
ATOM	161	N	GLY	24	11.118	45.024	35.519	1.00 24.95
ATOM ATOM	162 163	CA C	GLY	24 24	9.831 9.832	44.919	34.848 33.573	1.00 23.22
ATOM	164	ō	GLY	24	8.780	43.868	33.024	1.00 28.37
ATOM	165	11	HIS	25	11.000	43.691	33.071	1.00 20.89
MOTA	166	CA	нıs	25	11.042	42.840	31.877	1.00 19.30
ATOM	167	С	HIS	25	10.981	41.373	32.316	1.00 27.26
MOTA	168	0	HIS	25	11.898	40.850	32.951	1.00 25.47
ATOM	169	CB	HIS	25	12.268	43.060	30.958	1.00 24.20
MOTA	170	CG	HIS	25 25	12.313	44.382 45.514	30.218 30.758	1.00 33.04
ATOM ATOM	171 172		HIS	25 25	11.876	44.716	28.971	1.00 42.75
ATOM	173		HIS	25	12.801	46.497	29.867	1.00 39.14
ATOM	174		HIS	25	12.185	46.050	28.778	1.00 42.80
MOTA	175	N	LYS	26	9.872	40.728	32.028	1.00 25.90
ATOM	176	CA	LYS	26	9.675	39.355	32.446	1.00 26.27
ATOM	177	C	LYS	26	10.154	38.361	31.429	1.00 27.09
ATOM	178	0	LYS	26 26	10.027 8.230	38.576	30.232 32.863	1.00 25.75
ATOM ATOM	179 180	CB CG	LYS LYS	26	7.873	39.069 39.770	34.166	1.00 27.58
ATOM	181	CD	LYS	26	6.369	39.914	34.400	1.00 71.44
ATOM	182	CE	LYS	26	6.008	41.000	35.421	1.00 45.29
ATOM	183	N	PHE	2 <b>7</b>	10.703	37.250	31.910	1.00 22.04
MOTA	184	CA	PHE	27	11.164	36.236	30.978	1.00 13.72
MOTA	185	С	PHE	27	11.273	34.863	31.619	1.00 14.75
MOTA	186	0	PHE	27	11.293	34.722	32.842	1.00 15.94
ATOM	187	CB	PHE	27 27	12.495	36.638	30.287	1.00 21.58
MOTA MOTA	188 189	CG CD1	PHE	27	13.599 14.490	36.826 35.791	31.612	1.00 22.06
ATOM	190	CD2		27	13.722	33.029	32.005	1.00 17.55
ATOM	191	CE1		27	15.487	35.963	32.579	1.00 16.61
ATOM	192	CE2		27	14.747	33.234	32.931	1.00 19.78
ATOM	193	CI	PHE	27	15.621	37.127	33.234	1.00 13.83
ATOM	194	:1	SER	23	11.370	33.857	30.752	1.00 11.40

MOTA	195	CA	SER	28	11.492	32.479	31.186	1.00 15.59
MOTA	196	С	SER	28	12.579	31.749	30.379	1.00 15.96
MOTA	197	0	SER	28	12.699	31.933	29.167	1.00 18.99
ATOM	198	CB	SER	28	10.143	31.702	31.086	1.00 14.48
MOTA	199	OG	SER	28	9.510	31.678	32.353	1.00 31.95
MOTA	200	71	VAL	29	13.335	30.902	31.073	1.00 16.73
ATOM	201	CA	VAL	29.	14.361	30.093	30.435	1.00 14.06
MOTA	202	С	VAL	29	14.258	28.614	30.817	1.00 6.80
MOTA	203	0	VAL	29	14.058	28.266	31.987	1.00 10.85
ATOM	204	CB	VAL	29	15.768	30.570	30.839	1.00 17.96
ATOM	205	CG1	VAL	29	16.826	29.599	30.234	1.00 15.30
HOTA	206	CG2	VAL	29	15.989	32.001	30.357	1.00 16.37
ATOM	207	· 14	SER	30	14.462	27.781	29.824	1.00 11.31
MOTA	208	CA	SER	30	14.535	26.351	30.011	1.00 17.96
ATOM	209	С	SER	30	15.917	25.818	29.571	1.00 11.26
ATOM	210	0	SER	30	16.398	26.157	28.513	1.00 13.17
MOTA	211	CB	SER	30	13.471	25.603	29.202	1.00 19.91
MOTA	212	OG.	SER	30	12.249	25.667	29.882	1.00 48.74
MOTA	213	N	GLY	31	16.480	24.926	30.364	1.00 9.88
ATOM	214	CA	GLY	31	17.718	24.321	29.977	1.00 12.44
ATOM	215	С	GLY	31	17.737	22.816	30.249	1.00 13.16
MOTA	216	0	GLY	31	17.149	22.324	31.176	1.00 12.41
ATOM	217	11	GLU	32	18.459	22.112	29.433	1.00 13.44
MOTA	218	CA	GLU	32	18.622	20.670	29.570	1.00 13.73
MOTA	219	C	GLU	32	20.079	20.297	29.262	1.00 17.33
MOTA	220	0_	GLU	32	20.734	20.946	28.456	1.00 15.56 1.00 12.67
ATOM	221	СВ	GLU	32	17.761	19.893	28.543	
ATOM	222	CG	GLU	32	16.264	20.187	28.618	1.00 25.43
ATOM	223	CD	GLU	32	15.501	19.547	27.468	
ATOM	224		GLU	32	15.996	19.767	26.698 27.337	1.00 23.45
ATOM	225		GLU	32	14.292	20.022		1.00 15.36
ATOM	226	N	GLY	33 33	20.534 21.860	19.207 18.687	29.822 29.518	1.00 12.84
ATOM	227	CA	GLY	33	22.236	17.602	30.467	1.00 12.84
MOTA	228	C	GLY	33	22.236	16.919	31.011	1.00 13.56
MOTA	229	0	GLY	34	23.525	17.453	39.702	1.00 15.15
ATOM	230 231	N CA	GLU GLU	34	23.971	15.450	31.621	1.00 18.14
ATOM		C	GLU	34	25.220	15.874	32.367	1.00 16.26
MOTA	232 233	õ	GLU	34	25.926	17.760	31.944	1.00 18.67
ATOM	233	CB	GLU	34	24.180	15.114	30.927	1.00 22.53
ATOM ATOM	235	CG	GLU	34	24.948	15.261	29.624	1.00 33.78
MOTA	236	CD	GLU	34	24.879	14.020	28.796	1.00 55.15
ATOM	237	OE1		34	25.861	13.352	28.534	1.00 45.39
ATOM	238	OE2		34	23.653	13.719	28.430	1.00 56.26
ATOM	239	N	GLY	35	25.461	16.222	33.485	1.00 11.20
ATOM	240	CA	GLY	35	26.611	16.502	34.315	1.00 10.62
ATOM	241	c	GLY	35	27.293	15.192	34.662	1.00 19.92
ATOM	242	ō	GLY	35	26.650	14.161	34.750	1.00 16.69
ATOM	243	N	ASP	36	28.594	15.238	34.860	1.00 16.92
ATOM	244	CA	ASP	36	29.367	14.061	35.221	1.00 16.19
ATOM	245	c	ASP	36	30.396	14.505	36.233	1.00 13.94
ATOM	246	0	ASP	36	31.469	15.004	35.879	1.00 15.77
MOTA	247	CB	ASP	36	30.032	13.457	33.948	1.00 19.98
ATOM	248	CG	ASP	36	30.681	12.066	34.075	1.00 31.92
ATOM	249		ASP	36	31.236	11.519	33.141	1.00 30.97
MOTA	250	OD2	`ASP	36	30.587	11.515	35.248	1.00 25.32
ATOM	251	::	ALA	37	30.015	14.402	37.490	1.00 13.40
ATOM	252	CA	ALA	37	30.818	14.846	38.582	1.00 12.98
ATOM	253	С	ALA	37	32.181	14.145	38.637	1.00 21.94
MOTA	254	0	λLA	37	33.084	14.604	39.331	1.00 13.61
MOTA	255	C3	ALA	37	30.070	14.741	39.916	1.00 11.49
MOTA	256	::	THR	38	32.307	13.016 12.290	37.945	1.00 15.63
ATOM	2 <b>57</b>	CA	THR	38	33.581	12.290	37.943	1.00 19.54
MOTA	258	С	THR	3-8	34.705	13.114	37.335	1.00 25.61
MOTA	259	0	THR	38	35.850	13.069	37.775	1.00 17.89
MOTA	260	C3	THR	38	33.462	10.298	37.299	1.00 22.57
MOTA	261	OG 1	THR	38	32.543	12.146	33.067	1.00 29.86

FIG 5-5

								1.00 22.90
MOTA	262	CG2	THR	38	34.821	10.213	37.355	
ATOM	263	11	TYR	39	34.323	13.920	35.347	1.00 18.45
ATOM	264	CA	TYR	39	35.210	14.837	35.675	1.00 9.39
ATOM	265	С	TYR	39	34.874	16.291	35.991	1.00 14.41
						17.177	35.410	1.00 16.24
MOTA	266	0	TYR	39	35.454			
ATOM	267	CB	TYR	39	35.156	14.582	34.180	
ATOM	268	CG	TYR	39	35.426	13.137	33.929	1.00 28.73
ATOM	269	CD1	TYR	39	36.715	12.633	34.065	1.00 33.75
	270	CD2	TYR	39	34.392	12.249	33.642	1.00 39.19
ATOM							33.828	1.00 29.75
MOTA	271	CEl	TYR	39	36.982	11.276		
ATOM	272	CE2	TYR	39	34.635	10.885	33.435	1.00 45.41
MOTA	273	CZ	TYR	39	35.943	10.410	33.570	1.00 57.62
ATOM	274	ОН	TYR	39	36.199	9.070	33.364	1.00 70.77
		;;	GLY	40	33.935	16.525	36.929	1.00 9.94
ATOM	275						37.266	1.00 7.02
MOTA	276	CA	GLY	40	33.474	17.879		
ATOM	277	С	GLY	40	32.952	18.600	36.004	1.00 9.45
ATOM	278	0	GLY	40	33.068	19.830	35.829	1.00 12.63
ATOM	279	Ħ	LYS	41	32.380	17.823	35.092	1.00 5.44
	280	CA	LYS	41	31.954	18.335	33.842	1.00 6.63
ATOM					_		33.703	1.00 20.92
MOTA	281	С	LYS	41	30.414	18.554		
ATOM	282	0	LYS	41	29.617	17.693	34.085	1.00 12.94
ATOM	283	CB	LYS	41	32.360	17.357	32.827	1.00 8.27
ATOM	284	CG	LYS	41	32.099	17.771	31.419	1.00 13.19
ATOM	285	CD	LYS	41	32.521	16.644	30.481	1.00 20.20
						17.068	29.032	1.00 35.79
MOTA	286	CE	LYS	41	32.690			
MOTA	287	:1Z	LYS	41	33.113	15.954	23.147	1.00 47.55
ATOM	288	71	LEU	42	30.049	19.684	33.069	1.00 18.31
ATOM	289	CA	LEU	42	28.643	20.064	32.794	1.00 16.08
ATOM	290	c	LEU	42	28.456	20.422	31.330	1.00 14.23
			LEU	42	29.240	21.168	30.787	1.00 14.79
MOTA	291	0						1.00 13.22
MOTA	292	CB	LEU	42	28.223	21.300	33.621	
ATOM	293	CG	LEU	42	28.007	21.061	35.082	1.00 16.70
ATOM	294	CD1	LEU	42	27.894	22.406	35.782	1.00 13.79
ATOM	295		LEU	42	26.732	20.243	35.295	1.00 18.70
				43	27.395	19.914	30.672	1.00 8.04
ATOM	296	3	THR					
MOTA	297	CA	THR	43	27.103	20.275	29.282	1.00 4.67
MOTA	298	С	THR	43	25.636	20.656	29.186	1.00 17.23
MOTA	299	0	THR	43	24.811	19.818	29.442	1.00 14.35
MOTA	300	CB	THR	43	27.351	19.140	23.317	1.00 21.59
ATOM	301	0G1		43	28.692	18.743	23.415	1.00 42.74
							25.917	1.00 31.23
MOTA	302	CG2		43	27.073	19.675		
MOTA	303	::	LEU	44	25.327	21.934	13.830	1.00 11.83
MOTA	304	CA	LEU	44	23.944	22.409	22.847	1.00 13.81
MOTA	305	С	LEU	44	23.589	23.307	27.668	1.00 18.19
MOTA	306	0	LEU	44	24.416	23.989	27.107	1.00 13.86
							30.125	1.00 15.37
HOTA	307	CB	LEU	44	23.725	23.275		
ATOM	308	CG	LEU	44	23.369	22.584	31.456	1.00 24.69
ATOM	309	CD1	LEU.	44	21.869	22.381	31.601	1.00 23.20
ATOM	310	CD2	LEU	44	24.083	21.286	31.650	1.00 46.18
ATOM	311	N	LYS	45	22.294	23.331	27.339	1.00 10.29
	312	CA	LYS	45	21.752	24.224	25.358	1.00 11.94
ATOM			LIS					
ATOM	313	С	LYS	45	20.534	24.913	25.957	
ATOM	314	0	LYS	45	19.665	24.248	27.530	1.00 18.43
ATOM	315	CЗ	LYS	45	21.409	23.560	15.060	1.00 13.75
ATOM	316	CG	fiz	45	20.878	24.556	24.045	1.00 8.83
ATOH	317	CD	ĹYS	45	20.486	23.863	22.746	1.00 26.87
							22.740	1.00 16.58
MOTA	318	CE	LYS	45	19.574	24.688	11.842	1.00 16.32
ATOM	319	:12	LYS	45	19.318	24.024	20.555	1.00 18.33
ATOM	320	::	PHE	46	20.535	26.236	36.910	1.00 12.34
ATOM	321	CA	PHE	46	19.463	27.048	27.451	1.00 13.32
	322	c	PHE	46	18.759	27.718	26.343	1.00 18.26
ATOM							25.360	
ATOM	323	9	SHE	46	19.326	28.093	_ = . 360	1.00 16.83
ATOM	324	CB	SHE	46	19.934	33.101	12.473	1.00 15.29
ATOM	325	CC	PHE	÷6-	20.773	27.495	19.552	1.00 13.81
ATOM	326	co:		<b>÷6</b>	22.132	27.268	19.552	1.00 17.06
ATOM	327		2 PHE	46	20.209	27.121	13.774	1.00 8.24
								. 00 15 25
MOTA	328	CE.	1 SHE	<b>÷6</b>	22.924	26.693	30.331	1.00 15.95

ATOM	329	CE2		46	20.979	26.524	31.767	1.00 11.90
ATOM	330	CŽ	PHE	46	22.340	26.309	31.540	1.00 8.84
atom	331	N	ILE	47	17.440	27.845	26.498	1.00 13.24
HOTA	332	CA	ILE	47	16.588	28.453	25.479	1.00 18.02
ATOM	333	С	ILE	47	15.645	29.460	26.118	1.00 20.14
ATOM	334	0	ILE	47	15.039	29.162	27.148	1.00 17.67
ATOM	335	CB	ILE	47	15.737	27.386	24.801	1.00 22.67
ATOM	336	CG1	ILE	47	16.585	26.271	24.291	1.00 20.66
MOTA	337	CG2	ILE	47	15.024	28.002	23.641	1.00 33.79
ATOM	338	CD1	ILE	47	16.639	26.293	22.805	1.00 23.69
MOTA	339	N	CYS	48	15.564	30.653	25.561	1.00 14.68
ATOM	340	CA	CYS	48	14.681	31.635	26.170	1.00 16.93
MOTA	341	С	CYS	48	13.323	31.352	25.628	1.00 24.18
MOTA	342	0	CYS	48	13.122	31.513	24.453	1.00 20.63
MOTA	343	CB	CYS	48	15.063	33.116	25.885	1.00 16.85
MOTA	344	SG	CYS	48	13.913	34.268	26.712	1.00 22.06
MOTA	345	N	THR	49	12.424	30.871	26.484	1.00 27.31
MOTA	346	CA	THR	49	11.101	30.458	26.042	1.00 32.18
MOTA	347	С	THR	49	10.106	31.572	25.803	1.00.37.51
MOTA	348	0	THR	49	9.150	31.407	25.061	1.00 35.71
MOTA	349	CB	THR	49	10.537	29.417	26.972	1.00 23.66
ATOM	350	OG1	THR	49	10.387	29.989	28.258	1.00 30.10
ATOM	351	CG2	THR	49	11.512	28.226	27.022	1.00 29.98
MOTA	352	N	THR	50	10.314	32.693	26.447	1.00 32.34
atom	353	CA	THR	50	9.416	33.810	26.283	1.00 28.67
HOTA	354	С	THR	50	9.836	34.711	25.125	1.00 37.98
MOTA	355	0	THR	50	9.228	35.763	24.904	1.00 39.17
atom	356	СЗ	THR	50	9.251	34.611	27.589	1.00 36.23
atom	357	OG1	THR	50	10.512	34.980	28.118	1.00 35.37
HOTA	358	CG2	THR	50	8.507	33.773	28.602	1.00 27.78
ATOM	359	N	GLY	51	10.881	34.282	24.372	1.00 31.04
ATOM	360	CA	GLY	51	11.394	35.059	23.239	1.00 32.42
ATOM	361	C	GLY	51	12.865	35.542	23.427	1.00 48.45
MOTA	362	0	GLY	51	13.779	34.737	23.701	1.00 57.11
ATOM	363	И	LYS	52 53	13.087	36.862	23.282	1.00 36.08
ATOM	364	CA	LYS	52 52	14.416	37.460	23.415 24.861	1.00 29.65
ATOM	365	C	LYS		14.827	37.726		1.00 25.70
ATOM	366	0	LYS	52	14.140	39.420 38.714	25.620	
ATOM	367	C3	LYS	52 52	14.577 15.772	38.649	22.522 21.644	1.00 43.37
ATOM	368	CG	LYS LEU	53	15.772	37.190	25.250	1.00 19.22
ATOM	369	N CA	LEU	53	16.439	37.130	26.596	1.00 13.52
ATOM ATOM	370 371	C	LEU	£3	16.717	38.932	26.775	1.00 17.76
ATOM	372	Ö	LEU	53	17.392	39.539	25.973	1.00 21.59
MOTA	373	СВ	LEU	53	17.705	36.567	26.845	1.00 17.39
ATOM	374	CG	LEU	53	18.100	36.435	28.302	1.00 17.43
ATOM	375	CD1		53	17.048	35.621	29.053	1.00 20.12
ATOM	376	CD2		53	19.440	35.718	28.368	1.00 16.11
ATOM	377	N	PRO	54	16.197	39.525	27.817	1.00 16.69
ATOM	378	CA	PRO	54	16.324	40.962	28.092	1.00 18.60
ATOM	379	c	PRO	54	17.638	41.414	28.707	1.00 25.39
ATOM	380	Ö	PRO	54	17.865	42.609	28.861	1.00 18.88
ATOM	381	СВ	PRO	54	15.268	41.265	29.139	1.00 22.52
ATOM	382	CG	PRO	54	14.832	39.933	29.720	1.00 26.02
ATOM	383	CD	PRO	54	15.318	38.855	28.779	1.00 21.26
ATOM	384	ห	VAL	55	18.435	40.455	29.151	1.00 23.32
ATOM	385	CA	VAL	55	19.746	40.716	29.711	1.00 15.83
ATOM	386	C	VAL	55	20.688	39.868	28.973	1.00 19.38
ATOM	387	ŏ	VAL	5.5	20.268	39.035	28.219	1.00 20.34
ATOM	388	СB	VAL	5.5	19.814	40.409	31.147	1.00 17.67
ATOM	389	CG1		5.5	18.864	41.340	31.851	1.00 22.52
ATOM	390	CG2		5.5	19.402	38.959	31.397	1.00 19.11
ATOM	391	31	2RO	56	21.963	40.070	29.167	1.00 19.37
ATOM	392	CA	PRO	55	22.911	39.258	28.447	1.00 13.09
ATOM	393	c	PRO	3.5	23.059	37.834	29.038	1.00 5.83
ATOM	394	ō	PRO	3.5	23.067	37.631	30.254	1.00 12.35
ATOM	395	СВ	280	56	24.231	40.062	13.420	1.00 13.34

14/36

МОТА	396	CG	?RO	56	23.851	41.478	28.849	1.00 20.73
ATOM	397	CD	PRO	5 6	22.525	41.379	29.578	1.00 18.55
ATOM	398	t:	TRP	57	23.202	36.848	28.158	1.00 11.12
ATOM	399	CA	TRP	5.7	23.354	35.458	28.595	1.00 12.55
MOTA	400	С	TRP	57	24.411	35.239	29.700	1.00 14.13
MOTA	401	0	TRP	57	24.178	34.586	30.709	1.00 11.49
MOTA	402	CB	TRP	57	23.604	34.535	27.406	1.00 10.56
ATOM	403	ÇG	TRP	57	22.335	34.237	26.641	1.00 12.65
MOTA	404	CD1	TRP	57	21.999	34.714	25.426	1.00 16.24
MOTA	405	CD2	TRP	57	21.281	33.327	27.013	1.00 12.50
MOTA	406	NE1	TRP	57	20.784	34.200	25.018	1.00 14.25
MOTA	407	CE2	TRP	57	20.315	33.354	25.963	1.00 14.03
MOTA	408	CE3	TRP	57	21.052	32.521	28.129 26.007	1.00 14.36
MOTA	409	CZ2	TRP	57 5 <b>7</b>	19.148	32.583 31.767	28.170	1.00 14.23
MOTA	410	CZ3	TRP	57	19.887 18.945	31.818	27.128	1.00 10.01
ATOM	411	CH2	TRP	57 58	25.594	35.800	29.518	1.00 15.78
MOTA	412	N	PRO	58	26.629	35.616	30.503	1.00 9.53
ATOM	413	CA C	PRO PRO	58	26.241	36.010	31.878	1.00 9.71
ATOM	414 415	o	PRO	58	26.760	35.467	32.825	1.00 11.70
MOTA MOTA	416	СВ	PRO	58	27.833	36.441	30.040	1.00 10.83
MOTA	417	CG	PRO	58	27.597	36.748	28.582	1.00 18.50
ATOM	418	CD	PRO	58	26.137	36.432	28.278	1.00 15.82
ATOM	419	N	THR	59	25.336	36.977	32.021	1.00 7.54
MOTA	420	CA	THR	59	24.976	37.366	33.357	1.00 4.53
ATOM	421	c	THR	59	24.228	36.258	34.137	1.00 8.41
ATOM	422	Ō	THR	59	24.174	35.251	35.367	1.00 10.57
MOTA	423	CB	THR	59	24.197	38.691	33.384	1.00 15.54
ATOM	424	OG1	THR	59	22.295	38.480	32.844	1.00 15.51
ATOM	425	CG2	THR	59	24.917	39.731	32.542	1.00 15.76
MOTA	426	N	LEU	60	23.686	35.304	33.427	1.00 11.99 1.00 9.15
MOTA	427	CA	LEU	60	22.899	34.248	34.073	
MOTA	428	С	LEU	60	23.657	32.944	34.385	1.00 15.62 1.00 11.99
MOTA	429	0_	LEU	60	23.118	32.027	35.042 33.203	1.00 7.67
MOTA	430	CB	LEU	60 60	21.645	33.914 35.111	33.203	1.00 14.05
ATOM	431	CG	LEU	60 60	20.728 19.620	34.775	32.062	1.00 14.54
ATOM	432	CD1		60 60	20.142	35.456	34.394	1.00 10.67
ATOM	433	N	VAL	61	24.893	32.837	33.917	1.00 11.27
MOTA MOTA	434 435	CA	VAL	51	25.656	31.587	34.094	1.00 4.37
ATOM	436	C	VAL	61	25.678	31.013	35.496	1.00 6.02
ATOM	437	ō	VAL	61	25.355	29.805	35.743	1.00 10.75
ATOM	438	CB	VAL	61	27.050	31.643	33.406	1.00 7.14
MOTA	439	CG1		61	27.888	30.396	33.805	1.00 6.47
ATOM	440	CG2	VAL	51	26.890	31.745	31.876	1.00 6.63
ATOM	441	N	THR	52	26.053	31.843	36.442	1.00 7.02
ATOM	442	CA	THR	52	26.178	31.421	37.808	1.00 6.51
ATOM	443	С	THR	62	24.862	30.954	38.410	1.00 9.22
MOTA	444	0	THR	62	24.801	30.163	39.352	1.00 6.99
ATOM	445	CB	THR	52	26.816	32.520	38.660	1.00 16.97 1.00 12.00
MOTA	446	OGI		62	26.103 28.297	33.744 32.708	38.453 38.225	1.00 12.30
'ATOM	447	CG2		52 63		31.547	37.910	1.00 9.98
ATOM	448	N	THR	63	23.814 22.457	31.212	38.388	1.00 5.59
ATOM	449 450	CA C	THR THR	63	22.033	29.830	37.865	1.00 3.14
MOTA MOTA	451	Ö	THR	63	21.499	23.984	38.604	1.00 13.48
ATOM	452	CB	THR	53	21.458	32.312	37.925	1.00 11.14
MOTA	453	OGI		<b>63</b>	21.785	33.498	38.602	1.00 11.75
ATOM	454	CG		53	20.024	31.897	38.296	1.00 9.31
ATOM	455	N	PHE	54	22.250	29.620	36.583	1.00 10.19
ATOM	456	CA	PHE	5 <b>4</b>	21.895	28.371	35.995	1.00 3.00
MOTA	457	c	PHE	64	22.774	27.253	36.518	1.00 25.25
MOTA	458	ŏ	PHE	54	22.313	26.147	36.761	1.00 9.64
ATOM	459	СВ	PHE	54	22.114	23.438	34.513	1.00 £.28
MOTA	460	CG	PHE	54	21.233	29.357	33.750	1.00 10.95
ATOM	461	CD.		54	21.724	29.954	32.593	1.00 7.15
ATOM	462		2 PHE	54	19.899	23.563	34.106	1.00 14.43

ATOM	463	CEl	PHE	64	20.936	30.792	31.805	1.00 14.20
MOTA	464	CE2	PHE	64	19.077	30.375	33.317	1.00 13.95
MOTA	465	CZ	PHE	64	19.597	30.983	32.171	1.00 16.35
HETATM	<b>∔66</b>	31	CRO	6 <b>6</b>	24.077	27.513	36.610	1.00 11.86
HETATM	467	CG1	CRO	66	25.155	25.422	34.796	1.00 16.67
HETATM	468	OG1		66	26.679	27.129	35.461	1.00 14.22
HETATM	169	CB1		66	25.931	26.035	35.930	1.00 10.77
HETATM	470		CRO	66	25.011	26.478	37.078	1.00 7.34
HETATM	471	C1	CRO	66	25.718	26.991	38.253	1.00 17.70
HETATM	472	N2	CRO	66	26.975	27.732	38.216	1.00 9.21
HETATM	473	OH	CRO	66	32.894	30.804	36.971	1.00 13.84
HETATM	474	CD2		66	30.487	30.110	39.805	1.00 10.79
HETATM	475		CRO	66	31.614	30.563	39.085	1.00 10.01 1.00 9.48
HETATM	476	CZ	CRO	66	31.718	30.300	37.721 37.033	1.00 9.48
HETATM	477	CE1		66	30.707	29.546 29.103	37.742	1.00 17.44
HETATM	478	CD1		66 66	29.541 29.437	29.370	39.124	1.00 7.67
HETATM	479 480	CG2 CB2	CRO	66	28.329	28.822	39.960	1.00 10.75
HETATM HETATM	481		CRO	66	27.197	28.245	39.512	1.00 16.08
HETATM	482	C2	CRO	66	26.043	27.875	40.370	1.00 5.46
HETATM	483	02	CRO	66	26.022	27.962	41.566	1.00 13.20
HETATH	484	113	CRO	66	25.240	26.978	39.517	1.00 18.43
HETATM	485	CA3	CRO	66	23.840	26.511	39.734	1.00 10.40
HETATM	486	C3	CRO	66	23.413	25.550	40.817	1.00 11.96
HETATM	487	03	CRO	66	22.747	26.014	41.764	1.00100.00
ATOM	488	n	VAL	68	23.737	24.208	41.005	1.00 29.95
ATOM	489	CA	VAL	68	24.209	22.972	40.304	1.00 17.16
ATOM	490	С	VAL	68	25.692	22.550	40.734	1.00 14.88
MOTA	491	0	VAL	68	26.379	21.821	40.026	1.00 9.03
MOTA	492	CB	VAL	68	23.870	22.899	38.831	1.00 18.94
MOTA	493	CG1	VAL	68	24.685	22.088	37.942	1.00 17.17
MOTA	494		VAL	68	22.396	22.538	38.680	1.00 18.80
ATOM	495	51	GLN	69	26.129	22.965	41.914	1.00 11.04
MOTA	496	CA	GLN	69	27.465	22.764	42.394	1.00 15.00
ATOM	497	C	GLN	69	27.749	21.366	42.893	1.00 22.46
ATOM	498	0	GLN	69	28.876	21.025	43.154	1.00 15.84
ATOM	499	CE	GLN	69	27.929	23.852	43.414 42.615	1.00 14.13
MOTA	500	CG	GLN	69 69	28.202 28.216	25.174 25.385	43.520	1.00 17.01
ATOM	501 502	CD	GLN	69	27.433	26.476	44.448	1.00 18.94
MOTA MOTA	503		GLN GLN	69	29.151	27.300	43.241	1.00 8.52
MOTA	504	N	CYS	70	26.703	20.540	42.906	1.00 12.10
ATOM	505	CA	CYS	70	26.862	19.171	43.287	1.00 11.34
ATOM	506	c c	CYS	70	27.611	18.391	42.175	1.00 10.54
ATOM	507	ō	CYS	70	28.036	17.242	42.367	1.00 14.70
MOTA	508	СВ	CYS	70	25.476	13.584	43.596	1.00 14.52
ATOM	509	SG	CYS	70	24.325	19.012	42.251	1.00 15.61
ATOM	510	N	PHE	71	27.801	19.029	41.005	1.00 8.64
MOTA	511	CA	PHE	71	28.525	18.419	39.883	1.00 6.59
MOTA	512	С	PHE	71	30.041	18.754	39.876	1.00 16.43
MOTA	513	0	PHE	71	30.753	13.481	38.916	1.00 13.05
MOTA	514	CB	PHE	71	27.951	13.771	38.523	1.00 7.61
MOTA	515	CG	PHE	71	26.669	13.016	38.303	1.00 14.73
ATOM	516		PHE	71	26.693	15.642	38.050	1.00 10.34
ATOM	517		PHE	71	25.434	18.660	38.453 37.866	1.00 17.14
ATOM	518		PHE	71 71	25.506	15.931 17.961	38.300	
MOTA	519	CE2		71	24.238 24.282	15.598	37.990	1.00 20.92
MOTA	520	CZ	PHE			19.370	40.938	1.00 13.13
ATOM	521 522	:i CA	SER SER	72 72	30.500 31.889	19.715	41.075	1.00 13.13
ATOM A <i>TOM</i>	523	C	SER	72	32.689	3.446	41.357	1.00 14.56
ATOM	524	9	SER	72	32.256	7.566	42.122	1.00 10.90
MOTA	525	CB	SER	72	32.075	20.672	-2.257	1.00 3.65
MOTA	526	CG	SER	72	31.361	11.274	42.038	1.00 19.29
ATOM	527	;;	ÀRG	73	33.905	17.566 20.672 11.874 13.352	40.794	1.00 15.27
ATOM	528	CA	ARG	73	34.695	17.212	41.117	1.00 13.55
MOTA	529	C	ARG	73	35.414	17.212 17.426	42.443	1.00 19.96
						•	-	

16/36

3 50014	530	^	ARG	73	36.182	18.376	42.599	1.00 16.14
MOTA		0		73	35.694	16.817	40.013	1.00 16.80
HOTA	531	CB	ARG					1.00 20.13
MOTA	532	CC .	ARG	73	36.549	15.616	40.460	
MOTA	533	CD	ARG	73	37.489	15.093	39.381	1.00 28.47
ATOM	534	NE	ARG	73	38.743	15.85 <del>9</del>	39.260	1.00 25.48
ATOM	535	CZ	ARG	73	39.756	15.777	40.127	1.00 28.04
	536	NH1		73	39.688	15.004	41.195	1.00 28.76
ATOM				73	40.865	16.504	39.918	1.00 39.65
ATOM	537		ARG				43.424	1.00 12.05
MOTA	538	N	TYR	74	35.151	16.561		
ATOM	539	CA	TYR	74	35.861	16.659	44.690	1.00 11.57
ATOM	540	С	TYR	74	36.946	15.566	44.721	1.00 25.02
ATOM	541	ō	TYR	74	36.658	14.387	44.558	1.00 19.71
			TYR	74	34.978	16.528	45.934	1.00 15.51
MOTA	542	CB			34.395	17.850	46.402	1.00 16.59
ATOM	543	CG	TYR	74			45.631	1.00 14.44
ATOM	544		TYR	74	33.455	18.546	· ·	1.00 15.94
MOTA	545	CD2	TYR	74	34.799	18.399	47.618	
ATOM	546	CE1	TYR	74	32.901	19.756	46.059	1.00 7.99
ATOM	547	CE2	TYR	74	34.261	19.612	48.058	1.00 18.29
ATOM	548	CZ	TYR	74	33.294	20.276	47.298	1.00 13.87
	549	ОН	TYR	74	32.829	21.507	47.738	1.00 18.39
ATOM				75	38.181	15.947	44.902	1.00 19.20
ATOM	550	N	PRO		39.213	14.940	44.995	1.00 18.42
MOTA	551	CA	PRO	75				1.00 15.60
MOTA	552	С	PRO	75	38.958	13.993	46.175	
ATOM	553	0	PRO	75	38.373	14.361	47.174	1.00 11.99
ATOM	554	CB	PRO	75	40.514	15.681	45.195	1.00 18.31
ATOM	555	CG	PRO	75	40.242	17.158	44.863	1.00 24.81
		CD	PRO	75	38.742	17.306	44.694	1.00 15.41
ATOM	556			76	39.433	12.756	46.038	1.00 18.63
MOTA	557	N	ASP				47.062	1.00 16.19
atom	558	CA	ASP	76	39.269	11.770		1.00 15.92
ATOM	55 <del>9</del>	С	ASP	76	39.581	12.280	48.431	
ATOM	560	0	ASP	76	38.862	12.042	49.389	1.00 17.35
ATOM	561	CB	ASP	76	40.083	10.507	46.790	1.00 18.69
ATOM	562	CG	ASP	76	39.826	9.432	47.825	1.00 24.04
	563	OD1		76	40.523	9.268	48.817	1.00 29.72
ATOM				76	38.732	8.743	47.584	1.00 40.96
ATOM	564	OD2			40.647	12.984	48.561	1.00 18.79
atom	565	N	HIS	7 <b>7</b>				1.00 19.35
ATCM	566	CA	HIS	77	40.978	13.418	49.877	
MOTA	567	C	HIS	77	40.117	14.507	50.397	1.00 24.57
ATOM	568	0	∺IS	77	40.205	14.826	51.551	1.00 27.15
ATOM	569	C3	HIS	77	42.435	13.806	50.042	1.00 19.84
ATOM	570	CG	HI5	77	42.743	15.035	49.322	1.00 17.31
	571		. HIS	77	42.925	15.028	47.953	1.00 21.86
MOTA				77	42.925	15.295	49.774	1.00 18.70
ATOM	572	CD2		77	43.203	16.289	47.593	1.00 17.49
ATOM	573	CEI						1.00 18.11
MOTA	574	NE2		77	43.213	17.069	48.668	
ATOM	575	N	MSE	78	39.277	15.069	49.565	1.00 25.36
MOTA	576	CA	MSE	78	38.412	16.140	50.026	1.00 24.65
ATOM	577	С	MSE	78	36.920	15.774	50.066	1.00 26.47
	578	ō	MSE	78	36.070	16.636	50.260	1.00 28.16
ATOM	579	CB	MSE	78	38.596	17.331	49.121	1.00 26.38
ATOM				78	39.803	18.177	49.406	1.00 27.01
ATOM	580	CC	MSE					1.00 43.09
ATOM	581		MSE	78	39.987	19.608	48.117	
ATOM	582	CE	MSE	78	38.874	20.873	49.044	1.00 27.11
MOTA	583	11	LYS	79	36.606	14.509	49.856	1.00 18.68
ATOM	584	CA	LYS	79	35.216	14.061	<b>∔9.85</b> 3	1.00 21.54
ATOM	585	C	LYS	79	34.406	14.449	51.082	1.00 20.21
			LYS	79	33.186	14.652	51.025	1.00 21.08
MOTA	586	0				12.581	49.612	1.00 23.48
ATOM	587	CB	LYS	79 70	35.152	12.501		1.00 23.40
MOTA	588	CG	LYS	79	35.859	12.225	48.317	1.00 41.09
ATOM	589	CD	LYS	79	35.159	11.134	47.535	1.00 34.66
ATOM	590	CE	LYS	79	35.796		÷6.181	1.00 53.46
ATOM	591	:: 2	175	79	35.084	11.549	45.020	1.00 49.53
	592	::	ARG	30	35.069		52.213	1.00 19.77
ATOM				30	34.365		53.434	1.00 20.13
MOTA	593	CA	ARG				13.481	1.00 26.42
MOTA	594	С	ARG	30	33.898			
MOTA	595	0	ARG	30	33.251	15.717	54.467	1.00 23.51
ATOM	596	СЗ	ARG	30	35.155	14.549	54.700	1.00 24.58
-								

MOTA	597	CG	ARG	30	36.204	15.620	55.034	1.00 29.	71
ATOM	598	CD	ARG	30	36.964	15.344	56.335	1.00 61.	30
ATOM	599	NE	ARG	30	36.551	16.230	57.415	1.00 71.	
MOTA	600	CZ	ARG	50	37.398	16.882	58.192	1.00100.	
MOTA	601	NH1	ARG	٥٥	38.714	16.758	58.040	1.00100.	
MOTA	602	NH2	ARG	80	36.917	17.679	59.155	1.00 99.	
ATOM	603	N	HIS	31	34.275	17.121	52.473	1.00 18.	
atom	604	CA	HIS	<b>S1</b>	33.903	18.547	52.499	1.00 19.	
atom	605	С	HIS	91	32.841	18.883	51.486	1.00 18.	
ATOM	606	0	HIS	81	32.557	20.043	51.295	1.00 17.	
MOTA	607	CB	HIS	81	35.129	19.472	52.283	1.00 20.	
ATOH	608	CG	HIS	81	36.221	19.224	53.305	1.00 28.	
ATOM	609	ND1		81	36.127	19.701	54.618 53.202	1.00 29.	
MOTA	610			81	37.392	18.535 19.308	55.265	1.00 26.	
MOTA	611	CEI	HIS	81 81	37.218 37.991	18.603	54.452	1.00 28.	
ATOM	612	NE2	HIS ASP	82	32.298	17.843	50.841	1.00 12.	
MOTA MOTA	613 614	N CA	ASP	82	31.358	18.011	49.769	1.00 13.	
	615	C	ASP	82	29.922	18.148	50.259	1.00 24.	
ATOM ATOM	616	o	ASP	82	29.175	17.195	50.243	1.00 16.	
ATOM	617	СВ	ASP	92	31.480	16.917	48.730	1.00 12.	
MOTA	618	CG	ASP	82	30.642	17.209	47.518		92
MOTA	619	OD1	ASP	82	29.870	18.134	47.459	1.00 20.	31
ATOM	620	OD2	ASP	32	30.938	16.466	46.507	1.00 11.	12
MOTA	621	N	PHE	83	29.566	19.353	50.705	1.00 23.	
ATOM	622	CA	PHE	23	28.220	19.634	51.201	1.00 20.	
MOTA	623	С	PHE	83	27.154	19.333	50.168	1.00 20.	
MOTA	624	0	PHE	83	26.116	18.733	50.503	1.00 15.	
MOTA	625	CB	PHE	83	28.077	21.106	51.666	1.00 19.	
MOTA	626	CG	PHE	83	26.624	21.613	51.805	1.00 16.	
MOTA	627	CD1		<b>53</b>	25.946	21.498	53.021	1.00 17.	
MOTA	628	CD2		33	25.968	22.236	50.734	1.00 18.	
MOTA	629	CE1		83	24.635	21.960	53.156	1.00 24	
MOTA	630	CE2		83	24.650	22.690	50.840	1.00 19.	
MOTA	631	cz	PHE	83	24.001	22.575	52.068 48.921	1.00 20	
MOTA	632	N:	PHE	64 84	27.432	19.784 19.693	47.809	1.00 12	
MOTA	633 634	CA C	PHE	3 <b>4</b> 84	26.515 25.893	18.332	47.602	1.00.24	
MOTA MOTA	635	ò	PHE	84	24.674	18.200	47.534	1.00 21	
MOTA	636	СВ	PHE	34	27.085	20.265	46.513	1.00 13	
ATOM	637	CG	PHE	34	27.630	21.645	46.721	1.00 14	. 27
ATOM	638	CD1		34	29.001	21.845	46.890	1.00 15	. 17
ATOM	639	CD2		94	26.781	22.753	46.752		. 48
ATOM	640	CE1		34	29.520	23.129	47.073	1.00 14	. 63
ATOM	641	CE2	PHE	94	27.276	24.041	46.969		. 34
ATOM	642	CZ	PHE	34	28.650	24.221	47.137		.77
MOTA	643	N	LYS	35	26.738	17.330	47.482		.07
MOTA	644	CA	LYS	35	26.294	15.985	47.283		. 30
atom	645	С	LYS	35	25.657	15.371	48.547		. 43
MOTA	646	0	LYS	35	24.773	14.509	48.429	1.00 18 1.00 17	.40
MOTA	647	CB	LYS	35 35	27.434	15.089 15.372	46.757 45.323	1.00 17	
MOTA'	648	CG	LYS LYS	35 35	27.873 28.969	14.381	44.888	1.00 13	
ATOM	649 650	CE	LYS	35 35	29.766	14.819	43.662	1.00 10	
ATOM	651	NZ	LYS	35	30.319	16.185	43.773	1.00 12	
ATOM ATOM	652	21	SER	36	26.119	15.795	49.752	1.00 11	
ATOM	653	CA.	SER	36	25.610	15.267	50.998	1.00 12	
MOTA	654	c	SER	36	24.156	15.639	51.240	1.00 21	
MOTA	655	Õ	SER	36	23.452	14.979	52.013	1.00 19	.89
ATOM	656	СB	SER	36	26.448	15.661	52.208	1.00 15	
ATOM	557	OG	SER	36	26.308	17.042	52.495	1.00 22	.05
ATOM	658	:1	ALA	27	23.705	16.6 <del>9</del> 8	50.582	1.00 15	.09
MOTA	659	CA	ALA	<b>27</b>	22.333	17.138	50.762	1.00 19	
MOTA	560	С	ALA	<b>5</b> 7	21.337	16.399	49.870		.60
MOTA	661	0	ALA	€7	20.162	16.557	50.040	1.00 19	. 55
MOTA	562	C3	ALA	57	22.204	19.647	50.632		. 23
MOTA	663	::	MSE	38	21.935	15.536	48.976	1.00 14	.05

MOTA	564	CA	MSE	88	21.007	14.796	48.035		15.32
MOTA	665		MSE	88	20.496	13.448	48.579		21.48
MOTA	566		MSE	88	21.109	12.876	49.457		23.03
MOTA	667		MSE	88	21.848	14.593	46.791	1.00	15.92
ATOM	668		MSE	88	22.263	15.891	46.131		10.56
MOTA			MSE	88	20.737	16.894	45.394	1.00	31.99
MOTA	670		MSE	88	21.319	18.684	45.748	1.00	28.86
MOTA	671	N	PRO	89	19.363	12.930	48.084	1.00	14.78
ATOM	672	CA	PRO	89	18.552	13.475	47.008	1.00	14.80
ATOM	673	C	PRO	89	17.572	14.611	47.385		12.10
ATOM	674	ō	PRO	89 .	17.085	15.301	46.493	1.00	18.06
ATOM	675	CB	PRO	89	17.733	12.294	46.494	1.00	17.00
MOTA	676	CG	PRO	89	17.726	11.261	47.607	1.00	15.23
HOTA	677	CD	PRO	89	18.844	11.642	48.560	1.00	17.16
MOTA	678	N	GLU	90	17.278	14.795	48.695	1.00	14.63
ATOM	679	CA	GLU	90	16.348	15.838	49.157	1.00	20.68
ATOM	680	С	GLU	90	16.701	17.229	48.645	1.00	25.59
MOTA	681	ō	GLU	90	15.833	18.042	48.368		21.57
ATOM	682	СВ	GLU	90	16.031	15.816	50.682	1.00	22.21
ATOM	683	CG	GLU	90	15.782	14.403	51.228	1.00	37.59
MOTA	684	CD	GLU	90	17.071	13.641	51.447	1.00	83.49
MOTA	685	OE1	GLU	90	18.179	14.151	51.342	1.00	54.20
ATOM	686	OE2	GLU	90	16.875	12.373	51.749	1.00	64.65
MOTA	687	ห	GLY	91	17.977	17.509	48.510	1.00	21.39
MOTA	688	CA	GLY	91	18.394	18.769	47.906	1.00	17.77
ATCM	689	C	GLY	91	18.673	19.911	48.839	1.00	12.17
MOTA	690	0	GLY	91	18.769	19.764	50.055	1.00	16.81
ATOM	691	N	TYR	92	18.861	21.086	48.225	1.00	13.02
ATOM	692	CA	TYR	92	19.143	22.266	48.994	1.00	10.33
ATOM	693	С	TYR	92	18.575	23.478	48.347	1.00	9.87
MOTA	694	0	TYR	92	18.270	23.483	47.144	1.00	15.89
ATOM	695	CB	TYR	92	20.678	22.488	49.278	1.00	15.40
ATOM	696	CG	TYR	92	21.546	22.468	48.012	1.00	15.13
ATOM	697	CD1	TYR	92	21.620	23.576	47.166	1.00	14.75
MOTA	698	CD2	TYR	92	22.317	21.350	47.683	1.00	16.09
ATOM	699	CE1	TYR	92	22.404	23.561	46.006	1.00	6.50
MOTA	700	CE2	TYR	92	23.067	21.300	46.504	1.00	15.12
MOTA	701	CZ	TYR	92	23.156	22.424	45.683	1.00	18.13
MOTA	702	ОН	TYR	92	23.944	22.393	44.517	1.00	13.37
MOTA	703	N	VAL	93	18.447	24.504	49.189	1.00	11.93
MOTA	704	CA	VAL	63	18.025	25.822	48.778	1.00	14.74
ATOM	705	С	VAL	93	19.281	26.666	48.625	1.00	16.00
MOTA	706	0	VAL	93	20.172	26.625	49.451	1.00	15.16
MOTA	707	CB	VAL	93	17.073	26.480	49.791	1.00	23.45
MOTA	708	CG1	VAL	93	16.855	27.937	49.413	1.00	26.05
MOTA	709	CG2	VAL	93	15.716	25.764	49.771	1.00	22.90
atom	710	N	GLN	94	19.361	27.345	47.521	1.00	13.78
ATOM	711	CA	GLN	94	20.480	28.195	47.227	1.00	10.53
MOTA	712	С	GLN	94	19.948	29.583	46.998	1.00	12.23
MOTA	713	0	GLN	94	19.153	29.788	46.061		15.52
MOTA	714	CB	GLN	94	21.232	27.727	45.934	1.00	7.95
MOTA	715	CG	GLN	94	22.361	29.708	45.469	1.00	11.37
HOTA	716	CD	GLN	94	23.431	27.999	44.632	1.00	12.04
MOTA	717		GLN	94	23.805	26.879	44.946	1.00	13.60
MOTA	718	NE2		94	23.719	28.527	43.449	1.00	7.98
ATOM	719	N	,GTA	95	20.396	30.531	47.820	1.00	11.78
ATOM.	720	CA	GLU	95	19.974	31.899	47.643	1.00	13.47
MOTA	721	C	GLU	95	21.149	32.804	47.398	1.00	
MOTA	722	0	GLU	95 35	22.205	32.623	17.985		19.23
MOTA	723	CB	GLU	95 35	19.277	32.427	48.878	1.00	13.52
ATOM	724	CG	GLU	95	18.009	31.684	49.215	00	22.46 45.73
MOTA	725	CD	GLU	95 35	17.657	32.016	50.622		
HOTA	726		GLU	95 26	17.574	33.166	51.011		100.30
MOTA	727	OE2		5.5 5.5	17.764	30.987	51.423	00	61.33
MOTA	728	::	ARG	96	20.929	33.838	46.601	:.00	15.51
MOTA	729	CA	ARG	76	11.978	34.783	46.342	1.00	
MOTA	730	С	ARG	₹6	11.510	36.195	46.206	1.00	15.54

19/36

MOTA	731	0	ARG	96	20.389	36.488	45.806	1.00 15.01
MOTA	732	СЗ	ARG	96	22.582	34.463	44.967	1.00 16.19
ATOM	733	CG	ARG	96	23.495	33.247	44.929	1.00 17.61
ATOM	734	CD	ARG	96	24.615	33.453	43.908	1.00 9.06
ATOM	735	ΝE	ARG	96	25.411	32.277	43.766	1.00 9.88
ATOM	736	СZ	ARG	96	25.434	31.493	42.693	1.00 20.03
ATOM	737	H1	ARG	96	24.684	31.709	41.615	1.00 15.29
ATOM	738	::H2	ARG	96	26.236	30.430	42.714	1.00 11.03
ATOM	739	x	THR	97	22.470	37.068	46.344	1.00 13.39
MOTA	740 -	CA	THR	97	22.368	38.424	45.935	1.00 13.12
ATOM	741	С	THR	97	23.593	38.688	45.084	1.00 16.81
MOTA	742	0	THR	97	24.686	38.347	45.485	1.00 19.25
ATOM	743	CB	THR	97	22.282	39.442	47.066	1.00 26.27
ATOM	744	0 <b>G</b> 1	THR	97	21.225	39.101	47.945	1.00 31.43
ATOM	745	CG2	THR	97	22.038	40.804	46.445	1.00 15.90
ATOM	746	27	ILE	98	23.396	39.219	43.899	1.00 16.23
ATOM	747	CA	ILE	98	24.486	39.526	42.977	1.00 16.70
ATOM	748	С	ILE	98	24.533	41.017	42.686	1.00 21.10
ATOM	749	0	ILE	98	23.628	41.566	42.075	1.00 14.58
ATOM	750	CB	ILE	98	24.385	38.752	41.660	1.00 13.47
MOTA	751	CG1		98	24.480	37.236	41.890	1.00 16.09
ATOM	752	CG2	ILE	98	25.457	39.231	40.679	1.00 13.30
MOTA	753	CD1	ILE	98	23.875	36.431	40.738	1.00 13.93
ATOM	754	:1	PHE	99	25.613	41.678	43.110	1.00 14.86
MOTA	755	CA	PHE	99	25.719	43.098	42.896	1.00 12.44
ATOM	756	С	PHE	99	26.514	43.441	41.699	1.00 20.37
MOTA	757	0	PHE	99	27.696	43.164	41.700	1.00 20.07
MOTA	758	CB	PHE	99	26.401	43.770	44.084	1.00 15.96
ATOM	759	CG	PHE	99	25.638	43.624	45.356	1.00 21.41
ATOM	760		PHE	9 <b>9</b>	25.863	42.524	46.189	1.00 24.98
ATOM	761	CD2	PHE	99	24.698	44.585	45.743	1.00 22.94
ATOM	762	CEl	PHE	99	25.176	42.400	47.400	1.00 32.06
ATOM	763	CE2	PHE	99	23.992	44.469	46.946	1.00 24.26
ATOM	764	CZ	PHE	99	24.235	43.369	47.771	1.00 28.19
MOTA	765	X	PHE	100	25.906	44.085	40.704	1.00 12.53
MOTA	766	CA	PHE	100	26.679	44.522	39.554	1.00 8.75
MOTA	767	С	PHE	100	27.294	45.855	39.872	1.00 21.81
MOTA	768	O	PHE	100	26.599	45.775	40.308	1.00 20.31
ATOM	769	C3	PHE	100	25.927	44.572	38.226	1.00 5.94
MOTA	770	CG	PHE	100	25.537	43.183	37.764	1.00 12.75
ATOM	771	CD1		100	24.426	42.538	38.325	1.00 16.31
ATOM	772	CD2		100	26.317	42.484	36.843	1.00 15.27
MOTA	773	CEI		100	24.087	41.230	37.975	1.00 13.50
MOTA	774	CE2		100	25.965	41.192	36.435	1.00 21.25
MOTA	775	CZ	PHE	100	24.852	40.567	37.014	1.00 21.06
MOTA	776	!!	LYS	101	28.603	45.946	39.737	1.00 15.49
ATOM	777	CA	LYS	101	29.270	47.179	40.085	1.00 17.93
ATOM	778	C	LYS	101	28.732	48.349	39.287	1.00 13.71
ATOM	779	0	LYS	101		48.304	38.072	1.00 17.18
MOTA	780	CB	LYS	101 101	30.784	47.069	39.950	1.00 17.13
MOTA	781	CG	LYS	101	31.518	48.252	40.551 40.534	1.00 18.01
MOTA.	782	CD	LYS		33.036	48.060		1.00 26.70
MOTA	783	CE	LYS	101 102	33.797 28.353	49.116	41.332	
ATOM	784	::	ASP	102		49.403 50.618		1.00 18.09
ATOM	785	Cλ	ASP	102	27.805		39.368	1.00 23.08
ATOM	786	C	`ASP		26.559	50.356	33.549	1.00 25.42
ATOM	787	0	ASP	102 102	26.292	£1.061 £1.369	37.586	1.00 23.34
MOTA	788	C3	ASP		28.840	51.629	38.516	
MOTA	789	CG	ASP	102 102	30.109 31.206	51.233	39.296 38.931	1.00 57.01
ATOM	790			102	29.886	52.200	40.464	1.00 63.33
MOTA	791 792	::	ASP ASP	102	25.813	49.328	33.933	1.00 20.17
MOTA MOTA	793	 CA	ASP	103	24.602	43.949	33.233	1.00 20.17
ATOM	794	C	ASP	103	23.608	43.224	39.189	1.00 12.70
MOTA	795	0	ASP	103	23.749	43.431	+0.409	1.00 17.72
ATOM	796	СЗ	ASP	103	24.899	43.431	16.995	1.00 17.72
ATOM	797	CC	ASP	203	23.946	48.387	35.860	1.00 23.93
VIOW	, , ,	-0	r		23.540	40.55/		2.00 22.73

FIG 5-13

ATOM	798	OD1		103	24.238	48.274	34.688	1.00 19.05
MOTA	799	OD2		103	22.774	48.809	35.283	1.00 23.89
MOTA	300	31	GLY	104	22.612	47.542	38.646	1.00 20.17
MOTA	801	CA	GLY	104	21.598	46-900	39.498	1.00 20.22
MOTA	802	С	GLY	104	22.055	45.619	40.180	1.00 24.68
MOTA	803	0	GLY	104	23.202	45.211	40.085	1.00 18.06
MOTA	<b>SO4</b>	:1	ASN	105	21.125	44.967	40.872	1.00 15.71
ATOM	305	CA	ASN	105	21.425	43.703	41.510	1.00 8.89
ATOM	806	С	ASN	105	20.399	42.620	41.181	1.00 21.85
ATOM	807	0	ASN	105	19.255	42.911	40.824	1.00 15.17
MOTA	808	CB	ASN	105	21.605	43.840	43.001	1.00 8.58
ATOM	809	CG	ASN	105	20.359	44.366	43.697	1.00 43.57 1.00 36.67
MOTA	810		ASN	105	19.565	43.601	44.259	1.00 36.67 1.00 36.47
MOTA	811		ASN	105	20.178	45.674	43.659 41.328	1.00 16.80
MOTA	812	N	TYR	106	20.826 19.966	41.365 40.219	41.156	1.00 13.90
MOTA	813	CA	TYR	106 106	19.763	39.543	42.475	1.00 11.05
ATOM	814	C	TYR	106	20.678	39.404	43.281	1.00 13.86
ATOM	815 816	O CB	TYR TYR	106	20.547	39.128	40.246	1.00 15.88
ATOM ATOM	817	CG	TYR	106	20.619	39.398	38.793	1.00 15.57
MOTA	818	CD1		106	19.952	40.458	38.178	1.00 13.14
ATOM	819	CD2		106	21.373	38.524	38.006	1.00 13.35
ATOM	820	CE1	TYR	106	20.038	40.632	26.793	1.00 13.44
ATOM	821	CE2	TYR	106	21.481	38.692	36.628	1.00 10.87
ATOM	822	cz	TYR	106	20.814	39.751	36.025	1.00 15.93
MOTA	823	ОН	TYR	106	20.970	39.931	34.670	1.00 17.32
ATOM	824	H	LYS	107	18.538	39.115	42.709	1.00 12.39
ATOM	825	CA	LYS	107	18.194	38.349	43.897	1.00 11.51
ATOM	826	С	LY5	107	17.619	37.037	43.397	1.00 17.25
ATOM	827	0	LYS	107	16.704	37.010	42.562	1.00 13.14
MOTA	828	CB	LYS	107	17.217	39.063	44.823	1.00 14.82
ATOM	829	CG	LYS	107	17.860	39.631	46.060	1.00 40.71
ATOM	830	CD	LYS	107	18.528	40.974	45.793	1.00 43.48
ATOM	831	N	THR	108	18.205	35.951	43.835	1.00 14.95
MOTA	832	CA	THR	108	17.774	34.658	43.352	1.00 11.97
ATOM	833	С	THR	108	17.463	33.696	44.468	1.00 15.81
MOTA	234	0	THR	108	18.043	33.734	45.582	1.00 19.68
ATOM	835	C3	THR	108	18.847	34.034	42.410 43.137	1.00 23.81
ATOM	836	OGI		108	20.064 19.123	33.791 34.968	41.264	1.00 13.88
ATOM	837	CG2		108 109	16.560	32.804	44.154	1.00 13.57
ATOM	838	:: CA	ARG	109	16.212	31.751	45.048	1.00 12.56
ATOM	839 840	C	ARG	109	15.939	30.498	44.254	1.00 13.07
MOTA MOTA	841	0	ARG	109	15.239	30.509	43.249	1.00 12.52
ATOM	842	CB	ARG	109	15.069	32.100	45.959	1.00 17.32
MOTA	843	CG	ARG	109	14.767	30.995	46.932	1.00 17.92
ATOM	844	CD	λRG	109	13.400	31.160	47.610	1.00 19.99
MOTA	845	HE	ARG	109	12.821	29.854	47.883	1.00 36.05
ATOM	846	CZ	ARG	109	12.968	29.244	49.035	1.00 55.71
ATOM	847	NHI	ARG	109	13.630	29.815	50.046	1.00 44.11
ATOM	848	NH2	ARG	109	12.432	28.041	49.195	1.00 94.34
MOTA:	849	:1	ALA	110	16.577	29.414	44.635	1.00 13.26
MOTA	850	CA	ALA	110	16.377	28.207	43.870	1.00 12.68
MOTA	851	С	ALA	:10	16.346	26.979	44.734	1.00 13.15
MOTA	852	0	ALA	110	16.829	26.965	45.869	1.00 16.75
MOTA	853		, YTY	110	17.465	28.059	42.822	1.00 17.31
MOTA	354	::	GLU	111	15.770	25.939	44.175	1.00 15.39
MOTA	855	CA	GLU	111	15.741	24.655	44.823	1.00 15.24
MOTA	856	С	GLU	111	16.438	23.678	43.925	1.00 12.08
ATOM	857	0	GLU	111	16.086	23.545	42.771	1.00 15.70
MOTA	858	CB	GLU	111	14.303	24.123	44.993	1.00 19.20
ATOM	359 860	CO	GLU GLU	111 111	13.744 12.247	24.242	46.399 46.372	1.00 33.62
MOTA MOTA	861		L GIU		11.539	23.843	45.432	1.00 76.05
ATOM	262	UE:	672	111 111	11.742	24.956	47.380	1.00 54.87
ATOM	363	::	. GIS	112	17.438	22.965	44.457	1.00 10.78
ATOM	364	CA	∵AL	::2	18.063	31.978	43.631	1.00 10.98

FIG 5-14

ATOM	565	С	VAL	112	17.968	20.630	44.261	1.00 8.62
HOTA	866	0	VAL	112	18.271	20.438	45.432	1.00 15.63
MOTA	867	CB	VAL	112	19.428	22.358	43.012	1.00 22.75
ATOM	868	CG1	VAL	112	19.966	23.704	43.487	1.00 16.69
MOTA	869	CG2	VAL	112	20.452	21.232	43.078	1.00 18.47 1.00 14.67
MOTA	870	N	LYS	113	17.415	19.732	43.516	
MOTA	871	CA	LYS	113	17.175	18.421	44.045	
MOTA	872	С	LYS	113	16.822	17.485	42.931	1.00 7.11
MOTA	973	0	LYS	113	16.695	17.893	41.808	1.00 16.27 1.00 22.50
ATOM	874	CB	LYS	113	16.032	18.497	45.036	1.00 20.40
MOTA	875	CG	LYS	113	14.792	19.084	44.376	1.00 44.65
MOTA	876	CD	LYS	113	13.509	18.321	44.703	1.00 44.83
MOTA	877	CE	LYS	113	12.526	19.134	45.528	1.00100.00
ATOM	878	ΝZ	LYS	113	12.379	20.518	45.036 43.267	1.00 10.09
ATOM	879	N	PHE	114	16.683	16.208	43.207	1.00 11.41
MOTA	880	CA	PHE	114	16.325	15.175 14.975	42.181	1.00 14.18
MOTA	881	C	PHE	114	14.806 14.110	14.878	43.160	1.00 15.03
MOTA	882	0	PHE	114 114	16.866	13.838	42.838	1.00 12.89
ATOM	883	CB	PHE	114	18.231	13.536	42.338	1.00 16.80
MOTA	884	CG	PHE PHE	114	19.344	13.795	43.139	1.00 18.61
MOTA	885	CD1	PHE	114	18.403	13.009	41.056	1.00 19.50
ATOM	886		PHE	114	20.627	13.500	42.665	1.00 22.78
MOTA	887 888	CE2	PHE	114	19.673	12.708	40.572	1.00 25.36
atom atom	889	CZ	PHE	114	20.780	12.953	41.387	1.00 23.99
ATCM	890	::	GLU	115	14.354	14.819	40.966	1.00 15.29
ATOM	891	CA	GLU	115	12.978	14.473	40.642	1.00 11.40
ATOM	892	c	GLU	115	13.121	13.193	39.906	1.00 13.30
ATOM	893	ō	GLU	115	13.434	13.207	38.730	1.00 18.72
ATOM	894	CB	GLU	115	12.348	15.481	39.667	1.00 9.68
ATOM	895	CG	GLU	115	11.856	16.747	40.376	1.00 19.54
ATOM	896	CD	GLU	115	10.742	16.460	41.342	1.00 38.12
ATOM	897	OE1	GLU	115	10.181	15.395	41.431	1.00 34.84
MOTA	898	OE2	GLU	115	10.460	17.461	42.079	1.00 27.88
MOTA	899	N	GLY	116	13.005	12.087	40.585	1.00 14.51
ATOM	900	CA	GLY	116	13.225	10.861	39.869	1.00 15.91
ATOM	901	С	GLY	116	14.727	10.767	39.641	1.00 23.59 1.00 19.35
ATOM	902	0	GLY	116	15.516	10.922	40.570	1.00 20.25
MOTA	903	N	ASP	117	15.137	10.564	38.439	1.00 20.25 1.00 28.00
atom	904	CA	ASP	117	16.572	10.462	38.233 37.598	1.00 22.39
MOTA	905	С	ASP	117	17.237	11.677 11.672	37.265	1.00 21.33
ATOM	906	0	ASP	117	18.423	9.074	37.733	1.00 33.06
MOTA	907	CB	ASP	117 117	17.055 16.624	8.677	36.348	1.00 55.04
ATOM	908	CG	ASP ASP	117	16.230	9.468	35.495	1.00 59.57
ATOM	909	ODI		117	16.805	7.391	36.130	1.00 82.48
ATOM	910 911	OD2 N	THR	118	16.463	12.729	37.493	1.00 19.62
ATOM ATOM	912	CA	THR	118	16.889	13.981	36.910	1.00 18.21
MOTA	913	C	THR	118	17.186	14.988	37.976	1.00 18.92
ATOM	914	ō	THR	118	16.498	15.064	38.996	1.00 15.94
ATOM	915	СВ	THR	118	15.806	14.497	35.952	1.00 19.03
MOTA	916	OG 1		118	15.552	13.508	34.990	1.00 21.42
ATOM	917		THR	118	16.217	15.793	35.275	1.00 15.49
ATOM	918	N	LEU	:19	18.284	15.6Bl	37.805	1.00 13.66
ATOM	919	CA	LEU	119	18.679	16.706	38.759	1.00 13.50
MOTA	920	С	LEU	119	18.036	17.992	38.269	1.00 8.81
ATOM	921	0	LEU	119	18.194	18.368	37.091	1.00 12.49
ATOM	922	СВ	LEU	119	20.243	:6.815	33.839	1.00 12.25
ATOM	923	CG	LEU	:19	20.845	17.678	39.951	1.00 3.90
ATOM	924		L LEU	119	20.701	19.167	39.669	1.00 10.11
ATOM	925		2 LEU	119	20.366	17.311	41.333	1.00 7.36
MOTA	926	11	VAL	120	17.230	13.595	39.170	1.00 13.34
MOTA	927	CA	<b>VAL</b>	120	16.466	19.797	38.859	1.00 13.77
MOTA	928	C	VAL	:20	16.929	21.039	39.587	1.00 8.56
ATOM	929	0	"AL	130	17.135	21.039	40.752	1.00 13.32
ATOM	930	C3	VAL	120	14.939	19.566	39.082	1.00 17.50
ATOM	931	CC	1 TAL	120	14.133	20.790	33.642	1.00 17.32

FIG 5-15

MOTA	932	CG2	VAL	120	14.501	18.351	38.246	1.00 15.35
MOTA	933	:1	ASN	121	17.067	22.111	38.839	1.00 12.24
ATOM	934	CA	ASN	121	17.424	23.405	39.400	1.00 11.78
ATOM	935	C	ASN	121	16.301	24.382	39.060	1.00 11.18
	936	ō	ASN	121	16.195	24.802	37.934	1.00 11.09
MOTA				121	18.753	23.928	38.791	1.00 11.41
MOTA	937	CB	ASN		19.201	25.261	39.367	1.00 11.07
HOTA	938	CG	ASN	121			40.461	1.00 12.06
MOTA	939	OD1	ASN	121	18.773	25.654		
·ATOM	940	ND2	ASN	121	20.124	25.938	38.670	1.00 11.90
MOTA	941	N	ARG	122	15.470	24.706	40.029	1.00 13.69
ATOM	942	CA	ARG	122	14.348	25.610	39.825	1.00 12.99
ATOM	943	С	ARG	122	14.622	26.946	40.498	1.00 5.89
MOTA	944	õ	ARG	122	14.749	27.011	41.723	1.00 14.47
ATOM	945	CB	ARG	122	13.068	25.025	40.417	1.00 15.99
ATOM	946	CG	ARG	122	12.478	23.921	39.589	1.00 30.23
	947	CD	ARG	122	11.282	23.244	40.281	1.00 60.61
ATOM			ILE	123	14.663	27.992	39.680	1.00 11.46
ATOM	948	N		123	15.030	29.340	40.095	1.00 11.86
ATOM	949	CA	ILE		13.991	30.450	39.835	1.00 10.54
MOTA	950	C	ILE	123			38.765	1.00 12.83
MOTA	951	0	ILE	123	13.370	30.535		1.00 15.41
MOTA	952	CB	ILE	123	16.296	29.757	39.292	
MOTA	953	CC1	ILE	123	17.316	28.585	39.180	1.00 12.27
MOTA	954	CG2	ILE	123	16.944	30.993	39.918	1.00 14.01
MOTA	955	CD1	ILE	123	17.652	28.242	37.743	1.00 7.74
ATOM	956	N	GLU	124	13.953	31.358	40.793	1.00 11.36
MOTA	957	CA	GLU	124	13.189	32.572	40.700	1.00 15.20
ATOM	958	С	GLU	124	14.168	33.713	40.811	1.00 11.93
ATOM	959	ō	GLU	124	14.919	33.797	41.780	1.00 15.61
ATOM	960	СВ	GLU	124	12.028	32.677	41.751	1.00 19.74
ATOM	961	CG	GLU	124	12.387	33.337	43.089	1.00 72.94
ATOM	962	N	LEU	125	14.183	34.550	39.808	1.00 12.19
	963	CA	LEU	125	15.092	35.654	39.767	1.00 15.00
ATOM				125	14.420	37.011	39.722	1.00 19.35
YLOW	964	C	LEU	125	13.563	37.267	38.893	1.00 18.41
ATOM	965	0	LEU		15.976	35.533	38.510	1.00 14.29
ATOM	966	CB	LEU	125		36.683	38.375	1.00 17.55
MOTA	967	CG	LEU	125	17.003		37.849	1.00 13.46
ATOM	968	CD1		125	18.302	36.023 37.732	37.367	1.00 12.09
ATOM	969	CD2		125	16.511		40.554	1.00 12.73
ATOM	970	N	LYS	126	14.890	37.897	40.579	1.00 12.73 1.00 15.92
ATOM	971	CA	LYS	126	14.391	39.260		1.00 13.52
ATOM	972	С	LYS	126	15.563	40.276	40.445	1.00 19.86
ATOM	973	0	LYS	126	16.489	40.246	41.246	
MOTA	974	CB	LYS	126	13.611	39.487	41.877	1.00 17.31
ATOM	975	ÇG	LYS	126	12.853	40.786	41.923	1.00 33.94
ATOM	976	CD	LYS	126	11.356	40.601	41.675	1.00 60.87
MOTA	977	CE	LYS	126	10.652	41.929	41.521	1.00 52.70
ATOM	978	NZ	LYS	126	11.229	42.988	42.367	1.00 47.22
ATOM	979	N	GLY.	127	15.514	41.127	39.411	1.00 18.71
ATOM	980	CA	GLY	127	16.551	42.151	39.121	1.00 17.32
ATOM	981	С	GLY	127	16.012	43.572	39.272	1.00 25.32
· ATOM	982	0	GLY	127	14.981	43.908	38.693	1.00 20.14
ATOM	983	N	ILE	128	16.706	14.404	40.070	1.00 18.42
ATOM	984	CA	ILE	128	16.282	45.787	40.243	1.00 21.04
ATOH	985	c	ILE	128	17.405	46.789	40.196	1.00 25.93
HOTA	986	ō	ILE	128	18.562	46.496	40.429	1.00 19.37
	987	СВ	ILE	128	15.482	46.052	41.504	1.00 23.82
MOTA				128	16.408	45.888	42.701	1.00 23.86
MOTA	788	CGI			14.272	45.120	41.577	1.00 28.95
MOTA	989	CG2		128			44.013	1.00 29.89
MOTA	990	CDI		128	15.824	46.391		1.00 29.29
MOTA	991	11	ASP	129	15.999	48.002	39.918	1.00 20.25
ATOM	992	CA	ASP	129	17.861	49.124	29.882	1.00 12.53
MOTA	993	С	ASP	129	19.864	49.086	38.801	1.00 20.35
MOTA	994	0	ASP	129	19.949	49.632	38.953	1.00 24.22
MOTA	<del>99</del> 5	CB	ASP	129	13.498	49.407	41.253	1.00 20.57
MOTA	996	CG	ASP	129	17.545	50.077	42.226	1.00 43.70
NOTA	997	OD:	1 ASP	129	16.653	50.842	41.883	1.00 49.42
ATOM	398		2 ASP	129	17.770	49.740	<b>∔3.475</b>	1.00 32.07
<del></del>			-					

FIG 5-16

ATOM	999	::	PHE	130	18.510	48.493	37.693	1.00 16.40
ATOM	1000	CA	PHE	130	19.433	48.459	36.563	1.00 16.99
ATOM	1001	C	PHE	130	19.330	49.732	35.756	1.00 35.37
ATOM	1002	õ	PHE	130	18.242	50.318	35.623	1.00 27.34
ATOM	1003	СЗ	PHE	130	19.248	47.223	35.657	1.00 18.07
ATOM	1004	CG	PHE	130	19.809	45.980	36.312	1.00 19.10
ATOM	1005	CDI		130	19.021	45.210	37.171	1.00 16.15
ATOM	1006	CD2		130	21.126	45.572	36.073	1.00 19.17
ATOM	1007	CEI		130	19.536	44.074	37.801	1.00 23.37
ATOM	1008	CE2		130	21.665	44.445	36.703	1.00 21.11
MOTA	1009	cz	PHE	130	20.867	43.703	37.575	1.00 22.13
ATOM	1010	37	LYS	131	20.464	50.169	35.218	1.00 31.09
ATOM	1011	CA	LYS	131	20.477	51.371	34.400	1.00 27.52
MOTA	1012	c	LYS	131	20.105	51.045	32.992	1.00 25.57
ATOM	1013	o	LYS	131	20.695	50.169	32.343	1.00 22.97
ATOM	1014	СВ	LYS	131	21.796	52.109	34.438	1.00 32.64
ATOM	1015	CG	LYS	131	22.153	52.633	35.813	1.00 38.34
ATOM	1016	CD	LYS	131	23.646	52.886	35.975	1.00 75.76
ATOM	1017	11	GLU	132	19.116	51.751	32.509	1.00 26.88
ATOM	1018	CA	GLU	132	18.623	51.484	31.189	1.00 28.42
ATOM	1019	c	GLU	132	19.710	51.514	30.140	1.00 36.19
ATOM	1020	ō	GLU	132	19.617	50.862	29.101	1.00 39.24
MOTA	1021	C3	GLU	132	17.374	52.331	30.830	1.00 29.04
ATOM	1022	:1	ASP	133	20.752	52.254	30.438	1.00 40.08
ATOM	1023	CA	ASP	133	21.883	52.442	29.525	1.00 45.36
MOTA	1024	С	ASP	133	23.224	51.861	30.049	1.00 50.61
ATOM	1025	ō	ASP	133	24.299	52.243	29.572	1.00 52.14
HOTA	1026	СЗ	ASP	133	22.063	53.946	29.332	1.00 50.45
ATOM	1027	CG	ASP	133	22.109	54.642	30.670	1.00 87.10
MOTA	1028	ODl	ASP	133	21.408	54.314	31.624	1.00 91.27
ATOM	1029	OD2	ASP	133	23.047	55.552	3 <b>0.739</b>	1.00100.00
MOTA	1030	N	GLY	134	23.159	50.970	31.053	1.00 37.06
MOTA	1031	CA	GLY	134	24.349	50.376	31.639	1.00 30.22
MOTA	1032	С	GLY	134	24.845	49.228	30.803	1.00 23.10
MOTA	1033	0	GLY	134	24.360	48.990	29.685	1.00 19.23
MOTA	1034	N	ЛSИ	135	25.807	48.486	31.341	1.00 18.66
MOTA	1035	CA	ASN	135	25.339	47.370	30.563	1.00 18.03
MOTA	1036	С	ASN	135	25.372	46.199	30.406	1.00 15.75
ATOM	1037	0	λSN	135	25.485	45.430	29.461	1.00 16.03
ATOM	1038	CB	ASN	135	27.665	46.883	31.139	1.00 19.27
MOTA	1039	CG	ASN	135	28.743	47.943	31.108	1.00 20.99
MOTA	1040	OD1		135	28.969	48.595	30.078	1.00 25.69
MOTA	1041	:iD2		135	29.423	48.095	32.239	
MOTA	1042	11	ILE	136	24.444	46.052	31.362	1.00 18.14
ATOM	1043	CA	ILE	136	23.494	44.924	31.368 30.384	
ATOM	1044 1045	C	ILE	136 136	22.331 22.178	45.086	29.395	1.00 23.76
MOTA MOTA	1045	O CB	ILE	136	23.078	44.313	32.804	1.00 21.24
	1045		ILE	136	24.230	43.728	33.423	1.00 28.44
ATOM	1048		ILE	136	21.899	43.543	32.770	1.00 22.77
MOTA MOTA	1049		ILE	136	25.346	44.595	33.935	1.00 12.39
ATOM	1050	:1	LEU	137	21.543	46.117	30.640	1.00 18.21
ATOM	1051	CA	LEU		20.394	46.415	29.815	1.00 23.30
MOTA	1052	Ç	LEU	137	20.828	46.875	28.470	1.00 27.25
ATOM	1053	ŏ	LEU	137	20.181	46.619	27.488	1.00 27.00
ATOM	1054		, LEU		19.442	47.430	30.490	1.00 21.74
ATOM	1055	CG	LEU	137	18.828	46.852	31.762	1.00 22.56
MOTA	1056		LEU		17.856	47.837	32.415	1.00 22.27
MOTA	1057		LEU		18.113	45.554	31.424	1.00 37.52
ATON	1058	::	GLY		21.979	47.527	28.432	1.00 22.14
ATOM	:059	CA	GLY		22.510	48.C33	27.187	1.00 20.03
ATOM	1060	5	GLY	138	23.157	46.959	26.368	1.00 20.16
MOTA	1061	õ	GLY		23.600	47.202	25.264	1.00 22.44
ATOM	1062	::	HIS	139	23.246	45.756	25.903	1.00 19.27
MOTA	1063	CA	HIS	:39	23.859	44.655	26.148	1.00 20.24
NOTA	1064	:	HIS		25.301	44.929	25.616	1.00 20.13
ATOM	1065	5	HIS		25.605	44.745	24.439	1.00 17.97

24/36

FIG 5-17

MOTA	1066	CB	HIS	139	22.931	44.207	15.018	1.00 22.20
ATOM	1067	CG	HIS	139	21.708	43.551	25.550	1.00 25.52
ATOM	1068	ND1		139	21.666	42.182	25.785	1.00 25.67
ATOM	1069	CD2		139	20.525	44.092	25.927	1.00 28.09
	1070	CEI		139	20.474	41.918	25.275	1.00 27.50
HOTA				139	19.766	43.044	26.382	1.00 29.53
ATOM	1071	-	HIS	140	26.187	45.311	26.525	1.00 23.51
ATOM	1072	N	LYS		27.569	45.638	26.197	1.00 25.82
MOTA	1073	CA	LYS	140			26.560	1.00 26.28
MOTA	1074	С	LYS	140	28.600	44.537	26.391	1.00 22.29
MOTA	1075	0	LYS	140	29.824	44.730	26.911	1.00 27.56
MOTA	1076	CB	LYS	140	27.977	46.937		1.00 27.30
MOTA	1077	CG	LYS	140	27.269	48.217	26.445	
ATOM .	1078	CD	LYS	140	27.234	49.254	27.582	1.00 51.32
ATOM	1079	CΞ	LYS	140	26.924	50.696	27.169	1.00 47.92
ATOM	1080	NZ	LYS	140	27.112	51.663	28.284	1.00 73.76
ATOM	1081	N	LEU	141	28.116	43.403	27.115	1.00 19.33
ATOM	1082	CA	LEU	141	28.987	42.296	27.559	1.00 14.32
ATOM	1083	С	LEU	141	29.366	41.401	26.427	1.00 20.75
MOTA	1084	0	LEU	141	28.526	41.087	25.620	1.00 19.01
MOTA	1085	CB	LEU	141	28.313	41.488	28.676	1.00 12.53
ATOM	1086	CG	LEU	141	27.979	42.352	29.875	1.00 17.54
ATOM	1087		LEU	141	27.700	41.469	31.070	1.00 24.81
ATOM	1088		LEU	141	29.116	43.310	30.182	1.00 27.50
MOTA	1089	N	GLU	142	30.644	40.927	25.346	1.00 14.76
ATOM	1090	CA	GLU	142	31.040	40.059	25.311	1.00 13.43
ATOM	1091	C	GLU	142	30.462	38.691	25.641	1.00 15.69
	1092	ò	GLU	142	30.175	38.393	26.787	1.00 16.43
ATOM	1092	CB	GLU	142	32.558	39.866	25.204	1.00 14.73
ATOM	1094	CG	GLU	142	33.290	41.077	24.624	1.00 29.30
MOTA	1095	CD	GLU	142	34.787	41.003	24.825	1.00 56.32
ATOM		OE1		142	35.340	40.098	25.420	1.00 31.70
MOTA	1096	OE2		142	35.430	42.015	24.321	1.00 34.10
MOTA	1097			142	30.365	37.873	24.632	1.00 16.30
ATOM	1098	N.	TYR		29.837	36.542	24.764	1.00 20.04
ATOM	1099	CA	TYR	143	30.925	35.559	25.049	1.00 12.46
ATOM	1100	C	TYR	143		34.792	24.193	1.00 16.99
MOTA	1101	0	TYR	143	31.327			1.00 20.96
MOTA	1102	C3	TYR	143	29.035	35.113	23.498	
ATOM	1103	CG	TYR	143	28.187	34.857	23.674	
ATOM	1104	CD1		143	27.040	34.859	24.472	
atom	1105	CD2		143	28.512	33.684	22.986	
MOTA	1106	CE1		143	26.257	33.708	24.615	
ATOM	1107	CE2		143	27.735	32.530	23.104	1.00 16.58
MOTA	1108	CZ	TYR	143	26.603	32.551	23.914	1.00 17.35
MOTA	1109	ОН	TYR	143	25.861	31.432	24.035	1.00 23.40
ATOM	1110	N	ASN	144	31.392	35.597	26.251	1.00 12.40
ATOM	1111	CA	ASN	144	32.428	34.703	25.689	1.00 12.05
ATOM	1112	С	ASN	144	32.433	34.675	23.193	1.00 15.75
ATOM	1113	0	ASN	144	31.637	35.369	22.837	1.00 14.58
MOTA	1114	CB	ASN	144	33.823	35.038	25.068	1.00 18.45
MOTA	1115	CG	ASN	144	34.310	36.445	25.374	1.00 18.98
ATOM	:116	OD:	ASN	:44	34.150	36.951	27.488	1.00 20.34
ATOM	1117	ND2	N2A S	144	34.891	37.085	25.382	1.00 23.02
ATOM	1118	N	TYR	145	33.311	33.876	22.773	1.00 12.16
MOTA	1119	CA	TYR	145	33.343	33.765	30.195	1.00 10.63
ATOM	1120	С	TYR	145	34.765	33.458	30.730	1.00 14.58
ATOM	1121	0	TYR	145	35.510	32.751	30.090	1.00 19.83
MOTA	1122	CB	TYR	145	32.404	32.627	30.571	1.00 9.76
MOTA	1123	CG	TYR	145	31.698	32.916	31.826	1.00 11.86
MOTA	1124	CD:		145	30.515	33.658	31.208	1.00 9.04
ATOM	1125	CD		145	32.188	32.419	33.030	1.00 10.07
ATOM	1126	CE		145	29.860	33.948	32.999	1.00 2.36
ATOM	1127	CE		145	31.544	32.707	34.235	1.00 15.32
ATOM	1128	CZ.	TYR	145	30.375	33.469	34.206	1.00 11.69
•	1129	OH	TYR	145	29.730	23.735	33.376	1.00 15.23
ATOM	1130	::	ASN	146	35.086	23.931	31.933	1.00 15.36
ATOM	1131	CA	ASN	146	36.415	33.737	32.560	1.00 17.00
ATOM					36.426	32.618	33.539	1.00 17.68
MOTA	1132	С	ASN	146	30.420	22.516	7	00 17.00

FIG 5-18

MOTA	1133	0	ASN	146	35.395	32.043	33.848		14.71
MOTA	1134	CB	ASN	146	36.844	35.062	33.235		11.89
ATOM	1135	CG	ASN	146	37.013	36.147	32.215		35.45
MOTA	1136	OD1	ASN	146	37.533	35.890	31.105		31.63
ATOM	1137	ND2	ASN	146	36.547	37.349	32.553		19.74
MOTA	1138	N	SER	147	37.630	32.338	34.201		12.09
MOTA	1139	CA	SER	147	37.804	31.320	35.266 36.575	1.00	8.55
MOTA	1140	C	SER	147	37.769	31.999	36.671		16.56
ATOM	1141	0	SER	147 147	38.219 39.148	33.125 30.540	35.129	1.00	9.87
ATOM	1142 1143	CB OG	SER SER	147	39.212	29.980	33.828		33.20
MOTA MOTA	1144	N	HIS	148	37.195	31.365	37.583	1.00	5.53
ATOM	1145	CA	HIS	148	37.090	31.998	38.850	1.00	8.06
ATOM	1146	C	HIS	148	37.346	31.038	39.949		11.30
ATOM	1147	ō	HIS	148	37.328	29.844	39.754	1.00	16.87
ATOM	1148	CB	HIS	148	35.648	32.608	39.067	1.00	11.29
ATOM	1149	CG	HIS	148	35.215	33.554	37.972		10.84
ATOM	1150	ND1	HIS	148	34.548	33.121	36.836		12.77
MOTA	1151	CD2	HIS	148	35.403	34.887	37.851	1.00	8.82
MOTA	1152	CE1		148	34.389	34.178	36.060	1.00	8.84
MOTA	1153	NE2	HIS	148	34.882	35.242	36.647	1.00	8.82
MOTA	1154	N	ASN	149	37.534	31.579	41.125	1.00	
MOTA	1155	CA	ASN	149	37.626	30.805	42.345	1.00	13.35
ATOM	1156	C	ASN	149	36.409	31.157	43.205	1.00	14.47
MOTA	1157	0	ASN	149	36.099	32.320 31.093	43.327 43.184	1.00	12.67
MOTA MOTA	1158 1159	CB	NSA NSA	149 149	38.890 40.148	30.822	42.424	1.00	20.21
MOTA	1160	001		149	40.148	31.713	42.281	1.00	
ATOM	1161	ND2	ASN	149	40.210	29.641	41.818	1.00	16.44
ATOM	1162	N	VAL	150	35.773	30.144	43.741	1.00	14.65
ATOM	1163	CA	VAL	150	34.588	30.262	44.552	1.00	12.92
MOTA	1164	С	VAL	150	34.910	29.806	45.943	1.00	16.30
MOTA	1165	0	VAL	150	35.257	28.665	46.147	1.00	
ATOM	1166	CB	VAL	150	33.482	29.382	43.914	1.00	
MOTA	1167	CG1		150	32.252	29.297	44.765	1.00	
MOTA	1168	CG2		150	33.172	29.791	42.464	1.00	
MOTA	1169	N	TYR	151	34.796	30.716	46.900	1.00	17.64
ATOM	1170	CA	TYR	151	35.139	30.440	48.275	1.00	18.31
ATOM	1171	C	TYR	151 151	34.003	29.917 30.536	49.117 49.239	1.00	
ATOM ATOM	1172 1173	O CB	TYR TYR	151	32.963 35.793	31.681	48.920	1.00	
ATOM	1174	CG	TYR	151	37.025	32.033	48.141		25.B6
MOTA	1175	CD1		151	37.003	32.989	47.127		26.00
ATOM	1176	CD2		151	38.200	31.315	48.355		28.66
ATOM	1177	CE1		151	38.151	33.234	46.369	1.00	33.73
MOTA	1178	CE2		151	39.360	31.550	47.619	1.00	29.01
ATOM	1179	CZ	TYR	151	39.325	32.512	46.618	1.00	29.55
ATOM	1180	OH	TYR	151	40.449	32.737	45.877		38.69
ATOM	1181	N	ILE	152	34.250	28.791	49.753	1.00	17.71
ATOM	1182	CA	ILE	152	33.255	28.159	50.572		14.12
ATOM	1183	С	ILE	152	33.619	28.056	52.000		18.51
MOTA	1184	0	ILE	152	34.728	27.703 26.776	52.336		22.05
ATOM	1185 1186	CB CG1	ILE	152 152	32.979 32.431	26.875	50.060 48.638		16.66 11.30
ATOM ATOM	1187		ILE	152	32.431	26.078	51.021		17.96
ATOM	1188		, ITE	. 52	32.377	25.559	47.949		13.48
ATOM	1189	N	MSE	152 153	32.623	28.278	52.841		
ATOM	1190	CA	MSE	:53	32.789	28.162	54.269	1.00	22.61
ATOM	1191	c	MSE	153 153	31.534	27.648	54.916		27.31
ATOM	1192	ō	MSE	153	30.433	27.831	54.396	1.00	20.50
MOTA	:193	CB	MSE	153	33.145	29.490	54.855	1.00	19.11
ATOM	1194	CG	MSE	153	34.010	30.302	53.957	1.00	100.00
ATOM	1195	SE	MSE	153	34.060	32.117	54.524		100.00
MOTA	1196	CE	MSE	153 154	33.463	31.798	56.330	1.00	30.27
ATOM	1197	N	ALA	154	31.733	25.983	56.053	1.00	22.29
ATOM	1198	CA	ALA	154	30.669	25.389	56.796		22.56
ATOM	:199	С	ALA	154	29.820	27.401	57.552	1.00	29.00

FIG 5-19

ATOM	1200	0	ALA	154	30.274	28.457	57.960	1.00 27.02
ATOM	1201	CB	ALA	154	31.224	25.336	57.744	1.00 19.78
MOTA	1202	N	ASP	155	28.566 27.669	27.063 27.887	57.726 58.484	1.00 29.48
ATOM	1203	CA	ASP	155 1 <b>5</b> 5	26.976	27.019	59.511	1.00 44.51
MOTA MOTA	1204 1205	0	ASP ASP	155	25.898	26.492	59.274	1.00 39.55
ATOM	1206	СВ	ASP	155	26.659	28.617	57.597	1.00 31.70
ATOM	1207	CG	ASP	155	26.140	29.851	58.247	1.00 49.89
ATOM	1208	OD1	ASP	155	26.595	30.297	59.277	1.00 46.67
MOTA	1209	OD2	ASP	155	25.187	30.422	57.565	1.00 76.07
MOTA	1210	N	LYS	156	27.646	26.816	60.629	1.00 46.37 1.00 53.23
MOTA	1211	CA	LYS	156	27.116 25.750	25.954 26.369	61.654	1.00 65.62
MOTA	1212 1213	c o	LYS LYS	156 156	25.012	25.520	62.703	1.00 65.54
ATOM ATOM	1214	CB	LYS	156	28.147	25.612	62.725	1.00 59.51
ATOM	1215	N	GLN	157	25.398	27.655	62.138	1.00 68.32
ATOH	1216	CA	GLN	157	24.119	28.135	62.670	1.00 73.00
MOTA	1217	С	GLN	157	22.891	27.767	61.817 62.325	1.00 87.53
MOTA	1218	0	GLN	157	21.778 23.095	27.547 27.725	60.506	1.00 72.49
ATOM	1219	N	LYS LYS	158 158	22.040	27.725	59.593	1.00 66.19
MOTA MOTA	1220 1221	CA C	LYS	158	22.235	25.985	59.040	1.00 58.21
MOTA	1222	ō	LYS	158	21.447	25.524	58.226	1.00 59.85
HOTA	1223	N	ASN	159	23.303	25.294	59.502	1.00 40.00
MOTA	1224	CA	ASN	159	23.582	23.944	59.012	1.00 36.67
HOTA	1225	С	ASN	159	23.755	24.002	57.500 56.754	1.00 34.11
MOTA	1226	0_	ASN	159	23.223 22.431	23.167 22.952	59.367	1.00 46.42
ATOM	1227 1228	CB CG	asn Asn	159 159	22.842	21.485	59.428	1.00 80.46
ATOM ATOM	1229	OD1		159	23.850	21.121	60.054	1.00100.00
ATOM	1230		ASN	159	22.003	20.620	58.854	1.00 58.09
MOTA	1231	N	GLY	160	24.474	25.044	57.062	1.00 22.34
MOTA	1232	CA	GLY	160	24.686	25.247	55.663	1.00 17.58
MOTA	1233	С	GLY	160	26.055	25.791 25.664	55.433 56.271	1.00 26.75
ATOM	1234	N 0	GLY	160 161	26.960 26.200	26.395	54.277	1.00 23.23
MOTA MOTA	1235 1236	CA	ILE	161	27.442	26.975	53.909	1.00 16.45
MOTA	1237	c	ILE	161	27.200	28.354	53.395	1.00 15.77
MOTA	1238	0	ILE	161	26.118	28.680	52.962	1.00 15.95
MOTA	1239	CB	ILE	161	28.129	26.117	52.864	1.00 19.27
ATOM	1240	CG1		161	27.237 28.351	26.016 24.735	51.619 53.445	1.00 18.83
ATOM	1241	CG2		161 161	28.009	25.614	50.350	1.00 14.44
MOTA MOTA	1242 1243	CD1	LYS	162	28.226	29.169	53.471	1.00 17.56
ATOM	:244	CA	LYS	162	28.187	30.508	52.948	1.00 14.42
ATOM	1245	С	LYS	162	29.216	30.524	51.857	1.00 17.73
ATOM	1246	0	LYS	162	30.249	29.875	51.991	1.00 19.16
ATOM	1247	CB	LYS.	162 162	28.480 27.221	31.540	54.055 54.796	1.00 42.03
MOTA MOTA	1248 1249	CG CD	LYS LYS	162	27.493	32.787	56.039	1.00 70.42
. ATOM	1250	N	VAL	163	28.911	31.176	50.759	1.00 13.74
ATOM	1251	CA	VAL	163	29.798	31.201	49.629	1.00 11.95
ATOM	1252	С	VAL	163	29.928	32.610	49.103	1.00 19.30
MOTA	1253	0	VAL	163	28.944	33.318	48.983	1.00 19.84
MOTA	1254	CB	VAL	163	29.249	30.268 30.277	48.532	1.00 15.89
MOTA	1255	CG		163 163	30.105 29.029	23.852	47.261 49.077	1.00 15.86
MOTA MOTA	1256	CG2 N	VAL ASN	164	31.146	32.999	48.733	1.00 14.03
ATOM	1257 1258	CA	ASN	164	31.382	34.310	48.195	1.00 15.55
MOTA	1259	c	ASN	:64	32.396	34.271	47.050	1.00 20.08
ATOM	:260	0	ASN	154	33.268	23.386	46.988	1.00 23.49
MOTA	1261	CB	ASN	154	31.732	35.325	49.308	1.00 20.52
MOTA	1262	CG	ASN	: <del>6</del> 4	33.196	35.697	49.330	1.00100.00
MOTA	1263		l ASN	154 154	34.020 33.515	34.987 36.831	49.929 48.700	
ATOM	1264 1265	, עוה גי	2 ASN PHE	.54 155	33.515	35.207	46.109	1.00 17.37
nota Nota	1266		SHE	155	33.133	35.301	44.953	
ALVI.	-200					<del>-</del>		

АТОМ	1267	С	PHE	165	32.751	36.445	44.071		15.53
MOTA	1268	0	PHE	165	31.686	37.020	44.251		17.16
MOTA	1269	CB	PHE	165	33.207	33.960	14.187	1.00	12.36 14.35
ATOM	1270	CG	PHE	165	31.862	23.486	43.622 42.293		14.53
ATOM	1271	CD1	PHE	165	31.510	33.749 32.757	44.413		13.56
MOTA	1272	CD2	PHE	165	30.978 30.300	33.297	41.759		22.67
ATOM	1273 1274	CEI	PHE	165 165	29.774	32.282	43.893		15.78
ATOM	1275	CE2	PHE	165	29.426	32.572	42.573	1.00	16.20
ATOM	1276	CZ N	PHE	166	33.641	36.799	43.132	1.00	10.79
MOTA MOTA	1277	CA	LYS	166	33.417	37.864	42.162	1.00	10.74
ATOM	1278	C	LYS	166	33.603	37.344	40.774	1.00	15.95
MOTA	1279	ō	LYS	166	34.602	36.727	40.470	1.00	22.80
ATOM	1280	СВ	LYS	166	34.387	39.055	42.249	1.00	16.61
MOTA	1281	CG	LYS	166	34.573	39.688	43.573	1.00	18.11
MOTA	1282	CD	LYS	166	35.540	40.875	43.454		32.56
MOTA	1283	CE	LYS	166	35.272	41.966	44.476	1.00	48.19
ATOM	1284	NZ	LYS	166	34.823	41.435	45.782		85.21
MOTA	1285	:1	ILE	167	32.703	37.704	39.911	1.00	9.75
ATOM	1286	CA	ILE	167	32.768	37.340	38.558	1.00	9.35
ATOM	1287	С	ILE	167	33.203	38.542	37.823	1.00	14.36
ATOM	1288	0	ILE	167 167	32.811 31.379	39.640 36.929	38.170 38.005	1.00	12.16
MOTA MOTA	1289 1290	CB CG1	ILE	167	30.909	35.624	38.669	1.00	13.02
ATOM	1291	CG2	ILE	167	31.423	36.726	36.472	1.00	7.91
ATOM	1292	CD1	ILE	167	31.773	34.415	39.344		19.57
ATOM	1293	;;	ARG	168	34.005	38.299	36.815		12.19
MOTA	1294	CA	ARG	168	34.500	39.308	35.945	1.00	15.07
MOTA	1295	С	ARG	168	33.948	39.122	34.528	1.00	16.64
MOTA	1296	0	ARG	168	34.278	38.156	33.836		17.70
ATOM	1297	CB	ARG	168	36.024	39.287	35.944		16.54
MOTA	1298	CG	ARG	168	36.580	39.632	37.321		25.54
MOTA	1299	CD	ARG	168	37.894	38.910	37.601		63.52
ATOM	1300	HE	ARG	168	38.380	38.191	36.416	1.00	73.52
ATOM	1301	CZ	ARG	168	38.764 38.795	36.926 36.192	36.416 37.527		67.92 57.44
MOTA MOTA	1302 1303	NH2	ARG	168 158	39.192	36.375	35.271		59.15
ATOM	1304	11	HIS	169	33.090	40.064	34.098	1.00	14.88
ATOM	1305	CA	HIS	169	32.505	40.025	32.758	1.00	13.24
ATOM	1306	C	HIS	169	33.214	41.001	31.839	1.00	12.64
ATOM	1307	0	HIS	169	33.306	42.203	32.121	1.00	14.99
MOTA	1308	CЗ	HIS	169	30.970	40.374	32.760	1.00	10.46
ATOM	1309	CG	HIS	169	30.097	39.474	33.573	1.00	6.54
MOTA	1310	:ID1	HIS	169	29.724	33.246	33.111	1.00	12.63
MOTA	1311	CD2		169	29.474	39.695	34.764	1.00	10.21
MOTA	1312	CEI		169	28.892	37.718	34.031	1.00	10.53
ATOM	1313	NE2		169	28.734	33.566	25.063	1.00	11.84
ATOM	1314 1315	И СА	N2K N2K	170 170	33.691 34.349	40.513	30.737 29.812	1.00	15.87
MOTA MOTA	1316	C	ASN	170	33.356	42.224	29.067		25.06
ATOM	1317	ō	ASN	170	32.386	41.7C1	28.537		16.60
ATOM	1318	СB	ASN	170	35.110	40.550	28.755	1.00	
ATOM	1319	CG	ASN	170	36.245	39.717	29.312		18.70
ATOM	1320		ASN	170	36.702	38.752	28.684		48.29
MOTA	1321		ASN	170	36.695	40.073	30.480	1.00	19.13
MOTA	1322	::	ILE	171	33.662	43.527	28.947		18.75
MOTA	1323	CA	ILE	171	32.848	÷÷.460	28.168		16.74
MOTA	1324	С	ILE	171	33.459	44.638	25.791		19.51
MOTA	1325	0	ILE	171	34.643	44.596	25.642	1.00	21.06
ATOM	1326	C3	:LE	171 171	32.713	45.804	28.842	00	20.46
ATOM	1327	CG1		71	32.089	45.617	30.193		24.79
MOTA	1328	CG2	:::	171	31.852	46.727	27.997	1.00	19.03
ATOM	1329	CD1		171	32.630	46.599	31.229		41.65
ATOM	1330	::	GLU	172 172	32.632	44.818	25.804	00	16.54
NOTA	1331	ÇA	GLU	- /2	33.034	44.933	24.420		17.00
HOTA	1332	3	SLU	172	34.110	45.967	24.147 23.125		25.80 29.20
ATOM	لاتد	٠,	ال ما و	172	34.775	45.898	-3.143	2.00	25.20

MOTA	1334	CB	GLU	172	31.813	45.165	23.509	1.00 22.46
	1335		GLU	172	31.122	46.531	23.786	1.00 58.53
MOTA		CD	GLU	172	29.871	46.783	22.933	1.00100.00
MOTA	1336		-	172	29.415	45.970	22.156	1.00100.00
MOTA	1337		GLU	172	29.370	47.982	23.149	1.00100.00
ATOM	1338	OE2	GLU		34.277	46.934	25.034	1.00 24.41
ATOM	1339	И	ASP	173			24.852	1.00 25.03
MOTA	:340	CA	ASP	173	35.292	47.978	25.455	1.00 33.40
MOTA	1341	С	ASP	173	36.651	47.624	25.518	
MOTA	1342	0	ASP	173	37.561	48.451		1.00 30.42
MOTA	1343	CB	ASP	173	34.822	49.319	25.401	1.00 23.30
ATOM	1344	CG	ASP	173	34.743	49.358	26.912	1.00 32.47
MOTA	1345	OD1	ASP	173	34.406	50.355	27.513	1.00 37.58
ATOM	1346	OD2	ASP	173	34.949	48.196	27.504	1.00 49.22
ATOM	1347	N	GLY	174	36.766	46.410	25.956	1.00 23.87
MOTA	1348	CA	GLY	174	38.019	45.994	26.537	1.00 21.30
MOTA	1349	С	GLY	174	38.012	46.090	28.044	1.00 19.99
ATOM	1350	0	GLY	174	38.927	45.585	28.709	1.00 20.45
ATOM	1351	N	SER	175	36.972	46.767	28.598	1.00 13.88
ATOM	1352	CA	SER	175	36.898	46.931	30.034	1.00 8.70
ATOM	1353	С	SER	175	36.296	45.728	30.765	1.00 17.30
ATOM	1354	0	SER	175	36.136	44.655	30.175	1.00 18.77
MOTA	1355	CB	SER	175	36.288	48.235	30.450	1.00 14.07
ATOM	1356	OG.	SER	. 175	36.360	48.316	31.865	1.00 24.79
ATOM	1357	Ħ	VAL	176	35.963	45.912	32.051	1.00 13.74
ATOM	1358	CA	VAL	176	35.415	44.826	32.864	1.00 16.46
ATOM	1359	C	VAL	176	34.191	45.204	33.703	1.00 22.46
ATOM	1360	Ö	VAL	176	34.159	46.254	34.334	1.00 21.31
ATOM	1361	СB	VAL	176	36.477	44.285	33.818	1.00 24.43
MOTA	1362	CG1		176	35.847	43.344	34.827	1.00 27.45
	1363	CG2		176	37.532	43.536	33.035	1.00 25.65
MOTA MOTA	1364	N	GLN	177	33.234	44.269	33.787	1.00 15.47
	1365	CA	GLN	177	32.048	44.430	34.647	1.00 15.40
ATOM		C	GLN	177	32.102	43.457	35.813	1.00 10.60
ATOM	1366		GLN	177	32.027	42.243	35.634	1.00 13.65
MOTA	1367	0		177	30.709	44.283	33.872	1.00 15.57
MOTA	1368	C3	GLN	177	29.468	44.294	34.828	1.00 19.13
MOTA	1369	CG	GLN	177	29.103	45.673	35.361	1.00 14.91
ATOM	1370	CD	GLN GLN	177	28.759	46.588	34.574	1.00 20.17
MOTA	1371	CEI		177	29.128	45.821	35.690	1.00 17.28
ATOM	1372	NE2		178	32.227	43.993	37.018	1.00 8.17
MOTA	1373	N	LEU	178	32.313	43.180	38.181	1.00 16.66
MOTA	1374	CA	LEU		30.954	42.786	38.712	1.00 20.93
ATOM	1375	C	LEU	178	30.033	43.608	38.753	1.00 14.66
MOTA	1376	0	LEU	178		43.896	39.293	1.00 20.63
MOTA	1377	CB	LEU	178	33.089 34.286	43.110	39.815	1.00 39.28
MOTA	1378	CG	LEU	178			40.852	1.00 45.14
MOTA	1379	CD1		178	33.831	42.087		
MOTA	1380	CD2		178	35.018	42.426	38.648	-
MOTA	1381	51	ALA	179	30.869	41.550	39.171	1.00 16.72
MOTA	1382	CA	ALA	179	29.652	41.033	39.754	
ATOM	:383	C	ALA	179	29.932	40.277	41.040	1.00 15.70
ATOM	1384	0	ALA	179	30.337	39.119	41.028	1.00 15.91
MOTA'	1385	CB	ALA	179	28.853	40.197	38.731	1.00 14.08
MOTA	1386	::	ASP	180	29.694	40.946	42.155	1.00 8.88
ATOM	1387	CA	ASP	180	29.897	40.407	43.480	1.00 7.18
MOTA	1388	С	ASP	180	28.802	39.460	43.891	1.00 17.07
MOTA	1389	0	` ASP	180	27.651	39.844	43.987	1.00 18.22
MOTA	1390	CB	<b>ASP</b>	180	29.934	41.509	44.509	1.00 13.06
ATOM	1391	CG	ASP	180	31.285	41.902	44.935	1.00 46.28
ATOM	1392	001		180	31.981	41.206	45.655	1.00 60.46
MOTA	1393		ASP	190	31.574	43.121	44.560	1.00 46.61
MOTA	:394	::	HIS	181	29.173	38.242	44.197	1.00 14.51
MOTA	1395	CA	HIS	181	28.213	37.223	44.575	1.00 10.49
ATOM	1396	С	HIS	181	28.218	36.897	46.049	1.00 14.28
MOTA	1397	0	HIS	121	29.255	36.530	46.607	1.00 17.40
MOTA	:398	СЗ	HIS	181	28.450	35.915	43.769	1.00 9.89
MOTA	1399	CG	HIS	131	28.077	35.972	42.328	1.00 10.38
ATOM	1400		1 #15	181	28.606	36.926	41.455	1.00 12.24
	•							

MOTA	1401	CD2	HIS	181	27.279	35.146	41.606	1.00 10.42
MOTA	1402		HIS	181	28.093	36.678	40.269	1.00 9.97
MOTA	1403		HIS	181	27.314	35.594	40.316	1.00 9.38
ATOM	1404	N	TYR	182	27.029	36.897	46.668	1.00 10.40
MOTA	1405	CA	TYR	182	26.848	36.518	48.062	1.00 13.86
ATOM	1406	С	TYR	182	25.871	35.393	48.089	1.00 20.61
ATOM	1407		TYR	182	24.819	35.520	47.532	1.00 16.35
MOTA	1408	CB	TYR	182	26.359	37.664	48.934	1.00 21.12
ATOM	1409	CG	TYR	182	27.421	38.693	49.062	1.00 34.16
MOTA	1410	CD1	TYR	182	27.521	39.715	48.120	1.00 46.06
MOTA	1411	CD2	TYR	182	28.389	38.616	50.064	1.00 38.56
MOTA	1412	CE1	TYR	182	28.532	40.674	48.197	1.00 57.53
MOTA	1413	CE2	TYR	182	29.418	39.559	50.147	1.00 40.76
MOTA	1414	CZ	TYR	182	29.480	40.594	49.216	1.00 54.61
MOTA	1415	OH	TYR	182	30.461	41.534	49.308	1.00 61.92
ATOM	1416	N	GLN	183	26.246	34.277	48.686	1.00 17.63
MOTA	1417	CA	GLN	183	25.410	33.104	48.583	1.00 16.37
ATOM	1418	С	GLN	183	25.289	32.311	49.863	1.00 21.39 1.00 19.86
MOTA	1419	0	GLN	183	26.260	32.174	50.623	1.00 13.33
MOTA	1420	CB	GLN	183	25.984	32.219	47.422	1.00 17.38
ATOM	1421	CG	GLN	183	25.651	30.688	47.457 46.389	1.00 17.27
MOTA	1422	CD	GLN	183	26.411	29.884		1.00 13.80
MOTA	1423	OE1		183	26.975	30.454	45.456 46.473	1.00 13.00
MOTA	1424	NE2		183	26.361	28.553	50.055	1.00 19.74
ATOM	1425	N	GLN	184	24.080	31.739	51.168	1.00 16.55
MOTA	1426	CA	GLN	184	23.760	30.829	50.658	1.00 13.60
ATOM	1427	C	GLN	184	23.033	29.582 29.640	49.747	1.00 18.01
MOTA	1428	0_	GLN	184	22.219 22.949	31.444	52.330	1.00 20.11
MOTA	1429	CB	GLN	184	23.364	32.855	52.768	1.00 74.84
ATOM	1430	CG	GLN	184	22.312	33.517	53.657	1.00100.00
ATOM	1431	CD	GLN	184	21.159	33.054	53.752	1.00 97.99
MOTA	1432	OEl		184 184	22.689	34.625	54.286	1.00100.00
MOTA	1433	NE2		185	23.418	23.446	51.207	1.00 14.75
ATOM	1434	N	ASN	185	22.831	27.155	50.887	1.00 13.86
ATOM	1435	CA	asn Asn	185	22.421	26.463	52.166	1.00 16.06
MOTA	1436 1437	С 0	ASN	185	23.176	25.402	53.172	1.00 17.39
MOTA MOTA	1438	CB	ASN	185	23.761	26.212	50.119	1.00 15.20
ATOM	1439	CC	ASN	185	24.110	26.696	48.748	1.00 12.75
ATOM	1440	OD1		185	24.704	27.758	48.592	1.00 22.56
ATOM	1441	ND2		185	23.830	25.868	47.763	1.00 17.70
ATOM	1442	N	THR	186	21.227	25.941	52.139	1.00 18.01
ATOM	1443	CA	THR	186	20.707	25.227	53.288	1.00 17.40
MOTA	1444	С	THR	186	19.976	24.010	52.824	1.00 23.63
ATOM	1445	O	THR	186	19.389	23.991	51.730	1.00 24.57
ATOM	1446	CB	THR	186	19.856	26.100	54.206	1.00 28.82
ATOM	1447	OG 2	THR	186	18.874	26.752	53.446	1.00 35.65 1.00 28.86
MOTA	1448	CG2	THR	186	20.753	27.121	54.903	1.00 28.86
ATOM	1449	Ν.	PRO	187	20.101	22.951	53.620	1.00 22.40
MOTA	1450	CA	PRO	<u>:</u> 87	19.504	21.683	53.269	1.00 20.28
MOTA	1451	С	PRO	187	17.988	21.757	53.288	1.00 22.41
ATOM	1452	0	PRO	187	17.390	22.518	54.071	1.00 19.79
MOTA	1453	СВ	PRO	197	19.977	20.682	54.337	1.00 26.98
MOTA	1454	CG	PRO	187	20.840	21.449	55.338	1.00 22.04
MOTA	1455	CD	PRO	187	20.786	22.918	54.949	1.00 18.77
MOTA	1456	N	, ILE	188	17.382	20.957	52.453	. 00 20 12
HOTA	1457	CA	ILE	188	15.907	20.855	52.407	
MOTA	1458	С	ILE	188	15.470	19.766 19.966	53.389 54.202	1.00 31.58 1.00 38.58
ATOM	1459	0	ILE	188	14.596	20.574	50.991	1.00 21.52
MOTA	1460		ILE	188	15.385 15.555	21.775		1.00 16.10
MOTA	1461		1 ILE	198 188	13.555	20.141	50.981	1.00 28.85
ATOM	1462				15.139	21.471	48.660	
MOTA	1463			188 189	16.142	13.618		
MOTA	1464		GLY	139	15.833	_		
ATOM	1465			139	16.339			
NOTA	1466		GLY	139	17.016			
MOTA	1467	0	GLY	-37	_7.010			

								. 00 40 41
MOTA	1468	N	ASP	190	16.003	16.928	56.617	1.00 49.41
MOTA	1469	CA	ASP	190	16.392	17.047	58.021	1.00 55.01
MOTA	1470		ASP	190	17.556	16.115	58.338	1.00 56.15
ATOM	1471		ASP	190	18.083	16.100	59.463	1.00 58.30
ATOM	1472	CB	ASP	190	15.195	16.734	58.955	1.00 63.89
ATOM	1473	CG	ASP	190	14.592	15.365	58.686	1.00 99.67
ATOM	1474	001	ASP	190	14.599	14.466	59.514	1.00100.00
ATOM	1475	OD2	ASP	190	14.088	15.240	57.470	1.00100.00
ATOM	476	N	GLY	191	17.921	15.312	57.323	1.00 47.20
ATOM	1477	CA	GLY	191	19.015	14.347	57.419	1.00 44.96
	1478	C	GLY	191	20.359	15.044	57.587	1.00 34.43
ATOM	1479	ŏ	GLY	191	20.452	16.266	57.438	1.00 29.96
MOTA	1480	N	PRO	192	21.402	14.264	57.905	1.00 27.26
MOTA	1481	CA	PRO	192	22.737	14.834	58.100	1.00 24.01
MOTA MOTA	1482	C	PRO	192	23.444	15.274	56.787	1.00 20.55
MOTA	1483	0	PRO	192	23.323	14.648	55.740	1.00 23.84
ATOM	1484	CB	PRO	192	23.583	13.764	58.825	1.00 21.00
	1485	CG	PRO	192	22.739	12.501	58.915	1.00 27.49
ATOK	1486	CD	PRO	192	21.330	12.863	58.448	1.00 27.26
MOTA	1487	N	VAL	193	24.193	16.363	56.892	1.00 17.87
ATOM	1488	CA	VAL	193	24.964	16.902	55.792	1.00 19.51
MOTA			VAL	193	26.380	17.108	56.249	1.00 22.37
MOTA	1489	C	VAL	193	26.663	17.189	57.443	1.00 23.84
MOTA	1490	0	VAL	193	24.449	18.245	55.256	1.00 25.24
MOTA	1491	CB		193	23.059	18.118	54.632	1.00 21.90
MOTA	1492	CG1		193	24.497	19.322	56.346	1.00 24.81
MOTA	1493	CG2	LEU	194	27.253	17.241	55.277	1.00 19.04
HOTA	1494	N	LEU	194	28.654	17.438	55.516	1.00 20.29
ATOM	1495	CA	LEU	194	29.006	18.930	55.571	1.00 18.71
ATOM	1496	C	LEU	194	28.907	19.615	54.591	1.00 20.13
MOTA	1497	O CB	LEU	194	29.412	16.806	54.327	1.00 22.92
MOTA	1498	CG	LEU	194	29.994	15.423	54.542	1.00 30.60
MOTA	1499		LEU	194	29.227	14.642	55.595	1.00 35.19
YIOK	1500	CD2		194	30.048	14.672	53.211	1.00 25.61
ATOM	1501		LEU	195	29.453	19.430	56.713	1.00 17.39
MOTA	1502	N	LEU	195	29.881	20.808	56.785	1.00 18.83
MOTA	1503	CA	LEU	195	31.389	20.837	56.579	1.00 28.32
ATOM	1504 1505	C	LEU	195	32.161	20.152	57.281	1.00 21.98
MOTA		0	LEU	195	29.489	21.525	58.072	1.00 22.20
ATOM	1506	CB CG	LEU	195	28.055	21.349	58.444	1.00 26.40
ATOM	1507		LEU	195	27.937	21.508	59.941	1.00 31.99
ATOM ATOM	1508 1509	CD2		195	27.225	22.395	57.726	1.00 26.90
	1510	N N	PRO	196	31.789	21.610	55.597	1.00 21.58
MOTA	1511	CA	PRO	196	33.177	21.666	55.154	1.00 22.17
MOTA	1512	C	PRO	196	34.080	22.623	55.892	1.00 29.56
ATOM	1513	ŏ	PRO	196	33.635	23.588	55.490	1.00 29.04
ATOM ATOM	1514	СВ	PRO	196	33.054	22.265	53.752	1.00 22.77
ATOM	1515	CG	PRO	196	31.761	23.104	53.735	1.00 18.99
ATOM	1516	CD	PRO	196	30.910	22.567	54.861	1.00 16.42
ATOM	1517	N	ASP	197	35.379	22.410	55.716	1.00 22.95
ATOM	1518		ASP	197	36.364	23.370	56.134	1.00 19.71
ATOM	1519		ASP	197	36.556	24.295	54.931	1.00 24.74
ATOM	1520		ASP	197	36.251		53.800	1.00 24.88
MOTA	1521		ASP	197	37.711	22.730	56.446	1.00 22.28
ATOM	1522		ASP	197	37.690		57.687	1.00 43.93
ATOM	1523			197	36.912		58.608	1.00 53.47
ATOM	1524			197	38.634			1.00 31.58
ATOM	1525		ASN	198	37.062		55.168	1.00 19.74
ATOM	1526			198	37.254		54.118	
ATOM	1527		ASN	198	37.974			1.00 19.61
MOTA	1528		ASN	198	38.958			
MOTA	1529			198	39.013			
ATOM	1530			198	37.236	_		1.00 52.21
MOTA	1531			198	36.107			1.00 34.54
MOTA	1532			198	37.354			1.00 55.11
MOTA	1533		HIS	:99	37.462			
ATOM.	1534			:99	33.07			
ATO!	-227							

ATOM	1535	C	HIS	199	37.496	26.357	49.450	1.00 14.85
MOTA	1536	0	HIS	199	36.757	27.295	49.643	1.00 16.45
MOTA	1537	CB	HIS	199	37.988	24.103	50.471	1.00 16.53
MOTA	1538	CG	HIS	199	36.597	23.628	50.218	1.00 16.65
MOTA	1539		HIS	199	35.695	23.491	51.244	1.00 17.85
ATOM	1540	CD2		199	35.987	23.282	49.048	1.00 18.67
ATOM	1541	CEI	HIS	199	34.561	23.052	50.688	1.00 19.45
ATOM	1542	E2	HIS	199	34.716	22.905	49.364	1.00 18.74
	1543	3	TYR	200	37.879	25.998	48.247	1.00 12.56
ATOM		CA	TYR	200	37.334	26.689	47.100	1.00 14.01
MOTA	1544			200	37.207	25.824	45.870	1.00 15.57
MOTA	1545	C	TYR		37.793	24.751	45.768	1.00 20.20
MOTA	1546	0	TYR	200		28.011	46.779	1.00 19.79
MOTA	1547	CB	TYR	200	38.030	27.745	46.779	1.00 22.25
MOTA	1548	CG	TYR	200	39.382			1.00 22.23
MOTA	1549	CD1	TYR	200	39.543	27.526	44.835	1.00 25.73
MOTA	1550	CD2	TYR	200	40.473	27.605	47.057	1.00 25.75
ATOM	1551	CE1	TYR	200	40.800	27.222	44.317	
ATOM	1552	CE2	TYR	200	41.739	27.314	46.559	1.00 29.34
ATOM	1553	CZ	TYR	200	41.896	27.132	45.186	1.00 54.14
ATOM	1554	OH	TYR	200	43.153	26.820	44.703	1.00 62.66
ATOM	1555	7	LEU	201	36.393	26.309	44.946	1.00 15.07
ATOM	1556	CA	LEU	201	36.147	25.680	43.678	1.00 11.01
MOTA	1557	С	LEU	201	36.753	26.532	42.593	1.00 17.30
ATOM	1558	0	LEU	201	36.619	27.753	42.610	1.00 20.19
ATOM	1559	СЗ	LEU	201	34.628	25.518	43.354	1.00 10.69
ATOM	1560	CG	LEU	201	33.749	25.027	44.480	1.00 13.41
ATOM	1561	CD1	LEU	201	32.293	24.938	43.954	1.00 17.11
ATOM	1562	CD2		201	34.196	23.635	44.927	1.00 23.03
ATOM	1563	:1	SER	202	37.407	23.868	41.651	1.00 10.75
ATOM	1564	CA	SER	202	38.047	26.490	40.528	1.00 8.51
ATOM	1565	c	SER	202	37.222	26.189	39.294	1.00 11.56
ATOM	1566	o	SER	202	36.919	25.038	38.996	1.00 14.58
ATOM	1567	СB	SER	202	39.485	25.987	40.442	1.00 15.68
MOTA	1568	OG.	SER	202	40.067	26.353	39.228	1.00 36.44
ATOM	1569	:1	THR	203	36.798	27.241	38.601	1.00 12.36
MOTA	1570	CA	THR	203	35.879	27.067	37.499	1.00 15.60
MOTA	1571	c	THR	203	35.417	27.521	35.195	1.00 20.19
ATOM	1572	ō	THR	203	37.192	28.472	36.114	1.00 18.29
ATOM	1573	CB	THR	203	34.565	27.892	37.757	1.00 20.51
	1574	OG1		203	34.911	29.260	37.780	1.00 20.39
MOTA	1575	CG2		203	33.935	27.557	39.093	1.00 6.80
MOTA		::	GLN	204	35.913	26.883	35.164	1.00 10.30
ATOM	1576		GLN	204	36.173	27.271	33.807	1.00 14.85
ATOM	1577	CA		204	34.956	26.980	32.921	1.00 23.14
MOTA	1578	C	GLN		34.334	25.932	33.056	1.00 21.66
MOTA	1579	0	GLN	204				1.00 20.33
MOTA	1580	CB	GLN	204	37.475	26.696	33.237 32.518	1.00 40.16
ATOM	1581	CG	GLN	204	37.271	25.371 24.722	32.193	1.00 59.76
MOTA	1582	CD	GLN	204	38.588			1.00 41.80
ATOM	1583	OEI		204	39.011	24.716	31.035	
MOTA	1584	HE2		204	39.276	24.241 27.913	33.235	1.00 34.80
ATOM	1585	:1	SER	205	34.619		32.021	
ATOM	1586	CA	SER	205	33.447	27.762	31.172	1.00 14.60
ATOM	1587	C	SER	205	33.654	28.307	29.783	1.00 20.21
ATOM	1588	0	SER	205	34.282	29.337	29.581	1.00 17.82
MOTA	1589	СЗ	SER	205	32.197	23.445	31.758	1.00 11.88
MOTA	1590	QG.	'SER	205	32.121	28.406	33.177	1.00 15.45
ATOM	1591	::	ALA	206	33.065	27.630	28.827	1.00 13.00
MOTA	1592	CA	ALA	206	33.079	23.029	27.426	1.00 9.99
MOTA	1593	C	ALA	206	31.623	23.192	26.924	1.00 21.23
MOTA	1594	2	ALA	206	30.809	27.306	27.139	1.00 14.10
ATOM	1595	CВ	ALA	206	33.751	25.936	26.596	1.00 13.45
ATOM	1596	::	LEU	207	31.335	19.320	26.263	1.00 15.09
ATOH	:597	CA	LEU	207	30.036	23.617	25.706	1.00 12.07
HOTA	:598	\$	LEU	207	30.070	39.445	24.235	1.00 19.76
ATOM	:599	2	LEU	207	31.014	13.840	23.576	1.00 20.82
MOTA	1500	23	LEU	207	29.580	13.840 31.057	26.004	1.00 8.24
ATOM	1601	22	LEU	207	29.744	11.493	27.457	1.00 15.35
							• • • •	

FIG 5-25

MOTA	1602	CD1	LEU	207	28.955	32.790	27.707	1.00 13.73
	1603		LEU	207	29.268	30.406	28.400	1.00 18.79
MOTA						28.863	23.698	1.00 15.33
ATOM	1604	N	SER	208	29.011			
ATOM	1605	CA	SER	208	28.914	28.692	22.270	1.00 13.74
MOTA	1606	С	SER	208	27.449	28.852	21.794	1.00 20.16
					26.548	29.085	22.594	1.00 15.81
ATOM	1607	0	SER	208				1.00 17.82
ATOM	1608	CB	SER	208	29.495	27.367	21.822	
ATOM	1609	OG	SER	208	28.769	26.311	22.431	1.00 31.45
MOTA	1610	N	LYS	209	27.242	28.738	20.485	1.00 16.50
				209	25.907	28.828	19.906	1.00 18.02
MOTA	1611	CA	LYS					1.00 29.99
ATOM	1612	С	LYS	209	25.637	27.610	19.031	
ATOM	1613	0	LYS	209	26.578	27.004	18.502	1.00 32.55
ATOM	1614	CB	LYS	209	25.783	30.100	19.082	1.00 20.96
					24.746	31.055	19.606	1.00 34.50
ATOM	1615	CG	LYS	209				
ATOM	1616	CD	LYS	209	25.262	31.964	20.666	1.00 22.72
MOTA	1617	CE	LYS	209	24.370	33.159	20.896	1.00 18.96
	1618	NZ	LYS	209	23.565	33.067	22.116	1.00 27.39
ATOM						27.241	18.912	1.00 27.01
atom	1619	N	ASP	210	24.347			
MOTA	1620	CA	ASP	210	23.890	26.159	18.038	1.00 24.62
MOTA	1621	С	ASP	210	23.465	26.793	16.705	1.00 26.77
				210	22.468	27.514	16.605	1.00 23.00
ATOM	1622	0	ASP					1.00 24.43
ATOM	1623	CB	ASP	210	22.744	25.361	18.691	
ATOM	1624	CG	ASP	210	22.197	24.249	17.839	1.00 35.55
ATOM	1625		ASP	210	22.333	24.185	15.631	1.00 36.53
					21.499	23.400	18.535	1.00 45.51
MOTA	1626		ASP	210			15.708	1.00 30.25
ATOM	1627	N	PRO	211	24.306	26.618		
ATOM	1628	CA	PRO	211	24.120	27.224	14.397	1.00 30.30
MOTA	1629	C	PRO	211	22.733	26.982	13.770	1.00 39.72
				211	22.253	27.782	12.959	1.00 37.65
ATOM	1630	0_	PRO					1.00 29.99
MOTA	1631	CB	PRO	211	25.197	26.620	13.500	
ATOM	1632	CG	PRO	211	25.782	25.418	14.255	1.00 38.59
ATOM	1633	CD	PRO	211	25.158	25.405	15.647	1.00 35.05
				212	22.102	25.868	14.140	1.00 39.64
ATOH	1634	11	ASN					1.00 39.60
ATOM	1635	CA	ASN	212	20.808	25.515	13.592	
MOTA	1636	С	ASN	212	19.642	25.894	14.497	1.00 41.92
ATOM	1637	ō	ASN	212	18.485	25.518	14.263	1.00 42.30
					20.739	24.028	13.236	1.00 48.64
ATOM	1638	CB	ASN	212			-3.233	. 00 63 61
MOTA	1539	CG	ASN	212	21.883	23.678	12.230	1.00 53.51
ATOM	1640	N	GLU	213	19.947	26.675	15.520	1.00 27.84
ATOM	1641	CA	GLU	213	18.953	27.080	16.478	1.00 20.43
				213	18.485	28.527	16.241	1.00 29.95
MOTA	1642	С	GLU					
MOTA	1643	0	GLU	213	19.247	29.475	16.324	1.00 32.77
MOTA	1644	CB	GLU	213	19.535	26.878	17.894	1.00 16.45
MOTA	1645	CG	GLU	213	18.594	27.326	18.995	1.00 18.29
				213	17.229	26.703	18.853	1.00 38.01
MOTA	1646	CD	GLU					
MOTA	1647	OE 1	GLU	213	16.238	27.334	18.508	1.00 25.07
ATOM	1648	OE2	GLU	213	17.223	25.423	19.122	1.00 19.17
ATOM	1649	N	LYS	214	17.223	28.713	15.963	1.00 22.99
					16.721	30.081	15.726	1.00 22.84
MOTA	1650	CA	LYS	214				1.00 22.54
MOTA	1651	C	LYS	214	16.252	30.778	16.982	1.00 21.50
ATOM	1652	0	LYS	214	16.130	32.016	17.032	1.00 28.15
ATOM	1653	СВ	LYS	214	15.653	30.197	14.606	
				214	16.153	29.816	13.209	1.00 32.71
MOTA	1654	CG	LYS					
MOTA	1655	CD	LYS	214	16.752	30.979	12.431	1.00 55.31
MOTA	1656	N	ARG	215	15.947	30.028	13.014	1.00 14.52
ATOM	1657	CA	'ARG	215	15.518	30.726	13.014 19.209	1.00 15.58
					16.719		19.892	1.00 21.87
ATOM	1658	С	ARG	215				
ATOM	1659	0	ARG	215	17.848		19.572	1.00 26.69
MOTA	1660	C3	ARG	215	14.808	29.804	10.159	1.00 18.82
				115	13.660		19.475	1.00 23.30
ATOM	1661	CG	ARG					1.00 25.50
MOTA	1662	CD	ARG	215	13.220		20.205	1.00 15.45
ATOM	:663	NE	ARG	215	14.107	36.668	19.929	1.00 28.08
ATOM	1664	CZ	ARG	215	14.022	25.473	20.543	1.00 21.38
					13.074		21.455	1.00 23.92
MOTA	1665		1 ARG	215 215				
MOTA	1666	::H:	2 ARG	_15	14.893		10.225 10.330	1.00 20.46
MOTA	1667	::	ASP	116	16.466	32.275	10.330	1.00 16.72
ATOM	1668	CA	ASP	115	17.556	32.895	11.617	1.00 19.06
U T OI I	-000							


MOTA	1669	C	ASP	216	18.047	31.817	22.607	1.00 20.02
MOTA	1670	0	ASP	216	17.261	31.214	23.350	1.00 18.45
MOTA	1671	CB	ASP	216	17.066	34.169	22.383	1.00 21.33
MOTA	1672	CG	ASP	216	18.138	35.140	22.893	1.00 20.97
MOTA	1673	OD1	ASP	216	17.869	36.079	23.620	1.00 28.46
MOTA	1674	OD2	ASP	216	19.342	34.900	22.441	1.00 20.37
MOTA	1675	N	HIS	217	19.332	31.537	22.589 23.433	1.00 13.18
ATOM	1676	CA	HIS	217	19.813	30.482 30.614	23.723	1.00 21.35
ATOM	1677	C	HIS	217	21.313	31.471	23.163	1.00 15.03
MOTA	1678	0	HIS	217 217	22.014 19.587	29.168	22.690	1.00 13.03
ATOM	1679	CB	HIS HIS	217	20.525	29.025	21.542	1.00 15.49
MOTA	1680 1681	ND1		217	20.463	29.871	20.449	1.00 17.88
ATOM ATOM	1682		HIS	217	21.589	28.172	21.361	1.00 17.51
ATOM	1683	CE1		217	21.457	29.524	19.635	1.00 17.94
ATOM	1684	NE2	HIS	217	22.152	28.501	20.151	1.00 17.59
ATOM	1685	N	MSE	218	21.794	29.725	24.576	1.00 11.26
MOTA	1686	CA	MSE	218	23.186	29.642	24.887	1.00 11.49
MOTA	1687	С	MSE	218	23.560	28.198	25.094	1.00 24.15
ATOM	1688	0	MSE	218	22.822	27.446	25.751	1.00 20.70
MOTA	1689	CB	MSE	218	23.539	30.421	26.172	1.00 12.84
HOTA	1690	CG	MSE	218	24.809	30.004	26.907 28.434	1.00 12.59
ATOM	1691	SE	MSE	218	25.267 24.039	31.128 30.502	29.781	1.00 13.54
MOTA	1692	CE	MSE	218 219	24.727	27.824	24.558	1.00 15.62
MOTA	1693	N	VAL VAL	219	25.309	26.518	24.782	1.00 10.58
MOTA	1694 1695	CA C	VAL	219	26.473	26.689	25.753	1.00 16.54
atom Atom	1696	Ö	VAL	219	27.280	27.604	25.585	1.00 15.54
MOTA	1697	СВ	VAL	219	25.774	25.883	23.498	1.00 15.08
ATOM	1698	CG1		219	26.330	24.495	23.824	1.00 14.34
ATOM	1699	CG2	VAL	219	24.599	25.766	22.512	1.00 15.78
ATOM	1700	N	LEU	220	26.523	25.836	26.783	1.00 10.95
ATOH	1701	CA	LEU	220	27.490	25.939	27.850	1.00 11.01
MOTA	1702	C	LEU	220	28.206	24.643	28.184	1.00 21.26
MOTA	1703	0_	LEU	220	27.592	23.577	28.324	1.00 15.94
MOTA	1704	CB	LEU	220	26.807	26.545 26.578	29.100 30.402	1.00 21.10
MOTA	1705	CG	LEU	220 220	27.624 28.433	27.875	30.483	1.00 23.53
MOTA	1706	CD1	LEU	220	26.663	26.556	31.586	1.00 22.04
MOTA MOTA	1707 1708	N N	LEU	221	29.570	24.758	28.273	1.00 19.04
ATOM	1709	CA	LEU	221	30.498	23.666	28.697	1.00 13.22
ATOM	1710	C	LEU	221	31.309	24.178	29.887	1.00 10.73
MOTA	1711	0	LEU	221	31.846	25.267	29.857	1.00 12.98
ATOM	1712	CB	LEU	221	31.382	23.102	27.549	1.00 13.74
MOTA	1713	CG	LEU	221	32.580	22.257	28.045	1.00 18.64
ATOM	1714		LEU	221	32.149	20.868	28.496	1.00 17.38
MOTA	1715	CD		221	33.571	22.109	26.911	1.00 26.97
MOTA	1716	N	GLU	222	31.316	23.446	30.963	1.00 9.31 1.00 9.97
MOTA	1717	CA	GLÜ	222	31.936	23.929	32.144	
ATOM	1718	C	GLU	222 222	32.548 32.072	22.803	32.951 32.966	1.00 12.94
MOTA	1719 1720	O CB	GLU	222	30.836	24.762	32.896	1.00 12.14
MOTA '	1721	CG	GLU	222	31.092	25.119	34.364	1.00 13.88
ATOM	1722	CD	GLU	222	29.895	25.891	34.934	1.00 13.57
ATOM	1723		1 GLU	222	29.128	26.477	34.240	1.00 19.47
ATOM	1724		2 GLU	222	29.752	25.789	36.207	1.00 18.51
ATOM	1725	N	PHE	123	33.687		33.542	1.00 15.86
MOTA	1726	CA	PHE	223	34.476		34.373	1.00 9.34
MOTA	1727	С	PHE	223	34.711		35.722	1.00 11.08
MOTA	1728	0	PHE	223	35.028		35.828	1.00 19.86
MOTA	1729	CB	PHE	223	35.847		33.684	1.00 8.30
ATOM	1730	CG	PHE	223	35.703		32.431	1.00 10.50
ATOM	1731		1 PHE	223	35.570		32.469	1.00 13.56
ATOM	1732			223 223	35.750 35.481		31.184 31.287	1.00 11.32
ATOM	1733		1 PHE 2 PHE	-23	35.481 35.667		29.995	1.00 12.17
ATOM	1734 1735			123 123	35.521			1.00 10.27
MOTA	_/15		FRE	-60	125.		.5.030	1.30 10.07

FIG 5-27

MOTA	1736	N	VAL	224	34.54		22.081	36.765	1.00	9.28
MOTA	1737	CA	VAL	224	34.70		22.587	38.080		11.13
MOTA	1738	С	VAL	224	35.32		21.553	39.010		17.52 13.17
MOTA	1739	0	VAL	224	34.84		20.418	39.137 38.662		16.51
ATOM	1740	CB	VAL	224 224	33.37 33.62		23.078 23.736	40.022		13.90
ATOM	1741	CG1	VAL	224	32.67	14	24.048	37.697		13.85
ATOM	1742 1743	CG2 N	VAL THR	225	36.38		21.965	39.676		11.71
ATOM ATOM	1744	CA	THR	225	37.02		21.099	40.617	1.00	11.51
HOTA	1745	C	THR	225	37.36		21.798	41.927	1.00	14.76
ATOM	1746	ō	THR	225	37.70		23.002	41.962	1.00	16.64
ATOM	1747	CB	THR	225	38.16		20.279	40.014		20.33
ATOM	1748	OG1	THR	225	39.28		20.337	40.822	1.00	30.44
MOTA	1749	CG2	THR	225	38.46		20.722	38.631	1.00	10.89
MOTA	1750	N	ALA	226	37.22		21.065	43.011	1.00	7.89
MOTA	1751	CA	ALA	226	37.47		21.595	44.352	1.00	11.53
ATOM	1752	С	ALA	226	38.90		21.558 20.699	44.677 44.199	1.00	15.60
ATOM	1753	0	ALA	226	39.68 36.69		20.847	45.444	1.00	12.17
ATOM	1754	CB	ALA	226 227	39.39		22.490	45.479		13.95
MOTA	1755 1756	N CA	ALA ALA	227	40.7		22.550	45.871	1.00	19.64
atom Atom	1757	C	λLA	227	40.9		23.299	47.170	1.00	26.33
MOTA	1758	õ	ALA	227	40.0		23.715	47.840	1.00	25.39
ATOM	1759	CB	ALA	227	41.5	57	23.246	44.760	1.00	18.42
ATOM	1760	N	GLY	228	42.2	45	23.476	47.523		23.28
ATOM	1761	CA	GLY	228	42.6		24.292	48.658	1.00	21.61
MOTA	1762	С	GLY	228	42.8		23.562	49.939		32.93 32.53
MOTA	1763	0	GLY	228	42.9		24.201	51.009 49.842	1.00	33.59
MOTA	1764	N	ILE	229 229	42.8 43.0		22.231 21.375	50.998	1.00	31.81
ATOM	1765 1766	CA C	ILE	229	44.0		20.291	50.633		28.78
MOTA MOTA	1767	Ö	ILE	229	45.0		20.176	51.246		96.02
ATOM	1768	СЗ	===	229	41.6		20.772	51.519	1.00	35.70
ATOM	:769	CG1		229	40.8		21.807	52.325	1.00	30.66
ATOM	1770	CG2		229	41.9	90	19.549	52.392	1.00	33.37
ATOM	1771	CD1		229	39.3		21.715	52.092		38.74
ATOM	1772	0	HOH	301	27.5		12.735	38.010	1.00	15.09
MOTA	1773	0	HOH	302	23.9		34.589	37.331	1.00	
ATOM	1774	0	HOH	303 304	27.2 29.9		34.816 18.943	35.487 44.692	1.00	
MOTA	1775 1776	0	нон	305	30.9		21.886	49.900	1.00	
atom Atom	1777	ŏ	нон	306	20.0		31.196	43.592	1.00	
ATOM	1778	ō	нон	307	26.6		48.630	33.797	1.00	
ATOM	1779	ō	нон	308	22.3	129	33.239	41.399	1.00	
ATOM	1780	0	HOH	309	22.4	65	48.025	32.810	1.00	
MOTA	1781	0	нон	310	31.0		39.126	29.118	1.00	
ATOM	1782	0	нон	311	33.0		35.809	33.010	1.00	
ATOM	1783	0	нон	312	31.1		37.076	30.841	1.00	12.58 56.07
ATOM	1784	0	HOH	313 314	40.3 34.1		30.058 26.379	38.616 57.222	1.00	22.58
ATOM	1785 1786	0	нон	315	36.2		35.320	43.598	1.00	22.30
HOTA HOTA	1787	Ö	нон	316	33.8		29.786	34.671	1.00	12.21
ATOM	1865	ŏ	нон	317	42.3		20.166	43.534		25.67
ATOM	1788	o	HOH	318	10.2		28.684	30.403	1.00	43.55
ATOM	1789	0	нон	319	28.4		16.822	30.655	1.00	25.44
ATOM	:790	0	HOH	320	30.6		20.922	37.231		21.57
MOTA	1791	0	HOH	321	11.6		37.421	26.801	1.00	
MOTA	:792	0	нон	322	27.0		37.308	36.869	1.00	
MOTA	1793	0	HOH	323	33.1		14.524	43.070 53.352	1.00	
ATOM	1794 1795	0	HOH	324 325	37.9 32.0		14.036 49.100	37.02B	1.00	
MOTA	1795 1796	0	HOH	326	11.9		12.020	43.429	1.00	29.05
MOTA	1797	0	HOH	327	36.		29.941	31.666	1.00	
MOTA MOTA	1864	o	HOH	328	15.		26.513	15.694	1.00	
NOTA	1798	Ö	HOH	329	33.0		46.924	36.994	1.00	22.07
HOL	1363	ō	HOH		23.3		35.134	22.715	1.00	45.33
ATOM	1799	Ō	нон		33.		31.296	25.261	1.00	

FIG 5-28

		_		222	34 043	24 200	20 522	1 00 30 03
MOTA	1862	0	нон	332	34.942	24.780 12.919	29.532 54.611	1.00 38.93
MOTA	1800	0	нон	333	25.235 38.048	23.467	36.645	1.00 38.20
ATOM	1861	0	нон	3 <b>34</b> 3 <b>35</b>	12.284	43.511	38.338	1.00 37.73
MOTA MOTA	1801 1802	0	нон Нон	336	9.826	47.020	32.568	1.00 46.67
ATOM	1803	Ö	нон	337	7.671	41.532	29.806	1.00 40.88
ATOM	1804	Ö	нон	338	15.430	23.713	26.808	1.00 34.73
ATOM	1805	ŏ	нон	339	24.344	20.385	25.121	1.00 53.42
ATOM	1806	ŏ	нон	340	31.550	10.656	40.819	1.00 47.85
ATOM	1807	Õ	НОН	341	17.569	23.030	25.796	1.00 28.17
ATOM	1808	0	нон	342	19.174	38.552	23.965	1.00 45.54
ATOM	1809	0	HOH .	343	24.268	37.527	25.415	1.00 30.97
ATOM	1810	0	HOH	344	21.266	29.482	41.551	1.00 19.69
MOTA	1811	0	HOH	345	20.668	26.999	41.933	1.00 11.81
ATOM	1812	0	HOH	346	24.780	24.795	43.460	1.00 20.95
MOTA	1813	0	НОН	347	42.962	13.170	46.312	1.00 31.00 1.00 28.20
MOTA	1814	0	HOH	348	32.322	14.088	47.013 49.679	1.00 28.20
ATOM	1815 1816	0	нон Нон	3 <b>49</b> 3 <b>5</b> 0	31.708 22.408	13.186 35.801	50.514	1.00 40.71
MOTA MOTA	1817	0	нон	351	25.366	47.090	42.583	1.00 38.15
ATOM	1818	ŏ	нон	352	27.243	47.647	43.977	1.00 41.55
ATOM	1819	ō	HOH	3 <b>53</b>	29.868	45.076	42.906	1.00 29.32
ATOM	1820	ŏ	НОН	354	14.175	22.269	42.680	1.00 74.11
ATOM	1821	ō	нон	355	13.414	10.739	35.791	1.00 29.92
MOTA	1822	0	HOH	356	20.338	9.974	37.765	1.00 30.46
MOTA	1823	0	HOH	357	23.520	40.420	24.953	1.00 29.75
ATOM	1824	0	HOH	358	25.718	41.692	26.023	1.00 30.43
MOTA	1825	0	HOH	359	26.826	38.466	25.345	1.00 31.72
ATOM	1826	0	нон	360	37.768	42.373	25.123	1.00 41.53
ATOM	1827	0	HOH	361	40.078	42.268	25.852 22.083	1.00 37.12
ATOM ATOM	1828 1829	0	нон Нон	362 363	31.483 33.891	38.677 37.723	30.126	1.00 34.21
MOTA	1860	Ö	HOH	364	39.936	26.543	36.329	1.00 47.93
ATOM	1830	Ö	нон	365	36.631	34.210	41.636	1.00 62.74
ATOM	1831	ŏ	нон	366	37.038	29.783	52.197	1.00 40.07
ATOM	1832	ō	нон	367	37.289	37.407	40.231	1.00 37.59
ATOM	1833	0	нон	368	18.930	17.517	52.472	1.00 35.80
MOTA	1834	0	HOH	369	19.506	18.914	57.913	1.00 45.72
MOTA	1835	0	HOH	370	30.903	25.708	41.139	1.00 21.54
ATOM	1836	0	HOH	371	30.369	25.678	24.583	1.00 22.46
MOTA	1837	0	нон	372	21.000	33.705	20.826	1.00 26.00
ATOM	1838	0	нон	373 374	13.648	32.794 25.683	21.329 38.707	1.00 27.98
MOTA MOTA	1839 1859	0	нон Нон	375	29.735 33.670	24.419	60.503	1.00 50.04
ATOM	1840	ŏ	нон	376	30.034	11.047	37.420	1.00 43.28
ATOM	1841	ō	нон	377	8.662	35.846	35.068	1.00 51.94
ATOM	1842	ō	нон	378	10.847	36.466	39.503	1.00 42.32
ATOM	1843	0	нон	379	14.395	48.943	39.085	1.00 29.72
ATOM	1844	0	нон	380	36.676	11.660	40.172	1.00 39.81
MOTA	1845	0	HOH	381	35.968	7.212	34.763	1.00 58.66
ATOM	1846	0	нон	382	17.426	21.988	21.077	1.00 41.69
ATOM	1847	O	нон	383	29.837	22.623	39.378	1.00 32.82
ATOM	1848	0	HOH	384	23.855	29.386	55.164	1.00 55.00
MOTA MOTA	1849 1850	0	нон нон	385 386	17.408 27.900	35.360 49.720	47.495 42.448	1.00 61.61
ATOM	1851	0	нон	387	13.932	35.230	44.385	1.00 47.70
ATOM	1852	0	нон	388	12.650	23.021	43.288	1.00 49.86
ATOM	1853	ŏ	нон	389	16.974	42.367	43.435	1.00 34.38
ATOM	1854	ō	нон	390	37.335	42.653	28.295	1.00 64.46
MOTA	1855	Ó	нон	391	29.701	<b>49.856</b>	35.323	1.00 62.61
MOTA	1856	0	нон	392	27.267	50.835	33.976	1.00 66.60
MOTA	1857	0	нон	393	19.661	29.131	51.537	1.00 34.01
MOTA	1858	0	нон	394	29.412	17.505	59.089	1.00 51.78
TER								
END								

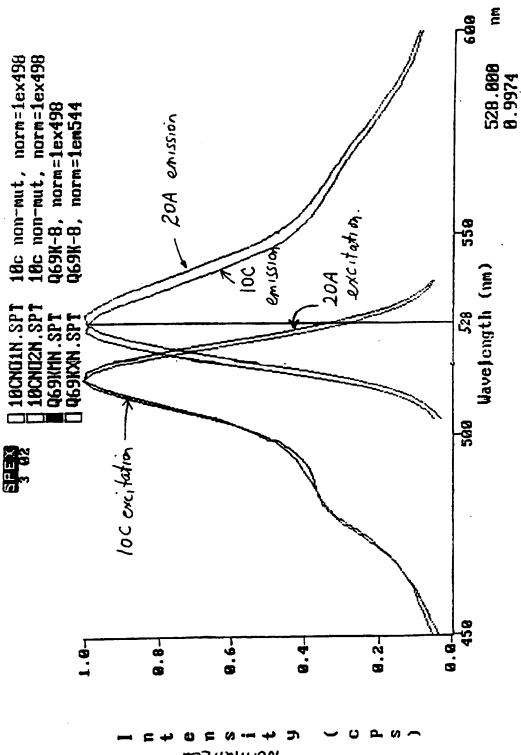

Normalized

International application No.
PCT/US97/14593

A. CLASSIFICATION OF SUBJECT MATTER  IPC(6): C07H 21/04; C07K 14/00, 16/00; C12N 1/20, 15/0  US CL: Please See Extra Sheet.  According to International Patent Classification (IPC) or to both							
B. FIELDS SEARCHED							
Minimum documentation searched (classification system follow	ed by classification symbols)						
U.S. : 435/252.3, 252.33, 325, 410, 320.1; 530/350, 387;							
Documentation searched other than minimum documentation to the	he extent that such documents are included	in the fields searched					
Electronic data base consulted during the international search ( Please See Extra Sheet.	name of data base and, where practicable,	search terms used)					
C. DOCUMENTS CONSIDERED TO BE RELEVANT							
Category* Citation of document, with indication, where	ppropriate, of the relevant passages	Relevant to claim No.					
X WO 96/23810 A1 (THE REGENTS CALIFORNIA) 08 August 1996, abs		46–53, 88, 92, 93, 94					
Y	$\mathbf{Y}$						
	X HEIM et al. Improved green fluorescence. Nature. Vol. 373, 23 88 February 1995, pages 663-664, see Figure on page 664.						
Y	.gato on page oo	2, 3, 10-16, 18- 26, 28-32, 34-37, 39-41					
X Further documents are listed in the continuation of Box	C. Soe patent family annex.						
<ul> <li>Special estagorise of cited documents:</li> <li>"A" document defining the general state of the est which is not considered</li> </ul>	"T" Inter document published after the inte date and not in conflict with the appli the principle or theory underlying the	ication but cited to understand					
to be of particular relevance  "B" cartier document published on or after the international filing date	"X" document of particular relevance; the considered novel or cannot be consider	claimed invention cannot be					
"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other	when the document is taken alone  "Y" document of perticular relevance: the	claimed invention connect has					
*O* document referring to an oral disclosure, use, exhibition or other manus.	special reason (as specified)  2 occurrent or particular reservance; use casimed investion commot be considered to investive stap when the document is combined with one or more other such documents, such combination						
*P* document published prior to the international filing date but later than the priority date claimed	"A" document member of the same putent						
Date of the actual completion of the international search 18 DECEMBER 1997	Date of mailing of the international sea	rch report					
Name and mailing address of the ISA/US Commissioner of Palents and Trademarks Box PCT	Authorized officer TWFO						
Washington, D.C. 20231	Telephona No. (703) 308-0196						

International application No. PCT/US97/14593

C (Continua	ation). DOCUMENTS CONSIDERED TO BE RELEVANT	
Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X  Y	HEIM et al. Wavelength mutations and posttranslational autoxidation of green fluorecent protein, Proc. Natl. Acad. Sci. USA. Vol. 91, December 1994, pages 12501-12504, see abstract	46, 48, 50, 52, 54, 56, 94
<b>Y</b> .	PEROZZO et al. X-ray diffraction and time-resolved fluorescend analysis of Aequorea green fluorescent crystals. Journal of Biological Chemistry. 05 June 1988, Vol. 263, No. 16, pages 7713-7716.	1-41, 46-57, 86- 94
X  Y	DELAGRAVE et al. Red-shifted excitation mutants of the gree fluorescent protein. Bio/Technology. February 1995, Vol. 13, pa 151-153, see Table 1 on page 152.	
Y	EHRIG et al. Green-fluorescent protein mutants with altered fluorescence excitation spectra. FEBS Letters. 1995, Vol. 367, pages 163-166, abstract.	1-41, 46-57, 89, 90
Y	WANG et al. Implication for bcd mRNA localization from spati distribution of exu protein in Drosophila oogenesis. Nature. 02 June 1994, Vol. 369, 400-403, see Figure 1.	ial 32-41, 54-57
P, Y	ORMO et al. Crystal structure of the Aquorea victoria green fluorescent protein. Science. 06 September 1996, Vol. 273, page 1392-1395, abstract.	1-41, 46-57, 86- 94
Р, Ү	YANG et al. The molecular structure of green fluorescent protein Nature Biotechnology. October 1996, Vol. 14, pages 1246-1251, abstract.	
Р, Ү	PALM et al. The structural basis for spectral variations in green fluorescent protein. Nature Struct. Biol. May 1997, Vol. 4, Number 5, pages 361-365.	1-41, 46-57, 86- 94
A	US 5,491,084 A (CHALFIE et al.) 13 February 1996, entire document.	1-41, 46-57, 86- 94
:		

International application No. PCT/US97/14593

Box   Observations where certain claims were found unsearchable (Continuation of Item 1 of first sheet)
This international report has not been established in respect of certain claims under Article 17(2)(a) for the following reasons:
1. Claims Nos.: because they relate to subject matter not required to be searched by this Authority, namely:
2. Claims Nos.:  because they relate to parts of the international application that do not comply with the prescribed requirements to such an extent that no meaningful international search can be carried out, specifically:
3. Claims Nos.:  because they are dependent claims and are not drafted in accordance with the second and third sentences of Rule 6.4(a).
Box II Observations where unity of invention is lacking (Continuation of item 2 of first sheet)
This International Searching Authority found multiple inventions in this international application, as follows:
Please See Extra Sheet.
1. As all required additional search fees were timely paid by the applicant, this international search report covers all searchable claims.
2. As all searchable claims could be searched without effort justifying an additional fee, this Authority did not invite payment of any additional fee.
3. X As only some of the required additional search fees were timely paid by the applicant, this international search report covers only those claims for which fees were paid, specifically claims Nos.:  1-41, 46-57, 86-94
4. No required additional search fees were timely paid by the applicant. Consequently, this international search report is restricted to the invention first mentioned in the claims; it is covered by claims Nos.:
Remark on Protest X The additional search fees were accompanied by the applicant's protest.
No protest accompanied the payment of additional search fees.

International application No. PCT/US97/14593

A. CLASSIFICATION OF SUBJECT MATTER:

US CL :

435/252.3, 252.33, 325, 410, 320.1; 530/350, 387; 536/23.1, 23.4

**B. FIELDS SEARCHED** 

Electronic data bases consulted (Name of data base and where practicable terms used):

APS, STN: Medline, Captus, Scisearch, Lifesci, Biosis, Embase, Wpids, Biotechds.

Search terms: aequorea and green fluorescent, T-203, Thr-203, T203, DNA, cDNA, sequence, s65t, t203h, s65t, t203y, s72a, f64l, s65g, 203y, s72a, s65g, v68l, t203y, t42x, v61x, t62x, v68x, q69x,n121x,y145x, v150x, f165x, i167x, q183x

BOX II. OBSERVATIONS WHERE UNITY OF INVENTION WAS LACKING This ISA found multiple inventions as follows:

This application contains the following inventions or groups of inventions which are not so linked as to form a single inventive concept under PCT Rule 13.1. In order for all inventions to be searched, the appropriate additional search fees must be paid.

Group I, claims 1-41, drawn to DNA coding for mutant fluorescent green protein having mutation at Thr-203, the fluorescent protein, antibody labeled with the fluorescent protein and the DNA coding for a fusion protein consisting of an antibody and the mutant fluorescent protein.

Group II, claims 42-45 and 58-61, drawn to a DNA probe labled with mutants fluorescent.

Group III, claims 46-57 and 86-94, drawn to drawn to DNA coding for mutants fluorescent green protein having mutation at an amino acid residue other than Thr-203, the fluorescent protein, antibody labeled with the fluorescent protein and the DNA coding for a fusion protein consisting of an antibody and the mutant fluorescent protein.

Group IV, claims 62-64 and 68-70, drawn to a method for engineering fluorescent protein.

Group V, claims 65-67, drawn to method of producing fluorescent resonance energy transfer.

Group VI, claims 71-74, drawn to a fluorescent protein crystal having the amino acid sequence SEQ ID NO: 2.

Group VII, claims 75-82, drawn to a computation method for the design of fluorescent protein.

Group VIII, claims 83-85, drawn to a storage device containing the atomic coordinate.

Group LX, claims 95-100, drawn to a method of identifying test chemicals.

The inventions listed as Groups I-IX do not relate to a single inventive concept under PCT Rule 13.1 because, under PCT Rule 13.2, they lack the same or corresponding special technical features for the following reasons:

Each of the above group has a special technical feature defined by the first claim in the Group. The following are the special technical feature for each Groups: (a) Group I is the nucleic acid coding for fluorescent protein having at least the mutation at Thr-203. (b) Group II is a fluorescent DNA probe labeled with a mutant fluorescent protein, (c)Group III is the nucleic acid sequence coding for mutant fluorescent protein having mutation at residues other than T-203. (d) Group IV is a method for the engineering of mutant fluorescent proteins, (e) Group V is a method for producing fluorescent resonance energy transfer, (f) Group VI is the protein crystal of the wild-type protein, (g) Group VII is the computation method to design mutants fluorescent protein with different fluorescent characteristics, (b) Group VIII is a storage device for data, and (i) Group IX is a method of identifying test chemicals.

Group I encompasses the nucleic acid coding for the mutant fluorescent protein, expression vector, recombinant host cell, the mutant proteins and a use for the DNA in making the fusion protein consisting of antibody and the fluorescent protein. Group II represent a second use for the mutant protein of Group I. Also, the special technical feature of Group I is different from that of Group III because the DNA of each Group codes for different sets of mutants that do not share common feature. The special technical feature for this Group I is distinct from those of Groups IV-IX.

The special technical feature of Group II, the fluorescent DNA probe is clearly different from those of Groups III-IX. The DNA probe of Group II represent a second use of the fluorescent protein of Group II.

International application No. PCT/US97/14593

The method of engineering fluorescent protein of Group IV is different from that of producing fluorescent resonance energy transfer of Group V because the resulting fluorescence is different in each case and vary in its characteristics. Similarly, the special technical features of each of Groups IV and V are different from those of the crystal of Group VI, the computation method of Group VIII, and the method of identifying chemicals of Group IX. Finally, the crystal of Group VI, the computation method of Group VIII, the storage device of Group VIII, and the method of identifying chemicals of Group IX are clearly unrelated to each other and there is no special technical feature that connects them together.